


PLEASE READ AND SAVE
THESE INSTRUCTIONS

INSTRUCTION MANUAL

Electric Foot Warmer

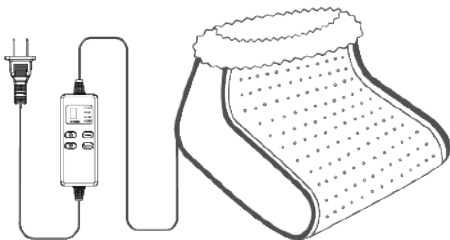

RoHS

MADE IN CHINA

IMPORTANT SAFETY INSTRUCTIONS

DO NOT DESTROY

DANGER! To avoid the risk of burns, electric shock, and fire, this product must be used in accordance with the following instructions.



Do not bleach



Do not iron



Do not fold the pad
Avoid sharp folds that may damage the product



Not to be used children (0 - 6 years)



Do not tumble dry

IMPORTANT SAFETY INSTRUCTIONS

DO NOT DESTROY



Do not dry-clean



Set 90°F (30°C) extra gentle washing program



Please read the instruction carefully before use



Class II electrical appliances



Avoid hot surfaces



PAP20
Dispose of packaging in an environmentally friendly manner

IMPORTANT SAFETY INSTRUCTIONS

DANGER! Non-observance of the following notes may result in personal injury or material damage (electric shock, skin burns, fire).

1. Read all the instructions.
2. Do not use while sleeping.
3. Never leave this appliance unattended, especially if children are around.
4. Do not use a heating pad for an infant.
5. Do not use a heating pad for animals.
6. This pad is not to be used on or by a disabled, sleeping or unconscious person, or a person with poor blood circulation unless carefully attended.
7. This heating pad must not be used by people who are not sensitive to heat, or by vulnerable people who may not react to overheating (e.g. diabetics, persons with skin changes related to illness, or persons with scarred skin in the area where the pad is to be applied, or after taking any pain-relieving medication or alcohol). Do not use on areas of insensitive skin.

IMPORTANT SAFETY INSTRUCTIONS

8. Check on skin under the pad frequently as burns can occur regardless of control setting.
9. Do not use in an oxygen-enriched environment.
10. Do not use pins or other metallic means to fasten this pad in place.
11. Never use the pad with liniments, ointments or salves, especially those containing heat-producing ingredients which could cause skin burns.
12. Never pull this pad with the supply cord.
13. Unplug when not in use.
14. Do not bend or pinch the cord, and do not use this pad if it has been dropped or damaged, or has malfunctions, a damaged power cord or plug. This pad has no user-serviceable parts. If for any reason this pad does not function satisfactorily, check the warranty for service information.
15. Do not tamper with this pad in anyway.
16. Loop the cord loosely when storing; tight wrapping may damage the cord and internal parts.


IMPORTANT SAFETY INSTRUCTIONS

- 17. This product is intended solely for non-commercial, non-industrial, and personal use. Do not use this product outdoors or for any other purpose.
- 18. Supporting controller during use.
- 19. Save these instructions.

OPERATING INSTRUCTIONS


How to use your heating pad?

The Control Connector must be connected to the Pad Connector before plugging the unit into an electrical outlet.

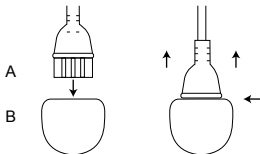

A. Control Connector
B. Pad Connector

NOTE: If the unit does not turn on when you press the power button, make sure the Control Connector button, make sure the Control Connector is securely plugged into the Pad Connector.

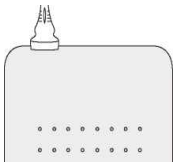


Body Side Fabric



Exterior Side Fabric

OPERATING INSTRUCTIONS

- 1.Plug the unit into an electrical outlet, the green indicator will show on.
- 2.To turn the heating pad on, press the "POWER" button.
- 3.Use the "Heat" button to select the desired heat setting: Heat Level 1-6. The higher heat level you choose, the warmer temperature you will gain.
- 4.Use the "TIMER" button to select the heating time: 1.0H, 2.0H, 4.0H, 6.0H.
- 5.To turn the heating pad off, press the "POWER" button again.

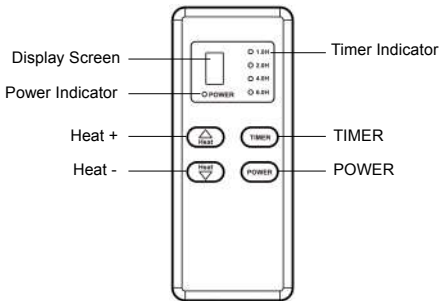

PRODUCT PARAMETERS

Voltage	100-120V(USA) ; 220-240V(UK,EUR)
Frequency	60 Hz(USA) ; 50 Hz(UK,EUR)
Material	Microplush, polyester
Heat Settings	6 Heat Levels
Time	1.0H, 2.0H, 4.0H, 6.0H
Cord	About 9ft

CARE AND MAINTENANCE

CAUTION: DO NOT DRY CLEAN! DO NOT BLEACH! DO NOT USE WRINGER! DO NOT IRON!

- Disconnect the power cord from the electrical outlet. Do not pull by the cord.
- Make sure the Pad Connector is disconnected from the Control Connector by pressing the tabs on the sides of the Control Connector while pulling the connectors apart.
- Machine or hand wash the Pad using a mild detergent, cool water, and a gentle/delicate cycle.
- Line dry, or lay the Pad flat to dry. Do not tumble dry.
- Make sure the product is completely dry before the next use.
- To store the product, loosely roll the Heating Pad and then store it in a cool, dry, dark place after each use. Avoid folding sharply the Heating Pad.

IMPORTANT SAFETY INSTRUCTIONS

DO NOT DESTROY

- **What temperature does this product reach?**

The maximum temperature does not exceed the Underwriters Laboratories temperature limit of about 140°F(60 °C) when set to "Heat Level 6". The temperature varies by heat setting. You may need to use a lower heat setting depending on the specific application and your personal preference.

- **How much does it take to warm up to each temperature selection?**

It takes about one minute to warm up, and the highest "Heat Level 6" takes about two minutes to warm up.

WARRANTY

The manufacturer does not bear any responsibilities when personal injury or property loss is caused due to the following reasons:

- Damage to the product due to improper use or misuse of the product
- Alteration, change, maintenance of the product with other equipment without following the instructions provided by the manufacturer.
- After verification, if the defects are due to improper operation during the transportation of the product.
- After verification, if the problem or dispute is caused by the quality specification or the performance of the parts and components that were produced by other manufacturer.
- In case of operation, repair or maintenance of the unit without following the manuals instructions or related regulations.
- If the damage is caused by natural calamities, improper environmental use or force majeure.

NO.:HP4555B 20220803