

Zeencare

7 Light Facial Wand USER MANUAL



Inclusions

- 1 Facial Wand
- 2 Charging Stand
- 3 Type-C USB Cord

7 LED Lights

Blue 460nm



Red 620nm

GREEN 520nm



Cyan 500nm

Yellow 620nm



Purple 405nm

White 390nm



Warmth Therapy



Microcurrent

Product Specifications:




90° Rotating Head

- 1 Rotating Head
- 2 Charging Indicator
- 3 Power Button
- 4 LED Light Mode
- 5 Charging Stand

Product Dimension:

Head: 5 in/12.7cm

Entire Wand: 17 in/ 43cm

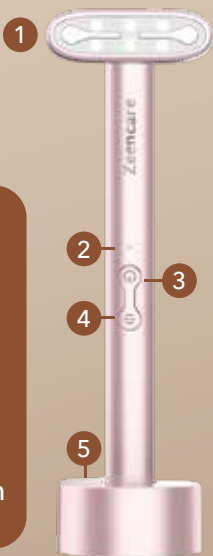

How To Use

1. Cleanse and moisturize skin.
2. Power on and select the mode.
3. Gently massage using the face wand.

Note: Use on hydrated or moisturized skin to avoid discomfort.

Charging Information

1. If the light is blinking, it's charging.
2. A steady light means the battery is full.
3. If the light flashes during use, charge it promptly to avoid shutdown due to low power.

Tips

- Charge the battery fully for the first use.
- Keep your skin moist while using the facial wand.
- Clean equipment with a dry cloth to remove residue and keep it dry.

Warning:

This product is not waterproof.
Don't use product while charging.