
## EMS foot massager instruction manual

1. Remove mat and main unit from plastic sleeve. Snap on the main unit to the mat.
2. Place your feet on the mat before switching on the device. Do not wear socks or pantyhoses.
3. Press the “+” button to turn on the main unit and select the mode by pressing the “M” button.
4. Adjust intensity level by pressing the “+” button. It is recommended that new modes are always started at intensity level 1 and increased as necessary.
5. Once intensity is selected, the device will do a 15 minute cycle and stop. Note that the higher the intensity, the higher the current consumption. The device can be switched off by pressing the “-” button once after the intensity level is set to zero.

The device can be used for approximately 160 minutes after one hour of charging.  
To charge the device, connect the micro USB cable to a suitable USB charger. A red LED will indicate that the device is charging and when the light goes off, it indicates it is fully charged. Do not use the device whilst charging
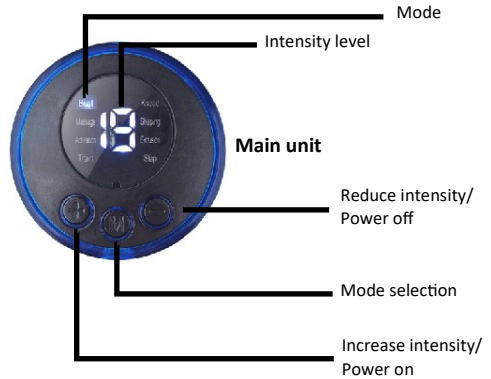
### What is EMS?

Electrical muscle stimulation, also known as neuromuscular electrical stimulation, uses electrical impulse to induce multiple muscle contractions. The central nervous system contracts muscles; EMS has been proved to be an active muscle strengthening tool as well as a circulation booster.

### What is the purpose of EMS?

- Muscle strengthening
- Improve blood circulation
- Preventing muscle atrophy
- Muscle relaxation
- Relieving soreness
- Muscle conditioning

### Warning !

The device should not be used by pregnant women, women who have just given birth or people with medical electronic devices such as pacemakers, etc.

### Troubleshooting

Problem	Possible reason	Solution
Not turning on	Battery has no charge	Charge the main unit
	Device used in low or high temperature environment	Use in ambient temperature between 15° and 35°C
No stimulation	Device switched off or intensity set too low	Press “+” button to turn on device or increase intensity
Stimulation is suddenly interrupted	Battery ran out of charge or the 15 minute cycle has ended.	Charge the main unit or press the “+” button to ensure the battery is charged.

### Technical specifications

Product name	EMS foot massager
Voltage	3.7V
Usage	A-F (8 modes) 19 levels
Frequency	1-100 Hz
Working time	15 minutes
Output	Maximum 35 mA
Battery capacity	150 mAh

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