Thank you for choosing our red light therapy panel. We appreciate your support and hope you'll enjoy this device.

Instructional Manual

Red-light - Therapy for Human & Animals !

Thank you for purchasing our light, we hope you can reap the benefits of red light therapy at home.

Light therapy
to relieve back & waist pain

Pain relief for neck

Relieve the shoulder arm
pain caused by Exercising

Red Light Therapy Pad

English ----- 2 - 12

Deutsch ----- 13 - 24

Français ----- 25 - 36

Italiano ----- 37 - 48

Español ----- 49 - 61

Instructional Manual

Red-light - Therapy for Human & Animals !

Thank you for purchasing our light, we hope you can reap the benefits of red light therapy at home.

Light therapy
to relieve back & waist pain

Pain relief for neck

Relieve the shoulder arm
pain caused by Exercising

Red Light Therapy Pad

Say hello to your red light!

Before starting your first session, please make sure to read this manual as it may explain important operating procedures.

Notes: Each LED inside with three chips at 1 pc 660nm and 2 pcs 850nm. Near Infrared 850nm is invisible to the human naked eyes.

This package includes

- 1* Red light therapy pad
- 1* Mobile pocket
- 1* Adaptor
- 1* Stretchy straps
- 1* Controller
- 1* User Manual

Specifications

Rating power	180W
Actual Power	60W
Product Size	800*300*4.95MM
LED QTY	360*0.5W
Net Weight	1.5KG
Wavelength	1*660nm,2*850nm
EMF	0 UT
Irradiance	0 Inch:>100mW/c m ²
Timer	5/10/15/20/25/30min
Material	Fabric
Lifespan	Over 50000 Hours

The red light device uses a combo of red 660nm and Infrared 850nm LEDs which is chosen to give the ideal wavelengths found in the scientific research of the last three decades.

660nm of red light

This light can be seen by the human eyes and targets skin cells.


850nm of near-infrared light

This light cannot be seen by the human eyes and targets deep tissue.

Red Light


Therapy Pad

360 Red Lights + 720 Infrared Lights , Unique pulse function , the effect is more obvious


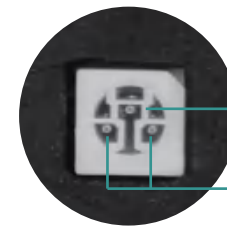
660nm Red Light

850nm Near Infrared Light

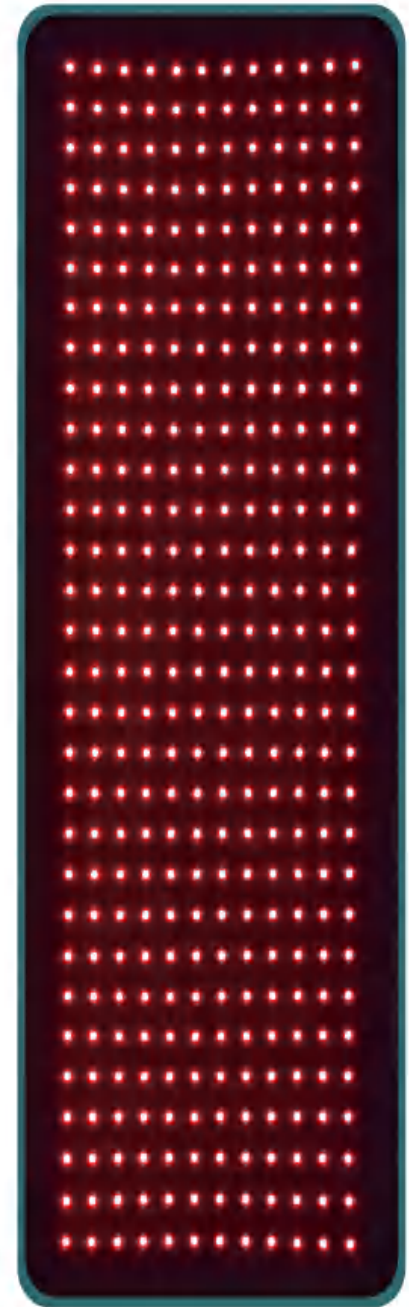CONTROL PANEL

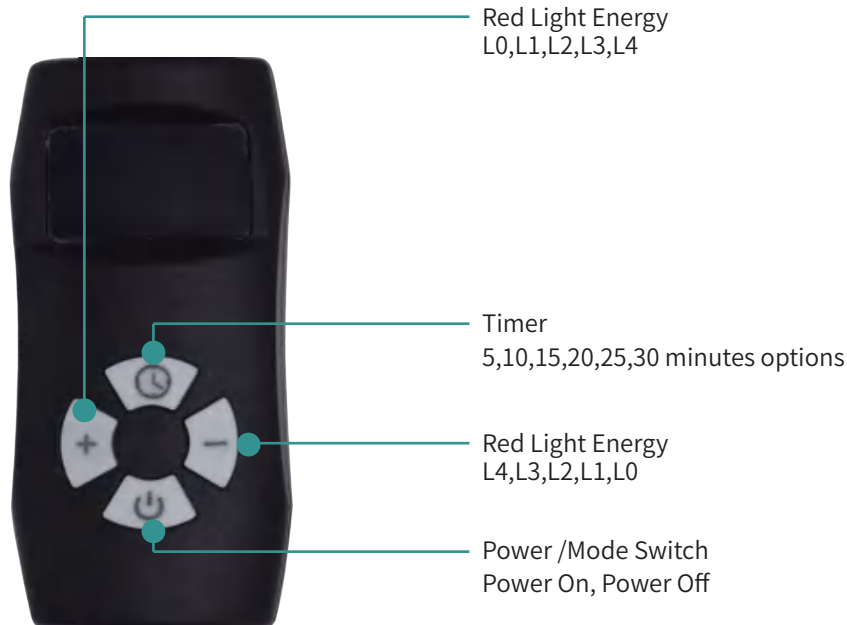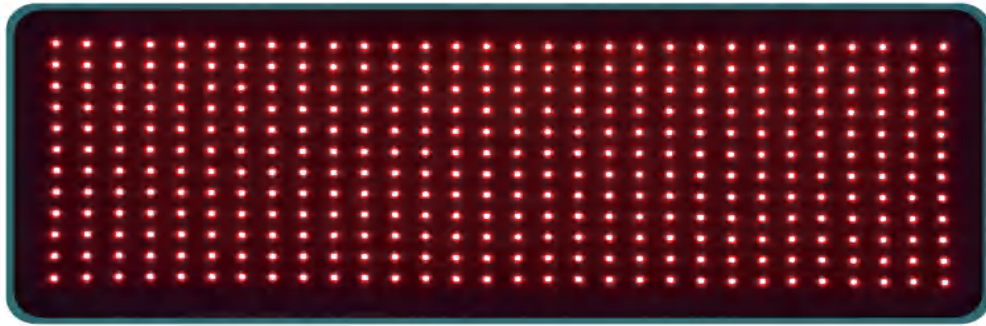Benefits backed by science:

Research shows that the wavelengths in the spectrum of 660-900nm show optimal biotical response.

- Muscle and Joints Pain
- Speed wound healing
- Enhanced blood circulate
- Weight loss
- Anti-inflammatory
- Decreases pain
- Agin sprts and wrinkle removal
- Increase collagen
- Hair Regrowth

Even though this technology doesn't have side effects, we suggest taking the following precautions:

- ▶ This device is not waterproof!
- ▶ Not intended for use by children and pregnant women.
- ▶ Do not exceed 30 Min. on one same treatment.

Safety instructions

Keep this manual in a safe place for future reference.

Before using this product, read all safety warnings and operating instructions. Failure to do so may result in burns, fire, electric shock, serious injury or death, or damage to the product or property.

Contraindications

- Do not use this product over open or bleeding wounds.
- Do not use this product if you are, suspect you may be, or are planning to become pregnant, or if you are, or plan to begin, breastfeeding.
- Do not use this product if you have cancer or suspect you may have cancerous lesions. Do not use over malignant or benign tumors.
- If you or anyone in your household has ever had symptoms related to epilepsy (seizures or loss of consciousness) when exposed to flashing lights, consult your physician before using this product.
- Consult your physician before using this product if you:
 - Are under the care of a physician for a pre-existing health condition.
 - Take steroid medication, use a topical steroid, or receive steroid injections.
 - Have recently had surgery, LASIK or PRK eye surgery, or skin filler injections, if you have breast implants, or if you have a tattoo in the treatment area.
 - Have a condition that may make you sensitive to light or sunlight, such as Lupus, some forms of eczema, light-induced migraines, albinism, and rosacea.
 - Take a medication that may make you sensitive to light or sunlight, including some antibiotics, chemotherapy drugs, cold and allergy medications, pain medications, and diuretics.
 - Use products that may cause photosensitivity, such as citrus essential oils, St. John's Wart, antibacterial soap, and coal tar and petroleum products.
- Are or may be sensitive to red light or near-infrared light.

Warning

To reduce the risk of burns, fire, electric shock, serious injury or death, or damage to the product or property:

- Only use the provided power adapter and USB power cable at the specified voltage to power the device. Using power adapters or cables other than the ones provided may result in fire or explosion.
- Do not use the provided power adapter or USB power cable to power other devices.
- Do not modify, tamper with, or damage the power adapter or USB cable.
- Never unplug the power adapter by pulling on the cord.
- Never use a modified or damaged power adapter or USB cable to power the device.
- Do not modify, tamper with, damage, or attempt to disassemble the device. This device has no user-serviceable parts.
- Always inspect the device, remote control, and power cables for damage, wear and tear, or fraying before each use. Do not use if damaged, malfunctioning, or if the device has been exposed to water.
- 10Hz pulsed-light therapy may affect users who are susceptible to photosensitive epilepsy or other photo-sensitivities. If you or anyone in your household experience symptoms during use, such as dizziness, blurred vision, eye or muscle twitching, loss of consciousness, disorientation, involuntary movements or convulsions, stop using the device immediately and consult your physician.
- Stop using the product immediately if you feel pain or discomfort.
- Keep the device, remote control, and power cables away from flammable or explosive items, chemicals, fire, and heat sources.
- Do not operate in the presence of aerosol (spray) products or corrosive gases, or where concentrated oxygen is being administered.
- Do not use in high temperatures or in direct sunlight. Do not operate under a blanket or pillow.
- Keep away from liquids. Do not use outdoors, in high-moisture environments or in any place where the device may become wet or damp. Do not immerse in water. If the device becomes wet, do not use.

- Do not touch the device, power adapter or USB cable with wet hands.
- Do not use the device over wet skin or clothing.
- Do not use outdoors. Do not use during lightning and thunderstorms.
- Do not use the power adapter or USB cable with an extension cord.
- Never leave the device unattended while running or plugged in. Unplug the device after use.
- Keep sharp objects away from the device, accessories, power cables.
- Do not operate the device while, sleeping, driving, smoking, or under the influence of drugs or alcohol.
- Do not use this product on sleeping or unconscious people or infants.
- For adult-use only. Keep away from children and pets.
- This device is not intended for use by adults with diminished mental, sensory, or physical capacity, unless they are closely supervised by a responsible adult.

Frequently asked questions

Q: WHAT IS RED LIGHT THERAPY?

A: Red light therapy (RLT), also called phototherapy, is a therapy that uses low wavelength red light to treat the at the surface and below the surface of the skin. Red light is absorbed by the body's cells, particularly the mitochondria, which use the red light to make energy, which stimulates cell regeneration and repair.

Red light therapy is a highly-effective, non-invasive, chemical-free therapy that can be used to safely and painlessly heal and rejuvenate many parts of the body.

Q: WHAT IS THE DIFFERENCE BETWEEN RED LIGHT AND NEAR INFRARED LIGHT THERAPY?

A: Both red and near-infrared light are wavelengths on the electromagnetic spectrum. Red light wavelengths, 600–700 nm (nanometers), are absorbed by the skin at a high rate. Near-infrared (NIR) wavelengths, 760–1400 nm (nanometers), are higher on the electromagnetic spectrum and not visible to the human eye. The human body easily absorbs near-infrared waves, up to 1.5 inches into muscle, fat, and joints.

Q: WHAT ARE THE BENEFITS OF RED LIGHT AND NEAR-INFRARED (NIR) LIGHT THERAPY?

A: The benefits of red light therapy and near-infrared (NIR) therapy are very similar, though red light therapy treats cells at the surface of the skin and near-infrared therapy treats cells up to 1.5 inches below the surface of the skin (your muscles, fat, joints, organs, etc.).

The benefits of red light and near-infrared (NIR) therapy include:

- Improve skin appearance, including sun damage, scars (including acne scars), burns, cellulite, and fine lines and wrinkles,
- Stimulate collagen and elastin production and cell regeneration, thus improving blood circulation to the skin and the overall appearance of the skin, including reduced cellulite.
- Energize the lymphatic system, aiding in detoxification.
- Lose inches off the waist, hips, or other treated areas.
- Decrease pain, inflammation, and stiffness in joints, muscles, and connective tissues.
- Reduce symptoms of chronic pain disorders such as tendinitis, arthritis, rheumatoid arthritis, and fibromyalgia.
- Boost immune system.
- Reduce symptoms of eczema and psoriasis, such as redness, itching, flaking, and swelling.
- Improve sleep quality.
- Reduce hormonal imbalances.
- Improve mental clarity and focus.
- Reduce cortisol and increase serotonin, reducing stress and depression.
- Promote faster recovery post-workout.

Q: HOW OFTEN CAN I USE IT?

A: We recommend 20 to 30 minute treatment sessions with your device for optimal results. We recommend using the device 3–5 times per week, depending on your wellness goals. Do not use continuously for more than 30 minutes, and wait 2–3 hours between 30-minute sessions.

Q: SHOULD I WEAR CLOTHING OVER THE TREATMENT AREA?

A: For red light therapy to work most effectively, it is best to wear the device directly against your skin. However, you may use the device over clean, light cotton clothing, if desired.


Q: HOW DOES RED LIGHT THERAPY HELP LOSE FAT AND INCHES?

A: Near-infrared (NIR) light can reach fat cells up to 1.5 inches below the surface of the skin. NIR wavelengths help breakdown stubborn fat cells, allowing them to enter the bloodstream and from there be eliminated via the body's normal waste-removal processes.

Application

Here show you some examples of how to use the red light device. You and your pets can enjoy the red light in different ways like shoulder, back, feet knee, etc.

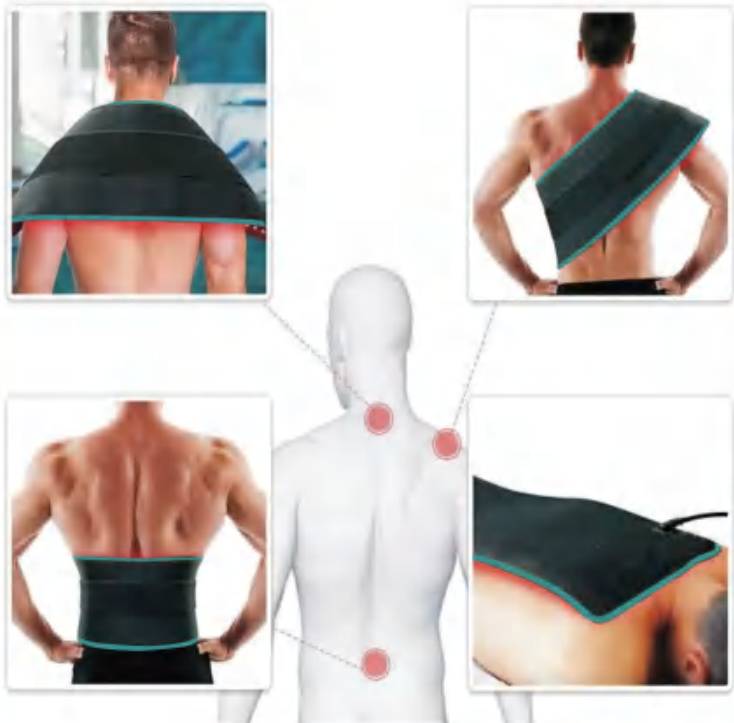LED THERAPY LIGHT PAD FOR JOINTS PAIN RELIEF

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