
Withings Thermo

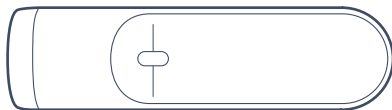QUICK INSTALLATION GUIDE

GUIDE D'INSTALLATION | INSTALLATIONSANLEITUNG | GUÍA DE INSTALACIÓN | GUIDA ALL'INSTALLAZIONE | GUIA DE INSTALAÇÃO

Read this manual before use

EN | FR | DE | ES | IT | PT

THANK YOU FOR CHOOSING WITHINGS THERMO

MERCI D'AVOIR CHOISI WITHINGS THERMO | DANKE, DASS SIE WITHINGS THERMO GEWÄHLT HABEN | GRACIAS POR ELIGIR WITHINGS THERMO
GRAZIE PER AVER SCELTO WITHINGS THERMO | OBRIGADO POR ESCOLHER O WITHINGS THERMO

Need help?



Besoin d'aide? | Brauchen Sie Hilfe? | ¿Necesitas ayuda? | Ti serve aiuto? | Precisa de ajuda?

➔ withings.com/support/thermo

This guide is available in additional languages at: www.withings.com/guides

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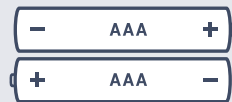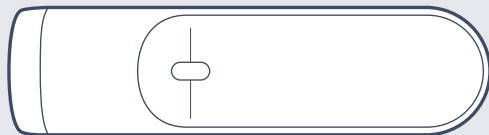How to install your Withings Thermo

Comment installer votre Withings Thermo | So installieren Sie Ihr Withings Thermo | Cómo instalar Withings Thermo |

1

Visit thermo.withings.com from your smartphone or tablet, and follow the instructions to install the Withings Thermo app.

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2

In the Thermo app, follow the instructions to install Thermo.

- FR** Dans l'app, suivez les instructions pour installer le Withings Thermo.
- DE** Folgen Sie der Anleitung in der App, um Ihr Withings Thermo Schläfenthermometer zu installieren.
- ES** Desde la aplicación, siga las instrucciones de instalación de Withings Thermo.
- IT** Nell'applicazione, seguire le istruzioni per installare Withings Thermo.
- PT** Siga as instruções no aplicativo para instalar o Withings Thermo.

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1

Press the button to turn the thermometer on.

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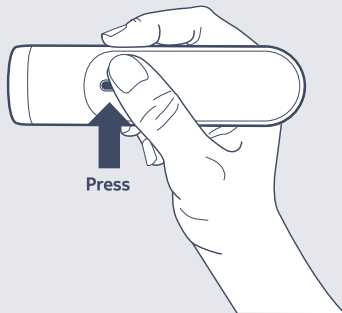
Clear any sweat or hair from the forehead before taking a measurement.

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2

Starting from the middle of the forehead, slowly scan straight across to the hairline, up to 1 cm above the skin.

- FR** En partant du milieu du front, balayez en ligne droite jusqu'aux cheveux.
- DE** Ausgehend von der Mitte der Stirn, parallel zum Haaransatz messen.
- ES** Explora a través de la frente, directamente desde el centro de la frente hasta el cabello.
- IT** Passare sulla fronte, a partire dal centro della fronte fino all'attaccatura dei capelli.
- PT** A partir do meio da testa, faça a leitura rente ao cabelo.

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3

Thermo will vibrate twice when the measurement is complete.

- FR Le thermomètre vibre pour indiquer la fin de la mesure.
- DE Das Thermometer vibriert, wenn die Messung abgeschlossen ist.
- ES El termómetro vibrará cuando la medición se haya completado.
- IT Il termometro vibrerà una volta che la misurazione sarà stata ultimata.
- PT O termômetro vibrará quando a medição estiver concluída.

4

To assign the temperature to a user, scroll through names on the touch sensitive zone (1). Press the button to assign the reading (2).

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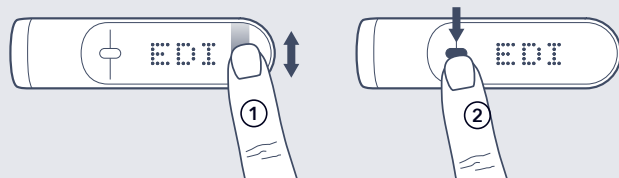
5

A colored LED indicates the fever level according to the user's age.

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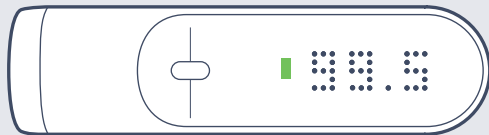
Person being measured must be indoors for 15 minutes prior to taking a reading. | Person being measured must have dry skin. If needed, dry with a towel prior to taking a reading. | Person being measured and the thermometer should be in the same ambient temperature for at least 15 minutes prior to taking a reading. | Person being measured should not have exercised within 15 minutes prior to taking a reading. | Don't take a measurement while or immediately after nursing a baby.

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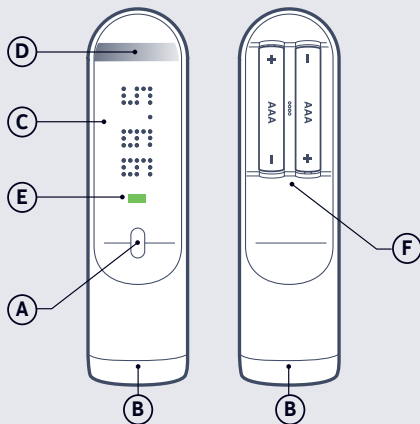
US Withings Two (2) Years Limited Warranty - Withings Thermo

Withings SA, 2 rue Maurice Hartmann, 92130 Issy-les-Moulineaux, ("Withings") warrants the Withings-branded hardware product ("Withings Product") against defects in materials and workmanship when used normally in accordance with Withings' published guidelines for a period of TWO (2) YEARS from the date of original retail purchase by the end-user purchaser («Warranty Period»). Withings' published guidelines include but are not limited to information contained in technical specifications, safety instructions or quick start guide. Withings does not warrant that the operation of the Withings Product will be uninterrupted or error-free. Withings is not responsible for damage arising from failure to follow instructions relating to the Withings Product's use.

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UK Withings Two (2) Years commercial guarantee - Withings Thermo

Withings SA, 2 rue Maurice Hartmann, 92130 Issy-les-Moulineaux, France («Withings») warrants the Withings-branded hardware product («Withings») against defects in materials and workmanship when used normally in accordance with Withings' published guidelines for a period of TWO (2) YEARS from the date of original retail purchase by the end-user purchaser ("Warranty Period"). Withings' published guidelines include but are not limited to information contained in technical specifications, safety instructions or quick start guide. Withings does not warrant that the operation of the Withings Product will be uninterrupted or error-free. Withings is not responsible for damage arising from failure to follow instructions relating to the Withings Product's use.

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EN Withings SA hereby declares that the device Withings Thermo is in conformity with the essential requirements and other relevant requirements of Directive 2014/53/EU. The declaration of conformity can be found at: withings.com/compliance

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The crossed out wheeled bin symbol on your product, battery and/or accessories means that the item is classified as electrical and electronic equipment

QIG_SCT01_42

withings

Withings Thermo | Smart temporal thermometer - Model: SCT01 - V1.0

© 2018 Withings SA. All rights reserved.