

The COSLUS logo is displayed in white, uppercase letters within a teal-colored circular shape that has a wavy, organic edge on its right side.

COSLUS

ACTIVATE YOUR

24-MONTH

WARRANTY FOR FREE



* Please contact customer service if you're unable to provide the order number

www.coslus.com/warranty



Win Free Accessories



Covered by 24-Month Warranty



Friendly Responsive Customer Service

Expire time: 7 days after the delivery date


 +1 855 594 3950 (9am-5pm Mon. - Fri. (EST))

 support@coslus.com


FLIP →

QUICK GUIDE

Scan QR code to watch the operation video. >>

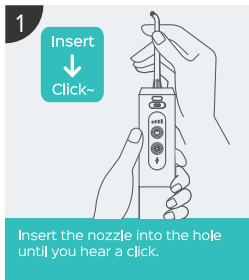


Insert the nozzle into the hole until you hear a click.

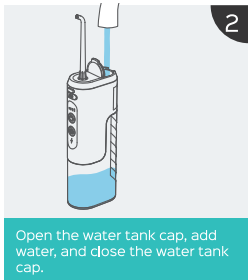

Open the water tank cap, add water, and close the water tank cap.




Adjust to the desired gear before use; short press the mode button to switch modes.

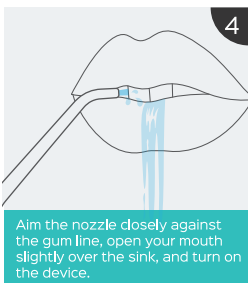

Aim the nozzle closely against the gum line, open your mouth slightly over the sink, and turn on the device.



After flossing for 1-2 min, remove the nozzle and store it by opening the fill cap.

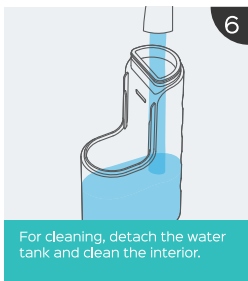

For cleaning, detach the water tank and clean the interior.