


User Guide of Portable Oral Irrigator

PACKAGE CONTENTS

- 1X Oral Irrigator (Model: HE-PCA028)
- 4X Nozzles
- 1X Charging Cable
- 1X Travel Pouch
- 1X User Guide

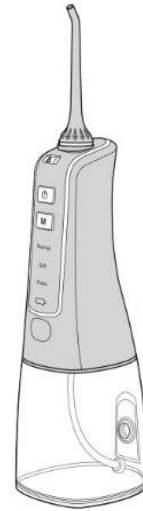

SPECIFICATIONS

Model	HE-PCA028
Input	DC 5V1A
Power	5W
Charging Time	Approx. 4 hours
Battery Life	2 weeks of lifetime between charges
Reservoir Capacity	300 ml/10.0 fl. oz
Battery Capacity	1400mAh
Working Cycle	2 minutes
Dimensions	8.5x7.2x22.6x24.9cm / 3.3x2.8x9.8in
Weight	249g / 8.8oz

Thank you for purchasing our product that has been tested to meet our stringent quality standards. This product is covered with our product and labor warranty for 12 months from the date of its original purchase. Please email us with your order ID and problem or issue, we will respond within 24 hours.

IMPORTANT

Safeguards

All information in this section is provided to instruct you on the correct and safe operation of this device so that you can avoid injuries to yourself and others as well as damage to property.

Dangers

- Toddlers, infants, people who are unable to operate should not be allowed to use this appliance. Please keep it out of reach of children.
- Never alter the appliance. Also, do not try to disassemble or repair it.
- Do not place or store the product where it can fall or be pulled into tub or sink.

Warnings

- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- The appliance is only to be used with the power supply unit provided with the appliance.
- The battery must be removed from the appliance before it is scrapped.
- The appliance must be disconnected from the supply mains when removing the battery.
- The battery is to be disposed of safely.
- Regarding how to remove the battery, please see the "Battery Maintenance and Disposal" section.
- Do not fill the tank with water warmer than 40 °C / 104 °F.
- If the cable is damaged, it shall be replaced by the manufacturer or its service agent or a similarly qualified person in order to avoid a hazard.
- Keep the cord away from heated surfaces.
- Please do not use it if the nozzle is lost.
- Do not use for any purpose other than oral cleaning.
- Keep nozzle away from children. The nozzle is small enough to be swallowed and can be lodged in the throat.
- People who experience tooth or gum pain may not be able to use this device and should consult their dentist.
- Do not press the nozzle too firmly against your teeth or gums. Doing so may injure your gum.

TROUBLESHOOTING

Problem	Cause	Solution
Water pressure too strong	Not used to the water pressure for the initial use	Select the Soft mode for the initial use and adjust to other modes after 2-3 uses
Water pressure too weak	Low voltage	Please fully charge the unit before use
Gum bleeding	1. Effective cleaning of the gingival sulcus during the initial use 2. Improper cleaning methods	Clean your teeth as instructed in our manual
No water or water leakage from the nozzle	Did not install the nozzle properly	Please reinstall the nozzle and make sure it is fit tightly into the body
Reservoir lid falls off	The reservoir lid is detachable	Install it as instructed in the "Adding Water" Section
Fail to charge the unit	1. Debris on the charging port 2. Battery malfunction	1. Clean the charging port regularly with Q-tip 2. If charging still doesn't start, consult our after-sales support team

CLEANING

1. Remove the nozzle and rinse it under tap water to clear all the debris.
2. Remove the reservoir and rinse it under tap water.
3. Brush the filter gently with soft brush.
4. Wipe the body with wet cloth.

Note:

- Make sure to unplug the unit before cleaning.
- Do not wash the unit in dishwasher.

