

FAQ Sheet for Tria Hair Removal Laser 4X

Powered by the same in-office Diode Laser Technology used by dermatologists, the Tria Hair Removal Laser 4X is the first and only FDA-cleared at-home laser to remove unwanted hair at a fraction of the in-office treatment.

Our advanced hair removal laser technology safely delivers over 3 times more hair-eliminating energy than any other at-home hair removal devices, making it the most powerful at-home hair removal technology. It's so safe and effective, you'll be amazed at what you can do all by yourself.

The Tria Hair Removal Laser 4X targets and permanently disables the hair follicle to stop hair growth on the body for permanent results, and even gentle enough to remove unwanted hair on the face for women.

SKINTONE INFORMATION

Recommended skin tone:

The Hair Removal Laser 4X is designed to work on naturally light brown to black body hair, on very pale white skintones to light brown skintones. Because the Hair Removal Laser 4X works with the dark pigment in the hair, lighter hair colors such as white, blond, red, or gray will not absorb enough of the laser's energy to disable the hair follicle.

Note: Hair color in one area can be darker or lighter than hair color in another area, so check the hair in the area you wish to treat.

The Hair Removal Laser 4X is safe and effective on most hair and skin tones.
Find out if Tria is right for you by selecting the skin tone and hair color in the area you wish to treat.

1 SKIN TONE: Select the skin tone of the face and/or body location you wish to treat.



2 HAIR COLOR: Select the natural hair color of the face and/or body location you wish to treat.



HOW TRIA 4X WORKS

Tria gets to the root of hair removal


Tria Hair Removal Lasers are the only at-home hair removal devices featuring diode lasers, which disable the hair follicle by targeting the melanin at the root without damaging the surrounding skin.



DEVICE COMPARISON

Tria is the most powerful at-home hair removal solution

Tria's advanced hair removal laser technology safely delivers over 3 times more hair-eliminating energy density (J/cm^2) than any other at-home hair removal devices. Our ultra high energy density, coupled with its precise laser wavelength selection, effectively targets the pigment in the hair follicle to permanently disable its ability to regrow hair. And, it's so safe and effective, you'll be amazed at what you can do all by yourself.

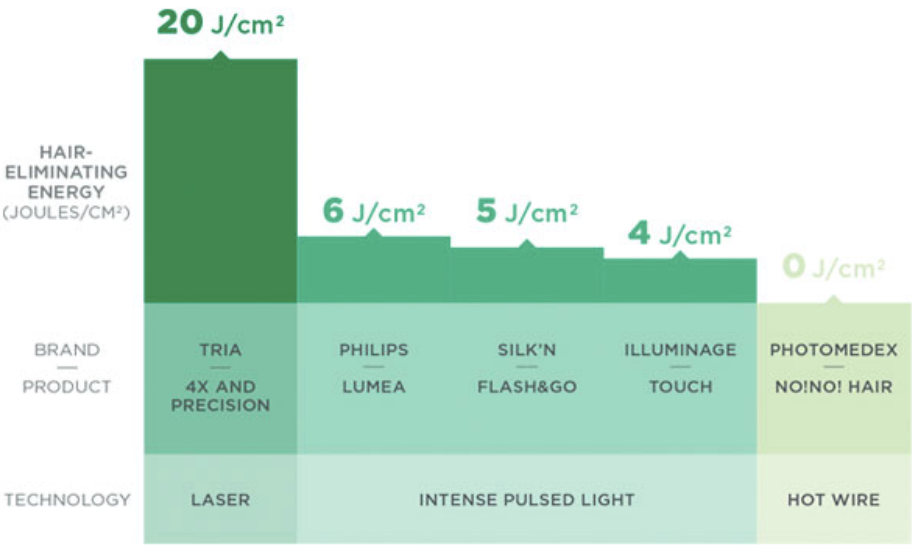

LASER TYPE: Diode laser

LASER CLASS: 1M

NOTE: This Class 1M laser device emits invisible laser radiation. Do not treat the eye, eyelids, or point directly at the eye. Do not view directly with optical instruments, such as microscopes and telescopes. This device is safe to eyes and skin if used in accordance with these Instructions for Use.

OUTPUT WAVELENGTH: 1440 nm

OUTPUT LASER ENERGY: 5-12 mJ/pulse

Tria Lasers vs. IPL Devices



More Energy

Tria's laser devices safely deliver over 3 times more hair-eliminating energy than any IPL or other at-home hair removal devices. This means each follicle that is within the treatment area will be exposed to more energy than within an IPL device.




Easier to Use

All Tria devices are powered by an internal battery while most IPL devices require a cord during use. Tria's cordless devices help to increase the range of places you can use them and creates a more comfortable experience.



Easy on the Eyes

Not only are Tria's lasers safe and effective for use on the face and body, but the light emission is not as intense compared to the bright light bursts from an IPL device.

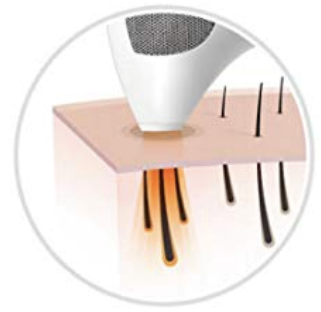

More Control and Precision

Tria's laser devices provide controlled precision on each treatment area, targeting each follicle at a much deeper and concentrated level than IPL for permanent results.

Affordable laser luxury...

Why pay \$1,000's for a professional treatment?

Professional laser hair removal is the "gold standard," but it's also expensive and can be a hassle. Tria Lasers feature the same core diode laser technology used in professional offices. While lasers used by professionals may offer higher peak power and/or greater speed of laser pulse, the Tria Laser 4X delivers excellent, permanent results – all on your schedule, in the comfort and the privacy of your home, at a fraction of the cost.



LASER 4X

[SHOP NOW](#)



PRECISION

[SHOP NOW](#)



IN-OFFICE VISITS

UPPER LIP
\$500

UNDERARMS
\$1,750

FOREARMS
\$1,750

BIKINI
\$2,500

UPPER LEGS
\$2,500

ADVANCED FEATURES

Tria Laser 4X features the same core diode laser technology preferred by dermatologists. While professional devices may offer higher peak power and/or greater speed of laser pulse, the Tria Laser 4X delivers excellent, permanent results in the comfort and privacy of your own home.



LED display

The 4X features a unique interactive LED display that guides you through each treatment, safely and effectively.



Powerful and effective

Specially designed so you can be well on your way to smoother more touchable skin.



Longer battery life

With a 30-minute battery life, Tria Laser 4X is ideal for covering multiple face and body parts during each treatment.

DIFFERENCES BETWEEN THE 4X AND PRECISION LASERS



Hair Removal Laser Precision is designed for small, sensitive areas. Small, slender and light, Tria Laser Precision is great for all your curved areas. With enough battery charge (about 15 minutes) to treat an area like both underarms, it's ideal to take along on a vacation.

- The treatment window is exactly the same size as the 4X
- There are three strength levels on the Precision as compared to the 4X's five levels, but the level three strength is exactly the same as 4X's level five.

QUICKSTART GUIDE ON HOW TO USE

Step 1: Prep your skin

Cleanse, shave and towel-dry for the most effective treatment.

Step 2: Unlock your device

Here's how to unlock the Tria Hair Removal Laser 4X:

1. Unplug the Tria Hair Removal Laser 4X from the battery charger.
2. Hold the device and press the power button to turn it on. The display window will show that it is locked.
3. The red light on the bottom of the device is the skin sensor. Hold it on the area you wish to treat to check

your skin tone. The skin sensor will automatically sense and check your skin.

If it's suited for your skin tone in the area you want to treat, the Tria Hair Removal Laser 4X will emit a ready tone and the lock indicator in the display window will disappear.

Step 3: Select a treatment level

Customize your comfort from 5 energy levels. Start on low, and work your way up to a higher setting.

Step 4: Choose an area and begin

Treat the desired body part every two weeks. Continue treatments as necessary for desired results.

Using the Tria Hair Removal Laser as directed for at least 8 treatments, most users will achieve a substantial visible reduction in hair regrowth, and additional treatments may result in additional reductions. Depending on skin tone, hair color, body part, treatment level, and other factors, individual results from laser hair removal may range from complete reduction to little or no reduction.

MOST FREQUENTLY ASKED QUESTIONS

Can I use the Tria Hair Removal Laser on my face?

Yes, the Tria Hair Removal Laser can be used on women's facial hair (cheek line downward). Darker, coarser facial hair is likely to show more pronounced reduction than lighter and finer hair known as “peach fuzz” on the upper lip and side of the face. Men should not use it on their face due to the density of facial hair.

How long does a battery charge last on the 4X?

Each battery charge on the Tria Laser 4X lasts about 30 minutes. This battery life allows users to stay focused and conscious of overlapping during treatment. A rule of thumb is that a single battery charge is enough to treat a larger area, such as the lower leg, or multiple small areas, such as the upper lip and underarms.

Why do I have to shave? How do I prep the skin for treatment?

We recommend that you prepare your skin before each Tria Hair Removal Laser 4X treatment for best results.

- Cleanse the treatment area thoroughly, removing all cosmetics, lotions, and creams from the area you wish to treat. (This will allow the laser treatment to easily penetrate the skin.)
- Shave the treatment area thoroughly using a new razor. This allows the laser treatment to be most effective because the laser works on the hair inside the follicle, under the surface of the skin.
- Wipe your skin with a cool, damp cloth to remove any shaving cream or soap, and then dry the area well.

Please note: Waxing and plucking are not advisable any time between treatments because both remove hair from the follicle. The hair inside the follicle needs to be present for the Hair Removal Laser Precision treatments to be effective.

Can I use the laser more often than every 2 weeks?

We recommend treating every 2 weeks for best results. This treatment protocol is optimal to target the hair follicle at the proper stage of its growth cycle. For permanent results, treat the area once every 2 weeks until desired result has been reached.

How long does it take before I will see any results?

After each treatment:

- The skin may experience mild redness that typically dissipates within 24 hours.
- Hair that was effectively treated begins to fall out. Hair that regrows may be finer and lighter. Continue your regular treatments.

After 3 months of treatment:

- With each treatment, additional follicles will be treated. Hair follicles that have been deactivated the Tria Hair Removal Laser 4X will no longer grow hair. Continue to use the Tria Laser until you achieve your desired results.

Can I use the Tria laser on my genital area?

No. Do not use the Hair Removal Laser 4X on genitals. These areas may have more sensitive skin, a darker skin tone and/or a greater hair density. Using the Hair Removal Laser 4X in these areas may cause serious skin injury. You may use the Hair Removal Laser 4X on your bikini line, which is the pubic hair visible beyond the boundaries of a swimsuit.

How old do you have to be to use the Tria laser?

The Tria Laser 4X should not be used by or on anyone under the age of 18. The use of the Hair Removal Laser 4X on or by children has not been studied and could result in serious injury.

Can I use the Tria laser on tattoos?

No. Do not use Tria Laser 4X on parts of your skin that are darker than your usual skin tone: tattoos, dark brown or black spots (such as large freckles, birthmarks or moles). Dark skin may absorb too much laser light, which can injure the skin.

What does that error beep noise mean while I'm treating?

Hearing a beep and buzz sound indicates that the laser window was lifted before the pulse was completed.

How can I prevent cracks on the housing of the device, and can I still use my device if it is cracked?

Cracks can happen over time with rapid heating and cooling of the device. There is no known way to prevent this at this time. If your device is cracked, it is safe to use as long as the cracks are not present near or on the treatment window.

Do you have to have the power level at 5 to get any results?

Treating at the highest level will deliver better results faster. You will may see results at the lower levels, but it will take much longer.

Will the hair removal laser work on women with PCOS or hormonal hair growth?

The HRL is not effective on women with PCOS since it can cause the follicle to regenerate. Instead, the laser will act as a means to manage the hair growth but will not stop it. Also, medications and birth control will also impact how effective the device is. Anything that pertains to hormones can affect results.

Do you get 90,000 pulses for the Precision and the 4X?

Both the Precision and 4X have 90,000 pulses and will hold 300 full charges.