Function Introduction

7 Color Beauty Device for Face & Neck is an over-the-counter hand-held device for the treatment of full-face wrinkles. It combines 7 colors LED light, vibration and heating to better solve wrinkles and other skin problems.

(A) Vibration

It contains 3 massage modes with different energy levels. The CLEAN Mode can clear the pores and keep your skin clean and glow; The EMS Mode can relax your skin, help your skin absorb skin care products and reduce fine lines and wrinkles; The HOT Mode can tighten the sagging skin and promote collagen production, thus rejuvenating your skin.

(B) Heating

45±5 °C constant temperature during the treatment promotes absorption of skincare products. The appropriate temperature also makes skin care more comfortable.

(C) 7 Color LED Light

RED LIGHT (630nm) Anti-aging & promote collagen production
BLUE LIGHT (415nm) Sooth inflammation & calm sensitive skin
GREEN LIGHT (525nm) Smooth wrinkles and fine lines
PURPLE LIGHT (390nm) Promote lymphatic drainage
YELLOW LIGHT (590nm) Reduce melanin production
CYAN LIGHT (490nm) Reduce skin tension and swollen capillaries
WHITE LIGHT (mix) Relax skin & accelerate growth and healing of skin cells

Using Steps

1. After washing your face, apply your favorite skin care products evenly to your face.
2. Press the power button for 2 seconds to turn on the device.
3. Press the button to select the light mode and massage mode according to your needs.
4. Place the massage head against your skin and massage from bottom to top.
5. Press the power button for 2 seconds to turn off the device after use.

Notes

1. This device cannot be used when charging.
2. The massage head needs to fit against the user's skin before it starts to vibrate
3. If there is no operation for one minute, the device will automatically shut down.

Using Tips

1. Massage from chin to cheek
2. Massage from cheek to outer eye corner
3. Massage from neck to chin
4. Massage the forehead from bottom to top


Safety Precautions

1. It is generally not advisable for people who have undergone facial plastic surgery, subcutaneous metal implant and those who are pregnant or within the period to use this product.
2. Wipe the massage head with tissues or cloth after use. Do not wash it with water.
3. Please not put it in the bathroom or other humid places, as it is easy to cause a short circuit.
4. Please not use corrosive solution (alcohol, banana oil, gasoline, nail polish remover and etc.) to wipe the device, which may lead to the failure of the product or the cracking, discolor of components and others.
5. Each mode can only be used once each day, and it is recommended that three modes are used in succession once each day.

Warranty Card		
Product Name	Intense Pulsed Light Wrinkles Reducing Instrument	Model No.
Distributor		
Name		
Address		
Tel	()	Mobilephone
Defective Phenomenon		
For repair purpose, please fill in detail. Main device: 1 year warranty		

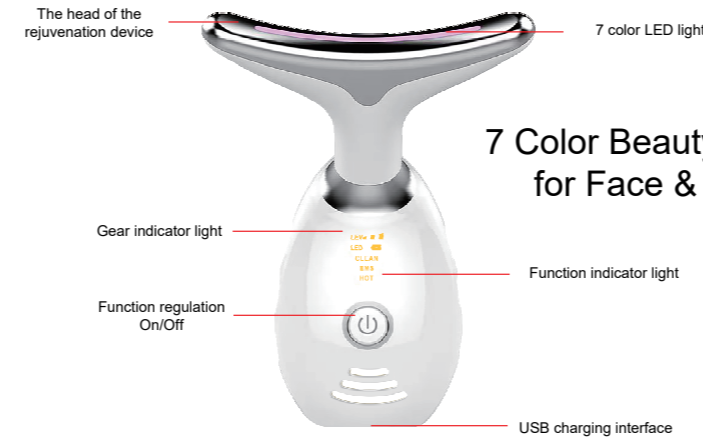Certificate	
Production Date	
Quality testing	

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