
Ultra Facial Lift
With Microcurrents & Light Therapy

EP-4008

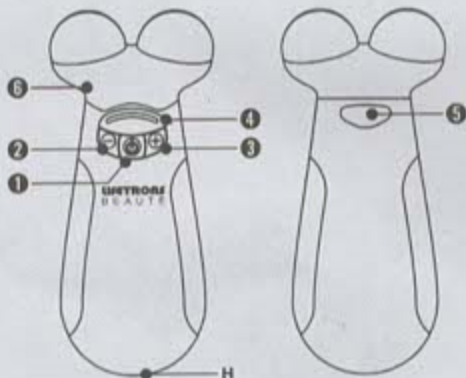Thank you for purchasing this Lifetons® product. This user guide outlines the operating procedures and safety precautions. Please read it carefully to make sure you are using your product properly and retain it for future reference.

How to Use

A. Changing the Treatment Head

1. Press the release button **1** on the back of the device to release the treatment head.
2. Select the treatment head **2** / **3** / **4** you would like to use and attach it to the body of the device.

Notes: The device will not vibrate when using the red photon light head. Do not change the treatment head when the power is on.

B. Microcurrent Facial Massage Treatment

1. Use the microcurrent facial massage treatment head to define and contour larger facial areas such as the jaw, jawline, cheeks and forehead.
2. Cleanse skin and be sure to remove makeup before use. Gel, serum, moisturizer or cream (not included) must be applied in combination with the steps below to ensure conductivity, proper operation and maximum benefits.
3. Attach the facial treatment head **1** to the body of the device. Press the power button **1** 1 time to turn on the device. Adjust the intensity level to your comfort by pressing the up (+) or down (-) button located on the front of the device. The default intensity level is level 0 (no microcurrents).
4. Gently press both spheres of the head onto your face while lifting and holding skin in place for 5 seconds or gently glide the device in a linear, upward motion. See drawing for more details.
5. The device will automatically shut down after 10 minutes of continuous use. Follow-up with your daily skincare routine after treatment.
6. For best results use this treatment for 10 minutes, 2-3 times per week.

C. Microcurrent Eye & Lip Massage Treatment

1. Use the eye and lip massage treatment head to target areas around the eyes and lips such as the lower lip, upper lip, smile lines, lower eyelid, upper eyelid and eyebrows.
2. Cleanse skin and be sure to remove makeup before use. Gel, serum, moisturizer or cream (not included) must be applied in combination with the steps below to ensure conductivity, proper operation and maximum benefits.
3. Attach the eye and lip treatment head **2** to the body of the device. Press the power button **1** 1 time to turn on the device. Adjust the intensity level to your comfort by pressing the up (+) or down (-) button located on the front of the device. The default intensity level is level 0 (no microcurrents).
4. Gently press both spheres of the head onto your face while lifting and holding skin in place for 5 seconds or gently glide the device in a linear, upward motion. See drawing for more details.
5. The device will automatically shut down after 10 minutes of continuous use. Follow-up with your daily skincare routine after treatment.
6. For best results use this treatment for 10 minutes, 2-3 times per week.

D. Red Photon Light Therapy

1. Use the red photon light therapy to target wrinkles, fine lines, dark spots and even out skin tone while helping to promote collagen production.
2. Cleanse skin and be sure to remove makeup before use.
3. Attach the red photon light head to the main body of the device.
4. Put on the eye protection goggles.
5. Turn on the device to start the treatment.
6. Place the lights in direct contact with the skin and hold for 5-8 seconds before moving to the next area. You can treat your entire face using this technique.
7. The device will automatically shut down after 10 minutes of continuous use. Follow-up with your daily skincare routine after treatment.
8. For best results use this treatment for 5-10 minutes, 2-3 times per week.

Notes: Always wear the eye protection goggles when using the photon light therapy. It is dangerous and harmful to your eyes to look directly into the lights.

E. Intensity Level Settings

There are 4 intensity levels for the microcurrent facial and eye treatment. We suggest starting with the first level and gradually increasing it over time depending on your skin's needs. Please note that the default level is level 0 (no current).

1. Press **1**: Power on and intensity level 0 begins (1 blue light appears)
2. Press **2**: 2 times: Switch to intensity level 1 (2 blue lights appear)
3. Press **2**: 3 times: Switch to intensity level 2 (3 blue lights appear)
4. Press **2**: 4 times: Switch to intensity level 3 (4 blue lights appear)
5. Press **2**: 5 times: Switch to intensity level 4 (5 blue lights appear)
6. Press **2**: to lower the intensity level

Notes: The device will automatically shut down after 10 minutes of continuous use. If you want to shorten the usage time, simply press the power button for 1 second. The device will beep 3 times before shutting down.

F. Recharging the Device

1. Connect part **1** (micro-USB jack) to H (micro-USB port).
2. Connect part **2** (USB jack) to a USB port with power supply (e.g., AC adaptor with USB power supply).

G. Battery Indicator

Low battery: All blue LED lights flash 5 times and the device will beep 5 times before shutting down (only when the device is in use)

When charging: Blue LED lights **1** appear to indicate the charging status

Fully charged battery: No lights appear

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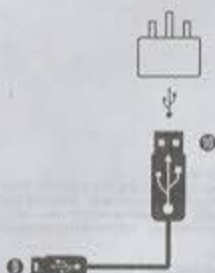
Safety Precautions

Please follow the rules outlined in this guide closely to ensure your safety and flawless operation:

- This product is only intended for use on the face.
- Discontinue using the product immediately if any abnormal situation, allergic reaction or discomfort occurs during use.
- This product is a skin care device and is not made for medical purposes. If you are under medical treatment, please consult with your doctor before use.
- This product is not suitable to use on body parts that are embedded with silicone, metal or other foreign objects. Do not use directly on the eyes.
- It is recommended to wait at least 6 months after any facial treatments such as a face lift or Botox injections before using this product. If you are unsure, consult with your doctor.
- Thoroughly clean the product after use.
- Do not use any kind of liquid to clean the product. Keep in mind that the surface of the product should only be wiped with a clean dry cloth.
- This product has a special lacquer - when cleaning, do not use any corrosive solution which may cause damage to the lacquer or appearance.
- Repair, removal or modification of this product is strictly prohibited. Improper handling may result in short circuits, product damage or even fire. Maintenance and warranty services will be void if unauthorized changes have been made to the product. Please contact an authorized dealer to obtain after-sales service if the product has a problem.
- Be cautious of any edges, uneven surfaces, metal parts, accessories or packaging to prevent possible injury or damage. Take note that this product should not be placed in your mouth.
- Do not use this product in locations with high temperatures or high humidity levels (e.g., bathroom). Ideally, the product should be used in a location free of excessive dust and not exposed to direct sunlight.
- Do not insert foreign objects or liquids into the product; this could damage the product.
- Refrain from dropping the product or bumping it against other objects. Avoid damage, collision and excessive pressure by always handling the product with care.
- Do not place this product near strong magnetic or electric fields - this may cause product failure to occur.
- Do not use this product in hazardous locations.
- The product may explode if disposed of in a fire.
- Do not use internally and avoid using near center of the throat to avoid thyroid stimulation.
- Do not use on open sores, wounds or infected areas. Never use internally (i.e., in mouth, nose, etc.).
- Do not use if pregnant, have a pacemaker, heart disease history or other medical condition.

Children

Keep out of reach of children. Do not allow children to play with your product or its accessories. They could hurt themselves or others and could accidentally damage the product. Your product and its accessories may contain small parts that could become detached and create a choking hazard.

Disposal of Electrical and Electronic Equipment

This symbol indicates that all electrical and electronic equipment included with this product should not be treated as household waste. Instead, it should be left at the appropriate collection point for the recycling of electronic equipment. By ensuring this product is disposed of correctly, you are helping to prevent negative consequences for the environment and human health. To recycle this product, please contact your local city office, environmental protection agency, household waste disposal service or where you purchased the product to learn more about appropriate recycling methods.

Battery

It is recommended that you completely use the initial charge up, then recharge the product at least two more times to ensure the best possible performance.

A rechargeable battery can be charged and discharged over one hundred times. However, it will eventually wear out - this is not a defect and corresponds to normal everyday use. When the operating or standby time is noticeably shorter, it may be a sign that the battery has reached the end of its life span. Please do not make repairs, modifications, or installations without contacting an authorized Lifetrons representative.

Limited Warranty

If within 12 months of purchase this product or any part thereof is proven to be defective by reason of faulty workmanship or defect in materials, we will repair or replace the same product free of charge of labor or materials on condition that:

- The original invoice or sales receipt with the date of purchase is presented with the defective unit.
- Lifetrons does not take any responsibility for lost information of any kind and will not reimburse you for any such loss.
- The limited warranty is only valid if the original proof of purchase for this product, issued by a Lifetrons authorized dealer, specifying the date of purchase, is presented with the product to be repaired or replaced. Lifetrons reserves the right to refuse warranty service if this information has been removed or changed after the original purchase of the product from the dealer.
- The warranty does not cover product failures which have been caused by use of accessories or other peripheral devices which are not Lifetrons branded original accessories intended for use with the product. The warranty does not cover any damage caused by misuse, abuse, neglect or accidental breakage.
- The warranty does not cover product failures caused by modifications, repairs or disassembly by any person not authorized by the manufacturer.
- Any defective unit or part shall become property of the manufacturer.

Additional Information

Register your Lifetrons Beauty product online at www.lifetrons.com.

Please visit our website www.lifetrons.com or contact our customer service center at customer cs@lifetrons.com to learn more about our new products and special offers.

This user guide was originally composed in English and was subsequently translated into other languages. The facility of subsequent translations cannot be guaranteed. In case of conflict between the English version and another language version, the English version takes precedence.

Product specifications are subject to change at the sole discretion and decision of Lifetrons Switzerland AG, Lifetrons International Ltd. or U-tac without prior notice. Any devices or accessories not listed under "What's Inside" are in reference to compatibility only and are not included. Product images are for reference purposes only; color may vary from the actual product.

Lifetrons Beauty, "Lifetrons Beauty Switzerland", "Lifetrons Beauty Professional Skincare", "Lifetrons Switzerland", "BEAUTÉ SWITZERLAND", "U-tac" are trademarks or registered trademarks owned by Lifetrons Switzerland AG, Lifetrons International Ltd. or U-tac. All other brands and product names listed herein are trademarks or registered trademarks owned by their respective holders and are hereby acknowledged and recognized as such. ©Copyright 2016 by Lifetrons Switzerland AG, Lifetrons International Ltd. and U-tac. All rights reserved.

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"Lifetrons BEAUTÉ", "Lifetrons BEAUTÉ SWITZERLAND", "Lifetrons BEAUTÉ PROFESSIONAL SKINCARE", "Lifetrons SWITZERLAND", "BEAUTÉ SWITZERLAND", "U-tac" are trademarks or registered trademarks owned by Lifetrons Switzerland AG, Lifetrons International Ltd. or U-tac. All other brands and product names listed herein are trademarks or registered trademarks owned by their respective holders and are hereby acknowledged and recognized as such. ©Copyright 2016 by Lifetrons Switzerland AG, Lifetrons International Ltd. and U-tac. All rights reserved.