



LED FACE MASK

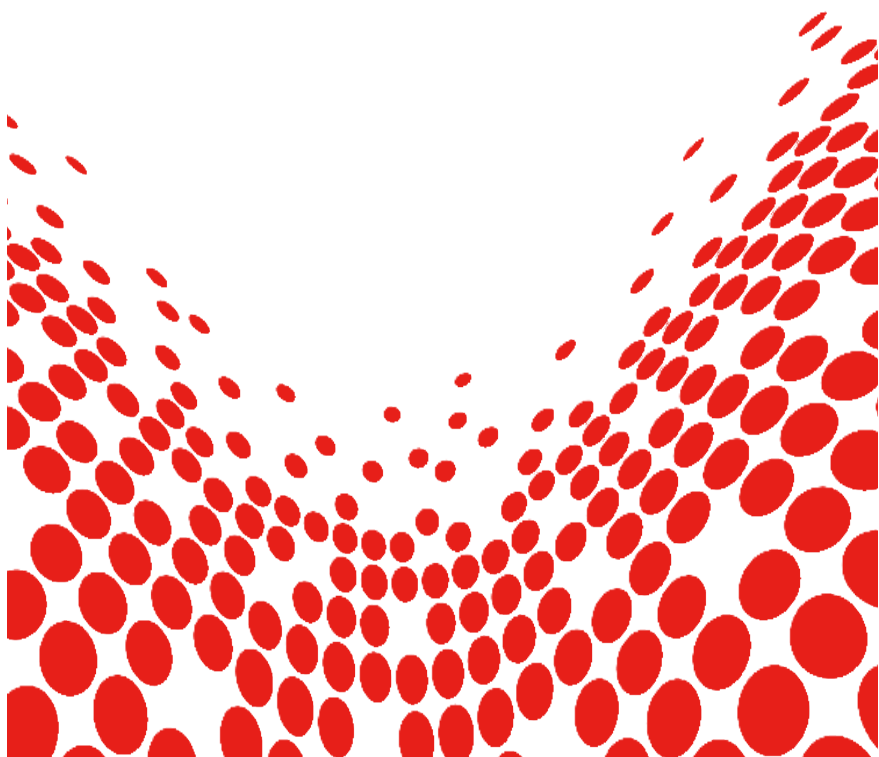

WHAT'S INSIDE?

Thank you for purchasing our LED face mask!

Our masks incorporate the optimal combination of 660nm red and 850nm near infrared wavelengths, which provides the most effective treatment. Give yourself or someone you love younger-looking clearer, brighter skin. Please read the user manual carefully before use, and check whether the product is damaged or not when you receive the product. If you have any questions or concerns, please feel free to contact us. We will be happy to assist you!

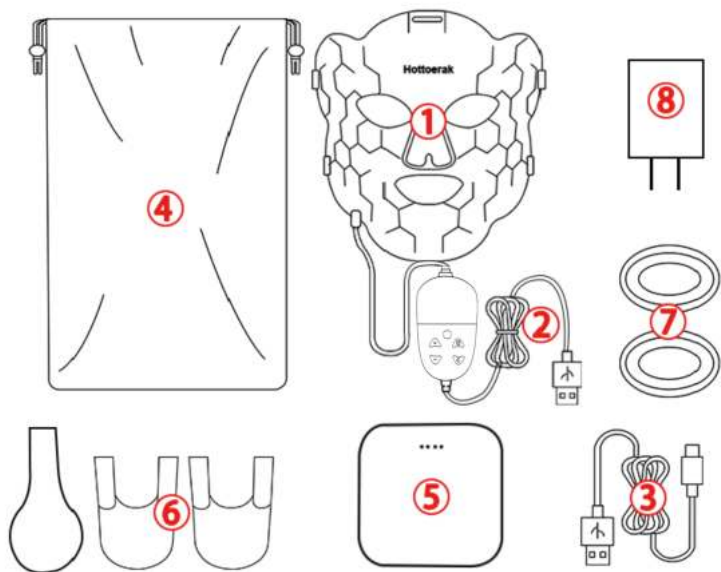

- ① 1x Light Therapy Mask
- ② 1x Controllers
- ③ 1x Charging Cable
- ④ 1x Carry Bag


- ⑤ 1x Power Bank
- ⑥ 1x Head Straps
- ⑦ 2x Pair of Eye Mask
- ⑧ 1x Adapter

PRODUCT SPECIFICATIONS

LEDs Wavelength	RED+NIR(660/850nm)
Mode	One Mode(Combines Wavelengths 660nm & 850nm)
Number of LEDs	74LEDs (3-in-1 Chips, 222Chips)
Working Time	2 Hours per Fully Charge
Timer	10/15/20 Minutes
Electrical Output	5V/2A
Battery Capacity	6000mAh
Operating Temperature	0°C to 40°C
Storage and Transport Temperature	-5°C to 55°C

WHAT IS LED LIGHT THERAPY?

660nm red and 850nm near infrared LED light therapy is a safe, non-invasive treatment for a variety of skin concerns. During LED therapy, the light waves generated by the LED device penetrate deep into the skin, triggering natural intracellular responses. The skin layer contains a lot of blood and water, which can absorb these light waves. This helps stimulate the production of collagen and elastin, which can help reduce the appearance of fine lines, wrinkles, and dark spots.(850nm near infrared is invisible light.)

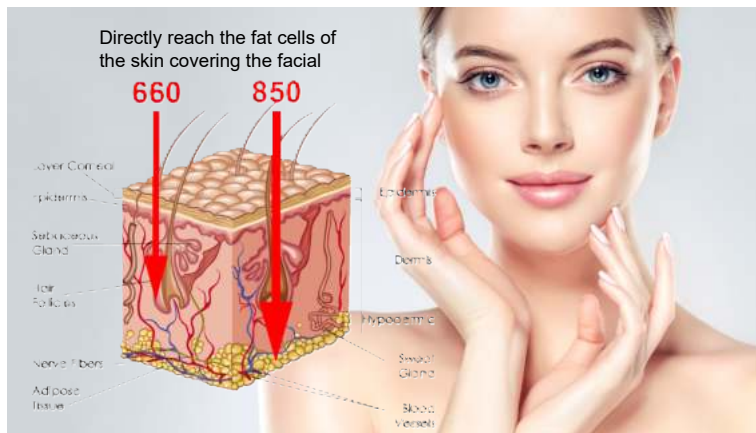

CONTROLLER INSTRUCTION




POWER BUTTON
Turn on/off



3 TIMING SETTINGS
10min/15min/20min options



BRIGHTNESS
"P1-P5"

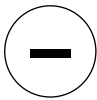

BRIGHTNESS
"P5-P1"



TIMING RESET

1. Press the power button to turn it on, use the default light intensity L1 for the first time, and use the default last memory function later.
2. To adjust the light intensity, use the plus or minus buttons provided. There are five levels of intensity to choose from.
3. Use the time button to select a timer setting of 10, 15, or 20 minutes. To exit timer mode, press the exit button.

ABOUT THE CHOICE OF TIMER USAGE

The duration of the red light therapy mask also affects its effectiveness for skin care. Here are the general effects of the product at different time intervals.

10 minutes: Helps to improve blood circulation and soothe skin.

15 minutes: This can promote collagen production and improve skin elasticity.

20 minutes: Helps to reduce inflammation and improve skin tone.

ABOUT THE CHOICE OF LEVEL USAGE

Generally higher intensity results in greater penetration of light energy into the skin, leading to better therapeutic outcomes. Our products are available in different strengths.

Level 1: Suitable for sensitive skin, it can relieve skin irritation, and redness, and improve skin tone.

Level 2: Suitable for daily skin maintenance, it can accelerate blood circulation, promote metabolism, and effectively reduce fine lines.


Level 3: Suitable for improving skin conditions, it can reduce acne, and blackheads, and improve skin elasticity and uneven skin tone.

Level 4: Suitable for treating deep skin problems, such as pigmentation and scars, and promoting skin cell regeneration and repair.

Level 5: Suitable for treating stubborn skin problems, such as deep wrinkles and pigmentation, with a stronger phototherapy effect that can promote skin cell renewal and repair.

RED AND NIR LIGHT & BENEFITS

Our masks use red and infrared LED light-emitting diodes, which produce uniform, cool, narrow-band light that penetrates deep into the epidermal layer of the skin to stimulate or activate skin cells. An LED mask that uses red and near-infrared light therapy has been scientifically proven to treat cosmetic indications including facial wrinkles.

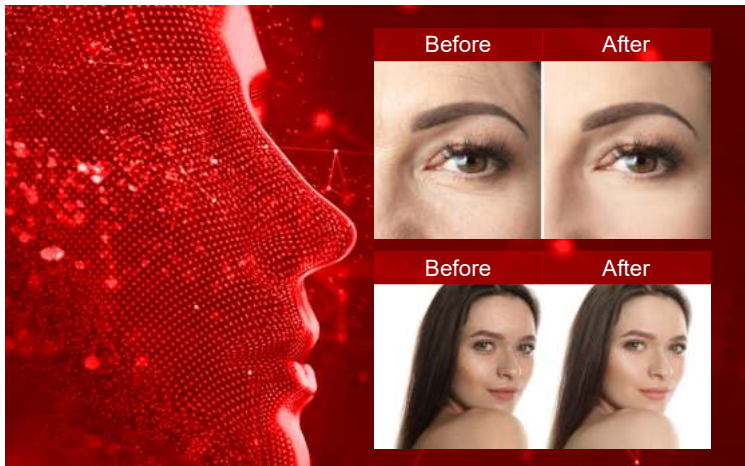

Red Light

620-750nm



Efficient: Whitening and rejuvenation, lighten dark spots, anti-wrinkle, firm skin, improve dark yellow, promote resolution of inflammation, improve the Immune system

Indications: Skin health and beauty, skin damage, acne, various inflammations, scald

Near-infrared
Lights

700-1000nm



Efficient: Repair skin, anti aging, anti-wrinkle, relieves allergy, symptoms, promote wound, healing, promote, metabolism, pain relief light, etc.

Indications: Joint pain, burn, sprain, diabetes, diffused redness visible capillaries dry, tight, flaky, or rough skin compromised, thin, or vulnerable skin, etc.

APPLICABLE & NON-APPLICABLE PEOPLE

Suitable for	Not suitable for
Patients with acne	Women who are pregnant
Psoriasis patients	Patients with epilepsy
Patients with sunburned skin	Patients with light sensitivity
Patients with red and swollen skin	Patients with cataracts
People with loose skin	Patients with glaucoma
People with itchy skin	Patients with acute skin diseases
People who need skin repair	Patients with eczema
People who need to produce collagen	Patients with dermatitis

+21.67% ↑
Skin elasticity

+23.66% ↑
Skin moisture content

+11.5% ↑
Smooth

+12.46% ↑
Gloss

-9.80% ↓
Nasolabial folds

-18.58% ↓
Canthus lines

Our mask allows every loved
beauty to have younger, clearer,
and brighter skin.


INITIAL USE VELCRO STRAPS INSTALLATION METHOD

Here are the steps:

Step 1: Place the face mask and the straps on a flat surface.

Step 2: Insert the end of the straps into the opening of the buckle and thread it through the buckle channel. Repeat the above steps for the remaining two sides, ensuring that the straps on both sides are correctly inserted into the buckles on the opposite side of the face mask.


Step 3: Press the strips firmly together to create a strong bond.

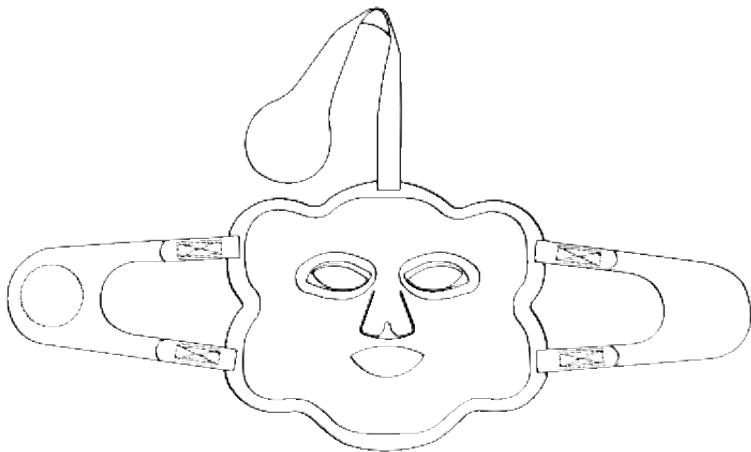

INITIAL USE EYE PAD INSTALLATION METHOD

Here are the steps:

- Step 1: Take out the face mask and the pair of eye pads from the packaging.
- Step 2: Flatten the face mask to make it ready for installation. (Make sure the inner side of the face mask is facing up.)
- Step 3: Insert the protruding end of each eye pad into the respective grooves on the eye area of the face mask, ensuring that the eye pads fit snugly to the upper part of the mask.
- Step 4: If you find that the eye pads are not correctly installed, readjust them until they are firmly fixed onto the face mask.
- Step 5: Finally, put on the face mask and adjust it for a comfortable fit on your face. The eye pads should gently touch your eyes without exerting excessive pressure.

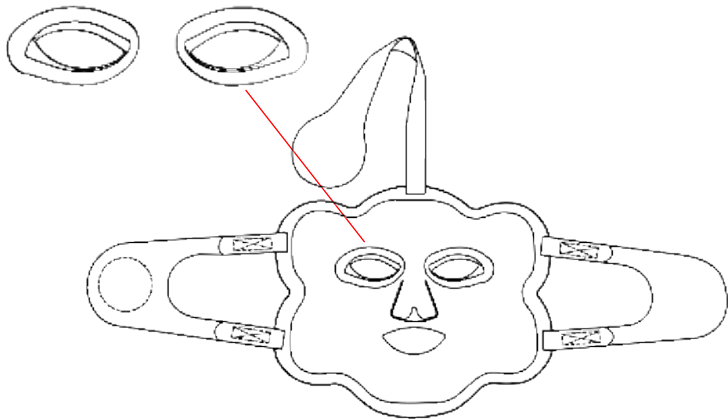


TO PROPERLY WEAR THE FACE MASK

Here are the steps:

Step 1: Place the mask gently over your face, ensuring it covers your nose and mouth completely. Take a moment to adjust the mask for a comfortable and snug fit, ensuring it sits comfortably on your face.

Step 2: Wrap the straps around the back of your head. Follow the step-by-step instructions in Strap Wearing Step to securely fasten the straps in the correct order.

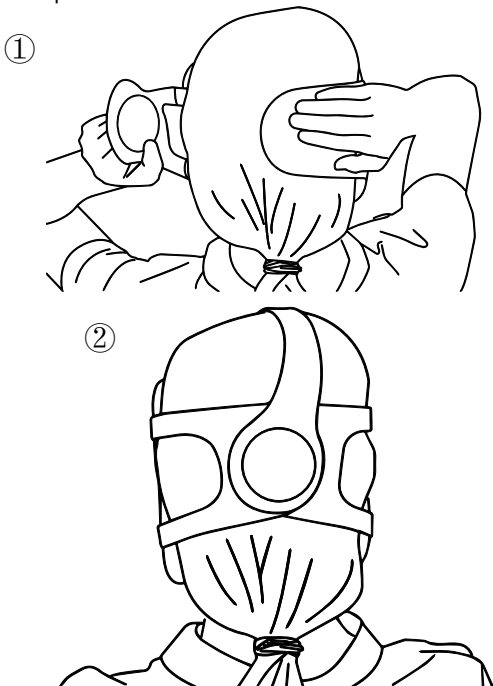

EASY-START SETUP GUIDE

01: GETTING READY TO USE THE RED LIGHT THERAPY MASK

Position the red light therapy mask in the area you wish to treat.

Utilize the velcro strap adjustment strap to fasten the red light therapy mask comfortably. It should be snug yet not overly tight.

02: CONNECTING THE CONTROL CABLE TO POWER

Connect the power adapter to the remote control cable and plug it into a wall outlet.

Or insert it into a fully charged power bank.

03: ACTIVATING THE RED LIGHT THERAPY MASK

Press and hold the On/Off button to power on the red light therapy mask.

Adjust the intensity level and time according to your preferences.


HOW TO USE AN LED FACE MASK?




1. Clean and dry the skin.



2. Adjust the eye pads and straps and place the mask on your face.

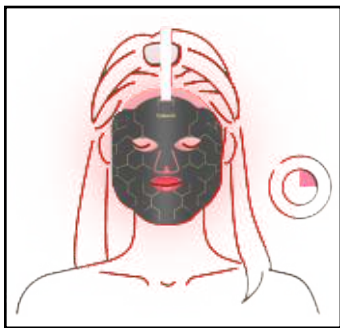

3. Press the button to turn on the mask.
Lie down or sit back and enjoy
the skin for 10-20 minutes.

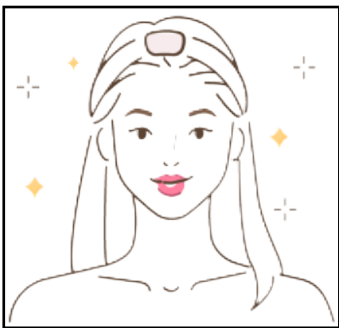

4. Follow your normal
skincare routine.

NOTE: Before using it the first time, insert the white velcro straps on the sides and top. Also, put in a set of eye pads.

ABOUT PRODUCT DETAILS

1. About The Peculiar Smell

Our products are made of food-grade silicone material, which is safe and non-toxic, It is normal to have a peculiar smell and will not affect your health, It is recommended that you open the package for 6-8 days after receiving the product before using it the smell will dissipate 95%

2. About The Nose And Eye Socket Have Indentations

It is normal to experience indentation when using eye pads. The indentation will disappear in 1-2 hours and will not leave a permanent mark. It is recommended that you lie down for more comfort.

3. About Wearing Comfort

We have three retention head straps and we recommend that you wear them all, this will soothe the hollows of your nose and jaw for a better fit against your skin.

4. About Eye Pads

Wearing eye pads allows you to do everything while you have your facial, and if you don't want to wear them, make sure to use them while lying down with your eyes closed.

5. About The Use Effect

Generally, after 4-6 weeks of regular use, the skin will become more delicate, and the effect of brightening and repairing will gradually appear. The most obvious effect is to promote the absorption of toner, which may also vary from person to person.

ADVANTAGES OF OUR MASK FACE

1. Made of medical-grade soft silicone, compared with stiff, hard plastic masks, it is more comfortable and fits any face shape, and the led light can penetrate and cover every area of the face.
2. 3-in-1 Chips technology and 660nm/850nm spectrum make facial skin more elastic. These two wavelengths work on the inner layers of the skin to rejuvenate it for a visibly radiant and healthier-looking complexion.
3. 10/15/20 minutes 3 timing modes. We usually recommend 15 minutes, choose longer after tolerance builds.
4. There is an adjustable velcro on the top of the head to avoid falling off. You can do yoga, read a book, or listen to music during light therapy.
5. Small design and equipped with a travel bag and power bank, you can take it with you to work or to a party.
6. Pay Attention: there will be some smell when you use it for the first time, and it will slowly dissipate later, which is a normal phenomenon and will not affect your use.
7. This mask offers the convenience of easy cleaning. Simply disconnect it from the power supply, gently wipe the surface with a water-based cloth or wet wipe, avoid harsh cleaning agents, let it air dry in a dry environment, and store it in a dust-free area using the provided case or pouch to ensure its well-maintained condition.

FREQUENTLY ASKED QUESTIONS

Will it harm my eyes?

Please do not look directly into the light (our devices come with eye pads to protect your eyes). Please wear eye pads when wearing eye masks to avoid glare from the light.

Can I use the device every day?

It is recommended to use for 15 minutes for each treatment with 3-5 treatments per week. Reduce to 3-4 treatments per week once results are visible. If you overuse the device, you may experience redness or dryness of the skin.

Can I put skin care products on my face before using the LED facial?

Yes. The red LED light stimulates new collagen production deep within the dermal layers of the skin. You can apply anti-aging serum and a hydrating sheet masque before using the red LED light therapy setting. **DO NOT** use heavy or greasy creams or lotions before treatments as these can cause a barrier between the skin and light penetration.

How long until I start to see visible results?

Consumer perception studies indicate that 95% of users showed effective improvement in skin tone and fine lines after 4~6 weeks.

What to expect from the LED facial treatment?

The LED Mask contains 222LEDs (3perLED bulb) pieces of 99% purity medical-grade LED light emitting diodes in Red, Near-infrared LED lights to smooth wrinkles, reduce pigmentation, and clear acne for a youthful complexion.

Is LED light therapy safe to use?

Yes. LED light therapy is generally safe and suitable for all skin types. We do however recommend you start slow with 10 minutes per session to test for individual tolerance.

CONTRAINDICATIONS

1. Stop using this facial mask if you encounter any prolonged skin redness or discomfort that lasts for more than 2 hours.
2. If a problem occurs for more than 2 hours you might be light sensitive or have another condition. We recommend you consult with a medical professional for advice.
3. Please do not use this product if you have lupus erythematosus, photosensitive eczema, or leukemia. Otherwise, it may cause severe skin reactions.
4. Please do not use this product if you are sensitive or dizzy to light. Otherwise, it may cause severe skin reactions.
5. It is recommended not to use this product if you have a migraine.
6. Please do not use this product if you have a history of seizures caused by light.

SAFETY WARNING

1. Do not place the mask where there is water or dripping water.
2. If the mask light is scratched or intruded on by liquid, please stop using it immediately.
3. Do not use the mask in the sauna or in the rain.
4. If the mask, controller, or wire is damaged, please do not use it again.
5. Do not pull the USB cable too hard to avoid the mask cable falling off.
6. Do not let children use the mask without supervision.

Thank you for choosing our products!

We are pleased to offer you the following:

- 60-day money-back guarantee
- 1-year free warranty



Made in China

