

Product introduction

All over body massager:

It is full body powerful vibrating massager that helps you relieve all of your body's tension and sore muscles after a long day work. It can be used on feet, legs, neck back & Shoulder.

Specification

1. Battery: Rechargeable lithium battery
2. Material: High-Grade silicone + ABS
3. Vibration Noise: Around 50dB
4. Vibrating mode: 8 different speeds and 20 vibration patterns
5. Water resistance
6. Size: 1.69*8.27 inch
7. Net weight: 6.35oz/180g

For any query contact:-shiatu@drphys.io

User Manual

Before using the product:


1. Change the massager each time before use by connecting it with USB charger cable to "DC" marked at the end of massager. We recommend to connect to the mobile phone charger (output voltage is 3-5V, current: 1-2A) or computer USB port (output voltage is 3-5V, current: 1-2A) for charging. (Wall Charger is not included)
2. Light blinks during charging & will glow white once the charging is done.
3. Clean the product after use with liquid soap and warm water. Gently pat dry with a clean towel.

How to Use:

1. Press the power button "⊙" to turn it off/on.
2. Short press the frequency conversion button "⊗" to switch between 20 patterns of the vibration, each time in turn.
3. Short press the speed change button "⊕" to switch between 8 kinds of speeds intensity, each time in turn.
4. Total 20 different vibration patterns and 8 different speeds intensities can be changed.

Note: Please charge the massager before first use.

The massager will not work if short of electricity.

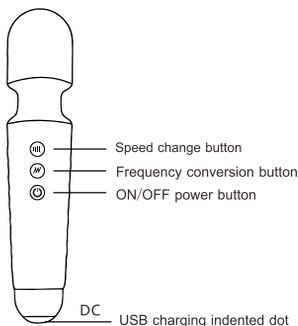

- Press the power button "⊙" to turn it off/on.
- Short press the frequency conversion button "⊗" to switch between 20 patterns of the vibration, each time in turn.
- Short press the speed change button "⊕" to switch between 8 kinds of speeds intensity, each time in turn.



SPORTS BODY WAND MASSAGER

