

REVITAGLOW

LIGHT THERAPY MASK

USER MANUAL





HELLO THERE.

Joel and I are athletes who suffered from severe ACL injuries a few months apart. We purchased multiple pieces of wellness equipment for our in-home rehab but were left disappointed and frustrated by wellness and recovery devices that made big promises but gave no support, no well-thought-out plan to follow, and no results. And that's how Lifepro came about, and why the Wellness Recovery collection was created — we've set out to change all of that. It's why every piece of equipment bearing our name not only undergoes rigorous scrutiny to make sure it delivers on power and efficacy, but also comes with built-in support from our 'guiding angels' and a built-in lifetime warranty so that you're never left hanging.

Welcome to Lifepro.
We're in this together now.

Abe & Joel

You don't need expensive trips to the spa to get smoother, softer, more supple skin. The Lifepro RevitaGlow Light Therapy Mask provides professional-level LED light therapy treatments from the comfort of your own home. Light therapy, also called phototherapy, is a highly-effective, non-invasive, chemical-free way to safely and painlessly repair and rejuvenate your skin. RevitaGlow provides three light therapy treatment options in one simple, hands-free solution:

- **Red (630–645 nm) light** improves blood circulation and boosts collagen and elastin production, reducing the appearance of fine lines, wrinkles, and scars. It tightens pores, smooths skin tone, and reduces symptoms of eczema, psoriasis, and rosacea including redness, itching, flaking, swelling, and broken capillaries.
- **Green (520–535 nm) light** boosts lymphatic drainage and flushes toxins, balancing oily skin and reducing hyperpigmentation.
- **Blue (470–485 nm) light** has anti-bacterial benefits and reduces inflammation. It eliminates the bacteria that causes acne, regulates over-active oil/sebum glands, and reduces pigmentation, sun damage, and skin scaling from psoriasis.

Use the RevitaGlow Light Therapy Mask three to four times a week for 15–20 minutes, and see results in as little as three weeks. The RevitaGlow is lightweight and portable. Take your beauty regimen to the next level and enjoy smoother, softer skin anytime, anywhere!

Redefine wellness and beauty as you experience the transformative power of the Lifepro RevitaGlow:

- **Treat acne:** kills acne-causing bacteria, regulates oil glands, and reduces acne scarring.
- **Reduce lines and wrinkles:** boosts blood flow and increases collagen and elastin production.
- **Detoxify skin:** improves lymphatic drainage and flushes toxins.
- **Improve symptoms of common skin conditions:** reduces itching, flaking, redness, swelling, and inflammation.

Check out the information in this manual for ways to get started using the RevitaGlow. Next, be sure to visit our website revitaglow.lifeprofitness.com for access to our **FREE library of videos** which demonstrate how to easily incorporate the RevitaGlow into your health and wellness routine.

WELLNESS RECOVERY REVITAGLOW LIGHT THERAPY MASK

GET STARTED

WHAT'S IN THE BOX

- Lifepro RevitaGlow Light Therapy Face Mask
- Controller
- USB Charging Cable
- Velcro Straps (2 Pairs)
- User Manual


FIRST STEPS

1. Remove the RevitaGlow and accessories from the shipping box and inspect each item for shipping damage.
2. Carefully read all safety warnings and operating instructions before using the RevitaGlow.
3. Go to revitaglow.lifeprofitness.com and register your product within 14 days of purchase to activate your lifetime warranty.
4. Complete the steps in the **EASY-START SETUP GUIDE** on page 9 to begin using the RevitaGlow.

CUSTOMER SUPPORT

If you have any questions about setting up the RevitaGlow, we're here to help. Call or email us and an experienced customer support representative will be happy to assist you: (800) 563-6604 or support@lifeprofitness.com.

PARTS OF THE REVITAGLOW

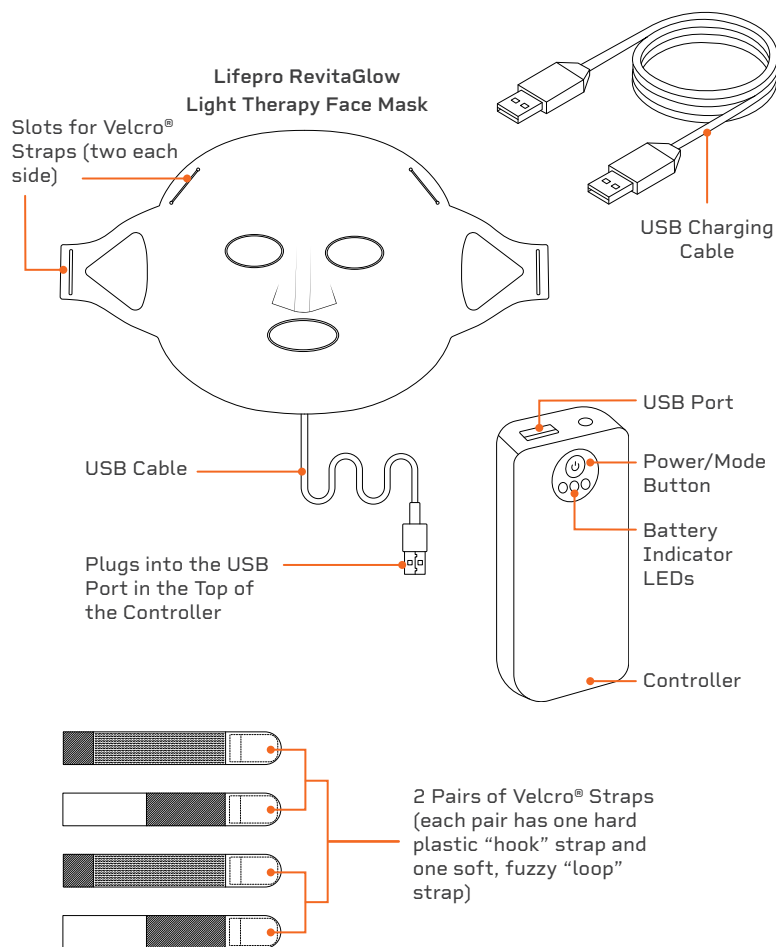

TIPS FOR USING THE REVITAGLOW

- Thoroughly wash and dry your face before using.
- We strongly recommend discontinuing any Retin-A use at least 5 days prior to using the RevitaGlow.
- Lie down on a flat comfortable surface, such as a bed or couch during treatment. We recommend spending the treatment time practicing blissful, uninterrupted relaxation and self-care, such as listening to music or meditating.
- Place the RevitaGlow in position on your face before turning it on.
- Check the treatment area frequently and stop immediately at the first sign of pain or discomfort.
- We recommend 15–20-minute treatment sessions, 3–4 times a week.
- We do not recommend using for more than 20 minutes per day.
- People with sensitive skin or skin conditions that can be exacerbated by heat, such as eczema and rosacea, may want to use the device only for short durations.
- Optionally, you may apply a light, soothing moisturizer before using.

⚠ CAUTION

- Do not look directly into the LEDs when the mask is turned on, as this may result in damage to the eyes.
- Do not use this product in combination with ointments, salves, liniments, balms, or other topical pain-relievers that contain ingredients that produce heat (such as capsaicin), as this may result in burns.
- Do not use sheet masks or harsh products such as acids, acne treatments, or retinols under the RevitaGlow.
- Do not wear metal jewelry or accessories, such as lip rings or nose rings, in the treatment area, as metal may heat and cause burns.
- Do not use directly on swollen or inflamed skin or open wounds.
- Keep out of reach of children and pets.

TESTING FOR SKIN SENSITIVITY

Before your first treatment, it is important to perform a skin sensitivity test to ensure you won't experience any adverse reactions to light therapy.

To do a sensitivity test:

1. Plug the face mask into the controller and place the mask on your inner arm, with the LEDs against your skin.
2. Press and hold the Power/Mode button on the controller for three seconds to turn on the RevitaGlow and allow the entire 10-minute treatment program to run.
3. If your skin feels uncomfortable or hot at any point during the test treatment, stop using the RevitaGlow immediately.
4. Wait at least 6 hours after the sensitivity test. If your skin becomes red and/or itchy, you may be sensitive light therapy and you should not continue to use the RevitaGlow.

NOTE: Please see a list of contraindications (people who should not use this device) on page 20.

FCC STATEMENT

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

NOTE: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.



EASY-START SETUP GUIDE

PREPARING YOUR FACE FOR TREATMENT

1

- Wash your face and carefully pat dry with a soft towel.
- Apply a light, soothing moisturizer to your face, if desired.

SETTING UP THE MASK

2

- Locate the two pairs of Velcro® straps.
- Attach the Velcro® straps to the mask following the instructions on page 10.
- Plug the USB cable exiting the bottom of the face mask into the standard USB port at the top of the controller.

BEGINNING YOUR TREATMENT SESSION

3

- Lie down on a flat comfortable surface, such as a bed or couch.
- Place the RevitaGlow in position on your face, with the LED-side against your skin.
- Fasten the Velcro® straps to secure the face mask to your face following the instructions on page 11.
- Press and hold the Power/Mode button on the controller for three seconds to turn on the RevitaGlow.
- Your light therapy treatment will begin immediately.
- Adjust the color mode as desired. See **OPERATING INSTRUCTIONS** on pages 12 and 13 for details.

NOTE: Charge the battery before first use if only one battery indicator LED is lit when the device is turned on for the first time. See **CHARGING THE BATTERY** on page 12 for instructions.



ATTACHING THE FACE MASK

ATTACHING THE VELCRO STRAPS TO THE MASK

1. Locate the two pairs of Velcro® straps.
 - a. Each pair has one “hook” strap (with stiff plastic hooks) and one “loop” strap (with soft, fuzzy loops).
 - b. One pair of straps will attach to the mask behind your head and the other will attach to the mask across the top of your head.

2. Attach the first pair of straps:

- a. Insert the square end of one strap through a slot in the side of the face mask and fold the end over to secure the strap.
- b. Insert the second strap through the other slot in the side of the face mask and secure it in the same manner.
- c. Make sure the two straps can be attached (the hooks and loops should face each other).

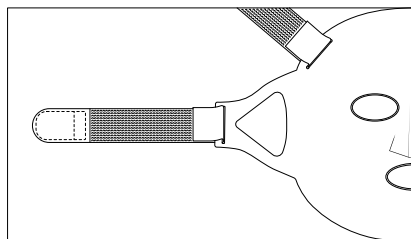
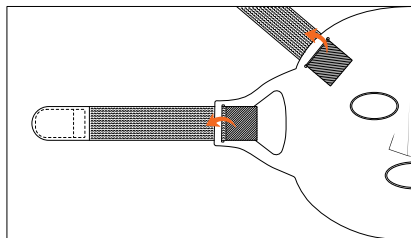


3. Attach the second pair of straps:

- a. Insert the square end of one strap through one slot in the top edge of the face mask and fold the end over to secure the strap.
- b. Insert the square end of the second strap through the second slot in the top edge and secure it in the same manner.

SECURING THE MASK TO YOUR FACE

1. Place the mask on your face with the LEDs facing your skin.
2. Fasten the straps at the sides of the mask behind your head.
3. Fasten the straps on the top edge of the mask across the top of your head.
4. Make sure the mask is snug to your face, but not too tight. Readjust the straps if necessary.

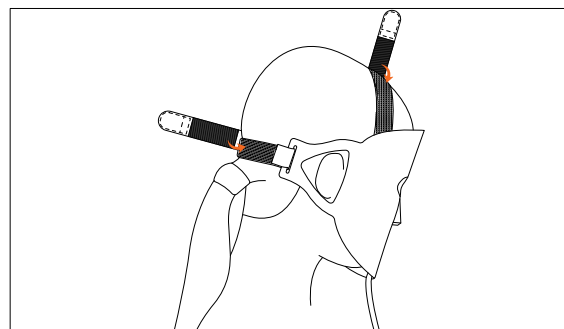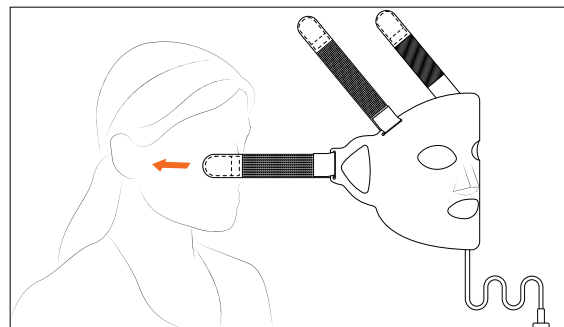

OPERATING INSTRUCTIONS

Please read the operating instructions before using the RevitaGlow Light Therapy Mask. Keep this manual for reference.

CHARGING THE BATTERY



- Do not use the RevitaGlow while it is charging. Always ensure it is powered off before charging.
- To charge, plug the USB charging cable into the USB port on the top of the controller. Then plug the other end of the charging cable into:
 - A USB port on a computer or laptop, OR
 - A 5V-1A or 5V-2A USB charging adapter. Then, plug the charging adapter into a wall outlet.

NOTE: While a 5V-2A charging adapter may be used, the device's maximum input current is 1A. This means that the device will not charge faster with a 2A adapter.

- Your RevitaGlow will typically take about 4 hours to fully charge.
- The battery indicator LEDs on the controller indicate the remaining battery:
 - 3 LEDs Lit = 75-100% battery
 - 1 LED Lit = less than 40% battery
 - 2 LEDs Lit = 40-74% battery
- We recommend charging the battery right away when only 1 battery indicator LED remains lit.
- The RevitaGlow will automatically power off when the battery is depleted.
- The average run-time on a full battery charge is 3-3.5 hours.
- For best results, only charge the battery at temperatures of 68-104° F.

USING THE CONTROLLER

1. **USB Port:** Plug the mask's USB connector into the USB port to begin using.
2. **Power/Mode Button:** Press and hold the button for three seconds to turn on / off the RevitaGlow. Tap the button to toggle through the three color therapy modes, red, green, and blue. The default color mode is red.
3. **Battery Indicator LEDs:** Three battery indicator LEDs indicate the remaining battery charge.

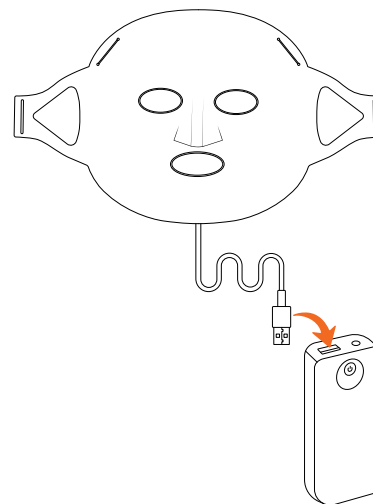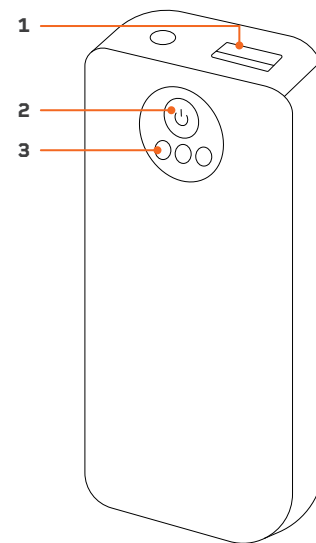

USING THE REVITAGLOW LIGHT THERAPY MASK

1. Prepare to use the RevitaGlow by following the instructions in the **EASY START SETUP GUIDE** on page 9.
2. Lie down on a flat comfortable flat surface, such as a bed or couch.
3. Place the RevitaGlow in position on your face, with the LEDs against your skin, and fasten the Velcro® straps to secure the face mask. See **SECURING THE MASK TO YOUR FACE** on page 11.
4. Press and hold the **Power/Mode button** on the controller for three seconds to turn on the RevitaGlow. Your 10-minute treatment session will begin immediately. The default treatment is red (630–645 nm) light therapy.
5. Tap the **Power/Mode button** to change the color therapy mode, if desired. Tap the button to toggle through the three color therapy modes, green (520–535 nm), blue (470–485 nm), and red (630–645 nm). See **CHOOSING A COLOR THERAPY MODE** on the next page.
6. The RevitaGlow will shut down automatically when the 10-minute treatment session completes. Press and hold the **Power/Mode button** for three seconds to stop the session at any time.
7. Always power off the RevitaGlow after each use.

DETERMINING TREATMENT FREQUENCY AND DURATION

ISSUE/ SEVERITY	ACUTE	CHRONIC	MAINTENANCE
Frequency <i>(How Often)</i>	1–2 times daily	1–2 times daily	3–4 times weekly
Duration <i>(How Long)</i>	10 minutes	10 minutes	10 minutes

NOTE: Do not use for more than 20 minutes per day.

CHOOSING A COLOR THERAPY MODE

The RevitaGlow Light Therapy Mask has three light therapy color options:

MODE	WAVELENGTH	BENEFITS
RED	630–645 nm	Anti-aging benefits including boosting collagen and elastin production. Improves blood flow. Reduces appearance of fine lines and wrinkles. Reduces symptoms of eczema, psoriasis, and rosacea, including redness, itching, flaking, and swelling.
GREEN	520–535 nm	Boosts lymphatic drainage and flushes toxins, balancing oily skin and reducing hyperpigmentation.
BLUE	470–485 nm	Has anti-bacterial benefits and reduces inflammation. Treat mild to moderate acne by eliminating the bacteria that causes acne, regulating over-active oil/sebum glands, and reducing the appearance of acne scars. Reduces pigmentation and sun damage. Reduces skin scaling from psoriasis.

TROUBLESHOOTING GUIDE

Troubleshoot simple issues you might experience with your RevitaGlow using the table below. Contact a customer support representative for additional assistance at: **(800) 563-6604** or **support@lifeprofitness.com**. Please do not attempt to repair the device on your own.

ISSUE	POSSIBLE REASON	HOW TO RESOLVE
DEVICE IS NOT WORKING	Device is not charged.	Charge the device.
	Device is not turned on.	Press and hold the Power/Mode button to turn on the device.
DEVICE STOPPED WORKING SUDDENLY	Device automatically shuts off after the 10-minute working time.	Turn the device on to begin a new treatment session.
TREATMENTS DO NOT APPEAR TO BE WORKING	The color mode selected is not optimal for the desired results.	Tap the Power/Mode button to choose a different color therapy mode from the three options: red (the default), green, or blue.
	Not enough treatments have been performed.	It may take up 3-6 weeks of treatments to see results.

FREQUENTLY ASKED QUESTIONS

Q: WHAT IS LIGHT THERAPY

A: Light therapy, also called phototherapy, is a highly-effective non-invasive, chemical-free therapy that can be used to safely and painlessly heal and rejuvenate many parts of the body.

Light therapy uses specific wavelengths, or colors of light to treat the skin, both at and below the surface. Different wavelengths of light have different treatment applications and benefits.

Q: WHAT ARE THE BENEFITS OF LIGHT THERAPY

A: The benefits of light therapy vary greatly depending on the wavelength of light used. The benefits of light therapy include:

- **Red (630–645 nm) light** improves blood circulation and boosts collagen and elastin production, reducing the appearance of fine lines, wrinkles, and scars. It tightens pores, smooths skin tone, and reduces symptoms of eczema, psoriasis, and rosacea including redness, itching, flaking, swelling, and broken capillaries.
- **Green (520–535 nm) light** boosts lymphatic drainage and flushes toxins, balancing oily skin and reducing hyperpigmentation.
- **Blue (470–485 nm) light** has anti-bacterial benefits and reduces inflammation. It eliminates the bacteria that causes acne, regulates over-active oil/sebum glands, reduces pigmentation and sun damage. Reduces skin scaling from psoriasis.

Q: HOW MANY LEDS DOES THE REVITAGLOW HAVE?

A: The mask has 52 triple-LED chip sets. Each chip set contains 1 red 630–645 nm LED, 1 blue 470–485 nm LED, and 1 green 520–535 nm LED.

Q: WHO SHOULD NOT USE THE REVITAGLOW (CONTRAINDICATIONS)?

A: Please see the full list of contraindications on page 20 of this manual. We recommend everyone perform a skin sensitivity test before using this device. See instructions on page 8.

Q: HOW OFTEN CAN I USE THE REVITAGLOW?

A: We recommend using the RevitaGlow for 15–20 minutes, three to four times a week for optimal results. We do not recommend using the device more than 20 minutes per day.

Q: WHAT IS THE BEST TREATMENT OPTION FOR ACNE?

A: Blue (470–485 nm) light is best for treating and preventing acne. It also works well to treat acne scarring.

Q: WHAT IS THE BEST TREATMENT OPTION FOR ANTI-AGING?

A: Red (630–645 nm) light is great for anti-aging. In addition, it boosts blood flow and production of collagen and elastin. This results in reduced appearance of fine lines and wrinkles, smoothed skin tone, and tighter pores.

Q: HOW DO I CLEAN THE REVITAGLOW

A: The mask is made from food-grade silicone, so it is easy to clean. Simply wipe down the inside of the mask with a soft, dry cloth after each use. Make sure to turn off the mask and allow it to cool down before cleaning. We recommend washing and drying your face before use to keep the mask clean. **WARNING!** Do not allow water or other liquids to come into contact with the device. Do not use abrasive, corrosive, solvent, or chemical cleaners or harsh detergents.

Q: HOW SHOULD I STORE THE REVITAGLOW

A: Always store the mask in its original box to prevent damage to the mask. Make sure to disconnect the controller before storing the mask in the box. Store the device in a cool, dry location. Do not store close to the ground, walls, or roof. Do not store other objects on top of the mask it could be crushed or damaged.

Q: HOW CAN I EXPERIENCE THE BENEFITS OF LIGHT THERAPY ON OTHER AREAS OF THE BODY?

Check out some of our other light therapy products at the end of this manual and on our website, lifeprofitness.com.

SAFETY INSTRUCTIONS

Keep this manual in a safe place for future reference.

Before using this product, read all safety warnings and operating instructions. Failure to do so may result in burns, fire, electric shock, serious injury or death, or damage to the product or property.

CONTRAINDICATIONS

- **Do not use this product if you:**
 - Have swollen or inflamed skin, or open or bleeding wounds in the treatment area.
 - Are, suspect you may be, or are planning to become pregnant; or if you are, or plan to begin breastfeeding.
 - Have cancer or suspect you may have cancerous lesions. Do not use over malignant or benign tumors.
- **Consult your physician before using this product if you:**
 - Are under the care of a physician for a pre-existing health condition.
 - Have ever had symptoms related to epilepsy (seizures or loss of consciousness) when exposed to flashing lights.
 - Take steroid medication, use a topical steroid, or receive steroid injections.
 - Have recently had surgery, LASIK or PRK eye surgery, or skin filler injections, if you have breast implants, or if you have a tattoo in the treatment area.
 - Have a condition that may make you sensitive to light or sunlight, such as Lupus, some forms of eczema, light-induced migraines, albinism, and rosacea.
 - Take a topical or oral medication that may cause photosensitivity (sensitivity to sunlight), including some antibiotics, antifungals, antidepressants, antipsychotics, oral contraceptives, psoralens, cardiac, diabetic, blood pressure, and cholesterol medications, chemotherapy drugs, pain medications, nonsteroidal anti-inflammatory drugs (NSAIDs), diuretics, cold and allergy medications, and retinoids (such as Retin-A and Retinol).
 - Use a product that may cause photosensitivity, such as alpha-hydroxy acids (like glycolic acid), beta-hydroxy acids (like salicylic acid), Retin-A and Retinol, melatonin, citrus essential oils, St. John's Wart, acne products, antibacterial soap, and coal tar and petroleum products.
 - Are or may be sensitive to red light or near-infrared light.

DISCLAIMER

- This product is not a medical device. This product and the information in this manual are not intended to diagnose, treat, cure, or prevent any disease. The health benefits suggested or implied in this manual, other product literature, and company website are not certified or endorsed by any regulatory authority or medical institute.
- The use of this product is entirely at the user's discretion. Please read all instructions and safety information carefully before using this product. We assume no responsibility and expressly disclaim any and all liability for personal injury, death, property damage or loss sustained through improper use, failure to follow the instructions and warnings in this manual, incorrect setup, or inadequate or incorrect maintenance or neglect of this product.
- It is the responsibility of the owner to ensure that all users of this product are adequately informed of all warnings and cautions.
- Use this product only as instructed in this manual.

⚠ WARNING

To reduce the risk of burns, fire, electric shock, serious injury or death, or damage to the product or property:

- Only use the provided charging cable at the specified voltage to charge the device. Using a charging cable other than the one provided may result in fire or explosion.
- Do not use the provided charging cable to charge other devices.
- Do not modify, tamper with, or damage the charging cable.
- Never unplug the charging cable by pulling on the cord.
- Never use a modified or damaged charging cable to charge the device.
- Do not modify, tamper with, damage, or attempt to disassemble the device. This device has no user-serviceable parts.
- Always inspect the device, charging cable, and controller for damage, wear and tear, or fraying before each use. Do not use if damaged, malfunctioning, or if the device has been exposed to water.
- Keep the device and charging cable away from flammable or explosive items, chemicals, fire, and heat sources, and direct sunlight at all times.
- Do not use or charge in high temperatures or in direct sunlight. Do not use or charge under a blanket or pillow. Excessive heat can cause battery explosion, fire, electrocution, serious injury, death, or damage to property.
- Do not operate or charge in the presence of aerosol (spray) products or corrosive gases, or where concentrated oxygen is being administered.

- Do not use the device while it is charging.
- Keep away from liquids. Do not use outdoors, in high-moisture environments or in any place where the device may become wet or damp.
- Do not immerse in water. If the device becomes wet, do not use.
- Do not plug in or unplug the charging cable with wet hands.
- Do not use the device over wet skin.
- Do not use during electrical storms.
- Keep sharp objects away from the device and charging cable.
- Do not use the charging cable with an extension cord.
- Never leave the device unattended while running or charging. Unplug the charging cable as soon as the device is fully charged. Do not overcharge.
- Stop using the product immediately if you feel pain or discomfort.
- Do not operate the device while sleeping, driving, smoking, or under the influence of drugs or alcohol.
- Do not use this product on sleeping or unconscious people or infants.
- For adult-use only. Keep away from children and pets.
- This device is not intended for use by adults with diminished mental, sensory, or physical capacity, unless they are closely supervised by a responsible adult.

⚠ CAUTION

To reduce the risk of minor or moderate injury:

- Do not look directly into the LED lights while the device is powered on as this may result in damage to the eyes.
- Ensure the mask is positioned properly on your face so your eyes are level with the eye openings to prevent light exposure to your eyes.
- Do not use this product in combination with ointments, salves, liniments, balms, or other topical pain-relievers that contain ingredients that produce heat (such as capsaicin), as this may result in burns.
- Do not use sheet masks or harsh products such as acids, acne treatments, or retinols under the mask.
- Do not wear metal jewelry or accessories in the treatment area, as metal may heat and cause burns.
- Keep away from long hair to avoid getting it caught in the Velcro®.

CLEANING AND STORAGE

⚠ WARNING! Do not allow water or other liquids to come into contact with the device. Do not immerse the device in water or other liquids.

- Use only over clean skin to keep the device clean.
- Turn off, unplug, and allow to completely cool before cleaning or storing.
- Wipe the device with a soft, dry cloth after use.
- Do not use abrasive, corrosive, solvent, or chemical cleaners or harsh detergents to clean the device.
- Store the device in a cool and dry indoor location.
- Do not store other objects on top of the device. Store in the original box to prevent damage.
- Do not attempt to disassemble or repair the device. This device contains no user-serviceable parts.

NOTICE

- Do not sit, stand, or place objects on the mask. Excessive force or weight applied to the device can damage it.
- Do not wrap the charging cable around the device to avoid damaging it.
- Keep away from sharp objects that could damage the device or LEDs.
- Do not insert objects into the USB ports.
- Unplug the device immediately in the event of a power failure.
- This device is intended for in-home use only. Do not use in any commercial, rental, institutional, or therapeutic setting.
- Dispose of device and all accessories according to state and federal regulations.
- Do not use for more than 20 minutes a day.



LIFETIME WARRANTY

Lifepro Fitness creates quality, durable exercise products and we stand by that quality with a lifetime warranty on all of our products. If the RevitaGlow ever breaks (and we doubt it will), we'll send you replacement parts and show you how to repair it. If the RevitaGlow cannot be repaired, we'll replace it—free of charge. Register the RevitaGlow at revitaglow.lifeprofitness.com to activate your warranty within 14 days of purchase.

PRODUCT
SPECIFICATIONS

UNIT SPECS

Size:	13.2" × 9.8" × 0.2"
Net Weight:	1.1 lb
Rated Power:	3 W
Total LEDs:	52 Triple-LED Chip Sets
1 Red:	630–645 nm (52 total)
1 Blue:	470–485 nm (52 total)
1 Green:	520–535 nm (52 total)
LED Life Expectancy:	50,000 hrs
Auto Shut-off Timer:	10 minutes
Modes:	3
Auto Shut-off Timer:	10 minutes
Operating Temperature:	68–104 °F
Material:	Food-grade Silicone
Certifications:	CE, ROHS

BATTERY

Rated Voltage:	3.7 V
Battery Type:	18650 Lithium Ion
Battery Capacity:	2000 mAh
Use-Time:	2–3.5 hours
Charging Time:	3–4 hours

CHARGING CABLE

Rated Input:	5V–1A
Rated Output:	3.7V–2A

WELLNESS RECOVERY
REVITAGLOW
LIGHT THERAPY MASK
MADE IN CHINA

FROM OUR WELLNESS RECOVERY COLLECTION

ALLEVARED

LIGHT THERAPY BELT

With the Lifepro AllevaRed Light Therapy Belt, you can experience the therapeutic benefits of red light therapy from the comfort of your own home. No need for expensive trips to the spa or doctor's office. Red light therapy provides a wide variety of benefits, including pain relief and quicker recovery after injury. Together, red (660nm) and near-infrared (850nm) light treat both the surface of the skin and deep below the surface, helping cells make additional energy to repair themselves, improving skin health, and providing anti-aging and fat loss benefits. Just a few treatment sessions a week can help you achieve improved wellness.

The easy-to-wear belt provides full coverage for your lower back and waist, and can be used on your upper back, arms, legs, and feet. Your pet can also take advantage of the therapeutic benefits of the AllevaRed. Help your pet find relief from chronic joint pain, lose inches, and reduce inflammation.

The AllevaRed offers customizable treatment options for your wellness recovery. Five intensity levels, three modes (red, near-infrared, and red + near-infrared), and an optional pulsed light feature allow you to choose the perfect settings for each therapy session.

ALLEVARED BENEFITS

- Reduce muscle and joint pain
- Reduce inflammation
- Rejuvenate skin
- Blast fat
- Boost circulation



BIOHEAL

RED LIGHT THERAPY PANEL

Use the Lifepro BioHeal Red Light Therapy Panel in conjunction with the RejuvaWrap add the benefits of red light therapy to your wellness program. With the BioHeal, you can experience the therapeutic benefits of red light therapy from the comfort of your own home—no need for expensive trips to the spa or doctor's office. Red light therapy provides a wide variety of benefits, including enhanced muscle growth and tone, pain relief, and accelerated workout recovery.

Together, red (660nm) and near-infrared (850nm) light treat both the surface of the skin and deep below the surface, helping cells make additional energy to repair themselves, improving skin health, and providing anti-aging and increased muscle tone and fat loss. The BioHeal also improves blood flow, sleep quality, mental clarity, and focus. The BioHeal comes with two pairs of protective eyeglasses, so you can safely treat your face or head.

The BioHeal therapy panel has 60 LEDs and multiple hanging options so that you may effectively target and treat any area of the body. It offers customizable treatment options for your wellness recovery, including pulsed red light therapy and adjustable treatment time from 1–30 minutes. Three modes, (1) red + near-infrared, (2) red light, and (3) near-infrared light, allow you to choose the perfect settings for each treatment session.



BIOHEAL BENEFITS

- Boost muscle growth and tone
- Speed up recovery
- Reduce joint pain and stiffness
- Boost blood flow
- Rejuvenate skin
- Lose fat and inches around the waist
- Reduce inflammation

BIOREMEDY

INFRARED SAUNA BLANKET

With the Lifepro BioRemedy Infrared Sauna Blanket, you can experience the myriad therapeutic benefits of infrared saunas from the comfort of your own home—no need for expensive trips to the sauna or spa. The nontoxic, sweat-resistant sauna blanket is heated by silicone heating wires, which produce infrared waves that penetrate up to an inch and a half into muscle and fat, warming your body from the inside out. There is no hot air, so it's less intense than a sauna and safer for sensitive skin. The BioRemedy is perfect for anyone who wants to improve their cardiovascular health, reduce pain from chronic conditions such as arthritis and rheumatoid arthritis, improve skin health, lose inches around the waist, and reduce stress. Just a few treatment sessions a week can help you achieve improved wellness.

The BioRemedy offers customizable treatment options for your wellness recovery. The temperature is adjustable from 113–176 °F and you may set a treatment time from 30–60 minutes. Use the BioRemedy 3 or 4 times a week for 30 minutes, and see results in as little as three weeks.

The BioRemedy is lightweight, folds up small, and comes with a carry bag, making it ideal for taking with you on the go. No other sauna is this portable and cost-effective.



BIOREMEDY BENEFITS

- Alleviate pain
- Increase metabolism
- Reduce inflammation
- Reduce joint pain and stiffness
- Boost mood and decrease stress
- Rejuvenate skin
- Lose inches around the waist



Access Lifepro TV, our library of free personal training videos on the web at revitaglow.lifeprofittness.com to learn how to get the most out of your Lifepro® RevitaGlow and achieve the results you want.



**YOU'VE
GOT
THIS**