---

## NASAL ASPIRATOR:

---

**To use:**

- Hold baby upright.
- While squeezing the aspirator bulb, gently place the tip of the aspirator into baby's nostril and tilt it to the side. Never point the tip directly in the center of the nostril. Be careful not to over-insert the aspirator.
- Slowly release the bulb allowing the suction to draw the mucus out of baby's nostril.
- Repeat with other nostril.
- Repeat as necessary until baby's experiences relief from congestion.
- To clean, put the aspirator under soapy, warm water. Squeeze and release the bulb several times to clean the inside of the bulb. Rinse with warm water. Squeeze out all excess water from the bulb.

---

## EMERGENCY CARD:

---

- Keep important phone numbers and medical information handy for quick reference.
- Keep in a safe place. Familiarize family members, baby-sitters and care givers with the location of the card and the information on it.

---

## MEDICINE SPOON:

---

**To use:** Clean before first use. Hold the Medicine Spoon upright and pour the correct dosage of liquid medicine into the handle, usin the dosage markings for reference. To help ensure accuracy, keep the medicine spoon upright and at eye level. Make sure your child is sitting up, place the spoon in the child's mouth and gently tilt it upward.

**To clean:** Wash in warm, soapy water. Air dry. Top rack dishwasher safe. Do not place in boiling water.

---

## DIGITAL THERMOMETER:

---

For complete instructions for use and care, please see separate instruction manual enclosed.

---

## STORAGE CASE:

---

- Keep items organized with the convenient storage case.
- To clean:** Surface wash with mild soap and warm water. Rinse thoroughly; air dry.

**14444, 14454**  
Please read the following instructions and warnings carefully. Keep this instruction manual for future reference. Adult assembly required.  
**For Technical Support Call 1-800-268-6237 or E-Mail at customerservice@summerinfant.com.**  
For use with children from birth & up.

---

## INSTRUCTION MANUAL

---

### Health & Grooming Kit

---

---

#### BRUSH & COMB:

---

Soft brush bristles are gentle on baby's scalp. Wide-toothed comb makes untangling hair easier. Handle design with soft grip texture is easy for toddlers to hold and use.

**To use:**

- Wash brush & comb before using.
- Untangle hair from the roots out. Keep a hand on the roots while brushing or combing to avoid pulling on child's scalp.
- Encourage your toddler to try brushing his or her own hair with your assistance (approximately 3 years old).

**To clean:** Wash with mild soap and warm water. Rinse thoroughly; air dry.

---

#### NAIL CLIPPERS:

---

This standard size nail clipper's handle rotates 360 degrees so you can safely trim nails from any angle. Soft over-molded handle is comfortable and helps you maintain a steady grip.

**Note:** Do not clip too close to baby's skin.

We recommend clipping baby's nails while he or she is sleeping.

**To clean:** Sterilize with isopropyl rubbing alcohol and rinse in warm water. Wipe dry.

---

#### 2 EMERY BOARDS:

---

- Gently glide across baby's fingernails to smooth them.

**CE** Digital thermometer manufactured by:  
0044 ONBO Electronic (Shenzhen) Co., Ltd.  
China

We **LOVE** our consumers!  
Get to know us at:  
**www.summerinfant.com**

Summer Infant stands behind all of its products. If you are not completely satisfied or have any questions, please contact our Consumer Relations Team at **1-800-268-6237** or **www.summerinfant.com/contact**  
(For Europe call +44 (0)144 250 5000 or customerserviceuk@summerinfant.com)

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**Summer Infant Europe, LTD.**  
1st Floor North Wing, Focus 31  
Cleveland Road  
Hemel Hempstead, Herts  
HP2 7BW UK  
0144 250 5000

01/13

# STANDARD DIGITAL THERMOMETER

14374, 14444, 14454, 14464, 14474

Read the instructions carefully before using this device.

### IMPORTANT NOTES:

Performance of this device may be degraded if it is operated outside of stated temperature humidity range, stored outside of stated temperature/humidity range, dropped or subjected to mechanical shock, and/or patient temperature is below ambient temperature.

- High, prolonged fever requires medical attention. Be sure to contact your physician.
- For parental use. Do not leave child unattended with thermometer.
- Avoid exercise and drinking hot or cold beverages before taking a reading.
- Clean the probe tip with isopropyl alcohol or mild soap and warm water. Clean before first use and after each use.
- Thermometer is water resistant.
- Store thermometer in protective case when not in use. Store at room temperature and out of direct sunlight.

### HOW TO TAKE TEMPERATURE:

Keep in mind that taking a temperature can be trying for both parent and child. When taking a reading be sure to remain calm and confident. The sound of your voice can be very soothing for an infant. Explaining what you are about to do and why can be helpful for a toddler or older child.

1. Remove thermometer from protective case. Be sure to clean thermometer before first use and after each use.
2. Push the ON-OFF button located on the large end of the thermometer. The display will read: "188.8  $\frac{9}{5}$ ". This is a function check indicating the unit is working properly. If a reading is in memory it will be displayed briefly. When the °F or °C begins to flash, the thermometer will begin measuring temperature when exposed to body heat.
3. The °F (or °C) symbol will continuously flash until completion of the temperature taking process. Please note that the "  $\nabla$  " indicates a low battery. An "ERR" message in the window indicates your unit is not functioning properly. Do not expect an accurate temperature reading if this low battery indicator is present.


### ORAL USE: (for children 4 years and older)

- Do not drink hot or cold fluids, exercise or perform other activities that will raise or lower temperature readings when compared to your child's normal, average temperature. The mouth should remain closed up to 5 minutes before attempting a reading.
- Place the probe tip well under the tongue as indicated by the "★" (see Figure A). Instruct children not to bite thermometer.
- The peak temperature should be reached in approximately 60 seconds. Opening the mouth or improper placement of the probe tip could result in a longer time for a reading.
- Follow numbered instructions below.

### UNDERARM (AXILLARY) USE:

- Wipe underarm with a dry towel.
- Place probe tip in armpit and bring arm down against the body. (Figure B) The probe should be pointing upward in the same direction as your body if possible.
- The time required for an accurate underarm temperature may be 90 seconds or longer depending on placement of the probe and closure of the arm against the body.
- Important Note: An underarm reading will be approximately 1 degree lower than an oral reading. To calculate the oral equivalent add 1°F or 0.5°C to the underarm temperature reading.
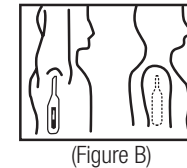
### RECTAL USE: (Infants and small children)

- Commonly used for babies, young children or when it is difficult to take an oral or underarm temperature.
- Lubricate probe tip with a water-soluble jelly or petroleum jelly.
- For easier insertion, position baby on his tummy, keeping his legs at a 45° angle to his back. It may be helpful to position the baby on your lap. Spread baby's buttocks with your thumb and forefinger. Gently insert probe a MAXIMUM of 1/2 in or 1.3 cm into baby's rectum and hold thermometer lightly in place. Do not force probe tip if you encounter resistance.
- The peak temperature should be reached in approximately 40 seconds.
- Follow numbered instructions below.
- Important Note: A rectal reading will be approximately 1 degree higher than an oral reading. To calculate the oral equivalent subtract 1°F or 0.5°C from the rectal temperature reading.
- Once you have used the thermometer to take a rectal reading, do not take any oral readings, for hygienic reasons.

### 4. When the °F or C° symbol stops flashing and the display reads a constant temperature, the peak temperature has been reached.

5. Read and record the temperature and time for reference. The reading will automatically be stored in the memory.
6. Push the on/off button to turn the unit off. This will prolong the battery life. If you do not turn the unit off, it will shut off automatically in approximately 10 minutes.

### BATTERY INFORMATION:

When the "  $\nabla$  " appears in the lower right of the display, the battery is exhausted. Dispose of the thermometer as the battery is non-replaceable. Battery Contains Mercury. Do Not Put in Trash. Recycle or Manage as Hazardous Waste. Properly dispose of product according to your local regulations concerning disposal of this product. The battery is well sealed inside the thermometer to reduce the risk of being swallowed by children. Keep away from small children and heat.

**Storage:** Room temperature 10 °C to 40 °C / 50 °F to 104 °F. Humidity: 15% to 95% (noncondensing)  
**Operating temperature:** -25 °C to 60°C / -13°F to 140°F. Humidity: 15% to 95% (noncondensing)

### WARRANTY

This Thermometer is guaranteed for a period of 24 months from date of purchase. This warranty does not cover batteries, or damage resulting from improper handling or use. Specific legal rights stipulated in your state may guide this warranty. If you have a legitimate claim under this warranty, please return the product to the place of purchase along with your sales receipt. If you have any questions or comments please contact our customer relations department.

Digital thermometer manufactured by: ONBO Electronic (Shenzhen) Co., Ltd. China

Type BF applied part

Batteries must be disposed of in accordance with the locally applicable regulations, not with domestic waste.

**Distributed by:**  
**Summer Infant, Inc.**  
1275 Park East Drive  
Woonsocket, RI 02895 USA  
1-800-268-6237  
© 2005 Summer Infant, Inc.

**Summer Infant Europe, LTD.**  
1st Floor North Wing, Focus 31  
Cleveland Road  
Hemel Hempstead, Herts  
HP2 7BW UK  
0144 250 5000

visit us at [www.summerinfant.com](http://www.summerinfant.com)

Please retain information for future reference. Colors and styles may vary. **MADE IN CHINA.**

FRONT

## RADIO AND TELEVISION INTERFERENCE

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to connect the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and the receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

You may also find helpful the following booklet, prepared by the FCC: "How to Identify and Resolve Radio-TV Interference Problems." This booklet is available from the U.S. Government Printing Office, Washington D.C. 20402.

Changes and Modifications not expressly approved by Summer Infant, Inc. can void your authority to operate this equipment under Federal Communications Commissions rules.

This device complies with part 15 of FCC rules. Operation is subject to the following two conditions; (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

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Termómetro digital fabricado por: ONBO Electronic (Shenzhen) Co., Ltd. China