
Health & Wellness

Quick Start Guide

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
English

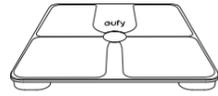
- **Not** intended for use by children under 16.*
- **Do not** use this device if you are pregnant.
- **Do not** use this device if you are wearing a pacemaker or any other internal medical devices.

* Considering the vast individual differences among children under 16, the measurement of body fat or other body composition measurements may be less accurate. Their weight-only measurement should remain accurate.

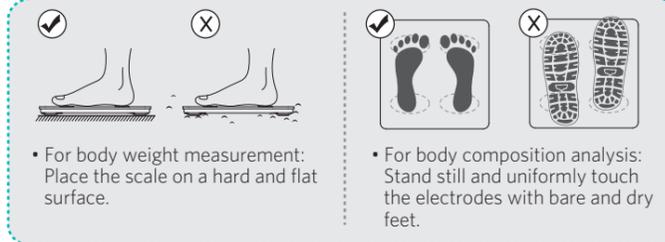
What's Needed

- 1 Compatible smartphone or tablet.
 - iOS 9.0 or above
 - Android 4.4 or above
- 2 Smart Scale P1 with batteries installed.

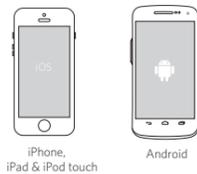iPhone, iPad & iPod touch


Android

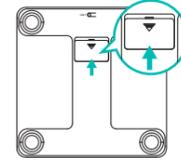01

Installing the EufyLife App and Pairing Devices

- 1 Download the EufyLife app from the App Store or Google Play.
- 2 Open the EufyLife app and then follow the on-screen instructions to pair.

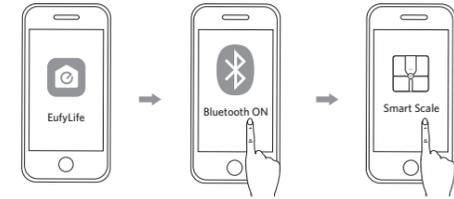Do not pair the scale from your smartphone's Bluetooth menu. Use the EufyLife app to pair the scale with your smartphone.

Using Your Scale

- 1 Place the scale on a hard and flat surface.
- 2 Open the EufyLife app and select Smart Scale on the device homepage. Then create a user ID if necessary.
- 3 Step barefoot onto the scale. An analysis of your body's composition will begin automatically.

Please read the owner's manual for more information


02

eufy

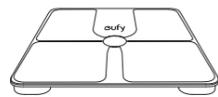For FAQs and more information, please visit:
www.eufylife.com

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