EMS-Halsmassagekissen mit Wärmefunktion
EMS Pulse Cervical Massage Pillow with Heat
Coussin Cervical Massant par EMS et Chauffant
Collo con Massaggio EMS e funzione Calore
Cojín de Masaje Cervical EMS con función de Calor
EMS-Pulsnekmassageapparaat met Warmte

BEDIENUNGSANLEITUNG – Bitte sorgfältig lesen!
INSTRUCTION MANUAL – Please read carefully!
MODE D'EMPLOI – A lire attentivement s.v.p.
ISTRUZIONI PER L'USO – Da leggere con attenzione!
INSTRUCCIONES DE MANEJO – ¡Por favor lea con cuidado!
HANDLEIDING INSTRUCTIE – Aandachtig lezen!

prorelax[®]

EMS Massagekissen **TRIO**

DEUTSCH	4
ENGLISH	12
FRANÇAIS	20
ITALIANO	28
ESPAÑOL	36
NEDERLANDS	44

Lieber *prorelax*® Kunde,

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prorelax[®]


Thank you for purchasing ***prorelax***[®]
EMS Pulse Cervical Massage Pillow with Heat

With normal care and proper treatment,
it will provide years of reliable service.

FEATURES	14
TECHNICAL DATA	14
OPERATION AND CONTROLLER INSTRUCTIONS	15
BATTERY CHARGING	16
PLEASE USE IN THE MANNER SHOWN BELOW...	16
SAFETY INSTRUCTIONS	17
CARE AND MAINTENANCE	18
REMOVABLE COVER	18
DISPOSAL	18
WARRANTY	19

FEATURES

- Electronic EMS pulse muscle Stimulator with **6 intensity levels**.
- 2 electrode plates plus 2pcs external electrode pads
- Optional Heat Therapy (2 heating levels)
- Built in rechargeable Lithium battery (USB cable charging)
- Adjustable drawcord with cordlock to fit neck size.
- Soft fleece cozy cloth, Removable and Washable

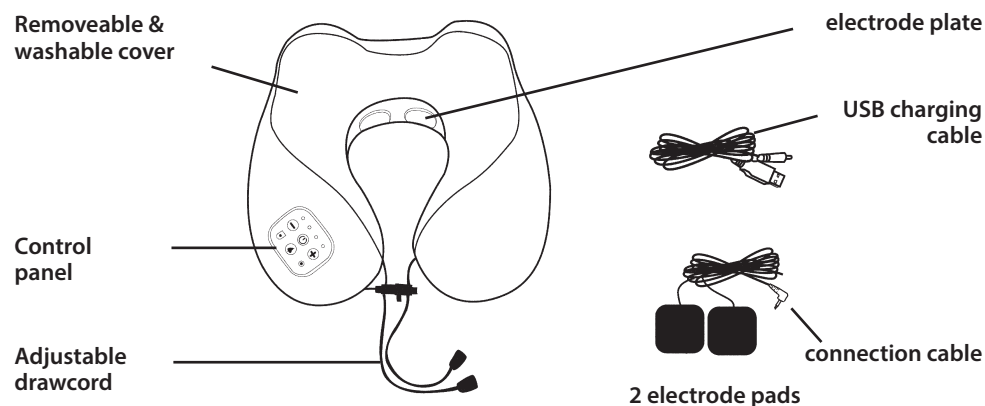TECHNICAL DATA

Measurements:	11.8*11*6.3inch
Weight:	1.38lbs
Battery:	3.7V,1000mAh.
Nominal power:	Max. 3 Watt
Battery Full Charging:	2.5-3 hours
Battery Max. Runtime:	3hours
Default Automatic Runtime:	15 minutes

OPERATION AND CONTROLLER INSTRUCTIONS

Indicator of vibration massage
GREEN → „ON“

EMS Indicator of Electrotherapy

GREEN → EMS - Mode 1
RED → EMS - Mode 2

- 1 Charging indicator +
- 2 Electrotherapy intensity indicator EMS

EMS
Increase intensity (+)
Decrease intensity (-)

Heat Indicator
YELLOW → Low
RED → High

Heat button
Press 1x → Low
Press 2x → High
Press 3x → Heat off

Port for external electrode pads

USB-charging-port

On-/off-button & program control (1-4):

Press 1x → Vibration massage + EMS-Mode 1
Press 2x → EMS-Mode 1
Press 3x → EMS-Mode 2
Press 4x → Vibration massage

Press and hold on 2-3 seconds to **turn off** the device.

EMS-Tip:

The electrical stimulation becomes more comfortable if you use an electrode gel or if you moisten the skin slightly.

1 Charging Indicator

While charging → Flash Red light
Fully charged → Green light


2 Electrotherapy Indicator EMS

The intensity of electrotherapy from low to high is displayed in **6 colours**:


- dark green → green → dark yellow → yellow → dark red → red +

BATTERY CHARGING


1. Connect the USB cable to the corresponding port on the device.
2. Connect the USB cable to any USB charger, it can be plugged in the portable power bank, mobile phone battery chargers, laptop and any USB port.
3. Normally it takes 2.5-3 hours to fully charge the battery.
4. While charging, the indicator light turns on flashing red color, when full charged, the indicator light will turn to green color.
5. If the battery is full, the device can work continuously for around 3 hours.
6. It's safe to use the device while charging.

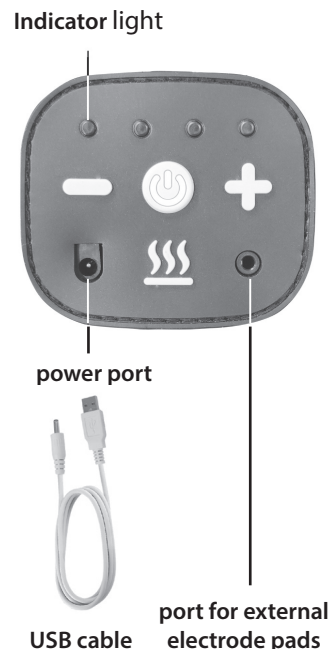PLEASE USE IN THE MANNER SHOWN BELOW...

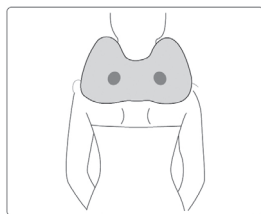NECK

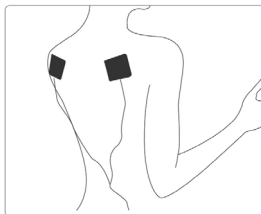BACK

- This device will switch off automatically after 15 minutes timer run out.
- Do not setup or use the massager in a bathroom or similiary wet / damp areas.
- Do not use this device while driving.

Likewise, you should allow your body some rest periods. To avoid straining your muscles too much, we recommend to not exceed a continuous massage of 15 Minutes.

SAFETY INSTRUCTIONS

Please read the following instructions carefully before putting your massage device into use to ensure trouble-free operation and optimum efficiency. Please retain these operating instructions for further use!

- The massage device conforms to recognized technical principles and the latest safety regulations.
- Do not wet, do not use pins.
- This items is NOT A TOY. Close supervision is necessary when this appliance is used by , on, or near children or disabled persons.
- This appliance should never be left unattended when plugged in.
- Any possible repairs may only be carried out by authorized expert staff. Improper use and unauthorized repairs are not permitted for safety reasons and lead to loss of warranty.
- Never touch the power plug with wet hands.
- Please avoid contact of the device with water, high temperatures and direct sunlight.
- Do not use any damaged cables, plugs or loose sockets.
- In case of malfunction, disconnect immediately from the mains.
- Do not use if you have skin disorders, open wounds, or swollen or inflamed areas.
- Misuse or incorrect use excludes any liability for damage.
- Do not use this device while driving.
- Do not use while you sleep.
- In order to avoid excessively stimulating the muscles and nerves, the recommended massage time should not be exceeded 15minutes at a time.
- Each massage – even a hand massage – must be refrained from during pregnancy or if one or more of the following complaints are present in the massage area: recent injuries, thrombotic diseases, all kinds of in flammations and swellings, and cancer. It is recommended to consult a doctor before having a massage for the treatment of ailments and illnesses.
- If you are dependent on electrical aids e.g. pacemakers, please be sure to consult your doctor for medical advice before taking a massage.
- The packaging material supplied should not be used as a toy. Failure to observe the above instructions may constitute misuse of the product and can cause serious injury or burns.

CARE AND MAINTENANCE

1. Do not wash in washing machine or immerse the whole product or the adapter in water.
2. Always switch off before cleaning and remove the adapter from the socket with clean, dry hands.
3. The device should be wiped with a slightly moist cloth and dried thoroughly. Do not wipe the product with alcohol or harsh chemicals such as thinner or bleach.
4. Keep away from heat, fire or exposure to direct sunlight.
5. Store the massager in a cool, dry place out of reach of children.
6. Avoid contact with sharp edges or pointed objects which might cut or puncture the fabric surface.

REMOVABLE COVER

- Unfasten the zipper and detach the cover
- Wash the cover in cold water with mild detergent.
Allow to air dry.

DISPOSAL

This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner.

EUROMEDICS GmbH, Hauptstrasse 169, D-59846 Sundern, Deutschland (hereinafter called "EUROMEDICS") provides a warranty for the Prorelax® EMS PULSE CERVICAL MASSAGE PILLOW WITH HEAT (hereinafter called the "Product") subject to the following preconditions and in the scope described below.

The warranty conditions described below do not affect the legally prescribed warranty obligations of the seller from the contract of purchase with the buyer. Moreover, the warranty applies irrespective of mandatory legal liability provisions, such as those according to the product liability law, or in cases of mala fide intent or gross negligence, or owing to injury to life and limb, or to health by EUROMEDICS or its agents.

EUROMEDICS guarantees the defect-free operability and the completeness of the product, especially that the product is free of material defects, manufacturing defects and construction defects. The decisive factor in this regard is the state of the art of the science and technology at the time of manufacturing. The product must have manifested the defect that caused the damage at that time itself.

The warranty period is 2 years from the date of purchase of a new, unused product by the purchaser.

The warranty applies worldwide, i.e. regardless of the country in which EUROMEDICS has marketed the product and the purchaser has purchased it.

This warranty only applies to products that the purchaser has procured as a consumer and exclusively for personal purposes in the framework of domestic use (i.e. not for products that were procured or used in the framework of a commercial or independent professional activity).

If the product is found to be deficient in operability in accordance with the following provisions during the warranty period, EUROMEDICS will provide a free replacement according to these warranty conditions.

If the purchaser wishes to file a warranty claim, he must first contact EUROMEDICS customer service:

EUROMEDICS KUNDENSERVICE

Biberweg 12
D-53842 Troisdorf
Tel: +49 (0)2208/9217999 Germany

E-mail: info@euromedics.de

The purchaser will then get further information on how to process the guarantee case, for example, where he can send the device free of cost and what documentation is required.

Filing a warranty claim can only come into consideration if the buyer can submit

- a copy of the invoice/proof of purchase and
- the original product

to EUROMEDICS or to a service centre nominated by EUROMEDICS upon demand.

Expressly excepted from this warranty are

- products that are used, cleaned, stored or maintained incorrectly and/or not in compliance with the conditions laid down in the operating manual, as well as products that were opened, repaired or modified by the purchaser or by a service centre not authorised by EUROMEDICS;
- products that were sold as second quality or used articles;
- consequential damages arising from a defect in the product (however, in this case, there may be entitlements from product liability or from other mandatory legal liability regulations).

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