

The background of the entire page is a solid blue color. It is covered with a repeating pattern of white line-art icons of various light bulbs. These icons include standard incandescent bulbs, CFL bulbs, and several different styles of LED bulbs, some with unique shapes like candelabra or globe forms. The icons are scattered across the page in a non-uniform but dense pattern.

Mebak Chic

USER MANUAL



www.mebakshop.com
amazon@mekabshop.com



TABLE OF CONTENT

WARNINGS	1
PRODUCT SPECIFICATION	2
PRODUCT DIAGRAM	2
OPERATING INSTRUCTIONS	3
CHARGING	4
MAINTAINANCE	4
WARRANTY	4
DISCLAIMER	5
QUESTIONS & ANSWER	5

Thank you for purchasing Mebak Percussive Massage Gun. This instruction manual contains important information about safety, usage and warranty. Use the product as described and keep this manual for future reference.

WARNINGS

- Do not expose the massage gun to sun or fire. Please keep it at room temperature.
- Do not immerse in water or place any objects into the ventilation holes.
- For adult use only. Children shall not use and play with the appliance. Do not use it if you got injured.
- Carefully examine the device before each use. Use only on dry, clean surface of the body over clothing by lightly pressing and moving for approximately 20 seconds per region.
- Bruising may occur regardless of the control setting or pressure applied. Use the device only on muscle soft tissue. Check treated areas frequently and immediately stop at the first sign of pain or discomfort.
- Do not use the Mebak massage device on your mouth, heart or near your genitals. Do not use the Mebak massage device above your Adam's Apple or spine, or any hard or bony area of the body.
- Do not use the Mebak Chic continuously for more than 1 hour. Please ensure the device rest for 5 minutes before re-using.
- Do not use any massage heads other than the ones provided by Mebak and only as instructed by Mebak. Use only the applicators that provide the best desired effect.
- Do not disassemble or tamper the device in any way, as this can increase the risk of electric shock, explosion, and exposure to hazardous battery chemicals.
- Never leave the Mebak operating or charging unattended.
- Keep the power off when storing and carrying.
- To protect the environment, please dispose of this product in accordance with local laws and regulations.
- Keep fingers, hair or any other body part away from the shaft and backside of the applicator as pinching may occur.

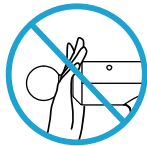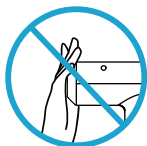

PLEASE DO NOT USE THE MEBAK, OR ANY VIBRATION DEVICE WITHOUT FIRST OBTAINING APPROVAL FROM YOUR DOCTOR IF ANY OF THE FOLLOWING APPLY:

Pregnancy, diabetes with complications such as neuropathy or retinal damage, wear of pacemakers, recent surgery, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylosis, recent joint replacements or IUD's, metal pins or plates or any concerns about your physical health. Frail individuals and children should be accompanied by an adult when using any vibration device.

These contraindications do not mean that you are not able to use a vibration or exercise device but we advise you to consult a doctor first.




PRODUCT SPECIFICATION

Adjustable Speed	1-6 Speed Levels
Amplitude	8 mm
Percussion Frequency	1200-2800 RPM/Min
Battery	2000mAh
DC Input	5V 2A
Item Net Weight	0.89kg (1.96lbs)
Item Size	17 * 6.2 *21 cm (6.7 *2.4 * 8.3 inch)


Package Include

- Mebak Chic Massage Gun × 1
- Massage Head × 8
- Portable Travel Case × 1
- Charging Cable × 1 (adapter not included)
- User Manual × 1

PRODUCT DIAGRAM

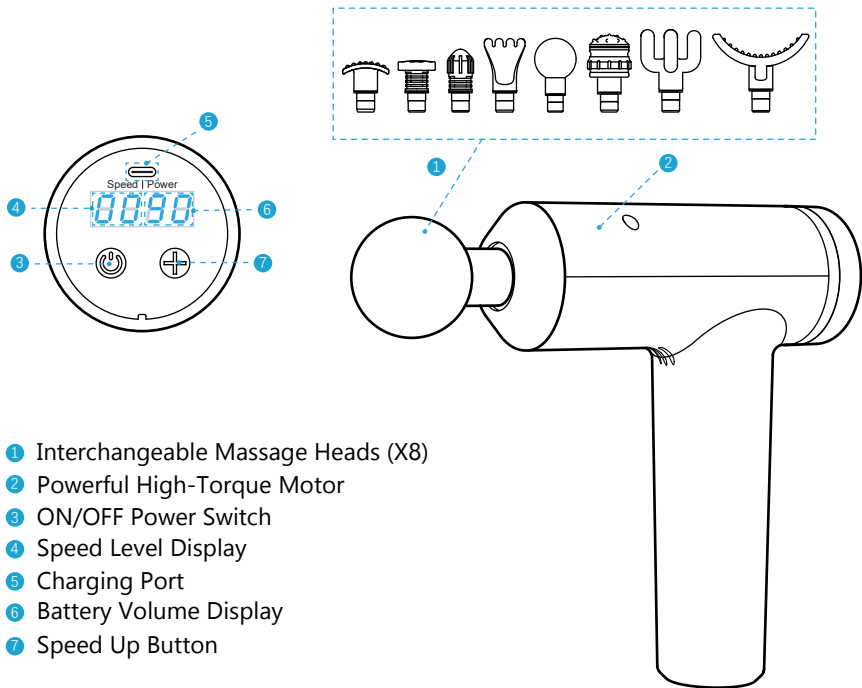

Massage Head

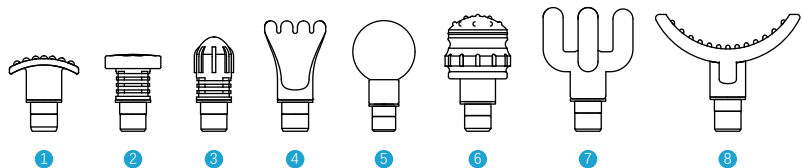

①	Bump	Best used to relax "flat area" muscle groups. Its multiple points are designed to help hypersensitive areas, such as the soles of the feet, with a less intense treatment.
②	Flat	Suitable for large and thick muscle groups, such as hip and leg area. Great for full body massage and relaxation. It's also effective for smaller and denser muscle groups.
③	Bullet	Used for pinpoint muscle treatment and smaller muscle groups like hands and feet, it allows you to dig deep into the tissues. Since it's used to massage specific areas, it shouldn't be used for longer than 15 seconds at a time.
④	Comb	Good for hard to reach areas, like shoulder blades or pelvis. Narrow, comb shape makes it great to deliver targeted relief.
⑤	Ball	Predominantly used for large muscle groups, like glutes, hamstrings, and quads. It covers a larger surface area, enabling you to target the large muscle groups effectively.
⑥	Cushion	Used mostly for tender or sensitive areas, but is also great for use all over the body. Being the softest and most flexible of all the massage heads, it is the safest to use close to bones or delicate areas. This also means that it could be very useful for non-flat parts of the body, as it adapts to the targeted surfaces better than any other head.
⑦	Claw	Comparing with normal U-shaped massage head, it has larger working area and more direct effect. At the same time, the radian of the four claws makes the force generated by them be more flexible, especially suitable for the small muscles of shoulder, neck and leg.
⑧	Crescent	Suitable for stimulating circulation and relaxing all muscle groups. Perfect for using on the calves and forearms.

OPERATION INSTRUCTIONS

Step ①

According to your needs to choose a suitable massage head, rotate the massage head firmly into the massage gun.



Step ②

- Turn on the device by pressing the Power Switch  for about 2 seconds.
- Shortly press the Speed Button  once to turn the device to speed level 1. Touch the " " again to enter the speed into a loop of 02-03-04-05-06-01. The corresponding speed will be displayed on the left side of the screen.

Step ③

- During the operation, you could press the Power Switch  at any speed level to set the device to "Stand-by (Speed 00)" status. Switch off the device after use by pressing the Power Switch  for about 2 seconds again.
- In consideration of your safety, the device will automatically turn off if there is no operation within 3 minutes in the Stand-by status.
- When the Power Display pulsates, please stop using the device and charge the device. When the display screen shows "00 00" or "E | 00", it means the battery has ran out and the device needs to be charged in time.

CHARGING

- Fully charge the device before first use. Please use the original charging cable and a 5V/2A output adapter to charge the device.
- When charging, the Power Display will pulsate to show the battery level and indicate the charging status. When the Power Display stops pulsating and stays at 99, it indicates that the battery is fully charged. To prolong the battery life, we recommend you to charge the device when the Power Display shows "20".
- The remaining battery volume will display in 99, 80, 60, 40, 20.
- The average run-time is 2-3 hours depending on the speed level and pressure applied during use. And the average charging time is about 4 hours. The charging time could be prolonged if the current output of the adapter is lower than 2A.
- Please be noticed that the device cannot work during charging.
- Please fully charge the device if it will not be used for a long time.

MAINTAINANCE

- Please do not use water to clean up the device, USB charging cable or the massage heads.
- Do not use gasoline, thinner, alkaline detergent, bleach, alcohol disinfectant when maintaining the device, as these may cause discoloration, deformation or cracking.
- Please pull out the USB charging cable when cleaning.
- Please do not plug in or pull out the charging cable with wet hands.
- If the surface is dirty, please soak a soft cloth in a neutral cleaner diluted with warm water. Wring it out and carefully wipe away the dirt.
- A small amount of grease may seep out near the massage head socket. Please don't worry. It is the grease of the internal part, not a defect of the device. You could wipe it away by a soft cloth.

WARRANTY

IF WITHIN SIX MONTHS FROM THE DATE OF PURCHASE, THIS PRODUCT FAILS DUE TO A DEFECT IN MATERIAL OR WORKMANSHIP, MEBAK WILL REPAIR OR REPLACE THE PRODUCT, OR NECESSARY COMPONENTS, FREE OF CHARGE. THIS WARRANTY EXCLUDES:

- (a) Damage caused by accident, abuse, mishandling, or transport;
- (b) Units subjected to unauthorized repair;
- (c) Units not used in accordance with Mebak instructions;
- (d) Damage exceeding the cost of the product;
- (e) Deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises, and
- (f) Failure to provide the dated proof of purchase.

Some states and countries do not allow a limitation of damages, so the foregoing limitation may not apply to you.

If you have any problems with the product, please contact us at this email address:

amazon@mebakshop.com

Or contact us through our website: **www.mebakshop.com**

This warranty guarantees specific legal rights, and other rights may vary from country to country and from state to state. This warranty service is available through the operating policies and procedures of Mebak and all Mebak international distributors. If you are a client, please contact the appropriate Mebak business entity, dealer or retailer, from whom you directly purchased the product for warranty and return authorization procedures. This warranty is initiated and executed by Mebak.

DISCLAIMER

To reduce the risk of injury, please read the instruction manual provided with the Mebak massage gun.

IMPORTANT SAFETY INSTRUCTIONS:

BEFORE USING THE DEVICE OR CHARGING THE LITHIUM-ION BATTERY, READ ALL INSTRUCTIONS AND CAUTIONARY MARKINGS IN THIS MANUAL, ON THE CHARGER, ON THE BATTERY, AND ON THE MESSAGE GUN.

QUESTIONS & ANSWERS

Questions		Reasons	Trouble shooting
Cannot be switched on	The display screen shows "00 00" or "E 00".	The battery has ran out.	Charge the device immediately.
	The display screen won't light up after charging.	Battery malfunction.	Please contact us for replacement.
	The display screen could light up but the device doesn't work.	The control panel is defective.	Please contact us for replacement.
Cannot be switched off		The device is in charging.	The device cannot be turned on/off or work during charging.
		Power Switch button press time is too short.	Keep pressing the Power Switch button for at least 2 seconds.
		The Power Switch button is defective.	Please contact us for replacement.
Unable to charge		The adapter or charging cable does not fit.	Please make sure that you are using the original charging cable and a 5V/2A adapter.
		The charging interface is defective.	Please use other 5V/2A output type-C charging cable to charge or contact us for replacement.
		Power cord is not plugged in tightly.	Make sure the power cord is plugged in place.
Stops working during massage		Touch the Power Switch button accidentally.	Turn on the device again.
		The pressure force is over stall force, the device turn off.	Turn on the device again.
		The device automatically turns off as it is in the Stand-by status for over 3 minutes.	Turn on the device again.
Abnormal noise		The massage head is not inserted in tightly.	Make sure the massage head is inserted in place.
		Shake the device and hear there are loose accessories inside.	Please contact us for replacement.

Mebak

Support E-mail: amazon@mebakshop.com



Name: Sea&Mew Accounting Ltd
Address: Electric Avenue Vision 25, London, Enfield EN3 7GD
Contact Person: Billy Han
Tel: +44 7399648608
E-mail: info@seamew.net