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AVITA CORPORATION
No.200.205 Provincial
Free way Song Ling Town,
Wu Jiang City,
Jiang Su, P.R.C. ZIP:215200

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4 in 1 infrared thermometer

Dear Customer

You have just bought this Terrillon infrared electronic thermometer and we thank you for your purchase. We hope you enjoy using the product. In order to be completely satisfied, we recommend you read these instructions carefully.

The infrared probe of the Thermo Color shows in a few seconds the heat emitted by the eardrum (this reflects the temperature of the nervous system located in the brain, which regulates body temperature) or emitted at the surface of the forehead. Being practical, it interprets the measurements in colour. It has 10 memories for monitoring temperature changes. In standby mode on its base stand, it measures the room temperature and displays the time.

READ BEFORE USE

Important notes


1. The temperature varies according to the part of the body where it is taken. Therefore, it is not necessary to compare the oral, underarm, rectal, ear or forehead temperature results.
2. The temperature also varies from one person to another, according to age and the time of day.

Advice


Determine your normal temperature range by taking your temperature when healthy. These readings can be used as a form of reference to compare with future readings when ill.

Is your temperature abnormally high or abnormally low compared to normal? Do not hesitate to consult a doctor in case of doubt.

PRODUCT DESCRIPTION

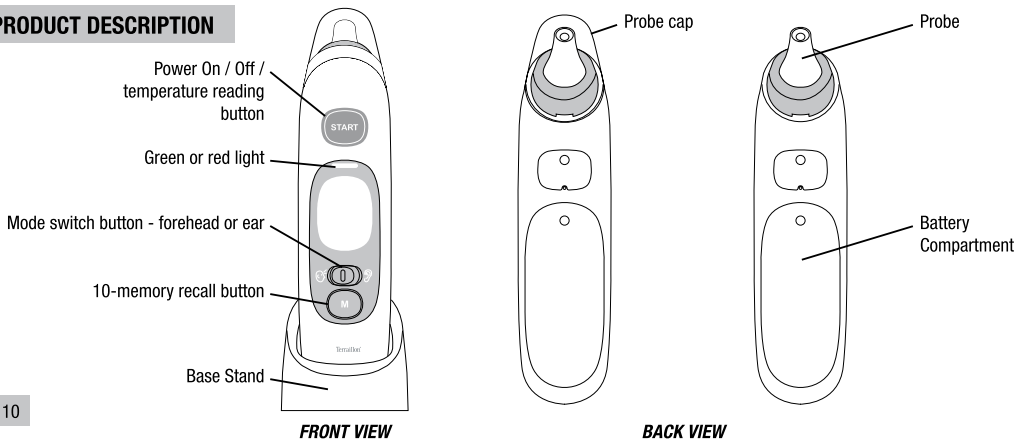VIEW INSIDE THE BATTERY COMPARTMENT

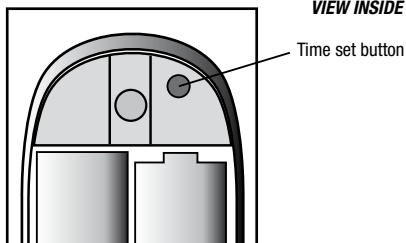FIRST USE

1. Inserting batteries

The unit works with 2 AAA batteries (not included).

Open the battery cover on the back of the Thermo Color using a screwdriver. Insert the batteries, matching the correct polarity (the polarities are shown on the inside of the battery compartment). If you do not intend to use this unit for a prolonged period of time, it is advisable to remove the batteries before storage.

2. Setting the clock

The first time that you use this thermometer or after having changed the batteries, you should set the clock.

- Open the battery compartment using a screwdriver.
- Press the time set button. The hour digits will flash.
- Press the Power On/Off «START» button on the thermometer in order to adjust them.
- Confirm by pressing the time set button. The minutes digits will flash.
- Press the Power On/Off «START» button on the thermometer in order to adjust them.
- Confirm by pressing the time set button.

If the time is set, the thermometer will go into standby mode after 30 seconds and the screen will display the room temperature and the time alternately.

If the time is not set, the thermometer will turn off automatically after around 30 seconds.

3. How to choose the temperature measurement unit (F° or C°). (p. 15)

TAKING A TEMPERATURE

Remove the cap from the probe to take a temperature, in forehead or ear mode.

1. Forehead temperature measurement







- Remove the cap from the probe.
- Switch to forehead mode (Fig. 1).
- Press the temperature reading «START» button.
- When first used to take a temperature, all the elements on the screen (Fig. 3) are displayed, followed by «--» (Fig. 4), a beep and «36,8°C» by default (Fig. 5).
- From its second use, all the elements on the screen are displayed (Fig. 3), followed by «--» (Fig. 4), a beep and the last temperature in the memory.
- Place the probe on the top of the temple.- Press and hold the temperature reading «START» button for approx. 2 seconds.
- You can release the «START» button and remove the Thermo Color from the temple. The temperature is displayed as well as a green or red colour depending on the result.


The forehead temperature may be affected by external climatic conditions so, in case of doubt, do not hesitate to take the ear temperature.

2. Auricular

- Remove the cap from the probe.
- Set the switch to ear mode (Fig. 2).
- Press the temperature reading «START» button.
- When first used to take a temperature, all the elements on the screen (Fig. 3) are displayed, followed by «--» (Fig. 4), a beep and «36,8°C» by default.
- From its second use, all the elements on the screen are displayed (Fig. 3), followed by «--» (Fig. 4), a beep and the last temperature in the memory (Fig. 6).
- Insert the probe into the ear.
- Press and hold the temperature reading «START» button.
- When the unit beeps, you can release the «START» button.
- Remove the probe from the ear. The temperature is displayed as well as a green or red colour depending on the result.

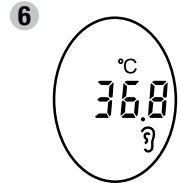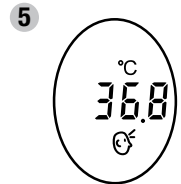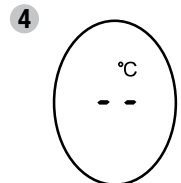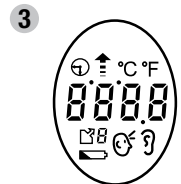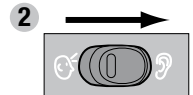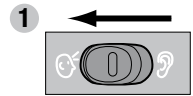Note: In order to obtain a precise measurement, ensure the correct positioning of the probe in the ear so that the probe tip senses the heat from the eardrum. In order to do this, you can pull the ear to straighten the ear canal and slightly rotate the probe inside. (It is not recommended to use the ear function of the thermometer on a baby under 12 months old.)

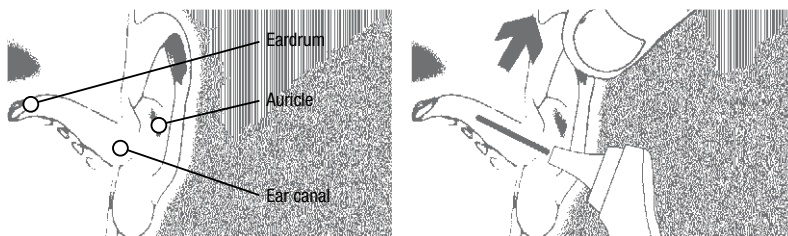Easy interpretation of the results

In both modes, if the temperature is lower or equal to 37.5°C, the green light comes on for around 5 seconds. If the temperature is higher than 37.5°C, the red light comes on.


- Temperature higher than 37.5°C - high temperature: the red light comes on
- Temperature lower or equal to 37.5°C - normal temperature: the green light comes on

MEMORY

To display the last memorised result, press the memory recall button «M». The last temperature reading will appear along with the memory number.

To display the previous memorised results, press the memory recall button «M» again and the memory numbers will be displayed.


MEMORY NUMBERS

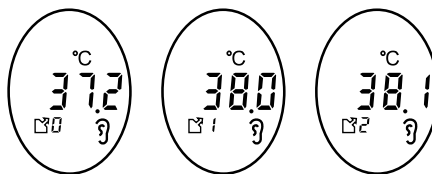REPLACING THE BATTERIES

When the symbol «» appears on the screen, this means that the batteries are nearing the end of their life.

When the screen displays the symbol «» and «Lo», this means that the battery is exhausted and should be replaced.

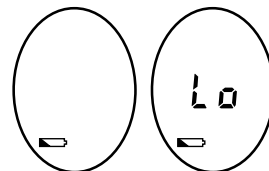To change the batteries of the Thermo Color:

- Open the battery compartment cover situated on the back of the unit using a screwdriver.
- Remove the exhausted batteries.
- Replace them with 2 AAA batteries, matching the correct polarity (the polarities are shown on the inside of the battery compartment).
- Replace the cover. The unit is ready to use.

ENVIRONMENTAL PROTECTION

Deposit exhausted batteries in a container designed for this purpose so they may be collected and recycled. Do not mix different battery types. Do not mix exhausted batteries with new batteries.

At the end of the product's serviceable life, take it to an electrical and electronic equipment recycling collection point.

DISPLAY SYMBOLS

Temperature measurement in progress

Ear mode

Forehead mode

Green light - temperature equal to or lower than 37.5°C

Red light - temperature higher than 37.6°C

Clock

Memory: result of the last temperature reading.

Degrees Celsius

Degrees Fahrenheit

This symbol appears when the battery begins to wear out. You can continue to use the unit but when you see the «Lo» symbol, the batteries should be changed because the unit can no longer function properly.

H_i The temperature is higher than 50°C

Lo The temperature is lower than 10°C

Err Temperature outside of operating conditions (16-40°C)

STOP

The thermometer automatically goes into standby mode after 30 seconds if it is not touched.

MAINTENANCE

Clean the thermometer probe with alcohol-soaked cotton wool before and after each use. Do not use abrasive substances. Do not immerse the thermometer in water or in any liquid, because only the probe is waterproof.

SAFETY PRECAUTIONS

1. Use of this thermometer is for monitoring purposes only and is not intended to provide diagnosis or treatment. Make sure that any abnormally high or low temperature readings are communicated to your doctor. Under no circumstances should you change your prescribed medicine doses.
2. Do not drop or break the thermometer.
3. Keep out of reach of children.
4. Keep thermometer away from direct sunlight and store well away from any source of heat in a clean, dry location.
5. Do not submerge the thermometer in water or any other liquid.
6. The Thermo Color only gives the forehead, ear and room temperature. Do not use it to take rectal, oral or underarm temperatures.
7. Avoid direct finger contact with the tip of the probe.
8. Do not take temperature within 30 minutes of eating, exercising or taking a bath.
9. Wait one minute between each temperature reading. Two consecutive temperature readings may give slightly different measurements. Should this occur, take the average of the two.

How to choose temperature units of measurement (Fahrenheit F° or Celsius C°)

When the unit is off, press the Power On/Off «START» button to turn it on. All the elements on the screen are displayed, followed by «--». Quickly press the Power On/Off «START» button (before the most recent measurement is displayed) in order to switch from one measuring unit to the other.

MANUFACTURER

AVITA CORPORATION
No. 200.205 Provincial
Free Way Song Ling Town,
Wu Jiang City,
Jiang Su, P.R.C. ZIP: 215200

Type BF

CE 0197

TECHNICAL DATA

Operating temperature range: 10°C - 50°C
Tolerance: Between 36 and 39°C: +/-0.2°C
Below 36°C: +/- 0.3°C
Above 39°C: +/- 0.3°C
Increments: 0.1°C or °F
Operating conditions: 16°C to 40°C (60.8°F~104°F), 15-95% RH
Storage conditions: -20 to 50°C 15-95% RH
Battery: 2 x AAA (not supplied)

WARRANTY

This product is guaranteed for 2 years against material and manufacturing defects. During this period, such defects will be repaired free of charge (proof of purchase must be presented if claim is under warranty). This warranty does not cover damage resulting from accidents, improper use or negligence. In the event of a claim, please contact the store where you bought the item first.

Applied Standard.

- IEC/EN 60601-1 (Safety)
- IEC/EN 60601-1-2 (EMC)

Keep these instructions in a safe place for future consultation.

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Terraillon SAS France & Headquarters
Terraillon UK Ltd
Terraillon Asia Pacific Ltd
Terraillon Corp USA