


PRESTIGE BRANDS INC
660 White Plains Rd. Tarrytown, New York 10591
914-524-6834 artdept@prestigebrands.com

DENTEK

PRODUCT NAME: ULTIMATE GUARD
COMPONENT: INSERT Back
SIZE: 1 CT
COUNTRY: USA
JOB #: DKUS038603
TECH DRAWING #:
UPC#:

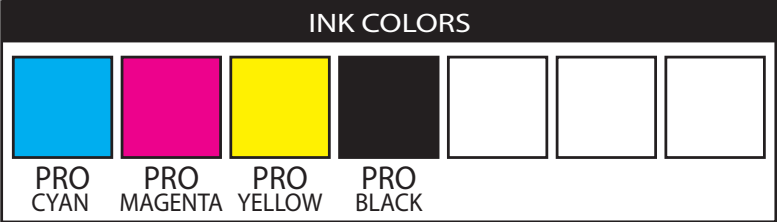

color mechanical	date	initials	comments
Graphics			
Marketing			
Packaging Dev			
Legal			
Regulatory			



CLEANING & CARE

After each use, brush your dental guard with toothpaste and/or mouthwash. Or, clean your guard weekly with Fresh Guard® Soak, specially formulated to clean dental devices. Never rinse or soak in hot water as it could alter the shape of your dental guard.

INDICATIONS FOR USE

The DenTek Ultimate® Dental Guard is indicated for use for protection against bruxism or nighttime teeth grinding. It is intended to reduce damage to the teeth and reduce the noise associated with bruxing or grinding.

WARNINGS

DO NOT USE:

- As a sports mouthpiece. The product is designed to absorb top-to-bottom pressure not shock absorption.
- If you are under 18 years of age.
- For more than three months without consulting your dentist.
- If you have loose teeth or teeth that can wiggle.
- If your dentist has told you that you have temporomandibular joint (TMJ) pain or clicking of the jaw.
- If you have any pain from teeth grinding, or other forms of tooth or jaw pain.
- If you have dental implants, or wear dentures, braces or other dental appliances.

ASK A DENTIST BEFORE USE IF YOU HAVE:

- Major health problems or serious breathing or respiratory conditions.
- Oral sores, bleeding gums or any gum disease.
- Cavities that have no fillings, or loose fillings or caps.
- Difficulty chewing or pain of the jaw, teeth or face.
- Two or more missing teeth.

STOP USE AND ASK A DENTIST IF:

- You develop soreness, bleeding gums or any other reaction inside your mouth.
- You develop jaw pain, teeth pain, ear pain, headache, neck stiffness, or joint clicking when using this product.
- The same symptoms persist after over a month of use.
- The product falls out of your mouth easily.
- The product disrupts breathing, feels uncomfortable or causes you to gag.
- You have a change in your bite that persists more than a few minutes after removing the product.

Wearers and users may experience choking, gagging, or swallowing while using this product. If this occurs, discontinue use and seek medical attention. See your dentist every six months.

Our dental guard only contains materials that have been proven to be safe and effective. They are BPA-free, latex-free and phthalate-free. The materials have been tested for biocompatibility and they do not contain any hazardous or toxic chemicals.

REFUNDS & REPLACEMENTS

REPLACING YOUR DENTAL GUARD

All dental guards should be replaced over time. Please check your dental guard monthly for wear. The life of your DenTek® dental guard for maximum protection depends on how often you grind and the force of your teeth grinding.

30 DAY MONEY-BACK GUARANTEE*

If you are not completely satisfied with the fit and comfort of the dental guard, we will gladly give you a full refund (not including taxes or postage costs) if you call us within 30 days of the purchase date.



6 Month Protection Guarantee*

From date of purchase, we guarantee that this dental guard will give you 6 months of teeth grinding protection. If not, we will send you one free dental guard as a replacement (one per purchase). For a refund or replacement details, call 1-800-433-6835. Please allow 3 to 4 weeks for delivery of your refund / replacement.

*Only one guarantee per 12-month period. This guarantee is only valid for products purchased in the U.S. by customers living in the U.S.

800-4DENTEK (433-6835)
www.DenTek.com



Mfd. for DenTek Oral Care, Inc., Tarrytown, NY 10591
A Prestige Consumer Healthcare company.
Patent Pending. Made in USA.
©2019 Trade dress is owned by DenTek Oral Care, Inc.
All Rights Reserved. DKUS038603