## FEATURES

- 1) 2 push buttons to control bubbles/heat and massage
- 2) full platform generates aerating bubbles
- 3) 3 attachments
  - dual node
  - single node
  - pinpoint
- 4) translucent splash cover
- 5) nonslip feet
- 6) line cord wrap

## RECOMMENDED MASSAGING TIMES

Generally, massage should not last longer than 10 to 15 minutes. Treatment may be applied once or twice daily as desired. Allow unit to cool before using again. As a person's metabolism varies, so does his or her response to massage. You will soon adjust your massage time to suit your needs and personal comfort. If you have any concerns regarding use, we recommend you consult your physician.

## LIMITED ONE YEAR WARRANTY (U.S. AND CANADA ONLY)

Conair will repair or replace (at our option) your unit free of charge for 12 months from the date of purchase if the appliance is defective in workmanship or materials.

To obtain service under this warranty, return the defective product to the service center listed below, together with your purchase receipt and \$12.00 for postage and handling. California residents need only provide proof of purchase and should call 1-800-3-CONAIR for shipping instructions. In the absence of a purchase receipt, the warranty period will be 12 months from the date of manufacture.

**ANY IMPLIED WARRANTIES, OBLIGATIONS OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED IN DURATION TO THE 12-MONTH DURATION OF THIS WRITTEN, LIMITED WARRANTY.** Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.

**IN NO EVENT SHALL CONAIR BE LIABLE FOR ANY SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS OR IMPLIED, WHATSOEVER.** Some states do not allow the exclusion or limitation of special, incidental or consequential damages, so the above limitation may not apply to you.

This warranty gives you specific legal rights, and you may also have other rights, which vary from state to state.

**Service Center**  
Conair Corporation  
Service Department  
7475 North Glen Harbor Blvd.  
Glendale, AZ 85307

Please register this product at: [www.conair.com/registration](http://www.conair.com/registration)

Visit us on the web at: [www.conair.com](http://www.conair.com)  
©2016 Conair Corporation  
East Windsor, NJ 08520/Glendale, AZ 85307

IB-10298A

FB30  
16PS016875

# CONAIR®

## relaxing spa foot bath

### INSTRUCTION BOOKLET FOR ALL FB30 MODELS

For your safety and continued enjoyment of this product, always read the instruction book carefully before using.

## IMPORTANT SAFETY INSTRUCTIONS

When using an electrical appliance, basic safety precautions should always be taken, including the following:

## READ ALL INSTRUCTIONS BEFORE USING DANGER

To reduce the risk of electric shock:


- 1) **Always unplug this appliance from the electrical outlet immediately after using and before cleaning.**
- 2) Do not reach for an appliance that has fallen into water. **Unplug immediately.**
- 3) Do not use while bathing or in a shower. Use only in a dry area.
- 4) Do not place or store appliance where it can fall or be pulled into a tub or sink. Do not place or drop into water or other liquids.

## WARNING

To reduce the risk of burns, fire, electric shock, or injury to persons:

- 1) **An appliance should never be left unattended when plugged in. unplug from outlet when not in use, and before putting on or taking off parts.**
- 2) This massager is not a toy!
- 3) Consult your doctor before using in case of pregnancy, illness, medical condition, or persistent pain.
- 4) Do not use on any unexplained pain or swollen muscles, or following a serious injury, before consulting your physician.
- 5) Do not use massager in close proximity to loose clothing or jewelry.
- 6) Keep long hair away from massager while in use.
- 7) Do not operate under a blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 8) This appliance should not be used by, on, or near children or individuals with certain disabilities.
- 9) Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by Conair.
- 10) Never operate this appliance if it has a damaged cord or plug, if it is not working properly, or if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination and repair.
- 11) Do not carry this appliance by the power supply cord or use the cord as a handle.
- 12) Keep the cord away from heated surfaces.
- 13) Never operate on a soft surface, such as a bed or couch, where air openings may be blocked.

- 14) Never drop or insert any object into the product.
- 15) Do not use outdoors.
- 16) Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- 17) Ensure that your hands are dry when operating the controls or removing the plug.
- 18) To disconnect, turn all controls to the OFF position, then remove plug from outlet.
- 19) This appliance is for household use only.
- 20) Do not stand on or in the foot bath. Use only while seated. Place feet in bath and remove feet from bath while in a seated position only.
- 21) Never use while sleeping or drowsy.
- 22) Keep out of the reach of children.
- 23) Use heated surface carefully. May cause serious burns. Do not use on sensitive skin or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
- 24) Burns may result from improper use.
- 25) To prevent scalding, never add boiling or steaming water to the foot bath.
- 26) Temperatures sufficiently high to cause burns may occur regardless of the control setting. Do not use on an infant, individuals with certain disabilities, or on a sleeping or unconscious person. Check the skin in contact with the heated area of the appliance frequently to reduce the risk of injury.

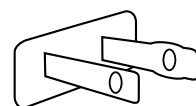
This unit has a two-wire, double-insulated polarized plug.

## OPERATING INSTRUCTIONS

To reduce the risk of electric shock, this appliance has a polarized plug (one blade is wider than the other). This plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician to install the proper outlet. Do not change the plug in any way.

## CAUTION

- 1) **Consult your doctor before using in case of illness or medical condition, including but not limited to pacemaker use, pregnancy, cancer, infection, fracture, or persistent pain.**
- 2) Do not use if you have any of the following conditions: diabetes, tuberculosis, benign or malignant tumors, phlebitis or thrombosis, hemorrhages, open or fresh wounds, ulcerated sores, varicose veins, poor circulation or bruised, discolored, burned, broken, swollen, or inflamed skin or heat sensitivity.
- 3) This product should not be used by any individual suffering from a condition that limits or altogether inhibits the user's capacity to feel or have sensation in any part of the body.
- 4) If you experience pain in a muscle or joint for a prolonged period of time, discontinue use and consult your doctor. Persistent pain could be a symptom of a more serious condition.
- 5) Use of this product should be pleasant and comfortable. Should pain or discomfort result, discontinue use and consult your physician.
- 6) If the user empties water contrary to instructions, a small amount of water may issue from the drain hole. This is normal, as it is a safety measure to prevent water ingress. Wipe floor with dry cloth after use to avoid damage.
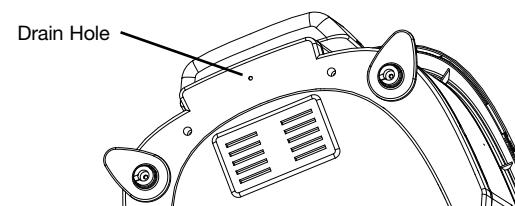
## USER/MAINTENANCE INSTRUCTIONS

- 1) Make sure foot bath is unplugged and is in the OFF position.
- 2) When transporting foot bath, position hands on opposite sides of foot bath for most secure grip.
- 3) Place foot bath on floor, fill foot bath with warm or cold water, and plug into wall outlet. Never plug in or unplug the foot bath while your feet are in water.
- 4) The heat setting is designed to maintain the temperature of warm water, not to heat water.
- 5) Sit in a chair and place your feet in the water. Never stand in foot bath.
- 6) Select foot bath feature by pushing ON/OFF button for bubbles and heat.
- 7) Choose one of the 3 attachments and place on post at the top of the foot bath. Select massaging by pushing ON/OFF for massaging. Place foot on attachment and apply pressure. Move your foot around for more pressure and treatment.

- 8) To clean the unit, make sure the massager is unplugged. Allow the unit to cool, and then wipe down with a damp cloth using soap and water only. Do not use strong chemical cleaners.
- 9) Fill foot bath with warm or cold water to the level indicator MAX. Do not fill above water-level indicator. Never overfill.
- 10) Pour water out of the foot bath from the side opposite the line cord only.
- 11) Do not hang unit from power cord.
- 12) Store unit in its box or in a clean, dry place. Wrap power cord loosely around tub feet.

## SAVE THESE INSTRUCTIONS

### SERVICING OF DOUBLE-INSULATED PRODUCTS

In a double-insulated product, two systems of insulation are provided instead of grounding. No grounding means is provided on a double-insulated product, nor should a means for grounding be added to the product.

Servicing a double-insulated product requires extreme care and knowledge of the system, and should be done only by qualified service personnel. Replacement parts for a double-insulated product must be identical to the parts they replace. A double-insulated product is marked with the words "DOUBLE INSULATION" or "DOUBLE INSULATED." This symbol (square within a square) may also be marked on the product.

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