

How to use the epilator.

Epilation head

The epilation head for legs/arms should be used only for legs/arms. The epilation head for underarms/bikini-line should be used only for underarms/bikini-line.


Areas suitable for epilation

Areas the epilation head for legs/arms can be used.



For models ES-EL8A, ES-EL7A, ES-EL9A

Areas the epilation head for underarms/bikini-line can be used. Before epilating the bikini-line, be sure to put on your swim wear or underwear etc. and check the areas you wish to epilate.

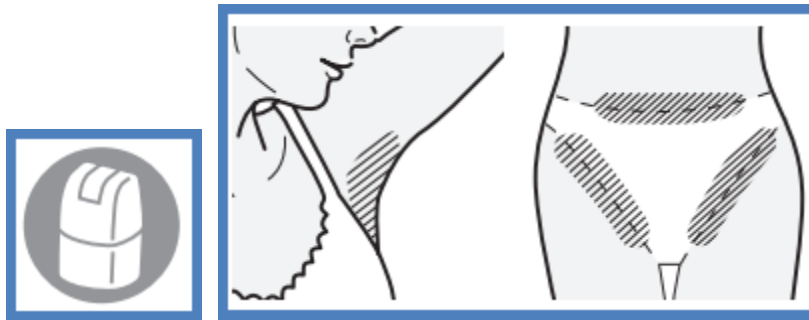

- Using the epilator on areas other than those marked areas may cause pain or skin trouble.

Areas not suitable for epilation

The epilator should not be used on the inner side of the upper arms and areas such as elbows or knees which are prone to sagging.



Wet Epilation

WET epilation (epilating after wetting the skin and epilator and then applying foam) makes the skin softer so it is gentler to your skin when epilating. For DRY epilation, begin at step 3. Before DRY epilation, wipe dry any moisture or sweat from the skin.

1. **Wet your skin.**
2. **Wet the discs and place a small quantity of liquid shower gel on the discs.**
 - Always use shower gel when wet epilating.
 - Do not use a shower gel with a scrub or with granular salt, this may cause a malfunction.
 - Rinse any shower gel out of your hands.

3. **Select the mode by pressing the power switch []**

- Foam is created. Foam makes the epilator slide better, so you can move it fast.
 - The appliance may not operate in an ambient temperature lower than approximately 5 °C (41 °F).
 - Apply the epilator at 90° to your skin and move it slowly so that it slides over your skin.
-

Epilating legs or arms


Upwards from the bottom of the leg. Inwards from the outside of the arm.

- For models ES-EL7A, ES-EL9A - Attach the fast cap for legs/arms or the gentle cap for legs/arms.



Epilating underarms or bikini line using ES-EL8A

Move in several directions as underarm/bikini line hair grows in different directions. Stretch the skin with your fingers, and move slowly from inside to outside going against the flow of hair.

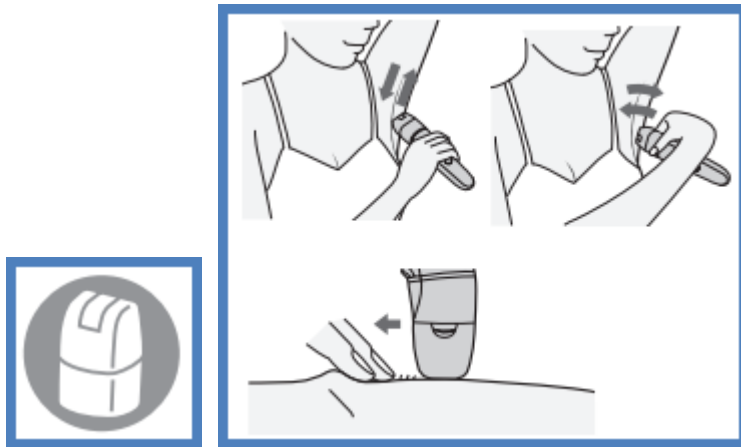

Epilating underarms or bikini line using ES-EL7A, ES-EL9A

- Attach the sensitive area cap for underarms/ bikini-line.
- Lightly apply the epilator at a 90° angle to your skin and be in close contact with your skin.
- Do not move the appliance back and forth while applying it to the skin. Remove it from the skin, change the direction, and then move.
- Efficiently stretch the skin with your fingers to prevent sagging of your skin, and move slowly going against the flow of hair.


How to charge the epilator.

For best performance, fully charge the appliance before use. You cannot operate the appliance while charging.

Charging time = Approx. 1 hour.

- The epilator can be used for approximately 30 minutes (40 minutes when using the gentle cap) after 1 hour of charging. (This will differ depending on the temperature.)
 - For models ES-ED70, ES-ED50 - 40 minutes when using the gentle cap.
 - For models ES-EL7A, ES-EL9A - 35 minutes when using the gentle cap for legs/arms.

Plug the AC adaptor into the epilator (a) and a household outlet (b).

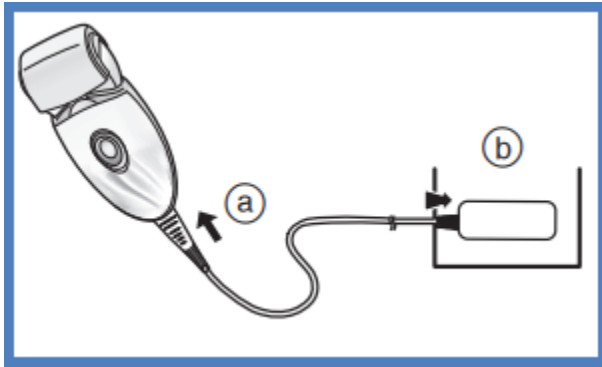

Model shown is ES-ED70.

- Recommended ambient temperature for charging is 15 - 35 °C (59 - 95 °F). The battery may not charge properly or not at all under extreme low or high temperatures.
- While charging, red light appears.
- After charging is completed, red light goes off.
- Abnormal charging, the light blinks twice every second.
- The charging time may decrease depending on the charging capacity.
- After charging is completed, turning the power switch ON with the epilator plugged in will cause the charging light to light up and then turn off after 5 seconds. This shows that the epilator is fully charged.
- To keep charging the battery will not affect battery performance.
- When charging the epilator for the first time or when it has not been in use for more than 6 months, the time it takes to charge the epilator may increase, or the charging light may not glow for a few minutes. It will eventually glow if kept connected.
- The AC adaptor is intended to be correctly oriented in a vertical or floor mount position.
- If the shape of the plug does not fit the power outlet, use an attachment plug adaptor of the proper configuration.

How to use the shaver head. - ES-EL7A, ES-EL9A

WET Shaving with soap lather makes the skin slippery for a closer shave. For DRY shaving, step 3 is unnecessary. Before DRY shaving, wipe dry any moisture or sweat from the skin.

1. **Attach the shaver head.**
 - Confirm the bikini comb is detached.
2. **Check that the outer foil is not deformed or damaged.**
3. **Wet your skin and put shower gel on your skin.**
 - Do not use a shower gel with a scrub or with granular salt. Also do not use shaving cream, skin cream or skin lotion. Doing so may clog the blade and cause a malfunction.

4. **Select the mode by pressing the power switch.**

- Press gently so that the whole blade is in close contact with the skin and the surface of the blade does not move up and down.


How to clean the epilator.

- Always clean the heads and caps after use to keep them in a hygienic state.
- Switch off and unplug the appliance before cleaning.
- Take care not to damage the skin protector, the frame, the fast cap and the gentle cap during cleaning.

Dry cleaning


Epilation head/shaver head

1. **Lift fast/gentle caps, the frame upward while holding the cap release ribs/the frame release ribs or remove the frame by pushing the frame release button.**

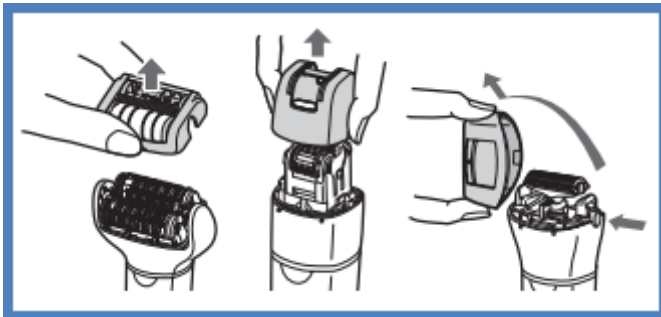

Model shown in photo is ES-EL3A.

2. **Clean fast/gentle caps, the epilation head and the shaver head with the cleaning brush.**

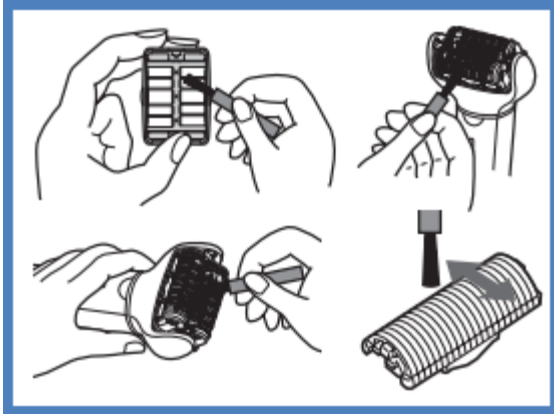


Model shown in photo is ES-EL3A.

Wet cleaning

Epilation head/Shaver head


1. **Wet the discs and blades, and then apply liquid hand soap.**
2. **Turn the switch on and then foam is created.**
3. **Wash the head with water to wash out hair.**
 - Do not use hot water.
 - If soap is not completely washed away, then a white deposit can accumulate and prevent the discs and the blades from moving smoothly.
 - Rinse the shaver head with the frame removed.
4. **Turn the switch off, wipe the appliance with a dry cloth and dry it well.**

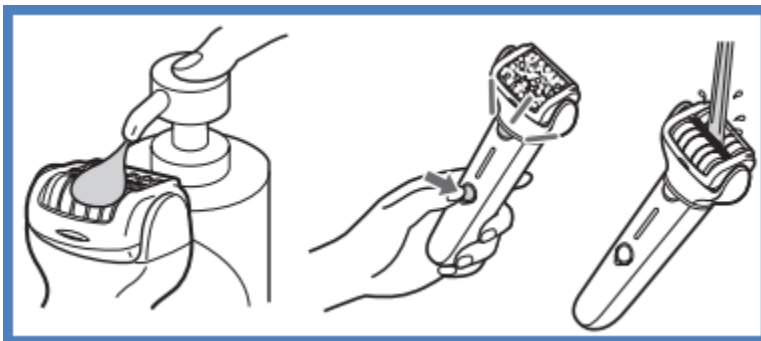

Model shown in photo is ES-EL3A.

Foot care head (ES-EL8A, ES-ED90, ES-EL7A/9A)

1. **Apply liquid hand soap to the central part of the file.**
2. **Pour a small amount of water on the file and turn the switch on.**
3. **Rinse any debris thoroughly with water or lukewarm water.**
4. **After turning the switch off, remove the frame and rinse inside.**
5. **Wipe with a dry cloth and dry in a well-ventilated location.**

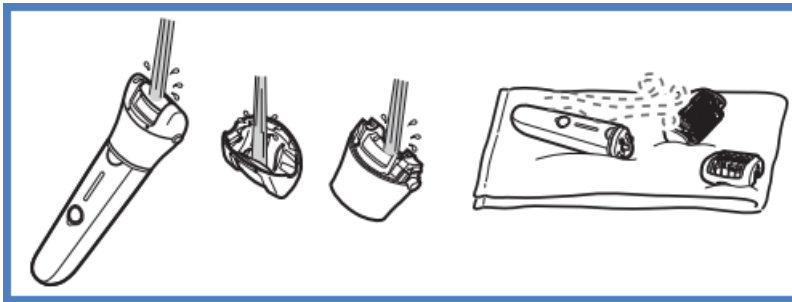


Model shown in photo is ES-EL3A.

How to replace the outer foil and inner blade.

We recommend replacing the outer foil every year and the inner blade every two years. Remove the outer foil of the shaver head only when replacing it.

Replacing the outer foil


1. **While pushing the outer foil gently with the fingers, use a fingernail to release the plastic panel (b) on the blade from the hooks (a) on the inside of the frame.**
2. **The new outer foil should be slightly bent and pushed in until it is hooked in the frame.**

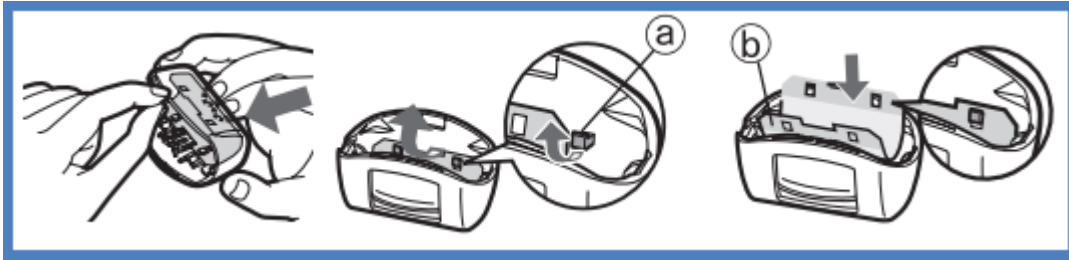

Model shown is ES-EL8A.

Replacing the inner blade (ES-EL8A)

1. Remove the inner blade by holding it at both ends and pulling straight away from the shaver.
2. Insert the inner blade by holding it at both ends and press downward.

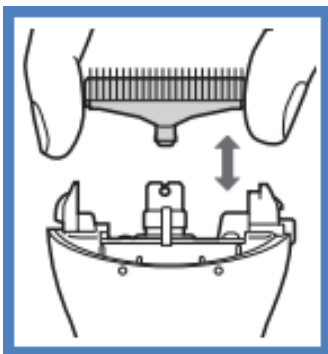


Model shown is ES-EL8A.

Replacing the file (ES-ED90/70/50, ES-EL7A/9A)

The file lifetime is approximately 1 year (assuming a single use every 2 weeks with approximately 10 minutes per use). We recommend replacing the file about once a year.

1. Press down with your finger and release the file.

2. Align the shaped part  with the mark side  and press down.

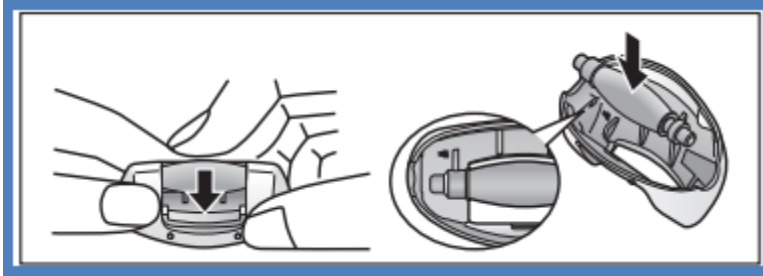

Model shown is ES-EL8A.

Replacing the sonic scrub brush (ES-EL9A)

Sonic scrub brush (Scrub brush)	once a year (assuming 3 minutes use per week)
---------------------------------	---