


RETAINER CLEANER CASE

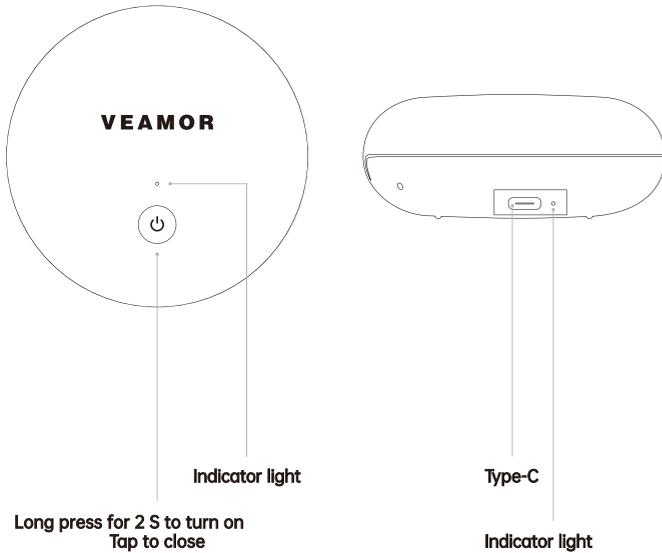

VEAMOR UV Retainer Case, Custom-Made Device for your Oral Health. VEAMOR has the function of UV sterilizing and storing retainers. In order to use it correctly and avoid damage, please refer to the following suggestions :

INSTRUCTION

First, open the lid of the box, put the retainer in and close the lid to check if it is completely closed. Press and hold the start button for 2 seconds and the white light will come on to indicate that UV sterilization is in progress. VEAMOR Will perform sterilization for up to 2 minutes. When the green light comes on and a set of sharp beeps sound continuously, the disinfection is complete

Note:

1. The device is not disinfected when it is turned off.
2. One minute after the disinfection is complete, the device will automatically shut down.
3. During the disinfection process, you can click the start button to turn off the machine.

CHARGING THE VEAMOR

Use the included Type-c charging cable to charge the VEAMOR. When the yellow light next to the Type-c port is on it indicates charging is in progress. The recommended charging time is 1-2 hours and the green light next to the Type-c port indicates that charging is complete. Can be used for about one month when fully charged (Depends on frequency and duration of use)

WARNING

1. VEAMOR retainer case IP64 waterproof. It's recommended to wipe the product case with wet paper. Please do not directly immerse the product in water. Avoid short circuit of the wire.
2. Please store under normal temperature, do not close to high temperature and humidity.
3. Please keep this product out of the reach of children.
4. Please keep your hands clean when using. When there is water on your hands, please wipe them clean. Avoid accidents.