
## Sonic Electric Toothbrush

### Operation Instructions

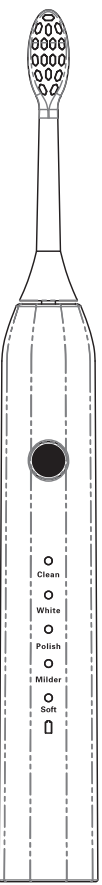Enjoy Your Life Every Day!

#### Operation Instructions (EN)

##### Important safety information

Before using the product, please read this important information carefully and keep it handy for future reference.

##### Danger

Keep the charger away from the water source. Do not place or store it above or near any container with water, such as bathtub, washbasin, pool, etc. Do not immerse the charger in water or any other liquid. After cleaning, make sure the charger is completely dry before connecting it to the power supply.

##### Warning

Do not replace the power cord. If the power cord is damaged, please discard the charger. You should always use the original-model charger for charging to avoid danger. If the product is damaged (brush head, toothbrush handle or charger), please stop using.

The charger of this product can not charge other brands of electric toothbrush or any product with inductive/wireless charging function, either.

If possible, avoid charging for more than 12 hours. Do not let children play with this product.

• Do not wash the brush head, toothbrush handle or charger in the dishwasher.

• If you have had oral cavity or gum surgery in the past 2 months, please consult your dentist before using this product.

• If the gum bleeds a lot after using this product or it will still bleed after 1 week, please consult your dentist.

• If you feel discomfort or pain when using this product, please also consult your dentist.

• For any questions related to pacemaker or other implantable device, please consult your doctor or the device manufacturer of the implantable device before use.

• If you have medical problems, please consult your doctor before using this product.


• This product is designed only for cleaning teeth and gum only. Do not use it for any other purpose.

• If you feel any discomfort or pain, please stop using this product and contact your dentist or head down hard to the metal axis until it will not go down anymore.

• Please stop using the brush head that deformed or with bent bristles. For bristles of intelligent and curing type, the bristles in the forward will become lighter with use when a brush head should be replaced, so it is recommended to replace the brush head every 3 months. If there are signs of wear, it should be replaced in advance. Do not use other brush heads that are not recommended by the manufacturer.

• If your toothpaste contains peroxide, sodium bicarbonate or other bicarbonates (a common ingredient in whitening toothpaste), be sure to thoroughly clean the brush head with soap and water after each use, which can prevent possible plaque buildup.

#### Parts

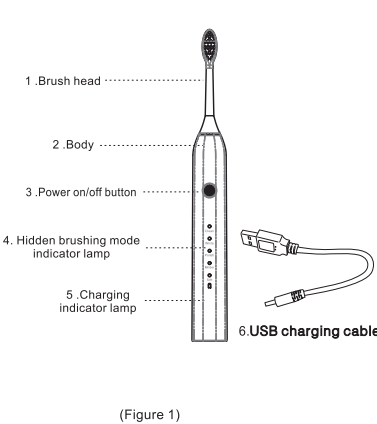(Figure 1)

##### Introduction

Intelligent sonic-wave electric toothbrush (Figure 1)

1. Brush head
2. Body
3. Power on/off button
4. Hidden brushing mode indicator lamp
5. Charging indicator lamp
6. USB charging cable

Note: Items in the packaging box may vary depending on the model purchased.

##### Product Instructions


|                          |                   |
|--------------------------|-------------------|
| Input voltage: DC 5V     | Rated power: 2W   |
| Product material: PC/ABS | Brush head: Nylon |

##### Preparation for use

Brush head installation steps:  
1. Align the brush heads so that the direction of the bristles is the same as the front of the toothbrush handle (Figure 2).

2. Press the brush head down hard to the metal axis until it will not go down anymore.

##### Instructions

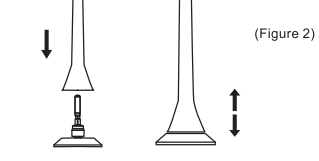
(Figure 2)

#### Charging the product

- Plug the charger into the power socket (Figure 3).
- Connect the body with the charging cable, and the charging indicator lamp will flash until it is fully charged when the indicator is on.
- The body is pre-charged, so you don't need to fully charge it before you use it for the first time.

##### USB charging cable

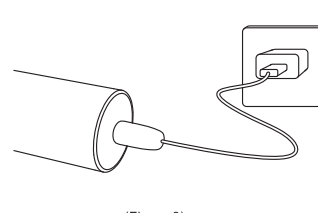(Figure 3)

##### Use this product

1. Insert the brush head carefully into the oscillating axis, and leave a certain gap between the brush head and the main body to reserve a working space for the normal work of the brush head.
2. To prevent toothpaste spilling, take a small amount of toothpaste on the brush head before use (generally speaking, sonic-wave electric toothbrush saves more toothpaste than manual toothbrush).
3. First soak the brush head and put it into the mouth, then start the toothbrush.
4. Move the toothbrush slowly with the bristles at 45° to the gum line.
5. Hold the toothbrush handle gently and move slowly in a circular motion. Do not use press it hard.
6. Presses your mouth with water after brushing.

This product can be used only with/for:  
Dental braces (when used with the dental braces, the brush head will wear faster).

Restoration of teeth (filling, crown, porcelain teeth).

##### Brushing mode

**Cleaning mode**  
Effectively remove periodontal plaque (about 30s intervals).  
Cleaning mode is 2 minutes, including four 30s intervals.

**Whitening mode**  
2 minutes to deep whitening teeth, effectively improving tooth surface color.

**Polishing mode**  
Can brighten and polish the front teeth, users can press and hold the power button to stop polishing in advance according to their own needs.

#### Milder mode (beginner mode)

It is suggested that this mode can be used for the beginner which can effectively clean teeth and gently massage gum to improve gum health.

**Soft mode**  
It is a super gentle mode for sensitive teeth and gums which lasts 2 minutes.

##### Toothbrush handle

1. Remove the brush head and wash the metal axis with warm water. Make sure to remove all remaining toothpaste.

**Caution:** Do not press sharp objects on the rubber sealing ring on the metal axis, otherwise it may cause damage.

2. Wipe the entire surface of the toothbrush handle with a wet cloth.

**Caution:** do not use rubbing alcohol of isopropyl, vinegar or bleach to clean the handle, otherwise it may cause discoloration of the handle.

##### Cleaning, service and storage

**Clean the brush head:** wash the brush head and bristles after each use, and dry them with natural air. Do not put the bristles, as the will deform them and cause them to fall off.

**Clean the body:** clean the handle surface with clean water, and scour it with soft cloth, and dry it with natural air. Do not press sharp objects on the rubber sealing ring on the metal axis, otherwise it may cause damage.

**Clean the charger:** clean the slots on the charger with a soft cloth. Be sure to unplug the power plug from the socket before cleaning the charger. Wipe the surface of the charger with a wet cloth and use it after it is completely dry.

**Replace brush head:** professional dentists recommend replacing the brush head every 2-3 months for the best results.

**Clean the bottom of the brush handle:** before charging, be sure to clean the bottom of the brush handle, wipe the bottom surface with a wet cloth, and then charge after the cleaning.

##### Storage

If you do not use this product for a long time, please fully charge the toothbrush handle, unplug the charger from the power socket, clean the two parts and store them in a cool and dry place away from direct sunlight.

##### Recovery

When discarding this product, please do not put it together with the general domestic waste. Instead, hand them to the recycling center designated by the government. This will be beneficial to the environment.

This product contains a built-in rechargeable battery, do not discard the rechargeable battery together with general domestic waste.

Follow your country's regulations on classifying and recycling electronic and electrical products and rechargeable batteries. Proper discarding of products will help to avoid adverse consequences on the environment and human health. Varying period of intelligent sonic-wave electric toothbrush.

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