

LX7 Touchscreen Hand Massager



SUPPORT@LUNIXINC.COM

Lunix Replacement Parts available at : www.lunixinc.com

Touchscreen Hand Massager

LX7



Lunix
you deserve to relax

Congratulations on your purchase of the
Lunix LX7 Electric Hand Massager!

This innovative hand massager squeezes away the stress and strain of daily life, to give you a robust massage experience unlike anything you've ever tried before!

Your satisfaction is our top priority and we are here to help with any issues!

We've started as a small American family business based in California with the CORE VALUE of enhancing each of our customer's life experience. Our team works hard every day to make the most thoughtful and the highest-quality products for those who matter the most:

OUR FAMILY & FRIENDS!

Table of Contents

Introduction	4
What's In The Box?	5
Terms of use	6
Warranty	7
Control Pad	8
Benefits	9
Before Using Your Hand Massager	10
First Step: Charging	10
Safety Precautions	11-13
Operations	14-17
How to Use the Finger Massager	18
Care, Maintenance & Cleaning Instructions	18
Troubleshooting	19
Lunix Wellness Products Collection	21-29

Please read the manual carefully
before using the Unix Hand Massager



Introduction

Thank you for purchasing the **LX7 Touchscreen Hand Massager**.

This manual contains the instructions you will need to operate and maintain your model so that you can enjoy it for years to come. However, please be advised that mishaps could happen in time and rest assured that LUNIX is a team of professionals who aim to provide the highest level of factory support possible.

We know you're excited about getting your hands on your new hand compression massager, but it's very important that you take some time to read through the User's Manual. This manual contains all the necessary setup and operating procedures that will allow you to use your new device. Also, be sure to read and follow the precautions and warnings in this manual and on any labels or tags attached to your model. They are there to educate you on how to operate your unit safely and also get maximum life and performance from it.

Thank you again for choosing LUNIX. We work hard every day to make sure that you receive the highest level of customer satisfaction possible. We truly want you to enjoy your new Hand Compression Massager!

You can find more information on our website at
www.lunixinc.com

If you have any questions about your model or its operation,
contact LUNIX Customer Support by email :
support@lunixinc.com

What's In The Box?



1x Hand Massager



1x Finger Massager



1x USB Charging Cable



1x Travel Bag



1x User Manual



PRODUCT NAME: Hand Massager
MODEL: LX7
POWER: 5W
WEIGHT: 2.5 lbs / 1.1kg
VOLTAGE: DC 3.7V
BATTERY: Lithium
DIMENSION:
11.8 x 6.3 x 4.7 in / 30 x 16 x 12 cm

Terms of Use

This product is surrendered by LUNIX to the purchaser with the understanding that the purchaser accepts the responsibility that using the LX7 Hand Massager with its accessories in a careless, improper, or unsafe manner can result in serious injury or death. Also, the purchaser assumes all liability resulting from any misuse, unsafe handling, failure to follow instructions, or any action that constitutes a violation of any applicable laws or regulations. LUNIX shall not be liable for personal injury, loss of property, or loss of life resulting from the use of this product under any circumstances, including intentional, reckless, negligent, or accidental behavior. LUNIX shall also not be liable for any special, indirect, incidental, or consequential damages arising out of the assembly, installation, or use of their products or any accessory or chemical required to use their products. By the act of operating & using the product, the user accepts all resulting liability and releases LUNIX of any and all liability associated with its use.

If you as the user do not accept liability of Usership, LUNIX requests that you do not use this product. Do not open any of the enclosed materials. Return the model to LUNIX but please be advised that we absolutely cannot accept an item for return or exchange after it has been used or is otherwise no longer in as-new condition.

All information contained in this User's Manual is subject to change without notice. For the most updated version of the User Manual, please visit our website.

LUNIX reserves the right to make changes and improvements to products without incurring any obligation to incorporate such improvements into products previously sold. If you have any questions about your model or its operation, contact Customer Support by E-mail: support@lunixinc.com

FCC Compliance

This device contains a module that complies with the limits for a Class B digital device as described in part 15 of the FCC rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

The limits for a Class B digital device are designed to provide reasonable protection against harmful interference in residential settings. This product generates, uses and can radiate radio frequency energy, and, if not operated in accordance with the instructions, may cause harmful interference to radio communications. The user is cautioned that changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

WARRANTY

To guarantee you have the best experience we offer a 1-year warranty covering defects in materials, workmanship and assembly on all our products*.

To ensure that you don't invalidate the terms of the warranty please read and understand this booklet.

*If a defect in materials, workmanship or assembly is found BEFORE the product has been used, we will exchange the entire product, if the product is returned complete and in perfect condition in its original packaging. **

IT IS THE PARENTS' OR GUARDIANS' RESPONSIBILITY TO ENSURE MINORS RECEIVE APPROPRIATE GUIDANCE AND/OR SUPERVISION WHEN USING THIS PRODUCT.

LUNIX DOES NOT operate a new for old warranty once the product has been used. Warranty is only valid if the items were purchased directly from LUNIX (purchases from third party sellers will not be covered under the warranty)

COVERED BY WARRANTY

- 1-year warranty for the Air Compression (excludes damages caused by the user, including but not limited to incorrect maintenance.)
- 6-month warranty for the Battery (excludes damages caused by the user, including but not limited to incorrect maintenance.)

THIS IS A COMPLEX HIGH PERFORMANCE HAND COMPRESSION MASSAGER AND MUST BE TREATED WITH CARE AND RESPECT.

LUNIX has made every effort to ensure your product is as strong and durable as possible. However, it is still possible that some mishaps could happen and we take full responsibility. If the air compression, the battery (or the programs) start to malfunction, please contact us at support@lunixinc.com and we will provide a replacement.


NOT COVERED BY WARRANTY

- All other parts not mentioned above
- Lack of maintenance, improper use or abuse
- Water damage or problems resulting from water intake

Control Pad



DESIGNED IN USA

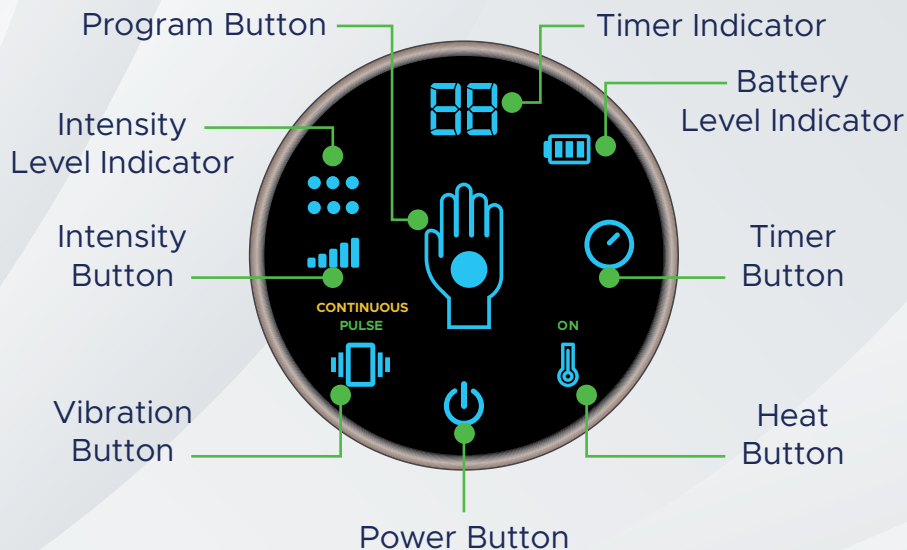

Heat



6 Levels



15'/20'/30' Session



2 Vibration
Modes



4 Programs



Shut Off
Protection


Benefits

With its ergonomic design, the Linux LX7 Full Hand Compression Massager is designed to give you the most comfortable massage. Based on the principles of human engineering, a unique design and the distribution of acupuncture points, the massager meets the demands of your hand according to its structure. It will be like having a professional masseur-physiotherapist at home.



Sit back and enjoy a good hand massage after a long day of work or before going to bed.

Use the LX7 Full Hand Compression Massager one or two sessions every day from 10 to 15 minutes for at least 7 days in a row in order to feel the benefits of the massage.

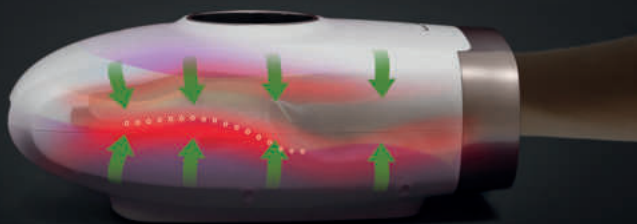

Before Using your Hand Massager

- After unpacking the hand massager, please check the integrity of the device and accessories and look for any damages.
- Before first use, charge the massager until the blue level indicator stops blinking (see below the section how to charge).
- Place the device on a flat, level surface.
- Insert one hand into the cloth-lined chamber.
- Press the Power Button, then the Program Button in order to start the massage.

NOTE: Please be sure to remove all sharp jewelry.

FIRST STEP: CHARGING

You should see a blue battery charge level indicator blinking while charging. Once the battery is fully charged, the level indicator stops blinking. In order to charge the massager, just connect the usb cable between the massager port and any usb plug (laptop, wall plug, etc.)

CHARGING TIME


5V / 1A: Around 7h

5V / 2A: Around 4h

WORKING TIME

Normal Use: Around 3h

With Heat & Vibration: Around 2h

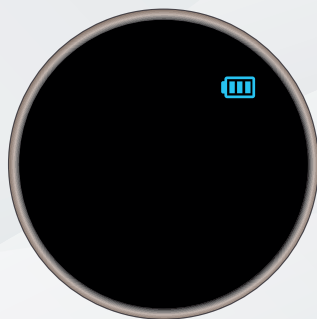

IMPORTANT:

Please fully charge the hand massager before the first use.

Keep in mind that you cannot use the massager while plugged-in and charging, due to safety reasons. If the blue light is not flashing while charging, please report to the troubleshooting page 19 of this manual.

Precautions



We want you to safely enjoy your new hand massager. Operate your model sensibly and with care, and it will be exciting, safe, and relaxing for you. Failure to operate your model in a safe and responsible manner may result in property damage and serious injury. The precautions outlined in this manual should be strictly followed to help ensure safe operation. You alone must see that the instructions are followed and the precautions are adhered to.

- 1.** Never charge the device with a USB plug or charging cable that exceeds 5V 2A as it could cause fire.
- 2.** Make sure the USB cable fits entirely into the port.
- 3.** This device should only be used with the provided USB cable.
- 4.** Do not let the device charge overnight or during the whole day.
- 5.** Make sure you avoid wooden surfaces while charging.
- 6.** Avoid placing any full glasses or liquid bottles next to the device as they could spill and damage the device.
- 7.** If the USB cable is damaged, do not charge the unit.
- 8.** Do not leave the product unattended when plugged in. Unplug from outlet when not in use.
- 9.** Do not operate after taking pain killers, sedatives or alcohol.
- 10.** Individuals with pacemakers should not use this product. Always consult a doctor before using it.
- 11.** Do not carry this product by its supply cord or use the cord as a handle.
- 12.** Keep cord away from heated surfaces.
- 13.** Do not use with wet hands or in humid places. Wetness and humidity can cause electric shock or mechanical malfunctioning.
- 14.** Do not attempt to disassemble or repair the unit, charging cord or screen. (doing so voids the warranty).
- 15.** Most importantly, use good common sense at all times.

16. Do not wear jewellery while using the massager.
17. If you have a medical condition like diabetes, inflammation, suffer from blood clots, recently had surgery and/or are pregnant, suffering from heart disease or osteoporosis, we recommend you consult your doctor before using the massager.
18. Never use the device more than 30 minutes in a row.
19. Never cover the massager while using it, this may result in overheating of the unit.
20. Never use the massager in a wet/humid room, in high temperature or next to corrosive gas.
21. Don't forget to turn off the device once you finish the massage, in order to protect the life battery.
22. Use this appliance only for its intended use.

WARNING:

Do not use the device to supplement or replace any medical treatment that you are currently receiving. Please consult your physician before use if you are experiencing any unexplained pain, already receiving medical treatment for circulation or vascular issues, or using at-home medical equipment.

Please contact your physician before using the LX7 Full Hand Compression Massager if you have physical health concerns or have any of the following conditions:

- You are pregnant.
- You have a pacemaker or other medical implant.
- You are suffering from heart or vascular disease.
- You have irregular blood pressure or are receiving treatment to normalize your blood pressure.
- You are suffering from osteoporosis.
- You are suffering from cancer.
- You have recently undergone surgery.
- You have concerns regarding areas of your health.

Never leave the appliance unattended, especially if children are present.

We recommend you limit the use of the massager to 15-minute sessions and for no more than an hour (60 minutes) consecutively.

Never use this product directly on swollen or inflamed areas of skin or open wounds. This product should NEVER be used by an individual suffering from any physical ailment that limits the user's mobility.

This product should NEVER be used on persons who have sensory deficiencies affecting their hands and appendages.

SAFETY PRECAUTIONS

This massager is capable of extreme performance. Responsibility, maturity, and common sense are required. This model must be treated with respect and caution to ensure your safety. Operating the hand massager in a careless, unsafe manner, without adequate care and preparation can result in serious injury or death.

The hand massager includes a lithium battery. Batteries have a severe risk of fire and injury if they are improperly handled, abused, or misused. The user must read and understand all included instructions, warnings, and precautions regarding the handling, care, and use of batteries.

- BEFORE you charge, ALWAYS confirm that the charger/plug settings exactly match the type (chemistry), specification, and configuration of the battery to be charged. DO NOT exceed the maximum manufacturer recommended charge rate which is 5V, 2A.
- ALWAYS charge batteries in a well-ventilated area.
- Charge the battery in a safe area away from flammable materials. Monitor the charging process and never leave batteries unattended while charging. Do not allow small children to charge the hand massager.
- DO NOT disassemble, crush, short circuit, or expose the battery to flame or other source of ignition.
- DO NOT leave the massager unattended while charging, discharging, or anytime the charger is ON with a battery connected. If there are any signs of a malfunction, unplug the power source and/or stop the charging process immediately.
- ALWAYS unplug the charger from the wall outlet and disconnect the cable when not in use.
- If the massager gets hot to the touch during the charging process (temperature greater than 110°F / 43°C), immediately disconnect the battery from the charger and discontinue charging.
- Keep a Class D fire extinguisher nearby in case of fire.
- DO NOT expose the unit to water or moisture.

WARNING: RISK OF INJURY!

We strongly recommend that people who have very sensitive hands use the lowest intensity.



WARNING! CAUTION! DANGER! FIRE HAZARD!



CHARGING OR DISCHARGING BATTERIES POSES A NUMBER OF SIGNIFICANT HAZARDS WHICH CAN INCREASE THE RISK OF FIRE, EXPLOSION, SERIOUS INJURY, AND PROPERTY DAMAGE IF NOT PERFORMED AS PER THE INSTRUCTIONS. BEFORE USE, READ AND FOLLOW ALL MANUFACTURER'S INSTRUCTIONS, WARNINGS AND PRECAUTIONS.

NEVER ALLOW CHILDREN UNDER THE AGE OF 14 TO USE THE MASSAGER WITHOUT THE SUPERVISION OF A RESPONSIBLE, KNOWLEDGEABLE ADULT. THIS IS NOT A TOY.


Operations

PLEASE WASH YOUR HANDS
EACH TIME BEFORE USING
THE HAND MASSAGER

Place your hand with your palm facing down into the hand massager before turning the device on. You will notice that your palm should fit perfectly inside. If you don't feel comfortable, you may adjust your hand until you find the perfect spot. The front opening provides a better comfort and a maximum ventilation.



POWER BUTTON

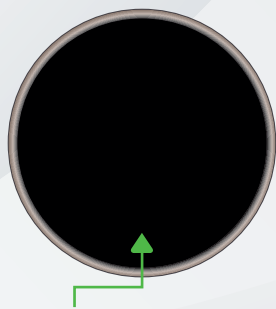


Press Power
Button Here

POWER ON: Long press the power button (located on the bottom of the black screen) for 2 seconds, the hand massager will switch on, and the power button will turn blue. Alternatively, short press the power button once, the power indicator will turn red, then long press the same button for 2 seconds, the hand massager will switch on, and the power button will turn blue. Then press the program button (the white hand), and the hand massager will start to work.

POWER OFF: Press the power button one time in order to shut down the device. The blue power button will change color to red and the airbags will progressively deflate.

After a few seconds the red power indicator will disappear.

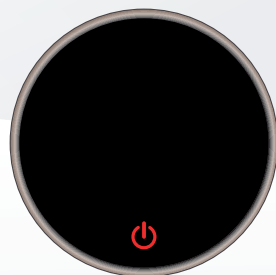




PROGRAM MODE

When turning the unit on, the display will show a white hand icon, signifying that the program mode is on standby.

Each program massages & releases pressure on different parts of the hand, palm & fingers depending on your hand shape. That means each person will feel a different and unique massage.

You can select your favorite program mode by pressing the program button (by pressing the white hand) and choose between 4 different programs.

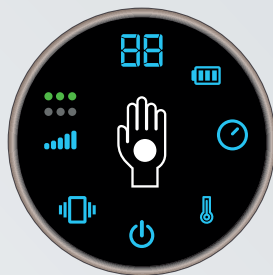

HAND MESSAGE MODES



Full Hand Massage



Fingers & Palm Massage



Wrist Massage



Heat Massage

This last program will apply heat to the hand coming from the top of the unit. The heat function provides a gentle warmth, around 104°F / 40° C, enough to relax & enhance the benefits of the massage. It may be subtle for people with sensitive hands.



HEAT BUTTON

Press the heat button to turn ON/OFF the temperature.

When the HEAT option is activated, ON in red will be displayed on the screen.

If you press the button a second time, the HEAT function will be activated



Please keep in mind that the heat function reaches the optimal temperature after 5/10 minutes of use. The heat function provides a gentle warmth, around 104°F /40°C - enough to relax & enhance the benefits of the massage.


If you want to feel more heat, we recommend you wear a pair of gloves. This will help your blood circulation.



VIBRATION BUTTON

Press the button to activate the vibration mode. When the option is activated the red light above the button will turn on.

The LX7 hand massager provides 2 different modes: a pulse mode and a continuous mode. Press the button several times to switch between the modes.

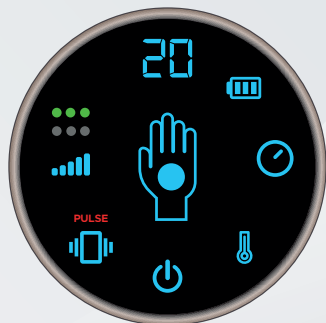

Pulse Mode



Continuous Mode



INTENSITY BUTTON

Use this function to increase or decrease the intensity level by pressing the button.

There are 6 different levels of intensity from the weakest 1 white light to the strongest 6 red lights.

We strongly recommend that people with very sensitive hands use the lowest intensity (one white light), as a higher intensity may be uncomfortable.

Wait for a few days in order to feel pain relief and to feel the full benefits of the massage.



TIMER BUTTON

When pressing the timer button, it allows you to choose between 15, 20, and 30-minute sessions.

Each time you press the button, it will change the time session. The hand massager automatically shuts off when the time has expired.



Lowest Intensity



Highest Intensity

How to Use THE FINGER MASSAGER

This Finger Massager is an easy and relaxing way to soothe your tired fingers after a long day. Just put your finger in the middle of the roller and scroll back and forth to enjoy a wonderfully comfortable massage.



Care, Maintenance & Cleaning Instructions



- Store the Hand Compression Massager in a cool and dry place when not in use.
- Do not place the unit into any type of liquid.
- Keep away from all solvents and harsh detergents.
- After use, wipe with a soft, damp cloth using non abrasive cleaners.

Troubleshooting



When operational problems occur, please refer to the chart below.

PROBLEM	POSSIBLE CAUSE	SOLUTION
I am unable to start the unit	The unit did not start properly OR battery is low	Press the Power Button for 2sec OR recharge the unit
The massage pressure is too low	Battery is low OR intensity level is too low	Recharge the unit OR increase intensity level
The massage pressure is too strong	You are using an Intensity level that is too high OR your hands are not correctly positioned inside the unit	Use I1 Intensity OR try to adjust your hand in the unit to feel comfortable
My unit suddenly stopped working	Power outage OR massage more than 15min OR long working time, the unit has been shut off due to overheating protection	The unit automatically shuts off after 15min of use. Use the unit after the motor has been naturally cooled
I hear a faint sound while using the massager	The air pump is working	Normal behavior
The airbags are not inflating	Battery is low	Recharge the unit. If the problem persists, please contact us at support@lunixinc.com
I cannot see the blue light flashing while charging OR the unit turns on for a few seconds, then all the lights flash and the unit powers down	The unit is not charging properly because the connection is not made properly between the usb port and usb cable OR the plug is higher than 5V 2A	The usb cable needs to fit entirely into the USB port OR the plug must not exceed 5V 2A
How do I know when the unit is fully charged?		The blue battery charge level indicator stops blinking once the unit is fully charged
Unit stopped charging	Unit not receiving the correct amount of voltage	Make sure the charger you used is 5V 2A and the charger (plug or laptop) must not exceed 5V
The display is constantly turning on and off making a beeping noise	PCB problem	Contact us at support@lunixinc.com for further assistance
Heat function not working	Unit not fully charged and the heat option takes some time to achieve the highest temperature	Wait at least 10-12min before feeling the heat and you can also wear a regular pair of thin gloves; thus, you can increase the heat- feeling
Unit changing intensity programs	PCB defective	Contact us at support@lunixinc.com for further assistance
Unit is hot	Excessive use	Power off the unit and let it cool before further use

If you have any questions, concerns or if you would like to order any spare parts or replacements. Please contact us at **support@lunixinc.com** and we will get back to you within 24 hours


Lunix Wellness Products Collection

LX3 Hand Massager - White



LX3 Hand Massager - Black

LX5 4Pcs Orthopedic Bed Wedge Pillow Set - Brown

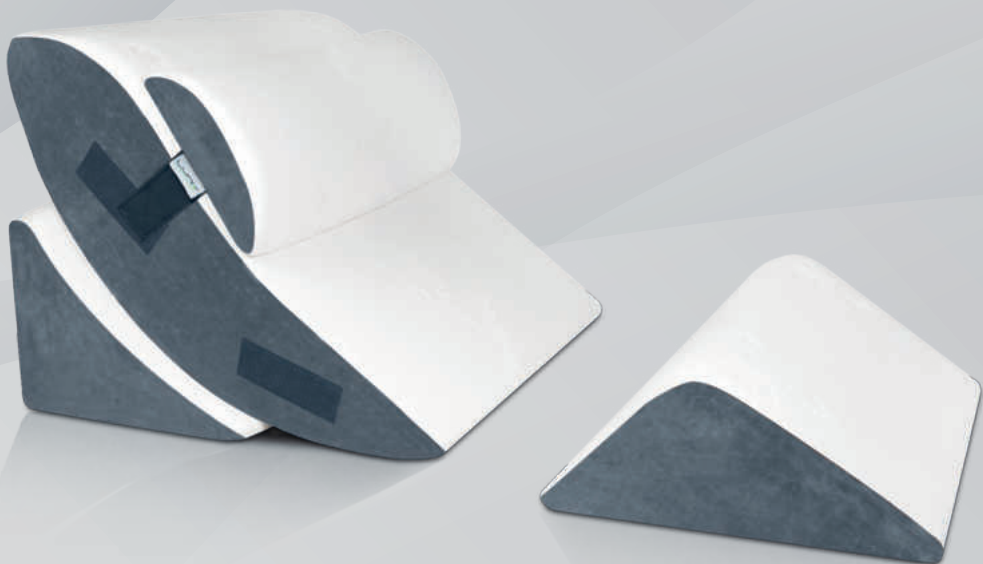

LX5 4Pcs Orthopedic Bed Wedge Pillow Set - Navy

LX6 3Pcs Orthopedic Bed Wedge Pillow Set - Brown



LX6 3Pcs Orthopedic Bed Wedge Pillow Set - Navy

LX7 Touchscreen Hand Massager - White



LX7 Touchscreen Hand Massager - Black

LX8 2 Layers Orthopedic Bed Wedge Pillow Set - Blue



LX8 2 Layers Orthopedic Bed Wedge Pillow Set - Gray

LX9 Extra Large Orthopedic Bed Wedge Pillow 7.5" - Gray



LX9 Extra Large Orthopedic Bed Wedge Pillow 12" - Gray

LX9 Extra Large Orthopedic Bed Wedge Pillow 7.5" - Navy



LX9 Extra Large Orthopedic Bed Wedge Pillow 12" - Navy

LX10 Cordless Full Leg Compression Massager - Green



LX10 Cordless Full Leg Compression Massager - Black

LX11 5Pcs Orthopedic Bed Wedge Pillow Set - Blue



LX11 5Pcs Orthopedic Bed Wedge Pillow Set - Gray

