

Oral-B®

VITALITY™

powered by **BRAUN**



Charger Type 3757
Handle Type 3709



문의사항이 있으시면 소비자상담실로
연락주시기 바랍니다.

소비자 상담실

(KR) 080-920-6000(수신자 부담)

If you should require further
assistance, please call our
Consumer Service line.

Internet:

www.oralb.com

www.oral-b.co.kr







www.braun.com

www.service.braun.com

Braun GmbH
Frankfurter Straße 145
61476 Kronberg
Germany/Allemagne

Charger Type 3757

Handle Type 3709

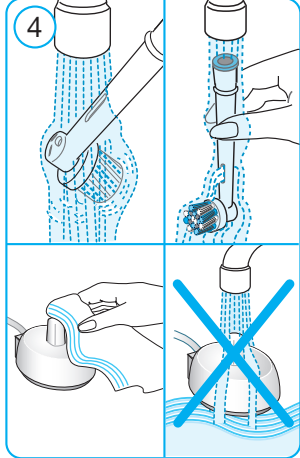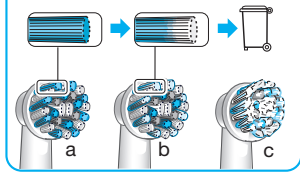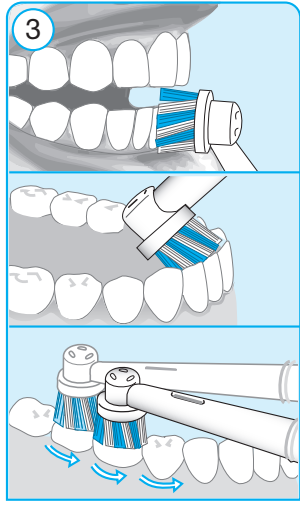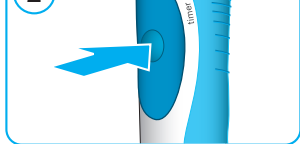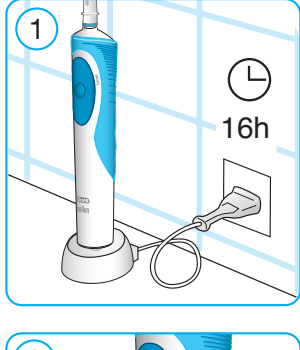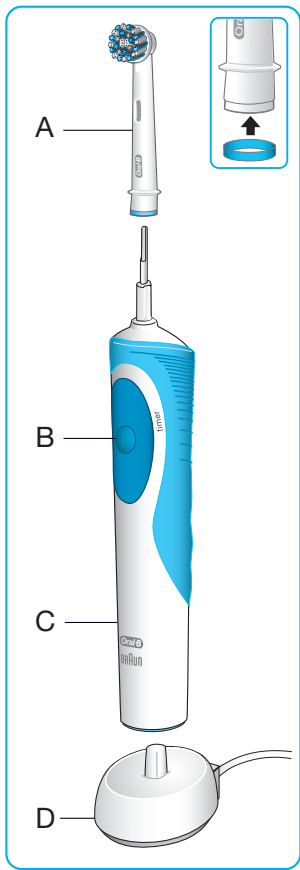

English

Welcome to Oral-B!

Before operating this toothbrush, please read these instructions and save this manual for future reference.

IMPORTANT

- Periodically check the cord for damage. If cord is damaged, take the charging unit to an Oral-B Service Centre. A damaged or non-functioning unit should no longer be used. Do not modify or repair the product. This may cause fire, electric shock or injury.
- Usage by children under age 3 is not recommended. Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge, if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Cleaning and maintenance shall not be performed by children.
- Children shall not play with the appliance.
- Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.

WARNING

- If the product is dropped, the brush head should be replaced before the next use even if no damage is visible.
- Do not place the charger in water or liquid or store where it can fall or be pulled into a tub or sink. Do not reach for it when fallen into water. Unplug immediately.
- Do not disassemble the product except when disposing of the battery. When taking out the battery for disposal of the unit, use caution not to short the positive (+) and negative (-) terminals.
- When unplugging, always hold the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.
- This toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- Replace brush (head) every 3 months or sooner if brush head becomes loose due to damage.

Description

A Brush head (depending on model/ interchangeable)

B On/off switch

C Handle

D Charging unit

Specification

For voltage specifications refer to the bottom of the charging unit.

Connecting and charging

Your toothbrush has a waterproof handle, is electrically safe and designed for use in the bathroom.

- Plug the charging unit (D) into an electrical outlet and place the toothbrush handle (C) on the charging unit (picture 1).
- A full charge takes about 16 hours and enables up to five days of regular brushing (twice a day, 2 minutes).

- For everyday use, the toothbrush handle can be stored on the plugged-in charging unit to maintain it at full power. Overcharging is impossible.

Battery maintenance

To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and discharge the handle by regular use at least every 6 months.

Brushing technique

- Wet brush head (A) and apply any kind of toothpaste. To avoid splashing, guide the brush head to your teeth before switching the appliance on (B).
- When using one of the rotating brush heads guide the brush head slowly from tooth to tooth, spending a few seconds on each tooth surface.
- When using the Oral-B «TriZone» brush head place the toothbrush bristles against the teeth at a slight angle towards the gumline. Apply light pressure and start brushing in back and forth motions, just like you would do with a manual toothbrush.
- With any brush head, first brush the outsides, then the insides and finally the chewing surfaces (picture 3). Do not press too hard, simply let the brush do all the work. Brush at least for 2 minutes. Do not press too hard, simply let the brush do all the work. You may also consult your dentist or dental hygienist about the right technique for you.
- Turn off the handle by pressing the on/off switch (B).

In the first days of using any electric toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or dental hygienist.

Timer (depending on model)

An incorporated timer memorizes the elapsed brushing time (depending on model).

«2-Minute» timer: A short stuttering sound indicates the end of the professionally recommended 2-minute brushing time.

«Professional» timer: A short stuttering sound at 30-second intervals reminds you to brush equally all four quadrants of your mouth. A long stuttering sound indicates the end of the professionally recommended 2-minute brushing time.

Brush heads

Oral-B offers you a variety of different Oral-B brush heads that fit your Oral-B toothbrush handle. Our **oscillating-rotating brush heads** can be used for tooth-by-tooth precision cleaning.

Oral-B CrossAction brush head



Our most advanced brush head. Angled bristles for a precise clean. Lifts and powers away plaque.

Oral-B FlossAction brush head



features micropulse bristles that get deep in hard-to-reach areas between your teeth, gently removing hidden plaque.

Oral-B Precision Clean brush head



cups each tooth with a curved bristle trim and reaches deep between teeth at the same time.

The Oral-B Sensitive brush head



features an extra soft bristle structure that is gentle on teeth and gums.

Oral-B 3DWhite brush head



features a dentist-inspired polishing cup to remove surface stains to naturally whiten teeth. We do not recommend using the Oral-B 3DWhite brush head with braces.

Our **dynamic «TriZone» brush head** can be used with all common brushing techniques.

The Oral-B TriZone brush head



features triple zone cleaning action for outstanding plaque removal, even in between teeth.

The Oral-B TriZone brush head



features triple zone cleaning action for outstanding plaque removal, even in between teeth.

Oral-B brush heads feature light blue INDICATOR® bristles to help you monitor brush head replacement need. With thorough brushing, twice a day for two minutes, the blue color will fade halfway approximately within 3 months, indicating the need to exchange your brush head. If the bristles splay before the color recedes, you may exert too much pressure on teeth and gums (picture 3).

We do not recommend using the Oral-B «FlossAction» or Oral-B «3D White» brush head with braces. You may use the Oral-B «Ortho» brush head, specifically designed to clean around braces and wires.

Cleaning

After use, rinse the brush head for several seconds under running water with the handle switched on. Switch off the handle and take off the brush head. Rinse both parts separately under running water and wipe them dry. From time to time, also clean the charging unit with a damp cloth (picture 4).

Subject to change without notice.

Environmental Notice

Product contains batteries and/or recyclable electric waste. For environment protection do not dispose in household waste, but for recycling take to electric waste collection points provided in your country.



Guarantee

We grant 1 year guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor. This guarantee does not cover: damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorized persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorized Oral-B Braun Service Centre.

Warranty replacement brush heads

Oral-B warranty will be voided if electric rechargeable handle damage is found to be attributed to the use of non-Oral-B replacement brush heads. Oral-B does not recommend the use of non Oral-B replacement brush heads.

- Oral-B has no control over the quality of non-Oral-B replacement brush heads. Therefore, we cannot ensure the cleaning performance of non-Oral-B replacement brush heads, as communicated with the electric rechargeable handle at time of initial purchase.
- Oral-B cannot ensure a good fit of non-Oral-B replacement brush heads.
- Oral-B cannot predict the long-term effect of non-Oral-B replacement brush heads on handle wear.

All Oral-B replacement brush heads carry the Oral-B logo and meet Oral-B's high quality standards. Oral-B does not sell replacement brush heads or handle parts under any other brand name.

Manufacturing date code:

To determine the year and month of manufacture, refer to the 4-digit production code located on the type plate of charger and handle. First two digits of the production code refer to the year of manufacture. Next two digits refer to the month of manufacture. Example: «1209» – means the product was manufactured in September of 2012.