KU2 Cosmetics

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
Product description

The IPL device for hair removal uses light pulses to treat hair follicles and prevent them from growing. The light pulses in a specific wave range cause the hair follicles to heat up. When the hair follicles reach a certain temperature, they enter the resting phase, which causes the hair to gradually stop growing and eventually fall out.


Main Machine

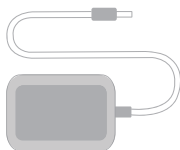
Adapter


Goggles
User Manual

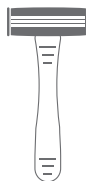
Razor

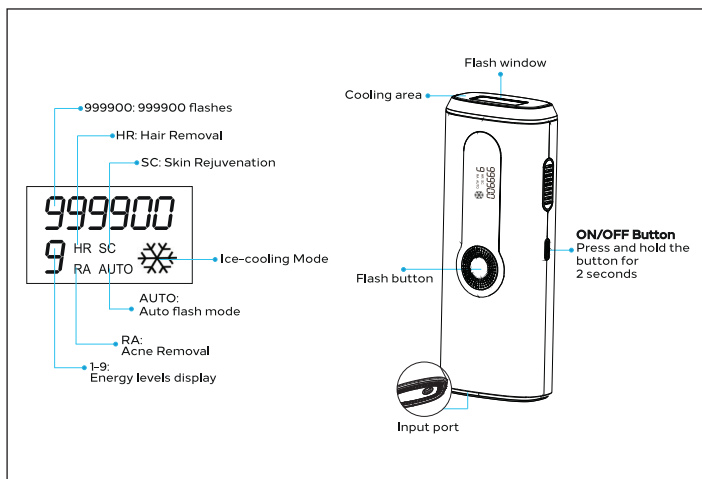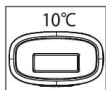The ice-cooling mode can reduce irritation and burns while in use.

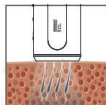
Light pulsing technology can effectively and safely remove unwanted body hair. In addition, light impulses can smooth the skin structure and improve acne scars.

How IPL technology works


IPL light pulse technology is highly effective at removing hair and leaving the skin silky smooth. The light pulse strikes the color pigment (melanin) located in the hair follicle. The light pulse absorbed by the pigment is converted into heat, causing the hair follicle to heat up briefly and destroy the hair or hair root.

Please note:


Hair regrowth is delayed or suppressed. The hair may only fall out after one to two weeks. Optimal hair removal cannot be achieved with a single treatment. To achieve long-term hair removal, you must perform the treatment at least every two weeks for the first few months.

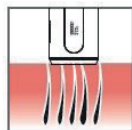Start (Weeks 1-4): Treat twice a week during the first 1-4 weeks.

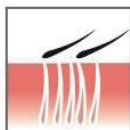In Between (Weeks 5-8): Weekly treatment during weeks 5 to 8.

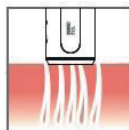Maintenance (After 8 weeks): Continue IPL treatment as needed (use every 1-2 months).

Technical data

Model name: IPL Supersilk KU20
Product name: IPL device
Rated power: 36 watts
Rated voltage: 100-240 V, 50/60 Hz
Adapter type: DC 12V, 3A
Light area: 3.3 m²
Technology used: IPL (Intense Pulsed Light)
Intensity density: 1.5-3.9 J/ cm²
Wavelength: 510 ~ 1100 nm
Lifetime of the lamp tube: 999999 flash times
Intensity levels: 9 levels
Operating humidity range: Relative humidity 30 %-60 %
Operating temperature 5 ~ 35 °C

General instructions

Read the operating instructions very carefully before using this appliance and keep them in a safe place, including the invoice and, if possible, the box with inner packaging. If you pass the appliance on to a third party, also pass on the operating instructions.

- Only use the appliance for its intended purpose at home. This device is not intended for commercial use.
- The adapter is only suitable for use in dry rooms.
- Ensure that the cord is not bent, clamped or comes into contact with heat sources.
- Do not use a damaged adapter.
- Do not repair the device yourself.
- The device and the power cord must be checked regularly for signs of damage.
- If damage is found, the device should no longer be used.
- Only use original equipment.
- For the safety of your children, do not leave any packaging parts (plastic bags, cardboard, polystyrene, etc.) within reach.

Special information for this appliance

This device is a powerful electrical appliance and should be used as such with special attention to the safety instructions. Read all warnings and safety instructions carefully before use. Follow them exactly during use! This device is not suitable for everyone. Follow the instructions below.

The device works best on darker hair types or hair that contains more melanin. Black and dark brown hair reacts best. Although brown and light brown hair may also respond, these usually require several sessions. Natural red hair may show no reaction. Natural white, gray or blonde hair usually does not react with this device.

Hair growth does not stop immediately after the first application. You will obtain a permanent and satisfactory result if you follow the instructions in this user manual.

Special safety instructions

WARNING: The appliance radiates strong visible and invisible optical radiation that can cause eye damage.

WARNING: Avoid direct exposure to the eyes.

CAUTION: Do not use the device on children.

CAUTION: Keep the device dry.

WARNING: Failure to follow the instructions may result in eye damage and possible loss of vision or skin damage.

WARNING OF ELECTRICAL SHOCK: Do not use this appliance near bathtubs, sinks or other receptacles containing water.

- Keep the power supply unit dry.
- Do not use any adaptor other than the one supplied.
- This device must not be used by children.
- Children must not play with the appliance.
- Appliances can be used by individuals with reduced physical, sensory or mental capabilities or lack of experience and / or knowledge if they are supervised or have been instructed in the safe use of the appliance and understand the hazards involved.
- Do not use the device in areas where the skin's health is compromised. For example, on areas of skin with burns, injuries, blisters, scars, cuts, open wounds, active skin diseases, recent tanning or sunburn and / or infections.
- This device is not a medical device. Medical advice should be sought for rashes, itchy skin, skin fungus and/or skin infections, skin bumps or skin spots.
- Prolonged or repeated use in one area may cause potential damage to the skin.

- The device and the light emission optics must be checked regularly for damage. Do not use a damaged device.
- Protect your eyes from the radiation! The device must not be used on the eyelids.
- Misuse of the device can lead to eye damage.
- Carry out a skin compatibility test before complete use.
- Refer to the table for the different hair and skin types. If you are unsure whether the device is suitable for you, consult your doctor or dermatologist.

Intended use

The device is used to remove unwanted body hair. The device is only suitable for people over 18 years of age. It is also used to reduce the regrowth of body hair. The device is designed for the following areas of application: Legs, arms, bikini line, armpits, chest and back. It is not suitable for the area of the nipples, anus or genitals (male or female).

Possible side effects

If the instructions in this manual are followed, side effects and complications associated with the use of this device are rare. However, any cosmetic procedure, including those intended for home use, carries a certain degree of risk. Therefore, it is important that you understand and accept the risks and complications associated with intense light pulse hair removal systems for home use.

Slight irritation of the skin

Although hair removal with intensive light pulses is generally very well tolerated, slight irritation of the skin can occur during treatment. The irritation of the skin is usually described as a slight stinging sensation on the treated skin area. The stinging sensation usually lasts as long as the flash of light is active, or for a few minutes afterwards. Anything more than this slight discomfort is not normal and means that you should either stop using the device as you cannot tolerate this hair removal method or that the light intensity has been set too high.

Reddening of the skin

Your skin may turn red immediately after using the device or within 24 hours of use. The redness usually fades within 24 hours. Consult your doctor if the redness does not disappear within 2 to 3 days.

Increased skin sensitivity

The skin of the treated area is more sensitive than usual, which may cause the skin to dry out or peel off.

Skin wounds and burns

In very rare cases, burns may occur following use of the device. The burn may take a few weeks to heal. In extremely rare cases, a visible scar will remain.

Scarring

Although very rare, permanent scarring may occur. When scarring occurs, it usually appears in the form of a flat and white lesion (hypotrophy) on the skin. However, it can also be large and red (hypertrophy) or large and extending beyond the edges of the injury (keloid).

Pigmentary changes

The device targets the hair shaft, in particular the pigmented cells in the hair follicle and the hair follicle itself. Nevertheless, there is a risk of hyperpigmentation (increased pigment or brown discoloration) or hypopigmentation (fading) of the surrounding skin. This risk of changes in skin pigmentation is higher in people with darker skin color. Usually, discoloration or pigmentation changes of the skin are temporary, and permanent hyper- or hypopigmentation rarely occurs.

Excessive redness and swelling

In rare cases, the treated skin may become very red and swollen. This occurs more frequently on sensitive areas of the body. The redness and swelling should subside within 2 to 7 days. Frequent cooling can help to reduce the swelling. Gentle cleansing of the skin is allowed, but you should avoid exposure to the sun.

Infection

Infection of the skin is extremely rare. However, there is a possible risk following a burn or skin wound caused by the device.

Hematoma

In very rare cases, bruising (hematomas) may occur after use of this device, which may last for 5 to 10 days. When the hematomas fade, a rusty-brown discoloration of the skin (hyperpigmentation) may occur.

Hair color and skin type table

Do not use the appliance on naturally dark skin or on tanned skin. Use on dark skin can result in burns, blisters and changes in skin color, for example.

Skin Tone \ Hair Color							
White/Grey		✗	✗	✗	✗	✗	✗
Light Blonde		✗	✗	✗	✗	✗	✗
Red		✗	✗	✗	✗	✗	✗
Dark Blonde Light Brown		✓	✓	✓	✓	✗	✗
Brown		✓	✓	✓	✓	✗	✗
Dark Brown		✓	✓	✓	✓	✗	✗
Black		✓	✓	✓	✓	✗	✗

Long-term hair removal

A typical hair growth cycle can last 18 - 24 months. To achieve long-term hair removal, several sessions with this device may be necessary during this time. Depending on hair color, the area of the body and how the device is used, the effectiveness of hair removal varies from person to person.

Avoid side effects:

- Do not use the device several times on the same area of skin during a session!
- Discontinue use immediately if your skin shows blisters or burns! Never use flammable liquids such as alcohol (including perfumes, disinfectants or other substances containing alcohol) or acetone to clean the skin before using the device.
- Shave the skin in the area where you intend to use the device. The hair should not be removed with wax or torn out in any way.
- The skin should be clean, dry and free of any antiperspirants or deodorants.

Carry out a skin tolerance test:

The purpose of the skin tolerance test is to determine how your skin will react to treatment with the device. Perform a skin tolerance test 48 hours prior to full treatment on each area of skin you wish to treat.


Avoid complications after application:

- Do not expose the skin to the sun after application. This also applies to tanning with artificial light (tanning bed).
- Do not use the device within 24 hours before or after a visit to the sauna or swimming pool.
- Cover moles and tattoos before you start using the device.
- Cover dark brown or black spots such as freckles, moles or warts before you start using the device.
- Never look directly into the flash of light! Always wear the safety goggles supplied when using the device.
- Do not use the device for any purpose other than its intended use.
- Always ensure that the skin contact surface is in full contact with the skin.
- Remove the device from the skin if either the skin contact surface or the skin feels hot to the touch.

Application Areas and Recommended Treatment Duration

For women: Application from head to toe, including upper lip, chin, armpits, legs, arms, and bikini area. For men: The ideal device for chest, back, arms, abdomen, and leg hair.

Note: The treatment times mentioned above for hair removal are based on various lab tests and trial runs. Individual treatment times may vary from person to person.

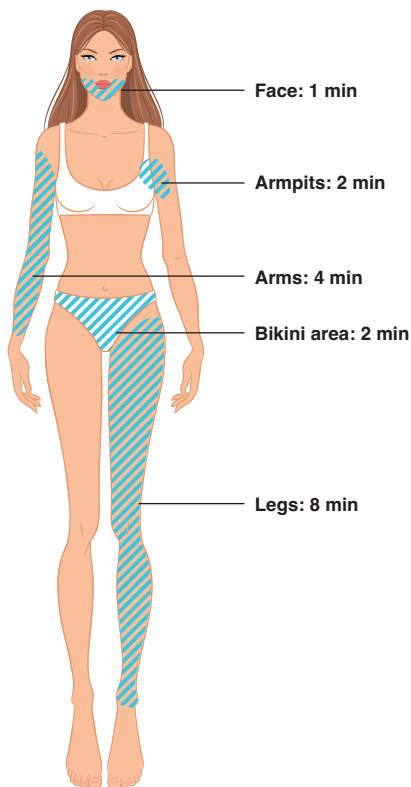Typical application schedule during a hair growth cycle:

Number of sessions	Time interval in weeks
1 - 4	2
5 - 7	4
8 +	if necessary, until the desired results have been achieved

Unpacking the device

1. remove the device from its packaging.
2. Remove all packaging materials such as foil, filling material, cable holders and cardboard packaging.
3. check the contents for completeness.
4. if the contents of the packaging are incomplete or damaged, do not operate the device. In this case, contact the manufacturer immediately.

Before use:

Carefully select the intensity level of the light flashes:

The intensity level refers to the intensity of the light flashes emitted during use, from the lowest (1) to the highest level (9). The level indicator on the display shows the intensity level set.


Always start the treatment at the lowest intensity level:

Only if you experience little or no discomfort during or after using the device at the lowest intensity level can you increase the intensity level by one level the next time you use it.

Usage


Shave the hair before each treatment to achieve optimal results.

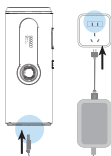
Then plug the device into a nearby power outlet.


Put on the protective glasses provided with the product or sunglasses to protect your eyes. (Note: Never look directly into the flash.)

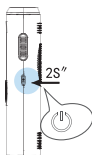
Press the „Power“ button for 2 seconds to turn on the device, and the energy indicator will light up. (When the device starts, the default setting is level 1.)

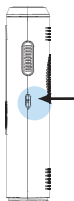
Briefly press the power button to set the energy level from 1 to 9.

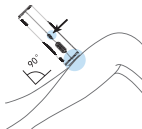Hold the device vertically close to the skin so that it is pressed against the skin. After scanning the skin, the display will switch to a constantly blinking hair count, indicating that the epilator is ready for use.

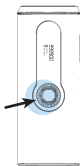Press the flash button to trigger the flash. Hold the button for 3 seconds to switch between automatic and manual modes. Press the flash button twice to switch to the HR/RA/SC mode.

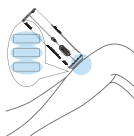Move the product over all areas of the skin where hair removal is desired while the device is running. Remember to gradually increase the light intensity only within the recommended range. Do not extend the flash interval duration at the same time.

After use

- Use a dry, clean cloth and, if necessary, a cleaner specially developed for electronic devices to gently wipe the surface of the device and in particular the surface where the device rests on the skin.
- After cleaning, keep the device in its original packaging. Store the device in a cool and dry place between sessions. Keep it away from water.

Problem handling

Problem	Cause	Solution
The device or the power supply unit heats up during use.	This is normal.	No measures are required.
The appliance makes loud noises when in use.	The built-in cooling fan of the appliance makes a noise similar to a hairdryer. This is normal.	No measures are required.
The device does not work/ cannot be switched on.	The device is not properly connected to the power supply/ power supply is faulty.	Check if the power supply unit is connected correctly. Check the power supply unit for defects. Contact customer service if necessary.
	The device is faulty.	Please contact our customer service.
No flash of light is emitted when the flash button is pressed.	Device does not rest properly on the skin.	Make sure that the flash light window is positioned correctly on the skin.
	The lamp is defective/ burnt out.	Please contact our customer service.
During application smoke can be observed.	The device is faulty.	Please contact our customer service.
Using the device causes pain.	Light intensity is too high.	Reduce the light intensity.
	Skin was not shaved thoroughly beforehand.	Make sure that the body area has been shaved before use.
	There is an intolerance.	Carry out a tolerance test before use. Do not use the device on injured/irritated skin.
No result can be determined after use.	The appliance has not been used often enough.	The appliance may need to be used more often at the start of use.
	The device is not suitable for the particular skin tone/ hair color.	Check whether the appliance is suitable for you before use. (see hair color and skin type table)

Warranty

In the event of a service / warranty claim, please contact our customer service department at support@ku-europe.de

Environmental protection

Do not dispose of the appliance in normal household waste at the end of its service life. Take it to an official collection point for recycling.
This will help to protect the environment.

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KU Europe GmbH