Massage Gun

Model: RP-GM171H

Thank you for your purchase of RENPHO Massage Gun. The massage gun offers a powerful percussion massage. While enjoying muscle and body relief, you should use it with care and safety in mind. To ensure safety and avoid any property damage, please read this manual carefully and use as indicated. Non-compliance may cause serious injuries or death and/or severe property damage.

Important Safety Instructions

When using an electric appliance, basic precautions should always be followed, including the following:

Read all instructions before using (this appliance).

1. This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
2. Children shall not play with the appliance.
3. Cleaning and user maintenance shall not be made by children without supervision.
4. The appliance is not intended for medical use. It is intended only for household to provide a luxurious massage.

DANGER - To reduce the risk of electric shock:

- 1) Always unplug the charger from the electrical outlet after charging.
- 2) Do not reach for the appliance that has fallen into water.
- 3) Do not use while bathing or in a shower.
- 4) Do not place or store appliance where it can fall or be pulled into a tub or sink. Do not place in or drop into water or other liquid.

WARNING - To reduce the risk of fire, electric shock, or injury to persons, please follow the instructions below:

- 1) The appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before cleaning, putting on or taking off parts.
- 2) Do not operate under a blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 3) Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.

- 4) Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- 5) Never operate this appliance if it is not working properly, has been damaged, or dropped into water. Return the appliance to a service center for examination and repair.
- 6) Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like.
- 7) Never drop or insert any object into any opening.
- 8) Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- 9) Do not keep the appliance near a magnetic field or where it can be affected by a magnetic field, for it can cause batteries to discharge.
- 10) Do not use the appliance before bed as it has a stimulating effect and can delay sleep.
- 11) Avoid touching the gap between the massage head and body with your finger or hair. Beware of pinching.
- 12) It is recommended to start the massage with the lowest intensity. Do not use for more than 30 minutes at a time. Excessive use could lead to slight bruise or appliance overheating that shortens its lifespan.
- 13) Please keep the appliance away from heat and avoid using it in places where the surrounding temperature is relatively hot.
- 14) Do not use the appliance with wet hands.
- 15) Only use the original massage head officially provided in the package. The use of a third-party massage head may cause damage to the appliance and void the warranty.
- 16) Do not drop, apply excessive force or heavy weight over the appliance to avoid malfunction and/or damages.
- 17) If a high temperature or foul odor occurs while using the appliance, stop using it immediately and unplug the appliance.
- 18) It is prohibited for non-professional personnel to repair the appliance or replace parts.
- 19) Use only massage heads that produce the best results.
- 20) Do not press the massage head vigorously or cover the appliance during a massage.
- 21) Do not use the appliance in close proximity to long hair, loose clothing or jewelry.
- 22) Do not massage on the same part for a long time in a row. If you feel unwell, please stop using it immediately.

Do not use the massage gun on

- 1) Sensitive skin, swollen or inflamed areas, in the area of poor circulation, where skin eruptions are present, or in the area of unexplained abdominal pain;
- 2) Head, face, cervical vertebrae, spine vertebral body;
- 3) Varicose veins;
- 4) Tumor site;
- 5) The throat area;
- 6) Hard parts such as the head or body bones;
- 7) Genital area (do not use this for sexual purposes. It is a powerful massaging appliance and may cause severe bodily harm);
- 8) An unconscious or sleeping person;
- 9) Any possible fracture site;
- 10) Positions close to the bony areas: the humerus, the back of the foot, the back of the hand and so on.

If you fall into any of the following user groups, please consult a doctor before use:

- 1) Blood stasis;
 - 2) Fractures;
 - 3) Heart disease, with a pacemaker, defibrillator, or metal contents in the body;
 - 4) Infants/minors/pregnant women;
 - 5) If you have experienced joint surgery within 90 days, do not use the appliance on the body within 3 inches of the surgical position.
- If you feel unwell or have any unusual discomfort, please stop using the appliance immediately.

About Battery

1. Warning Chemical Exposure: If the battery pack leaks, Do Not allow the liquid to contact with the skin or eyes. If contact has been made, wash the affected area with copious amounts of water and seek medical advice.
2. Do not throw the battery into fire or heat it. The used battery should be placed in a sealed plastic bag and safely disposed of in accordance with local environmental regulations.
3. Battery life can vary by application and/or operating conditions. If the appliance is not in use for a long period of time, please keep it powered off and charge it once a month to prolong its lifespan.

SAVE THESE INSTRUCTIONS

About the RENPHO Massage Gun

1.What's in the box

RENPHO Massage Gun x 1

Massage Head x 5

Type C Charging Cable x 1

Rubber Sleeve x 2

User Manual x 1

2.Product Introduction

Removable and replaceable massage head

TIP: None of the massage heads are to be used for your head. **DO NOT** use on your head.

Heating Head

Type-C charging port (Use 5V 1A/2A charger, not included)

When the massage gun reaches a certain amplitude, the massage head will automatically turn on the heating function.

- Built-in 600mAh lithium battery
- Charging time: about 90 minutes
- Battery life: about 90 minutes

* When the red light inside the heating head flashes, please charge it.

Ball head

Ball head: For large muscle groups like quads, glutes, extreme muscles, and joint soreness.

U-shaped head

U-shaped head: for spine, neck, and shoulder for sensitive muscles and soft tissues.

Air cushion head

Air cushion head: for sensitive body flesh relaxation.

Heating head

Heating head: for all parts of the body.

Bullet head

Bullet head: for joints, deep tissues, trigger points and small muscle areas like feet, and wrist.

Rubber Sleeve

Each massage head is pre-installed with a rubber sleeve, if it is damaged, replace it with a spare one.

Using the RENPHO Massage Gun

1. Important tips before use

1.1 Charging

- Please charge the massage gun before use.

DON'T use the massage gun when charging!

- To charge the massage gun, plug the type C cable into an adapter and into the unit's charging interface.
- The indicator light on the massage gun will **BLINK BLUE** during the charging process. The indicator light will be **SOLID BLUE** when charging is complete.
- For safety considerations, the massage gun will NOT work while charging. Additionally, the massage gun will automatically shut off after 10 minutes of continuous use. If you want longer massage, please turn it on then.

Laptop/Desktop

Adapter

RENPHO
Massage Gun

Note: Adapter is not included.

1.2 Battery Life

- This Massage gun has a built-in rechargeable battery. A full charge takes about 3 hours. Total battery life is up to 150 minutes after a full charge. Please fully charge it before first use.
- Battery life can vary by application and/or operating conditions. If the product is not in use for a long time, please keep it powered off and charge it once a month so that the battery will last longer.
- Please keep the massage gun away from heat and avoid using it in places where the surrounding temperature is relatively hot.


2. Preparation

- Please fully charge the massage gun for 3 hours before use.
- Select the appropriate massage head to connect with the host.

3. Operation Instruction

1. Connect the device with the massage head

Insert the massage head completely into the host when using, and plug out when stop using. Note: The installation of the U-shaped head should be positioned. When the buckle is inserted into the slot, it means the installation is in place. If the massage gun's head is not installed in place, there will be a risk of falling off while use.

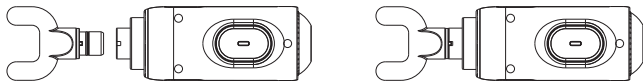2. Long press the  button on the panel for 2 seconds to power on and start the first speed, short press the  button on the panel, the speed will switch between 1-5 level. Long press the  button for 2 seconds to power it off.

*After 10 minutes of operation, the device automatically powers off.


3. Use the massage gun on any area of your body to relieve stiffness, soreness, and pain.

Shoulder Massage


Back Massage

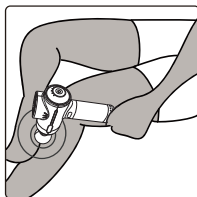
Foot Massage

Arm Massage

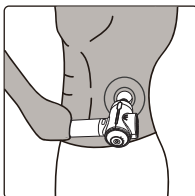Waist Massage

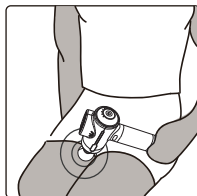Leg Massage

Product Care and Maintenance

- Place the massage gun in a safe, cool, and dry place when not in use.
- Always unplug the massage gun prior to cleaning.
- After use, wipe massage gun with a soft, damp cloth using non-abrasive cleaners. If necessary, clean with a sterilizer.
- **DO NOT** submerge the massage gun in water or any other liquid.
- Keep away from all solvents and harsh detergents.
- **DO NOT** attempt to repair this massage gun. There are no user-serviceable parts.
- **DO NOT** dispose the massage gun in a land fill. It contains a Li-ion battery and other electrical components, which is considered hazardous E-waste.

- 1) The appliance has a heated massage head. Persons insensitive to heat must be careful when using the appliance.
- 2) Use heated surfaces carefully. May cause serious burns. Do not use it on insensitive skin or on a person with poor blood circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
- 3) Temperatures sufficiently high to cause burns may occur regardless of the control setting. Do not use on an infant or invalid or on a sleeping or unconscious person. Check the skin in contact with the heated area of the appliance frequently to reduce the risk of blistering.

Troubleshooting

PROBLEM	CAUSE	SOLUTION
Massage gun doesn't turn on.	Power is off.	Long press the Power Button to turn on the massage gun.
	Battery is too low.	Charge the massage gun.
	Massage gun shuts-off after 10 minutes.	Long press the Power Button to turn on the massage gun for another 10 minutes.
Massage gun shuts off automatically.	The massage gun's overheat protection will switch off the massage gun automatically when too hot.	Wait for the massage gun to cool down before powering it back on.
	Massage gun shuts off after 10 minutes.	Long press the Power Button to turn on the massage gun for another 10 minutes.
Massage gun fails to run.	The massage head is to inserted too tight.	Plug out the massage head a little, then power it on.
If you encounter any other problems, contact: support-eu@renpho.com		

Specifications

Product Name: Massage Gun
Model: RP-GM171H
Input: 5V $\equiv$ 2A
Battery Voltage: 7.4V $\equiv$
Rated Power: 25W
Battery Capacity: 2500mAh
Amplitude: 10mm
RPM: Max 3200r/min
Auto Shut-off Timer: 10 minutes
Product Dimensions: 186*147*60mm
Product Weight: 680g

Warranty Policy

Your RENPHO product purchase is covered by a one-year limited manufacturer warranty from the date of delivery.

For warranty terms and conditions, please visit:

<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty. If you choose not to register your product, it will not diminish the product warranty.

Customer Service

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

 **Tel:** +44(0)7434-666088 (UK ONLY)
Monday-Friday 9:00AM-5:00PM(GMT)

 **Email:** support-eu@renpho.com(UK&EU)

***For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.**

Notice:

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

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EU Importer: Renpho GmbH
Manufacturer: Shenzhen Ruiyi Business Technology Co., Ltd.
UK REP Richers & Vision Ltd.
EU REP C&E Connection E-Commerce (DE) GmbH