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
ALDOM MESSAGE GUN

User Manual

Please read this manual carefully before using and keep this in a safe place for future reference

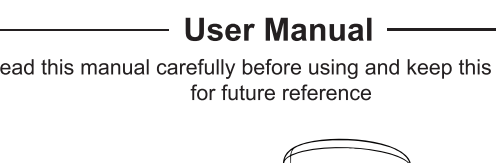Language Catalogue

English	P1-P7
Deutsche	P8-P14
Français	P15-P21
Italiano	P22-P28
Español	P29-P35
日本語	P36-P42

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HOW TO INSTALL MESSAGE HEAD

Insert the message head vertically.
To prevent, twist, other methods
the message head cannot be
directly and do not rotate.

Put it out

HOW TO USE

1. Press the power switch at the bottom of the handle, the touch screen will display the current power, the load is equal to eight times its fully charged capacity.
2. Press the power switch at the bottom of the handle, the touch screen LCD shows the power, press "+" or "-" to adjust speed according to your preference. The LCD shows the speed when you are using. The massage gun has 10 speed levels.
3. Massage your desired body part for 10 to 15 minutes.
4. When you finish the massage, press the power switch to turn off the massage power switch at the bottom of the handle. Turn it on it is a safe place.

MAINTENANCE

1. Cleanst start correctly: the battery may be used, please charge it if any other abnormal situation occurs, please contact us at sales@lodon.co.uk

DISCLAIMER

PLEASE READ PRIOR TO USE

WARNING & CAUTION:

- Do not use the machine while using.
- Do not use the machine or attempt to disassemble.
- Do not operate continuously for more than one hour. Allow to rest for 20 minutes before re-use.
- Do not charge unattended or leave charging overnight.
- Do not use the machine from power source before the charge is fully done.
- Do not use in a environment where the temperature is too high or environment is too hot.

• For adult use only. Do not use if injured. Consult your doctor before using the product.

2. SAFETY INSTRUCTIONS

- **⚠ DANGER - TO REDUCE RISKS OF ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, PLEASE READ AND FOLLOW THE FOLLOWING INSTRUCTIONS:**
- **Do not read while using.**
- **Do not use the device on the body while clothing is tightly pressing and moving.**
- **Do not use the body for approximately 20 seconds per person.**
- **Use the high power of the machine, it is recommended to use for 10-15 minutes at a time to avoid physical discomfort.**
- **Do not use device only on the body's soft tissue as desired without bruising or discomfort.**
- **Do not use on head or any other part of body area of the body.**
- **During or after use, regardless of control setting or pressure applied. Check treated areas for bruising or immediate injury. Do not use on the face or delicate skin.**
- **Check fingers, hair or any other body part away from the end and backside of the applicator or handle.**
- **Do not use the device on the face, neck, head, or back of the neck.**
- **Do not place any objects into the ventilation holes of the device.**
- **Do not immediately touch the device after the use.**
- **Do not get wet or misuse the device.**
- **Do not immerse device directly into water. The ventilation holes of the device.**
- **Turn the main power switch off with each use of the instrument.**
- **Do not tamper or alter the device.**
- **Please leave the device operating or charging unattended.**
- **NEVER USE TO REMOVE OR DISLOCATE, OR ANY VARIATION BEFORE FIRST OBTAINING PRIOR PERMISSION FROM YOUR DOCTOR IF ANY OF THE FOLLOWING ARE PRESENT:**
- Osteoporosis patients, Pregnancy, menstrual women with physical trauma or feel discomfort, skin diseases, blood diseases, heart diseases, high blood pressure, or any variation device. These contraindications do not mean that you are not able to use a vibration or exercise device but you should be aware of the contraindications.
- If you are a beginner, suggestion you use the machine from level one with the round massage head.

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