
INSTRUCTIONS FOR USE

Quick Start Guide



DO NOT use if you are pregnant.

DO NOT use the device if you are under the age of 18
or use it on anyone else under the age of 18.

DO NOT use if you have a pacemaker.

DO NOT use the device on or around the eyes, eyebrows, or eyelashes.


DO NOT use the device anywhere on the ears, nipples,
genitals or around the anus.

Frequency of treatment


- The device should be used 3 times a week for the first 2 weeks.
- For the following 2-3 weeks, use the device twice per week.
- After 5 weeks, use the device once a month.

* Depending on the hair condition, it is recommended to use the device more frequently at higher power levels if the hair is thick.




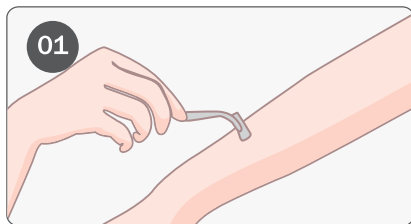


Ensure you shave the treatment area.
PLEASE DO NOT USE ANY OTHER FORM OF HAIR
REMOVAL (Waxing/Epilation etc).

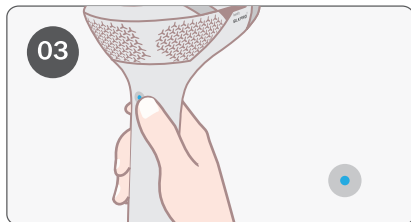

Plug in your Laser Hair removal device, and double
click the power button to start the device. The
device is available to be used immediately.

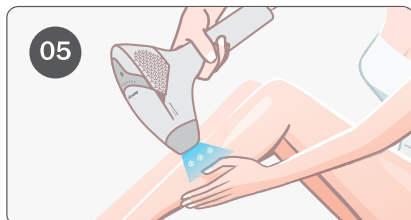


To commence treatment, press the device firmly
against the skin. Move the device slowly following
each individual flash of the laser across the treat-
ment area.

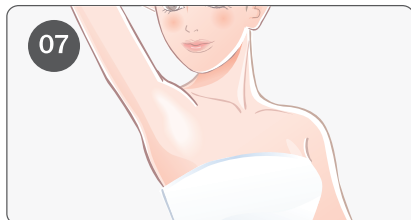


At the maximum power Level 6, a whole body
treatment only takes 15 minutes using the device.


* For reference only

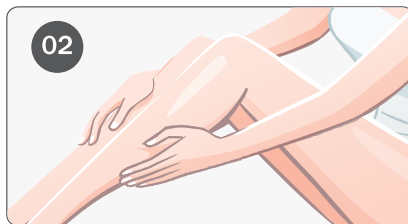


After shaving, keep the treatment area clean and
dry, and no residue of any creams or lotions
remain on the skin.

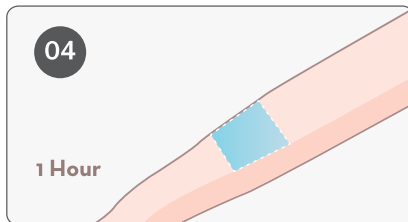

Before your first treatment, please conduct a
“patch test” on a small skin area and wait for an
hour to see if you have any adverse reaction to the
Laser Hair removal technology.

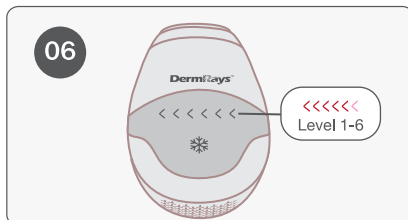

The Laser Hair removal device is automatically set
at the lowest (Level 1) for the gentlest treatment.
Press the power button once can increase one
level, with a maximum level of 6. A higher level of
power reduces the frequency of laser flashes,
which means slower treatments.

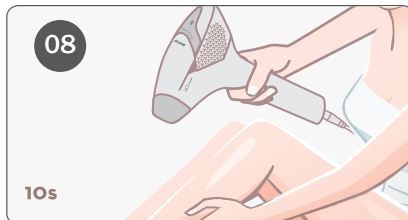

For safety purposes, if the device loses contact
with the skin for over 10 seconds, it will automati-
cally enter Standby mode. Simply press the power
button once to re-start your device.

* For more information, please refer to User Manual.