US	01-25
DE	26-50
ES	51-75
FR	76-100

US

✓STATEMENT

Thank you for purchasing this product, and in order to ensure proper use of this product, please read the user manual carefully before using it and keep this manual where you can refer to it at any time.

✓HISTORY OF LED LIGHT THERAPY

LED therapy began in 1903 when a scientist named Niels Finsen was awarded the Nobel Prize in Medicine for his use of light therapy. The field has expanded, and more clinics and the medical science industry are putting the effects of phototherapy to work with greater impact. Through extensive research and years of use of red light therapy and infrared therapy, it was found that not only make the tissues healed, but also make the condition of the skin surface was improved.

✓PRODUCT PRINCIPLE

LED red light therapy pad is through special LEDs, 660nm red light wavelength and 850nm infrared light wavelength, which has optical effects on organisms and produces important biological effects and phototherapy effects. Mitochondria in cells absorb red light the most, after the irradiation of red light, the catalase (CAT) activity of mitochondria increases, which can increase cell metabolism, increase glycogen content, increase protein synthesis and increase the decomposition of adenoglycine triphosphate (ATP), thereby strengthening the regeneration of cells, promoting the healing of ulcers, but also increasing the phagocytosis of white blood cells, improving the body's immune function, improving blood circulation, promoting local tissue metabolism, promoting hair growth, and promoting fracture healing, Accelerates the regeneration of damaged nerves, etc.

✓PRODUCT ADVANTAGES

1. The product can be used for phototherapy of the trunk and limbs (Which parts of the body the product can be used on is determined by the characteristics of the product)., as well as joint soft tissue health care and physiotherapy after injury, which can improve blood circulation and human metabolism to eliminate the pain of the injured part, so that the local damage of the human body can be repaired quickly, and the body's immune function can be effectively improved;
2. The technology of the product is very advanced, using LED red light module technology, combined with microelectronic heat compress, infrared and LED red light and other technologies in the performance of treatment, is a highly technological product, is a safe, effective, reliable, painless, non-toxic, natural body pain relief solution;
3. This product is a unique and high-end gift for your lover, parents, the elderly, friends, etc. on birthday, anniversary, Father's Day,

Mother's Day, Valentine's Day, Thanksgiving, Christmas and New Year.

☑SYMPTOMS ADAPTION

It is suitable for the adjuvant treatment of trunk and limb soft tissue injuries, acute and chronic shoulder, neck, back and leg pain. It has anti-inflammatory, anti-swelling, analgesic, improves microcirculation, promotes blood circulation, Anti-aging and cosmetic effect.

☑INTRODUCTION TO LIGHT SOURCES AND WAVELENGTHS

Color	wavelength	Function (different wavelengths, different skin depths)
Blue light	460nm	It is an extremely effective light source for activating cellular activity, and energy is available to many hair follicles and sebaceous glands. And can stimulate macrophages to release cytokines, thereby affecting the formation of inflammation, and play a certain healing and repair, anti-inflammatory, sedative, and have a certain effect on congestion and pain relief.
Red light	640nm	It has the characteristics of high purity, strong light source and uniform energy density. It plays an important role in skin care and recovery and is known as bioactive light. Red light can increase cell activity, promote cell metabolism, promote blood circulation, make the skin secrete a large amount of collagen and fibrous tissue to fill itself, accelerate blood circulation, increase skin elasticity, improve skin chlorosis and dullness, so as to achieve anti-aging, anti-aging oxidation and repair.
	660nm	
Infrared	850nm	It is invisible to the naked eye and penetrates deep into human tissues, which can expand local blood vessels, accelerate blood circulation, increase local arterial blood supply, help provide more energy and nutritional support, promote the healing of lesions, and relieve muscle tension, fatigue and strain.

☑ADAPTATIONS

1. Home; 2. Beauty salon; 3. Clinic; 4. Gym; 5. Sauna; 6. Physiotherapy center.

☑CONTRAINDICATION

1. Children under 13 years old and pregnant women should not use phototherapy products;
2. Cancer patients;
3. Tuberculosis patients;
4. Acute purulent inflammation;
5. Patients with bleeding;
6. Patients with heart disease and pacemakers;
7. Patients with arteriosclerosis;
8. Newborns and infants;
9. Patients allergic to light;
10. People who are deemed inappropriate to use these products by doctors.

WARNING

1. Do not irradiate the eyes;
2. Don't let children play, it's not a toy;
3. Do not contact with water to prevent short circuit leakage; Not waterproof, only for indoor use;
4. If other equipment such as wires or plugs is damaged, please do not use it, only use the equipment accessories provided;
5. This product contains SBR, if you are allergic to SBR, please do not use this product;
6. Do not use this product in sleeping, comatose or drunk people;
7. Do not wrinkle or fold the device, and need to roll it into a round shape when storing;
8. Do not use this equipment in open wounds, bleeding caused by increased heat and accelerated blood circulation;
9. When not in use, unplug the power supply to avoid tripping hazards;
10. If you feel unwell during use, discontinue use immediately and contact your doctor.


WARM PROMPT

1. After phototherapy, the skin will have a few traces of lamp beads, which usually disappear automatically within 30-60 minutes;
2. It is recommended to turn on infrared and red light at the same time;
3. The temperature is recommended to be gradually increased from 1 gear according to personal preferences and acceptance, and there are 3 gears;
4. Put the infrared pad flat, can not be cushioned on soft objects, the product has many FPC board and lamp beads composition, easy to make the therapy pad bend, damage the FPC board, resulting the therapy pad can not be used;
5. We recommend to enjoy the therapy with naked body to get better effect;
6. Turn on the therapy pad and warm it up for 2-3 minutes before use.

Controller Instruction

Power Connection Method

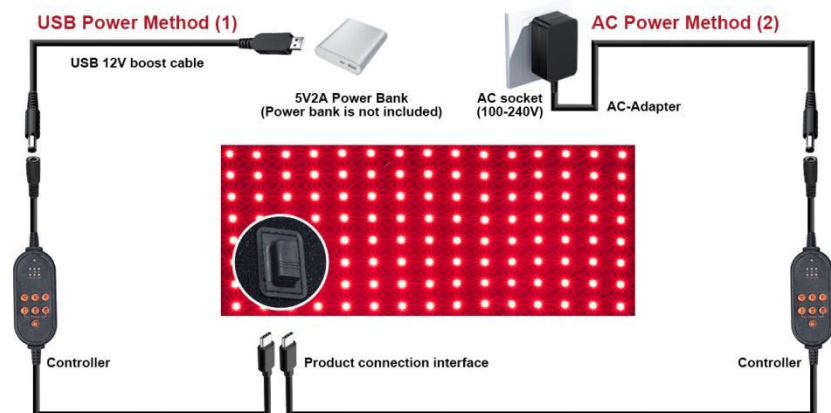INSTRUCTIONS FOR USE

For model numbers

CM-RED-M29; CM-RED-M39; CM-RED-M49; CM-RED-M59; CM-RED-M69;
CM-RED-M79; CM-RED-M89; CM-RED-C18

1. Press the power switch, the indicator light is on.
2. Wavelength selection: press the “R” key once, the LED power indicator will display green, and both red and infrared light will be on; press the “R” key twice, only the near-infrared light will be on (not visible to the naked eye, slightly glowing in dark light), the red light is off; press the “R” key 3 times, only the red light is on, and the near-infrared light is off; press the “R” key 4 times, all the LED lights are off. (Red light and infrared light can be selected according to individual needs. It is recommended to use red light and infrared light at the same time).
3. Temperature selection: the default LED power indicator will display green, and the heating temperature of the product is $39^{\circ}\text{C} \pm 2$; press the “C” key once, the LED power indicator will display blue, and the heating temperature of the product is $43^{\circ}\text{C} \pm 2$; press the “C” key twice, the LED power indicator shows red, and the heating temperature of the product is $47^{\circ}\text{C} \pm 2$. (It is recommended to choose the default temperature of $39^{\circ}\text{C} \pm 2$, and gradually increase it to a suitable temperature according to individual needs).
4. Pulse selection: Press the “R” key to set the wavelength mode, and then press the “P” key to start the corresponding wavelength pulse mode. (The pulse mode is selected according to individual needs).
5. Timing selection: press the timing button, press “I” for 10 minutes, the 10min indicator light is on, “II” is set for 20 minutes, the 20min indicator light is on, press “III” for 30 minutes, the 30min indicator light is on, the timer defaults to start for 30 minutes. Timing setting includes red light, infrared, temperature, pulse all modes. And product will auto shut off when the time is up.
6. Usage time: It is recommended to use it 1 to 2 times a day, no more than 30 minutes for each body part, and use it continuously for at least 4 weeks.
7. In the process of using red and infrared light, if you feel any discomfort, you should stop using it immediately.
8. The removable power bank pocket (only for some products) can be tied to the arm or pasted on the product in the direction of the Velcro, which is convenient for using the power bank during activities.
9. Please note: all plug connectors must be inserted tightly during use to prevent poor contact.
10. Do not fold during storage and use to avoid damage to the lamp beads.
11. There may be traces of light beads on the skin surface after use, which usually disappear automatically within 20-50 minutes.

US

CM-RED-M29

PRODUCT PARAMETERS

Model No	CM-RED-M29
Product Name	Red and infrared light foot pads
Voltage	AC 100-240V, DC 12V (can be connected to power bank or AC power)
Power	25W
LED QTY	150PCS
Wavelength	Red light 640nm*1; Red light 660nm*1; Infrared 880nm*1
Material	SBR
Size	760*130mm
Package contents	Foot pad*1; AC power adapter*1; Control switch*1; Elastic band*3; USB boost cable*1; Power bank * 1; User manual*1.
Function	Infrared, red light, pulse, heating, timing
Usage count	1-2 times a day, no more than 30 minutes per body part, at least 4 weeks.

US

CM-RED-M39

PRODUCT PARAMETERS

Model No	CM-RED-M39
Product Name	Red light & infrared light shoes
Voltage	AC 100-240V, DC 12V (can be connected to power bank or AC power)
Power	28W
LED QTY	158PCS
Wavelength	Red light 660nm*1; infrared light 850nm*2.
Material	SBR
Size	350*220mm
Package contents	Shoe*1 (not 1 pair); AC power adapter*1; Control switch*1; USB boost cable*1; Power bank bag*1; User manual*1.
Function	Infrared, red light, pulse, heating, timing
Usage count	1-2 times a day, no more than 30 minutes per body part, at least 4 weeks.

US

CM-RED-M49

US**CM-RED-M49**

AC CHARGING MODE


PRODUCT PARAMETERS

Model No	CM-RED-M49
Product Name	Red and infrared light full body pad
Voltage	AC 100-240V, working voltage 36V
Power	160W
LED QTY	960PCS
Wavelength	Red light 660nm*1; infrared light 850nm*2.
Material	SBR
Size	1700*400mm
Package contents	Mat*1; Control Switch*1; AC Power Adapter*1; Power Cord*1; User Manual*1
Function	Infrared, red light, pulse, heating, timing
Usage count	1-2 times a day, no more than 30 minutes each time, at least 4 weeks.

US

CM-RED-M59

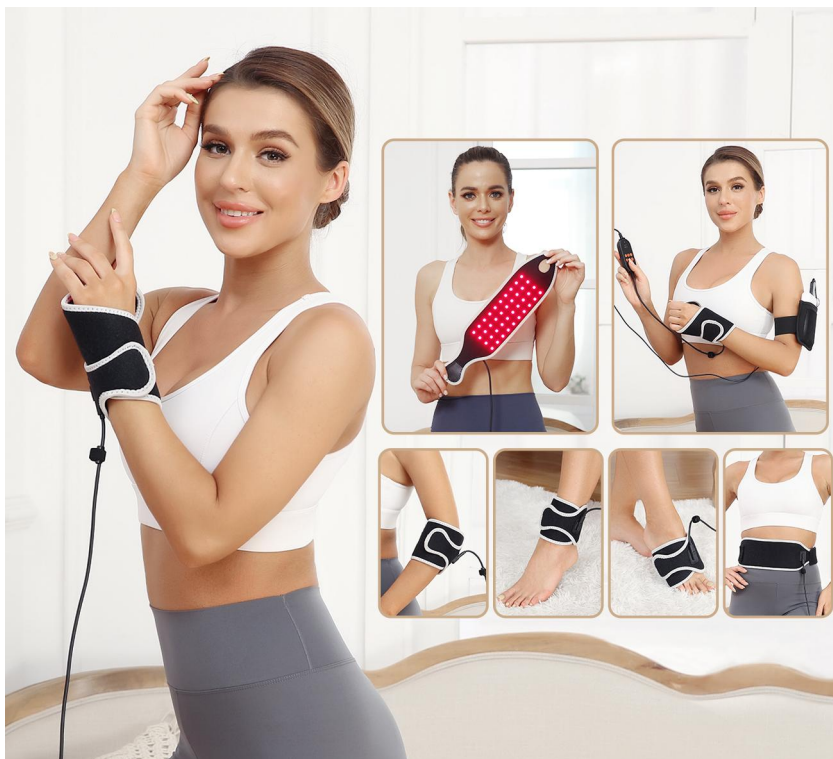PRODUCT PARAMETERS

Model No	CM-RED-M59
Product Name	Red light & infrared hand wrist band
Voltage	AC 100-240V, DC 12V (can be connected to power bank or AC power)
Power	9W
LED QTY	40PCS
Wavelength	Red light 660nm*1; infrared light 850nm*2.
Material	SBR
Size	100*400mm
Package contents	Hand Wrist band*1; AC power adapter*1; Control switch*1; USB boost cable*1;Elastic band*1;Power bank bag*1; User manual*1.
Function	Infrared, red light, pulse, heating, timing
Usage count	1-2 times a day, no more than 30 minutes per body part, at least 4 weeks.

US

CM-RED-M69

PRODUCT PARAMETERS

Model No	CM-RED-M69
Product Name	Red light & infrared foot ankle strap
Voltage	AC 100-240V, DC 12V (can be connected to power bank or AC power)
Power	15W
LED QTY	66PCS
Wavelength	Red light 660nm*1; infrared light 850nm*2.
Material	SBR
Size	320*340mm
Package contents	Foot ankle strap*1; AC power adapter*1; Control switch*1; USB boost cable*1; Power bank bag*1; User manual*1.
Function	Infrared, red light, pulse, heating, timing
Usage count	1-2 times a day, no more than 30 minutes per body part, at least 4 weeks.

US

CM-RED-M79

PRODUCT PARAMETERS

Model No	CM-RED-M79
Product Name	Red light & infrared gloves
Voltage	AC 100-240V, DC 12V (can be connected to power bank or AC power)
Power	19W
LED QTY	100PCS
Wavelength	Red light 660nm*1; infrared light 850nm*2.
Material	SBR
Size	170*320mm
Package contents	Glove*1(not 1 pair); AC power adapter*1; Control switch*1; USB boost cable*1; Elastic band*2; Power bank bag*1; User manual*1.
Function	Infrared, red light, pulse, heating, timing
Usage count	1-2 times a day, no more than 30 minutes per body part, at least 4 weeks.

US

CM-RED-M89

PRODUCT PARAMETERS

Model No	CM-RED-M89
Product Name	Red and infrared light belt
Voltage	AC 100-240V, DC 12V (can be connected to power bank or AC power)
Power	26W
LED QTY	150PCS
Wavelength	Red light 660nm*1; infrared light 850nm*2.
Material	SBR
Size	450*200mm
Package contents	Belt*1; AC power adapter*1; Control switch*1; USB boost cable*1; Elastic band*1; Power bank bag*1; User manual*1.
Function	Infrared, red light, pulse, heating, timing
Usage count	1-2 times a day, no more than 30 minutes per body part, at least 4 weeks.

US**CM-RED-C18****PRODUCT PARAMETERS**

Model No	CM-RED-c18
Product Name	Red and infrared light knee pad
Voltage	AC 100-240V, DC 12V (can be connected to power bank or AC power)
Power	16W
LED QTY	76PCS
Wavelength	Red light 660nm*1; Infrared 850nm*2
Material	SBR
Size	830*250mm
Package contents	Knee pad*1; AC power adapter*1; Control switch*1; Elastic band*1; USB boost cable*1; Power bank * 1;User manual*1.
Function	Infrared, red light, pulse, heating, timing
Usage count	1-2 times a day, no more than 30 minutes per body part, at least 4 weeks.

US

CM-RED-M19

PRODUCT PARAMETERS

Model No	CM-RED-M19
Product Name	Red and infrared light hat
Voltage	AC 100-240V, DC 12V (can be connected to power bank or AC power)
Power	25W
LED QTY	144PCS
Wavelength	Blue light 460nm*1; Red light 660nm*1; Infrared light 850nm*1
Material	SBR
Size	300*245mm
Package contents	Hat*1; AC Power Adapter*1; Control Switch*1; USB Boost Cable*1; Power Bank Bag*1; User Manual*1.
Function	Blue light, infrared light, red light, pulse, heating, timing
Usage count	1-2 times a day, no more than 30 minutes per body part, at least 12 weeks.

☑ INSTRUCTIONS FOR USE

1. Press the power switch, the indicator light is on.
2. Wavelength selection: press the "R" key once, the LED power indicator will display green, and the blue light will be all on; press the "R" key twice, the infrared and red lights will all light up; press the "R" key three times, the blue light, red light/infrared light turns on alternately; press the "R" key 4 times, all the LED lights are off. (Red light and infrared light can be selected according to individual needs. It is recommended to use red light and infrared light at the same time).
3. Temperature selection: the default LED power indicator will display green, and the heating temperature of the product is $39^{\circ}\text{C} \pm 2$; press the "C" key once, the LED power indicator will display blue, and the heating temperature of the product is $43^{\circ}\text{C} \pm 2$; press the "C" key twice, the LED power indicator shows red, and the heating temperature of the product is $47^{\circ}\text{C} \pm 2$. (It is recommended to choose the default temperature of $39^{\circ}\text{C} \pm 2$, and gradually increase it to a suitable temperature according to individual needs).
4. Pulse selection: Press the "R" key to set the wavelength mode, and then press the "P" key to start the corresponding wavelength pulse mode. (The pulse mode is selected according to individual needs).
5. Timing selection: press the timing button, press "I" for 10 minutes, the 10min indicator light is on, "II" is set for 20 minutes, the 20min indicator light is on, press "III" for 30 minutes, the 30min indicator light is on, the timer defaults to start for 30 minutes. Timing setting includes red light, infrared, temperature, pulse all modes. And product will auto shut off when the time is up.
6. Usage time: It is recommended to use it 1-2 times a day, no more than 30 minutes for each body part, and use it continuously for more than 4 weeks.
7. In the process of using red and infrared light, if you feel any discomfort, you should stop using it immediately.
8. The removable power bank pocket (only for some products) can be tied to the arm or pasted on the product in the direction of the Velcro, which is convenient for using the power bank during activities.
9. Please note: all plug connectors must be inserted tightly during use to prevent poor contact.
10. Do not fold during storage and use to avoid damage to the lamp beads.
11. There may be traces of light beads on the skin surface after use, which usually disappear automatically within 20-50 minutes.

US

CM-RED-C28

PRODUCT PARAMETERS

Model No	CM-RED-C28
Product Name	Red &Infrared Light Space Therapy Pad
Input voltage	110-240V (50-60Hz)
Output voltage	DC36V
Power	upper pad 270W, lower pad 280W
Current	7.8A*1,7.8A*1
LED QTY	1392PCS on the upper pad; 1440PCS on the lower pad
Wavelength	LED red light 660nm*1, LED infrared 850nm*2
Optical Density	1st gear: 120mw/cm ² ; 2nd gear: 135MW/cm ² ; 3rd gear: 150MW/cm ²
Temperature	109-131°F (43-55°C)
Timing range	0-60 minutes (Default time: 30 minutes)
Material	Therapy mat: SBR, Protective case: TPU + nylon cloth
Size	1700mm*600mm
Package contents	Therapy mat*2PCS (upper and lower);Protective case* 2PCS (up and down);remote control *1; Adapter*2; Goggles*1; Power cord *2; User manual *1; Storage bag*1; AAA battery*2
Function	Near infrared + red light + timing + pulse

US**CM-RED-C28**

✓HOW TO USE THE PROTECTIVE CASE

1. Place the protective case flat, one corner of the protective case with an opening, here is the power connection interface;
2. Make sure power supply interface of the therapy pad is coincide with the protective case opening before putting the therapy pad in the protective case;
3. Put the therapy pad completely flat into the protective case, and check whether the 4 corners are all aligned, and there can be no overlap part on the pad, so as not to damage the internal FPC board;
4. Zip the zipper after there are no overlap part on the therapy pad;
5. The protective case made in nylon, which is durable, wear-resistant, dirt-resistant, flat and easy to clean;
6. The transparent case is made in TPU material which is environmentally friendly, non-toxic, non-allergic, non-irritating, wear-resistant, good light transmission, and does not affect the penetration of light waves.

Protective Case Usage Steps

US

CM-RED-C28

Introduction of Remote Control

- ✓ A/B: A means control button for upper therapy pad.B means control button for Lower therapy pad
- ✓ Please install the 2*AAA battery before using the remote control, and pay attention to the positive and negative poles direction when installing.
- ✓ NOTE: Each LED beads contains 2pcs near infrared. Each therapy pad with 2 control buttons to turn on or off the 2pcs near infrared separately,or turn on them at the same time.

A/B	"A" on the remote control is the upper therapy pad control button; "B" is the lower therapy pad control button.
	Power Switch and Pulse Button Combination Power Switch: long press for 3 seconds to turn on the power supply, at the same time the button and the screen light on. On the contrary, long press for 3 seconds to turn off the power. Pulse Button: short press the button to start the pulse function. The button lights up after pressing, and the therapy pad starts the strobe function, and will see the light is flicking.
	The near infrared buttons, each LED beads has 2pcs 850nm infrared, a button to control a 850nm infrared alone, if you need turn on 2pcs 850nm at the same time, you need to turn on 2 buttons, the button lights on. The near infrared light is visible in the low light or dark light.
	The red light button, each LED beads has a 660nm red light, turn on the button, red light and the button light on at the same time, the therapy pad starts working.
	The power button, there are 3 gears, each pressing is a gear. The button lights up after pressing, the green indicator on the side of the display is the 1st gear, the blue indicator is the 2nd gears, and the red indicator is the 3rd gears. At the same time, you can obviously observe the brightness change, higher power, stronger light. According to personal preferences to choose different power.
	The time plus button, for timing the upper and lower therapy pads at the same time, add 5 minutes with each pressing. The button lights on after pressing and the display shows the time. Time is up, the device automatically stops working.
	The time minus button, for timing the upper and lower therapy pads at the same time, reduce 5 minutes with each pressing, the button lights on after pressing and the display shows the time. Time is up, the device automatically stops working.
1. Only one of the remote control mode and manual mode can be selected, and cannot be used at the same time. 2. Please install the 2*AAA battery before using the remote control, and pay attention to the positive and negative poles direction when installing.	

US**CM-RED-C28****☑MANUAL STANDBY MODE OPERATION METHOD**

1. Press the switch 1 time, the switch indicator is green: 1 group of red light + 2 groups of infrared (visible in the dark) work, the default low power.	
2. Press the switch 2 times, the switch indicator is green: 1 group of infrared work.	
3. Press the switch 3 times, the switch indicator is green: 1 group of red light work.	
4. Press the switch 4 times, the switch indicator is green: another group of infrared work.	
5. Press the switch 5 times, the switch indicator does not light up, and all LED lights stop working.	
6. Any of the above modes long press for 3 seconds, high power mode, with blue light indicator, and manual mode no mid-range power.	
7. Timing defaults to 30 minutes.	
8. Only one of the remote control mode and manual mode can be selected, and cannot be used at the same time.	

☑HOW TO USE THE PRODUCT

1. Put the lower therapy pad in the protective case, place it on hard objects, and never place it on the soft objects (such as soft Simmons, latex pad, sofa, etc.);
2. Plug in the power supply and check whether all interfaces are plugged in tightly;
3. Cover the upper therapy pad, which can be fixed according to personal habits;
4. The gap on the 2 sides of the upper and lower pads is specially designed for heat dissipation;
5. Wear the goggles;
6. Operate according to the instruction of the remote control;
7. Start enjoying your phototherapy moments;
8. It is recommended to use it once a day, 30-60 minutes each time, at least for more than 1 week.

US**CM-RED-C38**

PRODUCT PARAMETERS

Model No	CM-RED-C38
Product Name	Red and infrared light full body pad
Input voltage	110-240V (50-60Hz)
Output voltage	DC36V
Power	280W
Current	7.8A
LED QTY	1440PCS
Wavelength	LED red light 660nm*1, LED infrared 850nm*2
Optical Density	120mw/cm ² -150MW/cm ²
Temperature	118-132°F (48-56°C)
Timing range	0-30 minutes (Default time: 15 minutes)
Material	Therapy mat: SBR, Protective case: TPU + nylon cloth
Size	1700mm*600mm
Package contents	Therapy mat:*1PCS ;Protective case* 1PCS;remote control *1; Adapter*1; Goggles*1; Power cord *1; User manual *1; Storage bag*1; AAA battery*2
Function	Near infrared + red light + timing + pulse+Remote control

☑ HOW TO USE THE PROTECTIVE CASE

1. Place the protective case flat, one corner of the protective case with an opening, here is the power connection interface;
2. Make sure power supply interface of the therapy pad is coincide with the protective case opening before putting the therapy pad in the protective case;
3. Put the therapy pad completely flat into the protective case, and check whether the 4 corners are all aligned, and there can be no overlap part on the pad, so as not to damage the internal FPC board;
4. Zip the zipper after there are no overlap part on the therapy pad;
5. The protective case made in nylon, which is durable, wear-resistant, dirt-resistant, flat and easy to clean;
6. The transparent case is made in TPU material which is environmentally friendly, non-toxic, non-allergic, non-irritating, wear-resistant, good light transmission, and does not affect the penetration of light waves.

Protective Case Usage Steps

US

CM-RED-C38

✓MANUAL STANDBY MODE OPERATION METHOD

1. Press the switch 1 time, the switch indicator is green: 1 group of red light + 2 groups of infrared (visible in the dark) work, the default low power.	
2. Press the switch 2 times, the switch indicator is green: 1 group of infrared work.	
3. Press the switch 3 times, the switch indicator is green: 1 group of red light work.	
4. Press the switch 4 times, the switch indicator is green: another group of infrared work.	
5. Press the switch 5 times, the switch indicator does not light up, and all LED lights stop working.	
6. Under any of the above modes long press for 3 seconds, switch to the high temperature mode with blue light indicator. Manual mode without middle temperature mode (The default is low temperature mode with green light indicator).	
7. Timing defaults to 30 minutes.	
8. Only one of the remote control mode and manual mode can be selected, and cannot be used at the same time.	

Introduction of Remote Control

1. Only the remote control mode or manual mode can be selected at a time, and cannot be used 2 control modes at the same time.
2. Please install the 2*AAA battery before using the remote control, and pay attention to the positive and negative poles direction when installing.

US**CM-RED-C38****✓HOW TO USE THE REMOTE CONTROL**

	Power switch, long press for 3 seconds to turn on the power, while the button and display light up. Instead, long press for 3 seconds to turn off the power.
	The near infrared buttons, each LED beads has 2pcs 850nm infrared, a button to control a 850nm infrared alone, if you need turn on 2pcs 850nm at the same time, you need to turn on 2 buttons at the same time, the button lights on. The near infrared light is visible in the low light or dark light.
	The red light button, each LED beads has a 660nm red light, turn on the button, red light and the button light on at the same time, the therapy pad starts working.
	Temperature button, button light up after pressing. Level 1 (Low: 48°C/118°F), level 2 (Medium: 52°C/125°F), level 3 (High: 56°C/132°F). On the left side of the display, the green light means low gear, the blue light means medium gear, and the red light means high gear. At the same time, brightness changes can be clearly observed, higher temperatures and stronger light. Choose a suitable temperature according to your personal preference.
	The pulse function starts after short pressing the button, and the button P1, P2, P3 represents the flashing speed of low, medium and high. There are pulse indicator light on the right side of the display. Green light indicates P1 low gear. Blue light indicates P2 mid-gear. Red light indicates P3 high gear. According to personal preference, choose the suitable pulse gear. You will see the light flashing when the light therapy mat activates the pulse function.
	The time plus or minus button increases or decreases by 5 minutes each time it is pressed. When pressed, the button lights up and the display shows the time. When the time is up, the product will automatically stop working.
1. Only one of the remote control mode and manual mode can be selected, and cannot be used at the same time. 2. Please install the 2*AAA battery before using the remote control, and pay attention to the positive and negative poles direction when installing.	

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China·Dongguan Everest Information Technology Co., Ltd
Name: Sea&Mew Consulting GmbH