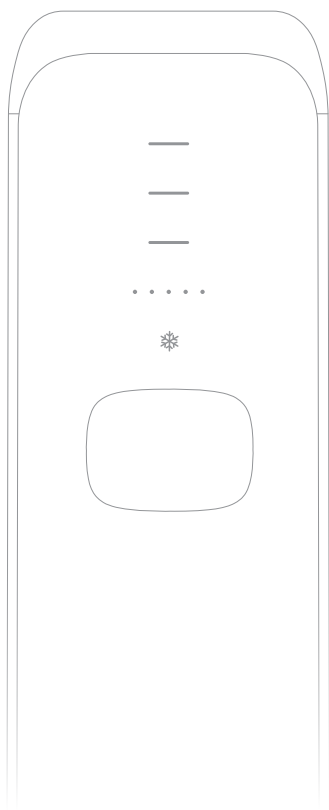
# Ulike

## Air3 IPL Hair Removal Device User Manual

Air Series

ENGLISH ————— 02-17

DEUTSCH ————— 20-35

FRANÇAIS ————— 38-53

ITALIANO ————— 56-70

ESPAÑOL ————— 72-87

Before using the device, please read the user manual carefully, including all contraindications, warnings and safety information, and keep them for future reference. Illustrations in this manual are for reference only.

**Thank you for choosing Ulike.**

**Founded in 2013, Ulike has focused on at-home IPL hair removal devices from day one.**

**We believe hair removal should be a safe, smooth, and pain-free experience.**

**And this device is going to be the hair removal system you like.**

Any questions or concerns about the products, please do not hesitate to contact us via site or email - we are here to assist you!

[de.ulike.com](http://de.ulike.com)

[uk.ulike.com](http://uk.ulike.com)

[de.support@ulikeglobal.com](mailto:de.support@ulikeglobal.com)

[uk.support@ulikeglobal.com](mailto:uk.support@ulikeglobal.com)

Shenzhen Ulike Smart Electronics Co., Ltd.  
810, Building 1, Xunmei Science and Technology Plaza, No. 8 Keyuan Road,  
Science Park Community, Yuehai Sub-District, Nanshan District, Shenzhen  
518000, Guangdong, P.R. China

## Safety Precautions

### ■ ATTENTION

- Please assemble and use the product precisely according to the instructions provided in this manual.
- Please read the User Manual carefully and keep this document for future reference.
- If you have any questions regarding the use of this product, please contact Ulike customer service for assistance.

⊗ **DANGER: Personal injury due to incorrect use**

ⓘ **WARNING: Danger due to incorrect use**

⚠ **CAUTION: Personal safety hazard or property loss due to incorrect use**

⊗ **DANGER**

- Do not use the device directly on your eyes; the strong light emitted by the product may harm the eyes.
- While the device is turned on and in use, do not directly touch the side ventilation outlet; the temperature of the outlet may scald skin.
- Do not insert any objects into the device; doing so may result in harm to you.
- No data regarding this product's use by children under 18 years of age is currently available. Children below 18 should not be allowed to touch or use it. Please keep this product away from children under 18 years of age.
- The product should not be used by people with physical, sensory or mental impairments, or lack of experience and knowledge of its use, unless supervised or directed by a person responsible for their safety.

ⓘ **WARNING**

- Please do not apply antiperspirant immediately after treating the armpit area;
- To avoid scalding and mechanical failure, please regularly clean the light window;
- Please do not place the light window near items that may easily cause demagnetization; otherwise, failure and damage may result.
- Please do not use the product when the ventilation outlet is blocked; otherwise, the device may cause smoke or fire.
- Please do not point the light window at black objects.
- To avoid blackening, discoloration, and breakage, please do not use the device towards clothes, hair, or other items.
- Please do not use the product on pets, as this may cause them harm.

- Eye (blindness) or skin damage may result if you don't follow the User Manual. Please carefully read and follow the User Manual.
- Please do not overuse this product or reuse on one area; doing so may harm your skin.
- Please protect your eyes. Do not use the product above the eyelids or near the eyes.
- Please test on your skin before use and check the area after 2 hours. You can use the product if no abnormal conditions (i.e. redness or blisters) occur in the irradiated area.
- Please wear safety glasses (provided) during the entire process. Keep the light window close to skin so that your eyes won't be irritated by strong light.
- After use, please wipe off the hair, dirt, and oily residue on the device to avoid mechanical failure and skin damage.
- Please set the intensity level according to your personal skin type and condition. Be sure to start with Gentle Mode the first time.
- It is normal if skin is slightly hot and red after use.
- The product cannot be left unattended while turned on, Please turn it off when not in use.
- If the operating environment temperature difference reaches more than 20 degrees (from low temperature to high temperature, or vice versa), please turn the device off and wait for about 2 hours before using it again.
- Please do not use the product near inflammable and explosive objects (i.e. sprayers, gas stations) which may cause fire, electric shock, et cetera.
- The product should be used at 5-30°C. Please check room temperature before use.
- Please avoid direct ultraviolet radiation on the irradiated area within 48 hours after using this product.
- Please do not expose the product to sunlight or UV light for an extended period of time; serious damage will result.
- If the built-in fan does not function after the product is turned on (there is no humming sound), please do not use the device; continued use in this case may result in smoke, fire, etc.
- It is not recommended to use the product for long stretches of time. Please turn the device off for 10 minutes after 1 hour of use.
- If the casing, lens, or power supply is damaged, please stop using the product immediately. Please do not attempt to repair the product without official Ulike authorization. The device must be repaired by professional personnel from the manufacturer's maintenance department or a similar department.
- This product is not suitable for everyone. Please do not use this product if you match any of the descriptions below.
- Do not use it in the bath, shower or other places near the water source: Please keep the product and the power adapter dry.

## ■ Unsuitable population and situation

- The product may cause discomfort for, or adverse effects on, people with allergies, skin diseases, allergic dermatitis, sensitive skin, etc. If any adverse effects occur during use, please stop using the product immediately.
- Please do not use the product on areas where skin integrity is damaged (in areas with burns, blisters, cuts, wounds, active skin diseases, or infections).
- If you are taking light-sensitive drugs, please seek medical advice before using this product; otherwise, using this device may damage your skin.
- The product is not suitable for pregnant, lactating, or menstruating women, or people with scars, diabetes, or fragile skin.
- People with severe heart disease, moderate to severe hypertension, malignant tumors, folliculitis, infectious skin disease, burned skin, and anticoagulant disorders are prohibited from using this product.
- Please do not use if you have abnormal collagen metabolism, including scars and poor wound healing.
- Please do not use it if you have vascular diseases such as varicosity and hemangiectasis in the area to be irradiated.
- Please do not use if your skin is sensitive to light and prone to rashes or allergic reactions.
- Please do not use it if your irradiated area is infected or you have herpes simplex, skin lesions, or hematoma.
- Please do not use it if you have immune response suppression diseases (such as HIV infection, lupus erythematosus, porphyria, etc.).
- Patients with hemorrhage syndromes and those who are using isotretinoin, alpha hydroxyl acid(AHA), beta hydroxyl acid(BHA) or azalea acid are prohibited from using this product.
- Please use this product carefully if you had chemical hair removal done within 30 days or recently applied skin whitening products or cosmetics with vitamin A.
- The product cannot be used with essential oil products. Please ensure that the skin in the irradiation area is clean and dry without applied lotions, lubricants, or other cosmetics.
- If the skin is not clean and dry, this may negatively affect your experience with the product and even damage the device.
- People who are taking aspirin drugs, immunosuppressive drugs, isotretinoin and other drugs, and recently took light-sensitive drugs such as tetracycline and sulfa, are prohibited from using this product.
- People recently injected with fillers and toxins are prohibited from using this product.

**■ Please do not use the product on any of the following body parts or any body parts with the following conditions**

1/ Oral cavity, genitals, mucosa, etc.	6/ Acnes	11/ Inside the navel
2/ Plastic surgery or other surgical procedure	7/ Spots and moles	12/ Large vein
3/ Throat	8/ Eyebrow, hair	13/ Tattoo, body painting
4/ Eyeball, eyelid, around the eyes	9/ Ear	
5/ Head	10/ Nipple, areola	

(Note: This list is not exhaustive. If you are not sure whether you can or cannot use this product, please consult a doctor)

**■ After using the product, you may experience the following side effects:**

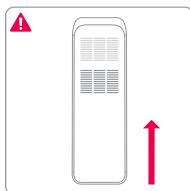
1/ Burns	3/ Swollen hair follicles and blisters	5/ Itchiness and pain
2/ Local redness and swelling	4/ Erythema and pigmentation	6/ Scarring

(Note: The list is not exhaustive. If you experience any discomfort, please consult Ulife customer service or a doctor)

**▲ NOTICE**

- 1.This product is not intended to treat diseases.
- 2.We advise against using this product to remove white or gray hair.
- 3.If you have naturally dark skin, please do not use this product.
- 4.Please do not expose yourself to sunlight before using this product.
- 5.On first use, please try this product in Gentle Mode on a small area of your arm to test your skin's reaction.
- 6.Make sure that the light window is close to the skin before beginning the hair removal process.
- 7.Do not use this product in dark conditions.
- 8.To avoid mechanical failure, please do not use the product in direct sunlight.
- 9.Please do not hit or violently shake this product.
- 10.Please turn off the device before unplugging it from an electrical outlet; doing otherwise may damage this product.

**Note:**When you finish using the product, please have the ventilation outlet facing upwards and the treatment window downwards for better heat dissipation and cooling.

## Product introduction

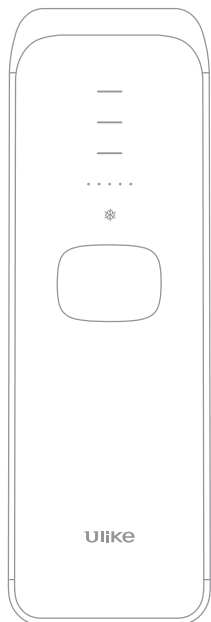
### ■ About the Ulike IPL Hair Removal Device

•The product, unlike concave-tip hot light hair removal devices, adopts the clinical sapphire contact cooling technology. Ulike's patented sapphire ice-cooling technology effectively cools the skin, reduces heat, avoids burns, and better protects the skin barrier. While removing hair, the sapphire flat tip compresses the skin, pressing down hair follicles to bring roots closer to the area with higher energy, thus improving the energy absorption efficiency of hair follicles. The sapphire flat tip is easy to clean and effectively transmits light to achieve better hair removal results.

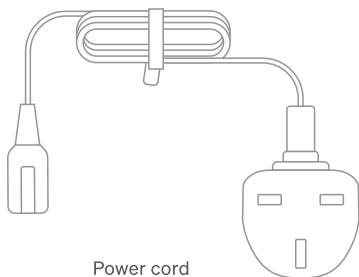
•This product's long lamp lifetime and versatility means that it does not need to be replaced and is suitable for multiple areas (such as the face, lips, underarms, bikini lines, arms, legs, etc.) to save expensive lamp head replacement costs. You don't need to buy accessories frequently or worry about additional fees. In addition, this product is designed to efficiently filter ineffective short wavelengths and harmful light with top-grade filter.

## ■ Product structure and accessories

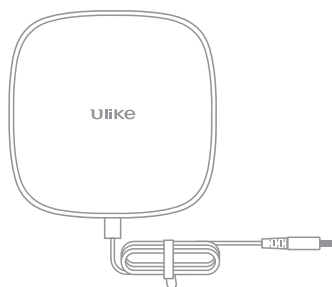
### Product composition

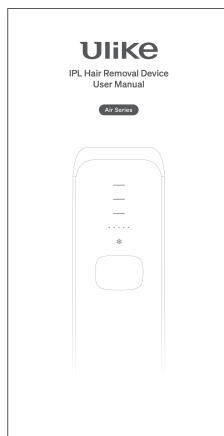
Host

Power cord

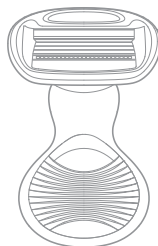
Power Adaptor

User manual

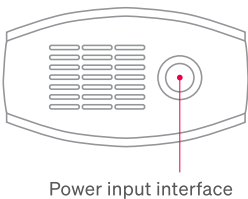
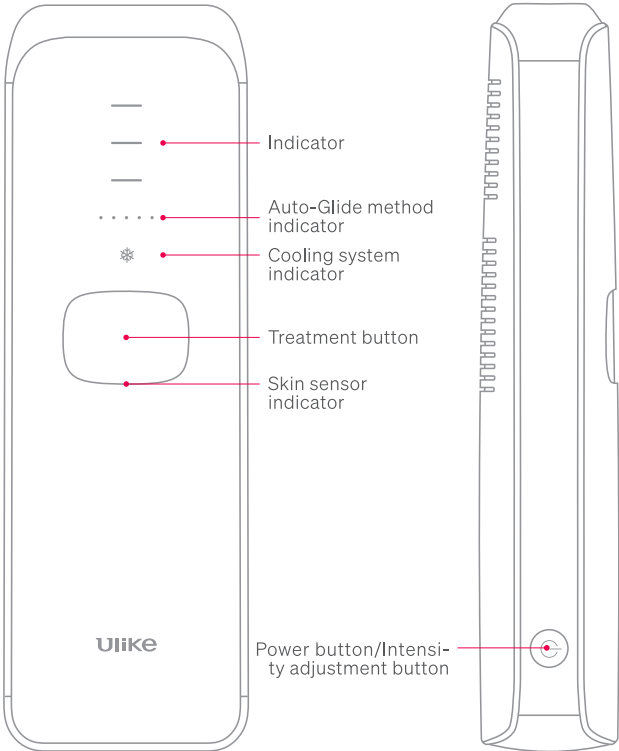
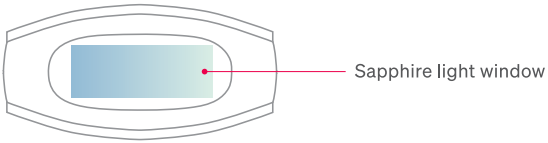
IPL Safety Glasses

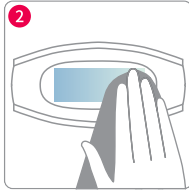
Razor

## ■ Product structure

## Before use

### ① Clean your skin thoroughly

Ensure that there is no hair or oily residue (perfume, lotion, deodorant, sun cream, etc.) left on your skin and it is completely dry.

### ② Make sure that the exposure port is free of dust, hair and dirt, and is not damaged or has no abnormal conditions

If there is dirt, please clean it up to avoid product failure, deterioration and skin problems.

### ③ Perform a skin test before using the depilator

Confirm your skin condition 2 hours after the skin test. If there is no abnormal signs (redness, blisters) at the irradiated site, you can use this product.

## ■ Set exposure intensity

This product offers 3 energy modes, which can meet the needs of different degrees of skin sensitivity and different body hair colors. The following are related instructions:

#### • Mode:

Soft Mode: Suitable for sensitive parts such as lips, people with dark skin and sensitive skin, etc.

Body Mode: Suitable for the body parts of people with thick hair in medium or low hardness, such as arms, legs, etc.

Strong Mode: Suitable for the body parts of people with thick and hard hair, such as arms, legs, etc. (Due to the higher power of Strong Mode, you may feel hotter than Soft Mode and Body Mode, but this is not product failure)

• Please refer to the table below to select a mode suitable for your skin nature and hair color:

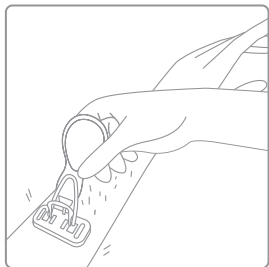
Skin color Hair color	White	Beige	Light brown	Medium brown	Dark brown	Brown black and darker
	White/Grey	×	×	×	×	×
Red	×	×	×	×	×	×
Light blond	×	×	×	×	×	×
Dark blond/light brown	regular mode power mode	regular mode power mode	regular mode power mode	soft mode regular mode	×	×
Brown	regular mode power mode	regular mode power mode	regular mode power mode	soft mode regular mode	×	×
Dark brown	regular mode power mode	regular mode power mode	regular mode power mode	soft mode regular mode	×	×
Black	regular mode power mode	regular mode power mode	regular mode power mode	soft mode regular mode	×	×

(Note: X means not applicable)

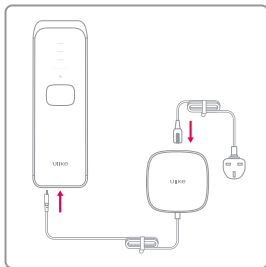
If you feel pain or discomfort during the use, please use a different mode (if you feel warm or slight burning without pain, it is normal), people with darker skin colors are prone to skin allergies, so it is recommended to use Soft Mode or Body Mode Please ensure you start with Soft Mode the first time.

## Usage

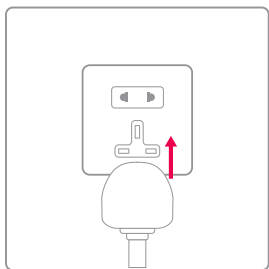
### ■ How to Use

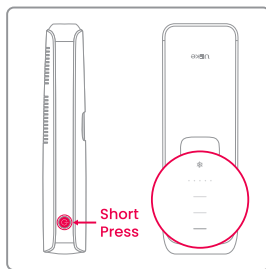
1. Use a scraper to remove hair from body areas in need of care. Please do not apply the scraper blade directly to the skin. Use with lather, or shave when taking a shower. The blade is sharp. Use it carefully to avoid scratches.

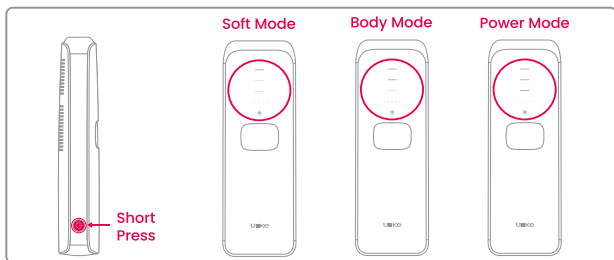
2. Connect the power chord to the adapter, plug the power chord into the power supply connection on the rear of the product

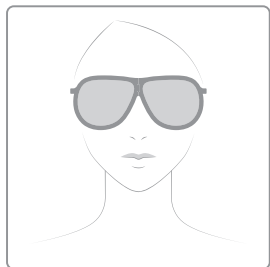
3. After connecting the power chord to the product's power supply connection, insert the plug into a power outlet.

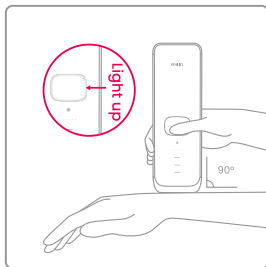
4. Short press Power On/Off on the side of the main unit to turn on the product, and the energy level and cooling function will be displayed when the product is turned on.

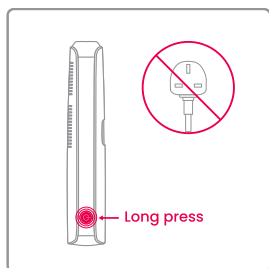
5. The product has three modes. Short press Power On/Off/Mode Switch on the side of the main unit to turn on the product and switch to the hair removal mode you need. The mode will change in the following order: Soft Mode - Body Mode - Power Mode - Body Mode - Soft Mode.

6. Wear protective glasses before removing hair to protect your eyes.

7. Place the light outlet on the area to be treated on the skin and fully fit the skin vertically. The white lamp of hair removal sensor lamp will light up, indicating that it has been fully contacted. Press the flash key and start the hair removal operation. If the white lamp is not light up, which indicating that lamp holder does not fully contact to the skin, and the device will not shine

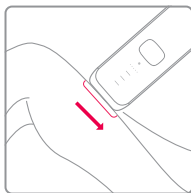
8. After finishing the hair removal process, long press the on/off button for 1.2s to turn the machine off (do not unplug it directly without turning off the machine to avoid damage.)

## ■ Treatment Methods

### •Auto-Glide Method (for large areas such as arms, legs):

Place the light outlet on the area to be treated and let it fully contact the skin. At this time, the hair removal sensor indicator light will light up, indicating that the lamp head fully contacts the skin. Long press the flashing button for 1.2 seconds and the continuous play indicator light will light up. Smart Continuous Play Mode will engage.

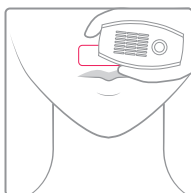
Please keep the light outlet contacting the skin. Flashes will be emitted at a certain interval. Please keep moving the product. Move it back and forth twice on the area to be treated to ensure that the area is covered evenly and avoid missing light.

Tip: You can adjust your movement speed according to the flashing speed, so that the product can cover the area more evenly. The flashing frequency is higher in Soft Mode and slightly lower in Body Mode. Furthermore, the higher the level, the slower the light will be emitted.

### •Spot Method (suitable for small areas such as fingers and upper lips):

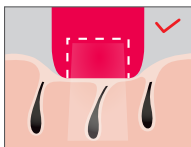
Place the light outlet on the area to be treated and let it fully contact the skin. The hair removal sensor indicator light will light up at this time, indicating that the lamp head fully contacts the skin. Short press the flashing button to start hair removal. Move the lamp head to the next area, wait for the white light of the hair removal sensor indicator light, and then press the flashing button again to repeat the above operation in order.

Tip: If you start to feel excessive heat on your skin during the hair removal process, pause for a few seconds to let the light outlet cool down again before using the product again

## ■ Tips for using Ulike Air Series

1. Due to continuous flashes, this product will heat up. In order to maintain an optimal cooling effect, we recommend you use the light window on the target area for about 30 seconds after hair removal before moving on to the next area.

2. For best results, press the light window tightly against your skin so that the lamp head is fully engaged. In this way the device can cool the skin in advance, reduce heat, and better protect the skin barrier.

## ■ Usage frequency

weeks 1-4  
3 times a week

from 5 weeks  
every half or one month

Use 3 times a week for weeks 1-4;

Use every half or one month from 5 weeks.

Note: Please use it multiple times according to the growth of the hair. Too thick hair can be appropriately increased in gear and frequency based on this, the specific situation varies from person to person.

## Care after using

Precautions to take after hair removal with this product

- We recommend you avoid having the treated area come into contact with water for 6 hours after hair removal. After that, you may moisturize the treatment area.
- Please do not use perfumes, facial masks, or cosmetics containing exfoliation beads on the treated area after hair removal.
- Please avoid exposure to ultraviolet light for 48 hours after hair removal. Apply sunscreen as necessary.
- Artificial tanning is prohibited while undergoing this hair removal process.
- If you require plastic surgery on the area being treated for hair removal, please consult a doctor.
- If you experience discomfort, please stop using the product and consult a doctor.

## ■ Cleaning and maintenance

For best results and to maximize the service life of the product, please clean the product before and after each use to avoid accumulation of hair and dust. Lack of proper maintenance of this product will negatively affect its functionality and effectiveness.

- Before cleaning the product, you must first turn it off, unplug it, and let the device cool down.
- Use a soft, clean, and dry cloth to wipe the mirror of the light window.
- Store the product in a dry and dust-free environment between -10 and 55 °C.
- Do not wash the product and its accessories under a faucet or in a dishwasher.

- Do not use steel balls, abrasive cleaners (toothpaste, etc.) or corrosive liquid (gasoline, motor oil, acetone, etc.) to clean the product.

**Note:** If water is insufficient to clean the product, you can use a few drops of high-concentration alcohol (lower than 95%) to wipe the parts.

## ■ Environmentally Safe Disposal

When disposing of this product, please do not include it with general household waste. Instead, please send it to a designated recycling center.

## FAQ

This section outlines the problems you may encounter when using this product. If you cannot resolve a problem with the solutions provided, please contact the Ulike customer service center in your country/region.

- If an abnormal situation occurs or a defect cannot be fixed, please contact the after-sales customer service department.
- Please keep your purchase invoice and warranty card. You will need to provide this documentation if the product fails.

This device complies with the standards

EN 55014-1

EN55014-2

EN 61000-3-2

EN61000-3-3

EN 62233

EN 60335-1

EN 60335-2-23

IEC 60335-2-113

IEC 62471

BS EN 60335-2-23

BS EN 60335-1

BS EN 62233

BS EN 55014-1

BS EN 55014-2

BS EN IEC 61000-3-2

BS EN 61000-3-3

## Faults and Solutions

Problem	Possible reason	Solution
The product will not turn on	The product is not connected to a power source.	Check that the product is plugged in to an electrical outlet
The product does not emit light	The temperature of the main unit head is too high. Or, the main unit head is not in full contact with skin	Wait for the device to cool down naturally before continuing. Check that the main unit is in close contact with skin—watch for the indicator light.
The working indicator light is not on	The product is defective or malfunctioning	Check if the power supply is working properly Repair or return the product immediately
There is an odor when the product is flashing	The light filtration mirror is contaminated by hair left on the area to be treated	Wipe off the light filtration mirror. Shave hair off area to be treated. Reduce product's intensity level
You experience unbearable pain during hair removal	The hair in the area to be treated has not been shaven	Before using this product, please shave off the hair in the area to be irradiated. If shaving causes skin irritation, cut the hair as short as possible and use a lower intensity mode
	The area to be treated is not suitable according to the guidelines presented above	Do not use this product on the inner labia, vagina, anus, nipples, and areola. Male users should not use this product on the face or on genitalia
	The intensity level is too high	Reduce the intensity level until you are comfortable
	The light window is damaged	Immediately contact Ulike customer service to arrange for repair or to return the product
After hair removal, the treated area is red	Slight redness is normal and harmless and will soon fade	Use aloe vera gel to soothe the reddened area or use the sapphire lamp head to cool down the skin
Unsatisfactory hair removal effect	You used an intensity level that is too low for your hair type	Select a higher intensity
	You didn't press the flashing button when moving to the next area	Press the flashing button when you move to the next area. There should be some overlap between different areas
	This device has no effect on your hair color	If your hair is light blond, gray, red, or white, irradiation will not be effective
	You did not use the product frequently enough	To successfully remove all hair, we recommend that you follow the provided irradiation regimen
	Your body's desired response to irradiation may take longer than expected	This product should be used continuously for at least 5 months in order to fully suppress hair regrowth
The product becomes hot and humid during use	The product has been on and in use for an extended period of time	Turn off the product and let it fully cool down before using it again

## Warranty regulations

- The warranty period of this product may vary according to the laws of the country in which it is sold.
- The warranty period of this product, the power cord, and included accessories is 2 year from the date of purchase.
- If the product has quality problems during the warranty period, Ulike will repair or replace it free of charge. However, damage caused by improper use is not covered by the warranty. In that case Ulike will replace the damaged part or the old device at the original price. The old device will be returned to Ulike.
- Please present your warranty card at the time you require service. The warranty period of the product is based on the date on which you bought it. If no date of purchase is available, the warranty period will be calculated based on the date of production.

**In the following cases, Ulike may not provide warranty service, refunds, or replacements.**

**Ulike is not responsible for reissuing the purchase invoice and warranty card. Please keep your purchase invoice and warranty card safe.**

- 1.Product damage caused by unpredictable events such as fire, earthquake, flood, etc.
- 2.Product failure and damage caused by negligence, misuse, or use under abnormal conditions.
- 3.Product damage caused by unauthorized disassembly and improper storage (such as rodent damage or liquid penetration).
- 4.Product damage caused by use of the product in inappropriate circumstances, including but not limited to excessive workload.
- 5.Product damage caused by other conditions that Ulike cannot control.
- 6.For damage associated with the function or use of this product (such as loss of business benefits, business interruption, other losses caused by the use of the product, costs of replacing equipment, devices or services related to using this product, customer claims, etc.), Ulike assumes no responsibility.

Depending on the type and extent of accidents, Ulike will bear corresponding responsibility according to relevant laws and regulations.

## Authorised Representative

Ulike Europe SAS  
4 rue Royale, 75008, Paris, France  
Contact: Raphael  
E-Mail: [contacteurope@ulikeglobal.com](mailto:contacteurope@ulikeglobal.com)

Name: Kingsmead Service Limited  
Add.: 19 mezzanine floor 19-21 Crawford Street London, England W1H 1PJ  
Contact: Peng Yan  
E-mail: [office@kingsmead-service.com](mailto:office@kingsmead-service.com)

## Basic parameters

<b>Air3 Series</b>	Specification	
	Product name	: IPL Hair Removal Device
	Product model	: UI06
	Product Size	: 179*58*34mm
	Rated Input	: 24V $\overline{\text{=}}$ 2.8A
	Power	: 67.2W
	Energy level	: 3 levels
	Spot size	: 3.3cm <sup>2</sup> $\pm$ 0.25cm <sup>2</sup>
	Spectrum	: 550-1200nm
	Pulse duration	: 1.12ms-7.44ms
	Pulse repetition rate	: 0.33Hz-2.5Hz
	Manufacturer	: Shenzhen Ulike Smart Electroics Co., Ltd
	Address	: 810, Building 1, Xunmei Science and Technology Plaza, No. 8 Keyuan Road, Science Park Community, Yuehai Sub-District, Nanshan District, Shenzhen 518000, Guangdong, P.R. China

Storage and transportation conditions		Operating conditions	
storage temperature	-10-55°C	Operating temperature	5-30°C
Humidity	$\leq$ 90%RH	Humidity	$\leq$ 80%RH
Relative pressure	50-106 kPa	Relative pressure	80-106kPa

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Skin color \ Hair color	Skin color						
	White	Beige	Light brown	Medium brown	Dark brown	Brown black and darker	
White/Grey	×	×	×	×	×	×	
Red	×	×	×	×	×	×	
Light blond	×	×	×	×	×	×	
Dark blond/light brown	regular mode power mode	regular mode power mode	regular mode power mode	soft mode regular mode	×	×	
Brown	regular mode power mode	regular mode power mode	regular mode power mode	soft mode regular mode	×	×	
Dark brown	regular mode power mode	regular mode power mode	regular mode power mode	soft mode regular mode	×	×	
Black	regular mode power mode	regular mode power mode	regular mode power mode	soft mode regular mode	×	×	

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Skin color \ Hair color	Skin color					
	White	Beige	Light brown	Medium brown	Dark brown	Brown black and darker
White/Grey	×	×	×	×	×	×
Red	×	×	×	×	×	×
Light blond	×	×	×	×	×	×
Dark blond/light brown	regular mode power mode	regular mode power mode	regular mode power mode	soft mode regular mode	×	×
Brown	regular mode power mode	regular mode power mode	regular mode power mode	soft mode regular mode	×	×
Dark brown	regular mode power mode	regular mode power mode	regular mode power mode	soft mode regular mode	×	×
Black	regular mode power mode	regular mode power mode	regular mode power mode	soft mode regular mode	×	×

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Skin color Hair color	Skin color					
	White	Beige	Light brown	Medium brown	Dark brown	Brown black and darker
White/Grey	×	×	×	×	×	×
Red	×	×	×	×	×	×
Light blond	×	×	×	×	×	×
Dark blond/light brown	regular mode power mode	regular mode power mode	regular mode power mode	soft mode regular mode	×	×
Brown	regular mode power mode	regular mode power mode	regular mode power mode	soft mode regular mode	×	×
Dark brown	regular mode power mode	regular mode power mode	regular mode power mode	soft mode regular mode	×	×
Black	regular mode power mode	regular mode power mode	regular mode power mode	soft mode regular mode	×	×

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