INVOSPA

YOUR PERSONAL MASSAGE EXPERT

Multipurpose KNEADING MASSAGER

Contact info: support@invospa.com

INVOSPA

Introduction

Thank you for purchasing the InvoSpa Multipurpose KNEADING MASSAGER.

This manual contains the instructions you will need to operate and maintain your model so that you can enjoy it for years to come. However, please be advised that mishaps could happen in time, and rest assured that InvoSpa is a team of professionals who aim to provide the highest level of support possible.

We know you're excited about getting started with your new massager, but it's very important that you take some time to read through the User's Manual. This manual contains all the necessary setup and operating procedures that will allow you to use your new device. Also, be sure to read and follow the precautions and warnings in this manual and on any labels or tags attached to your model. They are there to educate you on how to operate your unit safely and also get maximum life and performance from it.

Thank you again for choosing InvoSpa. We work hard every day to make sure you receive the highest level of customer satisfaction possible. We truly want you to enjoy your new multipurpose Kneading Massager!

Find more information on our website:
www.invospa.com.

Scan the QR code with your phone
camera to get assistance.

For premium support & lifetime warranty
Register now!

INVOSPA

What's in the Box

InvoSpa Kneading Massager

AC and Car
Adapters

User Manual

Product Name: InvoSpa Kneading Massager
Model: JC-668 (Blue), JC-779 (Dark Gray), TE-336 (White),
TE-339 (Black and Gray), TE-886 (Light Gray),
TE-889 (Black), ZC-3006 (Green), ZC-3007 (Yellow),
ZC-3008 (Red)

Power Consumption: 24W

Voltage: 12V 2A

Box Size: 40*19.5*17 cm \ 15.75 x 7.68 x 6.69 in

Product Weight: 2.031 kg \ 4.478 lbs

Terms of Use

This product is surrendered by InvoSpa to the purchaser with the understanding that the purchaser accepts the responsibility that using this model and the enclosed accessories in a careless, improper, or unsafe manner can result in serious injury or death. Also, the purchaser assumes all liability resulting from any misuse, unsafe handling, failure to follow instructions, or any action that constitutes a violation of any applicable laws or regulations. InvoSpa shall not be liable for personal injury, loss of property, or loss of life resulting from the use of this product under any circumstances, including intentional, reckless, negligent, or accidental behavior. InvoSpa shall also not be liable for any special, indirect, incidental, or consequential damages arising out of the assembly, installation, or use of their products or any accessory or chemical required to use their products. By the act of using the product, the user accepts all resulting liability and releases InvoSpa of any and all liability associated with its use. If you as the user do not accept liability of usership, InvoSpa requests that you do not use this product. Do not open any of the enclosed materials and return the unit to InvoSpa.

All information contained in this User's Manual is subject to change without notice.

For the latest version of the User Manual, please visit our website: www.invospa.com. InvoSpa reserves the right to make changes and improvements to products without incurring any obligation to incorporate such improvements into products previously sold. If you have any questions about your product or its operation, contact InvoSpa Customer Support at support@invospa.com.

Warranty

LIMITED WARRANTY

Our InvoSpa return policy allows returns for a full refund within 30 days from the delivery of the order, no questions asked, just by providing proof of purchase. After 30 days, you can still return the product for replacement or repairs up to 12 months (US) / 24 months (EU/UK) ONLY IF the unit is defective.

Please note that our InvoSpa products are regulated by the U.S. Consumer Product Safety Commission.

COVERED BY WARRANTY

To guarantee that you have the best experience, we offer a one-year warranty covering defects in materials, workmanship, and assembly on all our products.

*To ensure that you don't invalidate the terms of the warranty, please read and understand this booklet.

IT IS THE PARENTS' OR GUARDIANS' RESPONSIBILITY TO ENSURE MINORS RECEIVE APPROPRIATE GUIDANCE AND/OR SUPERVISION WHEN USING THIS PRODUCT. THIS IS A COMPLEX HIGH-PERFORMANCE MASSAGER AND MUST BE TREATED WITH CARE AND RESPECT.

InvoSpa has made every effort to ensure that your product is as strong and durable as possible. However, it is still possible that some mishaps may occur, and we take full responsibility. If the product starts to malfunction, please contact us at support@invospa.com and we will assist right away.

NOT COVERED BY WARRANTY

- All other parts not mentioned above
- Lack of maintenance, improper use or abuse
- Water damage or problems resulting from water intake

INVOSPA

Control Panel

Power ON/OFF Button

Massage Direction Button

Speed/Intensity Button


Heat ON/OFF Button

MASSAGER FEATURES


1. Shiatsu-type Massage
2. 8 Kneading Rollers
3. Heat
4. 2 Massage Directions
5. For neck, shoulders, back, waist, arms, thighs and legs.
6. Auto shutdown after 15 minutes
7. Home AC adapter and car adapter included.

Benefits

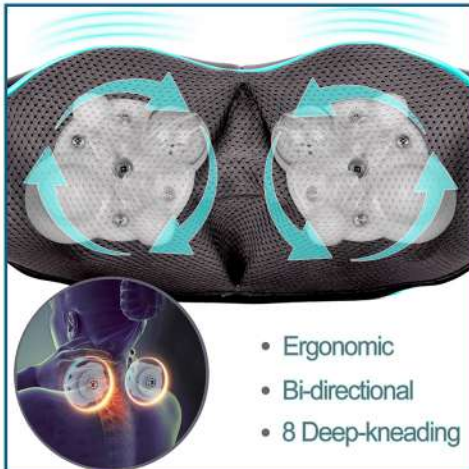8 MASSAGE ROLLER BALLS – This shoulder massager comes with 4 big nodes and 4 small nodes, which provide deep tissue massages on your neck, shoulders, upper back, lower back, waist, foot, thighs, calves, legs, feet and arms.

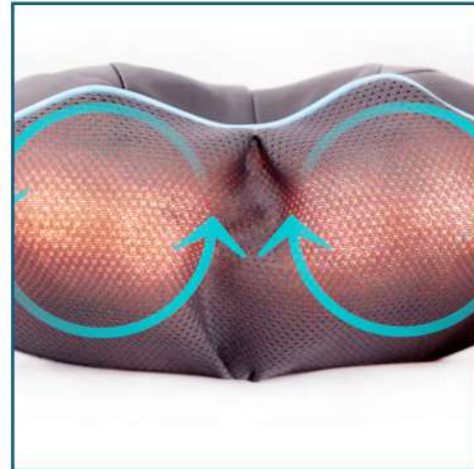HEAT FUNCTION – By turning on the heat function (optional), you can enhance your massage session.

DESIGNED FOR YOU – Massaging For Upper and Lower Body.

ADJUSTABLE INTENSITY – This neck massager pad has 3 speed strength levels, which allow you to get the right and appropriate amount of pressure on your muscles.

PORTABLE AND DURABLE – This neck and back massager includes a wall power adapter and a car charger, so you can have an electric self-massage in your home, office, or car when you are traveling. This Shiatsu Kneading massager pillow is the perfect gift for your loved ones!

Safety Instructions


1. The massager has a built-in timer and it will automatically turn off after 15 minutes. You can push the power button to start another 15 minute session.
NOTE: To avoid possible overheating of the unit, it is recommended NOT to use the massager for more the 30 consecutive minutes.
2. Avoid using the massager under direct sunlight or in any high temperature setting.
3. Do not use the massager in the bathroom, bathtub or any setting with high levels of moisture.
4. Disconnect the power when the massager is not being used.
5. Always disconnect the power before cleaning the massager.
6. Do not clean the massager with detergent or any sort of liquid.
7. Wipe the massager with dry or damp cloth to clean it.
8. Do not use the massager if any part of it (including power cord) is damaged.
9. Do not let children use the massager.
10. Stop using the massager if you are experiencing discomfort or pain.
11. Do not use the massager if your body is wet.
12. Do not use the massager on sensitive skin.
13. Do not use the massager on joints and bony parts of the body.
14. Do not use the massager on your abdominal region if you are pregnant.

UPPER BACK/NECK


MID BACK

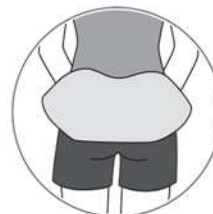LOWER BACK

THIGH / LEG

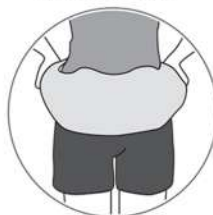WAIST

SHOULDERS / ARMS

Safety Precautions

SAFETY PRECAUTIONS: This massager is capable of extreme performance. Responsibility, maturity, and common sense are required. This model must be treated with respect and caution to ensure your safety. Operating this massager in a careless, unsafe manner, without adequate care and preparation, can result in serious injury or death.

- ALWAYS confirm that the charger/plug settings exactly match the type (chemistry), specification, and configuration of the device to be charged. DO NOT exceed the maximum manufacturer recommended charge rate, which is DC 12V, 2A.
- DO NOT disassemble, crush, short circuit, or expose the device to flame or other sources of ignition.
- ALWAYS unplug the charger from the wall outlet and disconnect the cable when not in use. If there are any signs of a malfunction, unplug the power source.
- If the massager gets too hot to the touch during operation (temperature greater than 100-104°F / 38-40°C), immediately disconnect it from the power source.
- Keep a Class D fire extinguisher nearby in case of fire.
- DO NOT expose the unit to water or moisture.

WARNING: Do not use the device if you have any of the following conditions:

- You are pregnant.
- You have a pacemaker or other medical implant.
- You are suffering from heart or vascular disease.
- You have irregular blood pressure or are receiving treatment to normalize your blood pressure.
- You are suffering from osteoporosis.
- You are suffering from cancer.
- You have recently undergone surgery.
- You have concerns regarding areas of your health.

Never leave the appliance unattended, especially if children are present.

We recommend you limit the use of the massager to 15-minute sessions and for no more than 30 minutes consecutively.

Never use this product directly on swollen or inflamed areas of skin or open wounds.

This product should NEVER be used by an individual suffering from any physical ailment that limits the user's mobility.

WARNING

We strongly recommend that people who are very sensitive and/or have sensitive limbs use the lowest intensity.

Note: The feeling of pressure when using the massager is perfectly normal and not a sign of injury or malfunction. We recommend using the intensity/speed levels gradually, going from the lowest to higher levels.

Care, Maintenance, and Cleaning Instructions

- Store the massager in a cool dry place when not in use.
- Do not store the unit in the bathroom or in a humid setting with high levels of moisture.
- Do not place or immerse the unit into any type of liquid.
- Do not place the unit near a fireplace, heater, stove, oven, furnace, or other sources of ignition.
- Keep away from all solvents and harsh detergents.
- After use, wipe with a soft, damp cloth using non-abrasive cleaners.

Troubleshooting

PROBLEM	CAUSE	SOLUTION
Massager is not turning on	Power plug is loose	Reconnect the plug to the outlet.
	Adapter is not plugged in correctly	Unplug the adapter and plug it in again.
	Power is off	Push the on button to turn the massager on.
Massager is turning off on its own	Faulty adapter	Connect the massager to the car plug to verify massager functionality.
	No power	Try a different power outlet or check the car power sockets.
	Power plug is loose	Follow instructions above and check the connection, adapter, and sockets.
	The 15-minute timer went off	The massager has a timer and will turn off on its own after 15 minutes. You can push the Power button to start another 15 minute session.

NOTE: IT IS NOT RECOMMENDED TO USE THE MASSAGER FOR MORE THAN 30 CONSECUTIVE MINUTES.

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INVOSPA

INVOSPA Wellness Products Collection

Shiatsu Back, Shoulder, and Neck
Massager (Dark Gray)

Shiatsu Back, Shoulder, and Neck
Massager (White Fabric)

35 Shiatsu Back, Shoulder, and Neck
Massager (Black and Gray Fabric)

Shiatsu Back, Shoulder, and Neck
Massager (Light Gray Fabric)

INVOSPA

INVOSPA Wellness Products Collection

Shiatsu Back, Shoulder, and Neck
Massager (Black Fabric)

Shiatsu Back, Shoulder, and Neck
Massager (Green)

Shiatsu Back, Shoulder, and Neck
Massager (Yellow)

Shiatsu Back, Shoulder, and Neck
Massager (Red)

INVOSPA

INVOSPA Wellness Products Collection

**Electric Throw Heated
Blanket Gray - 50"x60"**

**Electric Throw Heated
Blanket Blue - 50"x60"**

**Electric Throw Heated
Blanket Gray - 72"x84"**

**Electric Throw Heated
Blanket Blue - 72"x84"**

INVOSPA

INVOSPA Wellness Products Collection

Shiatsu Foot Massager

Muscle Massage Gun for Athletes

Leg Massager for Circulation

INVOSPA

YOUR PERSONAL MASSAGE EXPERT

JC-668

JC-779

TE-336

TE-339

TE-886

TE-889

ZC-3006

ZC-3007

ZC-3008

Shaft Innovations HK Limited

HK Address
Room 409, Beverley Commercial Centre,
87-105 Chatham Road South, Tsim Sha Tsui,
Kowloon, Hong Kong

EU Address
ORCA - Calle Gil-Vernet 54/55, Les Tapies 1
Hospitalet De L'infant, Tarragona, Spain 43890

Email: support@invospa.com
Web: www.invospa.com

Design by InvoSpa in Austin, TX
Made in China

Model Number:

JC-668 (Blue), ZC-3006 (Green),
JC-779 (Dark Gray), ZC-3007 (Yellow),
TE-336 (White), ZC-3008 (Red)

TE-339 (Black and Gray),

TE-886 (Light Gray),

TE-889 (Black)

Voltage: DC 12V 2A

Power: 24W

FCC UK CA CE RoHS

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

The above declaration relates to both the massager device and the DC wall adapter.