Coslus Group

CONTACT

✉ support@coslus.com 🌐 www.coslus.com

☎ +1(855) 594-3950 (9am-5pm Mon. - Fri. (EST))

Made in China


COSLUS

3+8 DIY Levels

Oral Irrigator
Munddusche

C30

USER MANUAL
Benutzerhandbuch

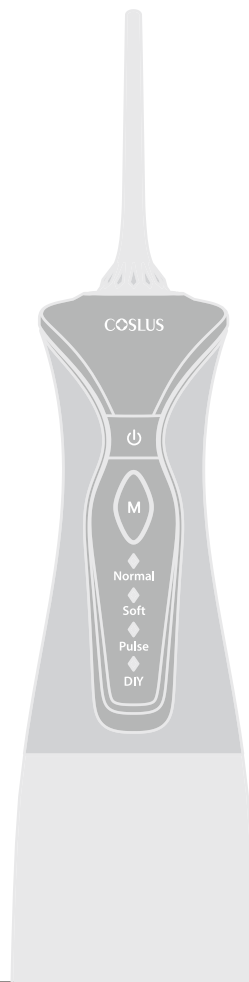Contents

Important safeguards	01-03
Product description	04
Product accessories and specifications	05-06
Operation steps	06-09
Cleaning and maintenance	09
Environmental protection	10
Troubleshooting	11

In case of any typographical errors or misunderstandings in the contents of the manual, we reserve the right of final interpretation. In case of product upgrade or information update, it will be incorporated into the new version of the manual without prior notice. If case of any change in the appearance and color of the product, the physical objects in the package shall prevail.

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

PLEASE READ ALL INSTRUCTIONS BEFORE USING THE PRODUCT.

DANGERS-To reduce the risk of electrocution:

1. Do not place or store product where it can fall or be pulled into a tub, shower stall or sink.
2. Do not place in or drop into water or other liquid.
3. Do not reach for a product that has fallen into water. Unplug immediately.
4. Do not handle charger with wet hands when charging.

WARNING-To reduce the risk of burns, electrocution, fire, or injury to persons.

1. A product should never be left unattended when plugged in.
2. Close supervision is necessary when this product is used by, on, or near children or invalids.
3. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
4. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the product to a service center for examination and repair.
5. Keep the cord away from heated surfaces.
6. Never block the air openings of the product or place it on a soft surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint, hair, and the like.
7. Never use while sleeping or drowsy.
8. Never drop or insert any object into any opening or hose.
9. Unplug this product before filling. Fill (reservoir) with water only unless otherwise specified by manufacturer. Do not overfill (or specify filling instructions).
10. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
11. Children should be supervised to ensure that they do not play with the appliance.
12. Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
13. This product needs to be charged with the charger specified by the manufacturer, do not use the charger not specified or other charger.
14. Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger.

15. The appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
16. This appliance contains batteries that are non-replaceable.
17. For household use only. DO NOT use the product outdoor.

Danger: To reduce electric shock and other hazards, the following issues should be noted

1. Before or during the first use of the product, please check whether the charger and charging cable are damaged. In case of any damage, please contact the after-sales service for replacement before using;
2. Do not put the product into fire or heat the product;
3. Do not immerse this product in water to avoid functional failure of the product;
4. Please do not use wet hands to plug and unplug the charger, when unplugging the charger, please be sure to hold the charger plug rather than the charging cable;
5. Please do not damage, destroy or modify the charging equipment, do not forcefully stretch, bend or twist the charging cable, and do not press heavy objects on the charging cable or let the charging cable get stuck between objects;
6. Do not drop or insert foreign matter into any opening of the product;
7. Do not place the product, the charger and the charging cable on the heated surface;
8. Do not place the product in a place with high temperature or humidity and avoid direct sunlight;
9. Before cleaning and maintenance, please ensure that the product is uncharged and powered off.

Warning: To reduce the risk of burns, electric shock, fire, and personal injury, please remember the following

1. Do not operate the product in a way other than the Manual;
2. Do not use the product without installing the nozzle;
3. People (including children) who do not have full capacity and are lack of relevant experience and knowledge must use, clean and maintain the product under the supervision of guardians;
4. Please place the product out of the reach of children and babies, and do not use it as a toy;
5. Patients with periodontal disease or those who have had oral surgery in the past 2 months, please consult a doctor before using the product;
6. Only nozzles and accessories donated, sold and recommended with products are allowed to be used;
7. Do not disassemble the product or change the parts of the product to avoid damage or danger;
8. Do not fill hot water above 40°C into the water tank and clean your mouth to avoid scalding;
9. Do not use the product for other purposes than oral cleaning;

10. Do not press the nozzle hard against gums or teeth to avoid injury;
11. Do not use and operate the product in outdoor places where aerosol (spray) or oxygen supply is used;
12. Please remove all dental ornaments when using the product;
13. Do not use iodine, salt water or concentrated essential oil insoluble in water in products. Using these liquids will reduce product performance and shorten its service life;
14. Do not fill pure mouthwash into the product. If mouthwash is required, please dilute it proportionally before filling.

SAVE THESE INSTRUCTIONS

FCC STATEMENT


This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Caution: Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:


- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Product description

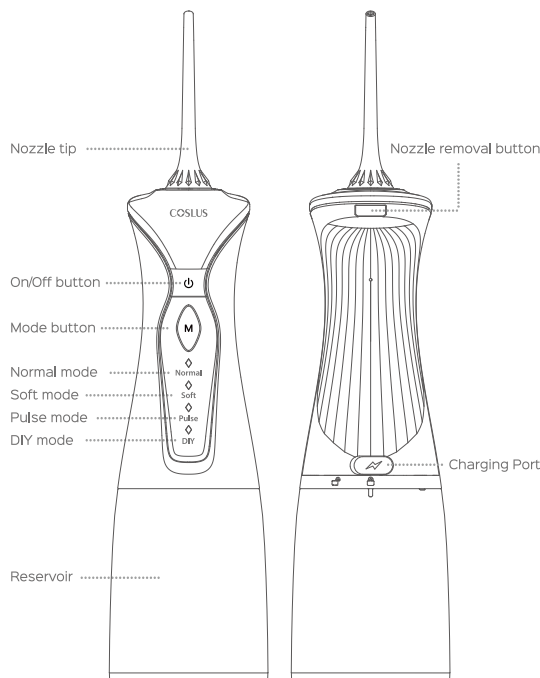Product accessories and specifications

Product fittings

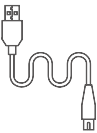			
Standard Nozzle x 2	Orthodontic Nozzle x 1	Tongue Scraper x 1	Charging Cable x 1


You can buy the fittings from the seller

Product specification

Product name		Oral irrigator
Oral irrigator	Product model	C30 (FC5360)
	Rated voltage	3.7V ===
	Waterproof grade	IPX7
	Charging time	≤4h
	Input parameters	5V === 1A max 2A
	Adapter ratings	Input: 100-240V~50/60Hz 0.3A
		Output: 5V === 1A 5.0W
Working time		40-120min (depend on mode)

Note: Charger plug is not included.


Nozzle description

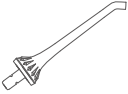Replace every 6 months		
Standard Nozzle		Daily cleaning
Orthodontic Nozzle		For washing braces, bridge and residue
Tongue Scraper		Clean coating and bacteria effectively

Operation steps


Start work: Very important! Ensure the device is fully charged before using it for the first time by connecting the USB cable to the recharging socket and turning on the power. This device cannot be used while being charged.

Step 1: Fill water into reservoir

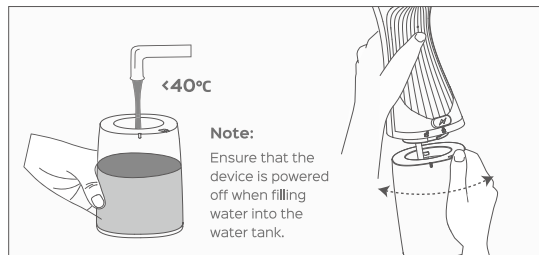Rotate the water tank clockwise and pull to remove the water tank, and then fill water into the water tank; After filling water, align the water tank with the buckle point, and rotate counterclockwise to install it back.


Step 2: Install the nozzle

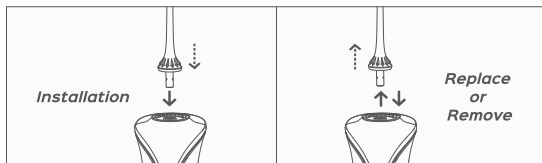Insert the nozzle into the top hole of the main unit and when you hear "click" sound, the installation is successful; Rotating the nozzle can adjust the direction of water outflow; To replace the nozzle, press and hold the nozzle removal button, and then remove the nozzle.

Step 3: Select the mode

There are 4 modes for the product. You can select the required mode by tapping the mode button when the device is powered off, and the functions of each mode are as follows:

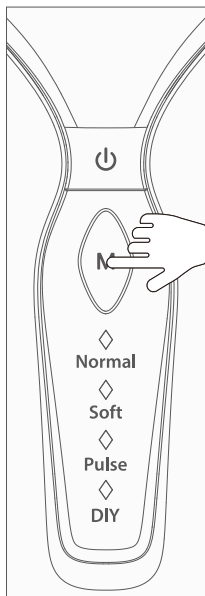1. Standard modes: press the On/Off button, the product enters the working mode; short press the mode button to change between 4 modes: normal - soft - pulse - DIY, cycle mode in turn. In any working mode, the machine will be powered off if pressing the ON/OFF button.


2. DIY mode: if there is no suitable pressure under standard modes, you can enter DIY mode to choose 8 water pressure gears.

In the working state, press the mode button to select the DIY mode, and long press the mode button for about 3 seconds. Until the DIY indicator flash quickly, you can change water pressure by short pressing. If you have selected the required gears, loosen the mode button. After 3 seconds, it will exit from choosing gears, and the indicator will stop flashing.

3. Long press the mode button to display the current power level, which is indicated by the number of lights on. When charging, the lights will flash in sequence. After fully charged, all lights are on. If the power is too low during use, the indicator light will flash to remind the user to charge. When the indicator light flashes 3 times, the device automatically shuts down.
4. Working time of about 2min; memory mode: when power on, the product will enter the working mode before last shutdown.
5. When in use, it is erected vertically to ensure that the equipment is in the best working state. After running in any mode for 120s, the device is automatically powered off.

Step 4: Preparation before use

Before use, please bend over the sink slightly, extend the nozzle into your mouth, aim at your teeth, gums or slit between the teeth, rotate the nozzle to adjust the direction of water outflow, ensure that the water flow is perpendicular to the flushing position, and open your mouth slightly to facilitate the flushing water to flow out of your mouth.

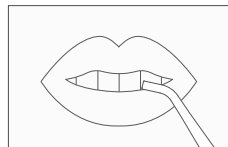Note:

When using the oral irrigator for the first time, please start from the lowest gear, and gradually adjust the water pressure after the oral cavity adapts; When using the product, itching or slight gum bleeding may occur, which is a normal phenomenon and usually stops within 1-2 weeks.


Step 5: Using


After pressing the On/Off button, the device is powered on. If you need to pause during use, please press the On/Off button lightly. If you continue to use it, please press the On/Off button again. The product has a timed power-off function in 120s. Before the automatic power-off, water in the tank will be used up. If you need to continue to use, please first lightly press the On/Off button to turn off the machine, add water and then start the equipment to use. After the use, please turn off the device.

Notice

To obtain a better cleaning effect, please start from the molars (posterior teeth) and start cleaning in the direction of incisors. When cleaning, move the nozzle 90° along the gum line and stay between the teeth for a short time. When using the product, the lips are slightly closed to prevent the water from flowing sideways, leaving a certain gap to facilitate the water in the mouth to flow freely to the wash basin.

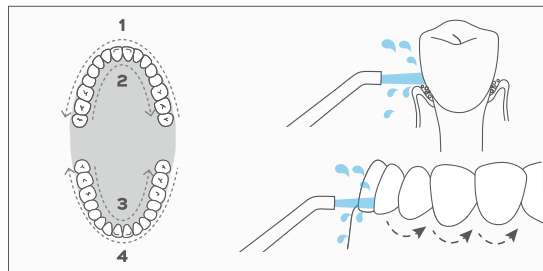Step 6: End

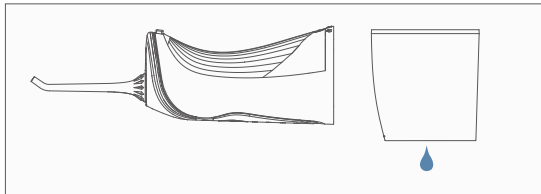After the use, open the water tank, pour out residual liquid in the water tank, open the equipment to drain the residual liquid in the equipment until no more water is sprayed, and then power off the equipment.

Step 7: Charging

The indicator light is white, and the light is always on when in use. When the light flashes 3 times and then turns off, it means the product is in low voltage. Press the On/Off button when the voltage is too low, and the product can not work. The indicator light flashes 3 times, indicating that the product needs to be charged in time. When charging, the indicator light flashes; after fully charged, indicator light is always on. Always on for 20 minutes and then goes out.

Cleaning and maintenance

Clean the main unit:

1. Please use water or neutral detergent to remove the product, do not use corrosive or abrasive detergent, do not use hot water over 40°C to clean any part of the machine to avoid damage to the product and reduce the performance of the product.
2. Do not soak the oral irrigator in water. If necessary, wipe the surface of the device with a soft cloth dipped in a mild non-corrosive detergent.

Clean the water tank:

1. Grab the product body with your left hand, grab the water tank with your right hand, rotate it in the counterclockwise direction to remove the water tank.
2. Please rinse the inside of the water tank with clear water of normal temperature, and do not clean it with any detergent.
3. After cleaning, align the water tank with the buckle point and then rotate clockwise to put it back.

Removal of hard water precipitates:

1. Hard water precipitates may accumulate on the device, depending on the mineral content in water. If it is not removed in time, the product performance will be reduced.
2. Align the nozzle with the sink, turn on the machine for running until the water is used up. Then fill the water tank with clean warm water and repeat the above operation to carry out cleaning. Cleaning operation should be repeated every 1-3 months to ensure performance.


Environmental protection

Battery maintenance

- The battery must be fully charged before use for the first time;
- If you are going to store it for a long time, the battery must be fully charged before storage;
- To prolong the service life of the battery, please charge it immediately in case of low power;

Battery recycling

The product is equipped with recyclable rechargeable battery. At the end of the product life and before discarding the product, you should take out the battery and recycle or dispose of it according to the local waste management regulations.

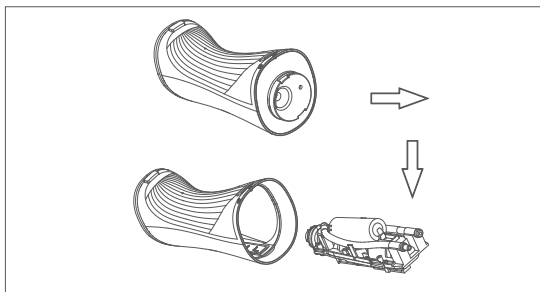Recovery steps: When the equipment is powered off, remove the nozzle and the water tank, open the bottom cover of the main unit with tools such as screwdriver and remove the built-in battery for recycling.

(Note: Do not short circuit the positive and negative terminals of the removed battery. Stick tape on the lead terminals for insulation)

Problem	Cause	Solution
Leakage between nozzle and host	The sealing ring in the nozzle fixing seat is damaged, or the nozzle is not installed in place	Contact the seller for repair, or try to install the nozzle in place
Leakage at pump inlet	Sealing ring is damaged or there are impurities at inlet	Contact the seller for repair
Leakage of water tank	Sealing ring is damaged	Replace the sealing ring or contact the seller for repair
Low pressure	Use under low voltage	Use when fully charged
Device cannot be started	Use under low voltage	Fully charged or contact the seller for repair
No water coming out	The product is used for the first time or has not been used for a long time	Fill proper amount of water to the nozzle installation port and use it after installing the nozzle

If the problems cannot be solved according to the above methods, please consult the selling store or go to the designated after-sales repair point for repair.

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