Coslus Group

CONTACT

✉ support@coslus.com 🌐 www.coslus.com

☎ +1(855) 594-3950 (9am-5pm Mon. - Fri. (EST))

CE FC IC

RoHS


Made in China

V2.0

COSLUS

E40

Oral Irrigator User Manual

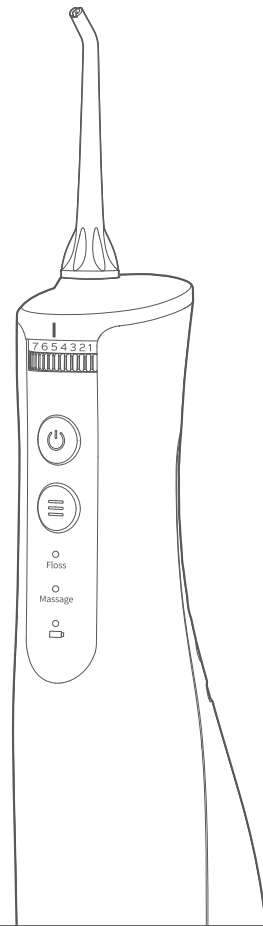Model: E40

Oral Irrigator

Munddusche

Irrigateur buccal

Irrigador bucal

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User Manual (EN)	01
Bedienungsanleitung (DE)	13
Manuel d'utilisation (FR)	25
Manual del usuario (ES)	37
Manuale dell'utente (IT)	49

IMPORTANT SAFETY CAUTIONS

All information in this part is provided to instruct you to use the appliance in a correct and safe way so that you can avoid injuries to yourself and others as well as damage to property.


- The appliance can be used by children aged 8 years and above. Person with reduced physical, sensory or mental capabilities or lack of experience and knowledge are not allowed to use unless they have guardian supervision or instruction concerning safe using.
- Children shall not play with the appliance. Cleaning and maintenance shall not be made by children without supervision.
- The appliance is only to be used with the power supply unit provided with the appliance.
- Never alter the appliance. Also, do not try to disassemble or repair it.
- The battery must be removed from the appliance before it is scrapped.
- Unplug the appliance from power supply when removing the battery.
- The battery should be disposed of safely.
- Regarding how to remove the battery, please refer to disposal instruction.
- Do not fill the tank with water warmer than 40°C.
- If the external flexible cable or cord is damaged, please ask the manufacturer, service agent or qualified professional person for help to avoid hazard.
- Keep the cable or cord away from heated surfaces.
- Please do not use it if you lost the nozzle.
- Do not submerge the device in water.
- Do not use for any purpose other than oral cleaning.
- Keep nozzle away from small children. The nozzle is small enough to be swallowed and can be lodged in the throat.
- Certain types of mouthwash may damage the appliance by causing the housing or water reservoir to crack. Therefore, do not fill water tank with mouthwash.
- People who is suffering from tooth or gum pain may not be able to use the appliance and should consult their dentist before use.
- Do not press the nozzle too firmly against your teeth or gums, or it may cause injury.

Product specifications

Model	E40(WFP14)
Charge Input	5V $\mp$ 2A (Max)
Voltage	3.7V
Power	5W
Battery capacity	1400mAh
Noise	<72dB
Waterproof rating	IPX7
Water pressure range	30-150PSI
Storage temperature	0-40°C
Humidity	20%-85%
Atmospheric pressure	70-106kPa
Water tank capacity	300ml


Package list

- 1 Main unit x 1
- 2 Standard jet tip x 2
- 3 Orthodontic tip x 1
- 4 Periodontal tip x 1
- 5 Tongue scraper x 1
- 6 User manual x 1
- 7 USB cable x 1
(adapter is not included)

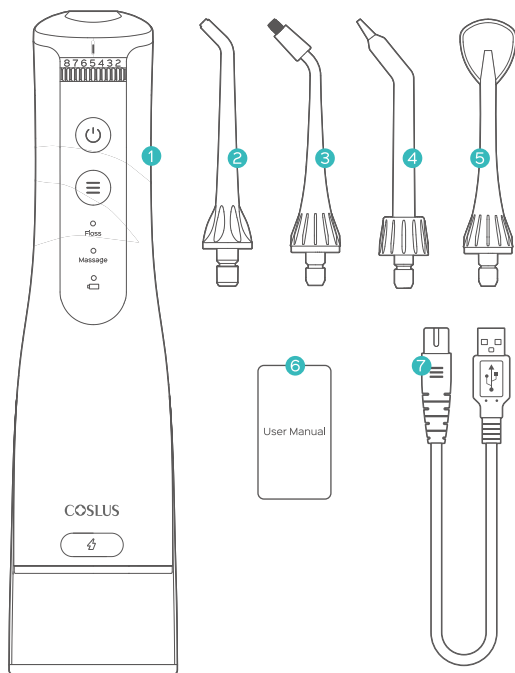Product Overview

- a. Speed Control Knob
- b. On/Off button
- c. Mode button
- d. Mode indicator
 - Floss
 - Massage
- e. Charging base
- f. Charging base cover
- g. Nozzle
- h. Rotatable ring
- i. Nozzle suction tube
- j. Filter
- k. Water tank
- l. Water tank cap
- m. Nozzle release button
- n. Battery indicator

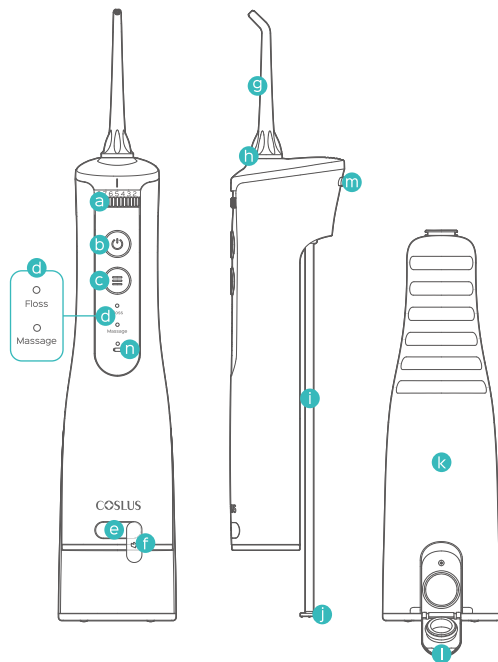Note:


There may be some liquid remaining inside the main appliance before use (tank, water suction tube, etc), this is the remaining distilled water for performance (water flow) test, so there are no health concerns.

Nozzle introduction

Standard Jet Tip: daily clean teeth, gum line, fissure, remove food residue, stain, plaque, etc.


Standard jet tip
x 2

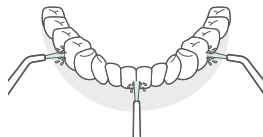Orthodontic Tip: clean bridge and teeth with braces, suitable for those who wear braces or have dental implant.


Orthodontic tip
x 1

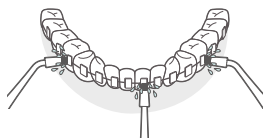Periodontal Tip: clean gingival sulcus, especially suitable for those who are suffering from periodontal diseases. When cleaning periodontal pocket, please use the lowest water pressure. Direct water flow along the gum line with nozzle tip at 45° and turn it on.

Periodontal tip
x 1

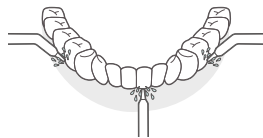Tongue Scraper: clean tongue surface, remove bad breath, food residue and thick tongue fur. When cleaning tongue, please use the soft mode. Lean forward, put scraper on your tongue and turn it on, then move gently from back to front. You can adjust water pressure according to your need.

Tongue scraper
x 1


Note: It is recommended to replace cleaning tip every 6 months for better cleaning experience.

Use instruction

How to charge


1. Take charging base cover down, securely plug USB cable into charging port of the appliance. (Make sure no liquid remaining inside the port before plugging.)
2. Connect the charging adaptor securely with the power outlet.
3. Battery indicator blinks during the oral irrigator is charging. When fully charged, it lights up for 1min and then turns off.
4. Cover charging base again after charging is finished.

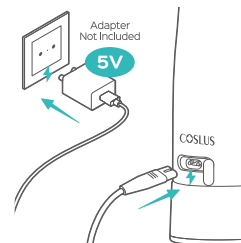
Note: When the battery indicator lights up or blinks, it means the oral irrigator needs charging. A full charge takes about 4 hours.

How to operate

1. Attach nozzle onto the main unit.

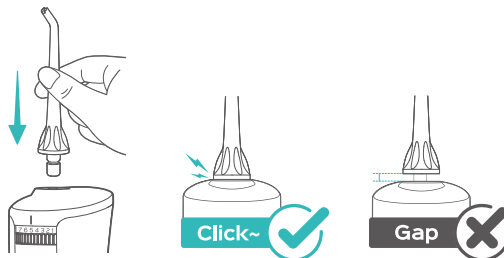2. Open water tank cap to fill the tank with water(not exceed 40°C), then close the cap. Or slide water tank down off the main unit to fill with water(not exceed 40°C), then attach the tank onto the main unit.

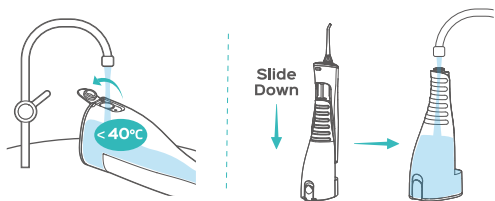Note: Steps below are very important. Improper operations may cause water leaking or spilling during use.

3. Press "Mode" button to select desired mode. Turn speed control knob to adjust water pressure. During irrigating, you can change the working mode and water pressure.

Mode introduction:

Floss: for normal deaning.

Massage: for cleaning and massaging gums.


Speed control knob introduction (Factory default pressure 3 levels): Turn speed control knob to the right to increase water pressure. Turn speed control knob to the left to reduce water pressure.

Pressure introduction:


Level 1-3: Low water pressure for beginners, people with sensitive gums, and children.

Level 4-7: Medium water pressure for most people skilled in using oral irrigators.

Level 8-10: High water pressure for people with experience in tooth cleaning and strong gum tolerance.


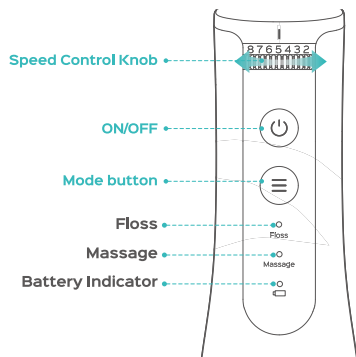
4. Upright hold the main unit with nozzle tip point toward your teeth.
5. Keep your mouth open slightly so that water can flow out.

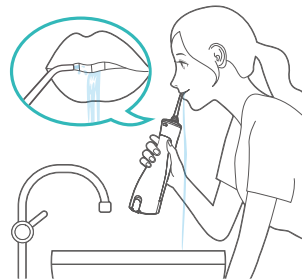
6. Press "On/Off" button to start irrigating.


- 1). Direct water flow at 90° to the teeth and gums
- 2). Direct water flow slowly along the teeth
- 3). Direct water flow along the gum line with nozzle tip at 90°
- 4). Clean between the teeth and around any braces, crowns or bridgework thoroughly

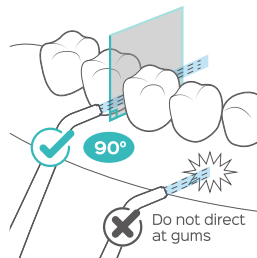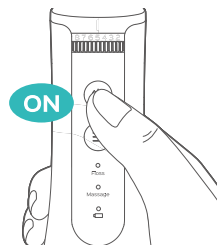7. Press "On/Off" button to turn it off, then press nozzle release button to replace other nozzle you need.

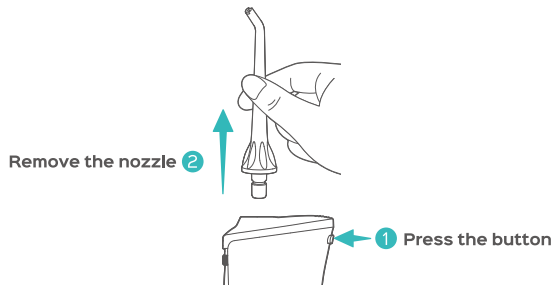8. After irrigating is finished, press "On/Off" button again to turn off.

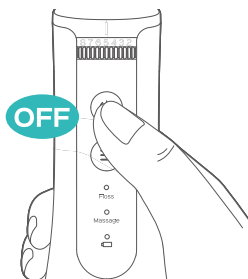Maintenance instructions

After using

1. Open water tank cap and dispose of water inside the tank.
2. Tilt and shake the main unit to drain the water from the speed control knob.
3. Press "On/Off" button to turn on and drain water inside the appliance.
4. Wipe the appliance with dry cloth.

Note: To prevent bacteria breeding, make sure there is no remaining water inside. If you're not going to use the appliance for a long time, be sure to wipe and dry it before storing.

Maintenance

Use only water or neutral cleaning agent to clean the appliance. Never use caustic or abrasive agents (e.g., vinegar or lime scale removers), since they may cause damage to the appliance. Do not use water over 40°C. Keep the appliance away from high temperature and avoid direct sunlight.

• Main unit

1. Wipe the dirt by cloth with soapy water.
2. Do not submerge the main unit in water to clean.
3. Dry after washing with water.

• Water tank cap

Wipe off any grime and dirt with a cloth.

• Nozzle

1. Wash with water and wipe with soft cloth.
2. Replace the nozzle every 6 months.

• Water suction tube

1. Wash with water and wipe with soft cloth.
2. Do not bend, pull or twist the water suction tube.

• Tank


1. Wash with water and wipe with soft cloth.
2. Dry the tank when not using for a long period.

Troubleshooting

Problem	Cause	Solution
Water does not come out	Negative pressure caused by not using for a long time.	Solution 1: Full fill the tank with water, invert the device and shake it to eliminate negative pressure.
		Solution 2: Install the nozzle tip on the main unit, suck out the air from nozzle tip with mouth to eliminate negative pressure.
	You've just purchased the appliance or it has been left unused for a period.	Pull out the nozzle, then flush nozzle hole on the body.
It does not work	You've just purchased the appliance or it has been left unused for more than 3 months.	Charge the appliance for at least 4 hours.
It only works for a few minutes even after charging	The battery has reached the end of its service life.	Contact service center for support.
	Charging time was insufficient.	Charge the unit for at least 4 hours.
Water pressure is very low	There is dirt stuck in the nozzle.	Replace the nozzle.
	The nozzle is deformed.	
	There is dirt stuck in the filter.	Clean the filter.



Disposal

The rechargeable battery inside your appliance is easy to remove for recycling. At the end of its service life and prior to disposal, please remove the rechargeable battery following the instructions below.

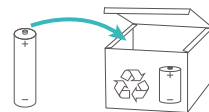Note: During below operation, be sure to protect your eyes, hands, fingers, and the surface on which you work.

1. Remove 2pcs rubber feet on the bottom of the appliance.
2. Use screwdriver counterclockwise to loosen bottom cap.
4. Open top cover and put inner part out of the main unit, then unscrew all screws of the inner part.

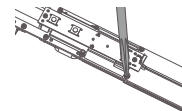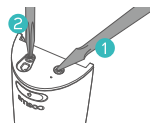3. Slide water tank down off the main unit and unscrew two screws beside nozzle suction tube.
5. Cut off battery wire and pull out batteries from the inner.

Warm Tip: If you have any question about recycling, please contact your local waste management office.

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