CS10G

User Manual

EN

Smart Body Fat Scale

1- 6

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FR

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
Precautions

Read carefully before use:


-  Do not use the body fat scale if you have a pacemaker or any other internal medical device in your body or are pregnant. When you take a measurement, the scale passes a very low, harmless electrical current through your body. The current may cause your internal medical devices to malfunction.
-  Do not use the scale on a wet surface. This could lead to slipping.
-  Dispose of the scale properly following your local regulations.
-  Do not jump on the scale. This scale has a weight limit of 396 lbs.

Scale Overview

GE CS10G Smart Body Fat Scale measures weight, BMI, body fat %, and other body compositions using bioelectrical impedance analysis (BIA) technology.

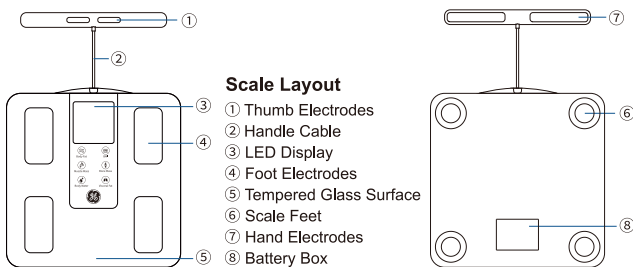LED Display

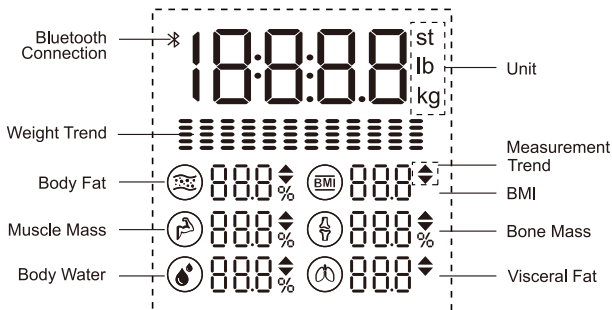Scale Overview

Model:	CS10G
Dimensions:	12.2" x 13.2" x 1.3" / 310 x 336 x 33.8 mm
Weight:	4.5 lb / 2.05 kg
Weighing Range:	0.4 - 396 lb / 0.2 - 180 kg
Power:	3*1.5V AAA battery

Get Started

It needs you to create a profile in the app and open it when weighing at first use to get more data displayed on the screen.

1. Download the app and create an account

- (1) Download the app. Go to the App Store for iOS or Google Play Store for Android, and search for "Fit Profile", or scan the QR code below.
- (2) Sign up and complete your profile.

Scan QR code to download the app

Sign up

2. Setup your scale

If you use an Android phone, please open "Location". And if iPhone, please agree with the app "Fit Profile" to use your phone's Bluetooth first.

- (1) Place the scale on a hard flat floor surface.
- (2) Open the "Fit Profile" app and turn on the Bluetooth.
- (3) Turn on the scale by briefly stepping on and off it.



Your scale will be paired with your smartphone automatically within seconds. If not,

- (1) Tap "+ Pair the Device"; Or click the icon at the top-right-hand corner and tap "+ Add device"
- (2) Turn on the scale by briefly stepping on and off it. On the app, tap the smart scale icon, then select "CS10G" and click "Connect."

Your scale has now been successfully paired with your smartphone and is ready to start taking your measurements.

3. Take measurements


- (1) Open the "Fit Profile" app, and then keep the app open on the "Home" screen.
- (2) Place the scale on a hard flat floor surface, do not place it on carpet or any other soft surfaces.

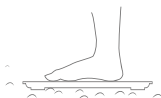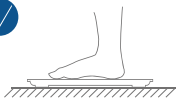
- (3) Turn on the scale by briefly stepping on and off it. Then, wait for the scale to set to 0.
- (4) Please stand on the scale barefoot and ensure your hands and feet are dry.

Guidance on Correct Measuring Posture


a. Step on the correct corresponding electrode area with both feet.

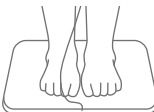b. Grip the handle: Position your thumbs on both oval electrode pads and your other fingers on both slender electrode pads.

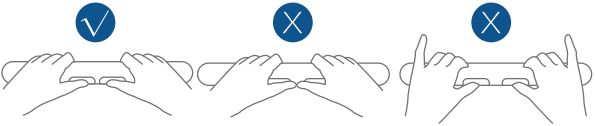c. Stand upright with arms naturally hanging down at an angle from the body, ensuring arms do not touch the waist.

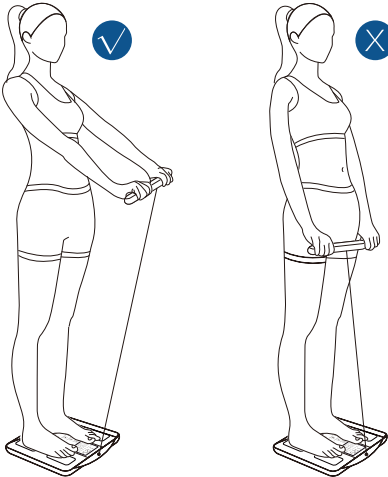(5) Wait for the weight to stabilize; blinking weight does not mean the measurement is complete, you need to wait for other body composition to be measured.

Multiple User Profiles

Separate Fit Profile Accounts:

Users can create or log in to their individual Fit Profile accounts and pair the smart scale with their personal accounts. Weight data for each user remains private and can be synced with third-party apps.

Shared Fit Profile Account:

Users with a Fit Profile account can create multiple user profiles within the same account. All users can log in to the shared account on their own smartphones without needing others to log out.

To add more users and share your account, click the icon at the top left corner on the Home screen in the app.

Note: Users sharing a Fit Profile account can access each other's data recorded under this account. Only the main user's data can be synced with third-party apps.

Switching Units

By default, the scale uses pounds (lb) as the unit of measurement, but it also supports kilograms (kg).

To switch the unit of measurement on the scale, first connect it to the Fit Profile app, then go to the app "Me" page -> "Settings" -> "Unit".

Recommendations

To have accurate, consistent, and comparable measurement results, please follow the recommendations below.

1. DO' s

- Ensure the personal information in your profile is correct.
- Place the scale on a hard flat floor surface when you weigh.
- Make sure the scale displays a 0 before stepping on to measure. While barefoot, stand still on the scale until the measurement is completed.
- Weigh yourself under similar conditions and around the same time every day.
- Empty your bladder before weighing.

2. DON' Ts

- Do not place the scale on a carpet or a soft surface when you weigh.
- Do not measure with socks or shoes on.
- Do not measure when your body is wet. For example, immediately after a shower, bath, sauna, or strenuous exercise.
- Do not measure after drinking a lot of water or alcohol.
- Do not measure if you have a cold or fever.

3. This smart scale is intended for home use only.

- For people ages 18 years and above, it provides full body composition measurements.
- For children under 10, it only provides weight and BMI measurements.
- For children between the ages of 10 and 17, it provides full body composition measurements. (Please note that the body mass-related items like body fat%, etc., are only for your reference.)

4. This smart scale is not intended for use by pregnant women.

Troubleshooting

Symptom	Cause	Solution
CAL	Calibrating	/
Lo	Low power	Please change and use new 1.5V AAA battery.
EEr	Overload	Use the scale within 396lbs.
EEr!	Unstable measurement	Place the scale on the flat, hard surface and measure again.
----	Weight detected after stepping off	Place the scale on the flat, hard surface and let the scale calibrate.
----	/	Please return the handle to the scale when this icon is scrolling.
No Display	Auto-off, or no power	Auto-off: Step on then off to turn on the scale; No power: Please change and use new 1.5V AAA battery.
Display Weight Only	Unable to determine body composition	No user profile: Create app profile and weigh once with Bluetooth. If a profile exists, open the app and weigh once. All others: Need to stand on the scale barefoot/ Not for children under 10.
Unable to Sync with App	Lost connection, or scale turns off	Reconnect the scale: For iPhone, authorize Bluetooth; for Android, enable location. Click the icon on the top-right corner and tap "My Device" > "+ Add device" on the app.

Warranty

This product is warranted by the manufacturer against defects in materials and workmanship for one (1) year from the original date of purchase. Proof of purchase is required.

This warranty does not cover accidental damage, mistreatment, or unauthorized repair. Detailed warranty policy: www.fit-profile.com/warranty

If you have any issues or questions, please contact Global eFit Inc.

Contact : Global eFit inc.

After sales service - US: cs@fit-profile.com

1-866-786-5988 (tollfree US & Canada; M-F 9:00am-5:00pm Central Time)

After sales service - Europe: cs-eu@fit-profile.com

FCC Compliance Statements

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, under Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Caution: Changes or modifications not expressly approved by the party responsible for compliance could void the user authority to operate the equipment.

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

Global eFit inc.
17304 Preston Road, Suite 800, Dallas, TX 75252, USA
www.fit-profile.com
After sales service - US: cs@fit-profile.com
1-866-786-5988
(tollfree US & Canada; M-F 9:00am-5:00pm Central Time)
After sales service - Europe: cs-eu@fit-profile.com

FCC ID: 2ANDX-CSX00

GE is a registered trademark of General Electric Company.
Manufactured under trademark license.

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Kundenservice: Global eFit Inc.