
T40

Oral Irrigator

Munddusche

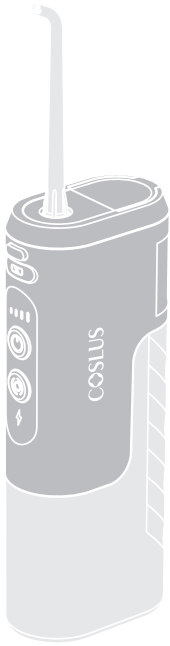USER MANUAL

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WhatsApp contact

Coslus Group

www.coslus.com

Tel: +1(855) 594-3950
(8am-5pm Mon. - Fri. (EST))

E-mail: support@coslus.com

Made in China

EN

1. Nozzle
2. Nozzle release button
3. Charging port
4. Mode indicator
5. On/off button
6. Mode button
7. Power indicator
8. Water tank
9. Water inlet & nozzle storage

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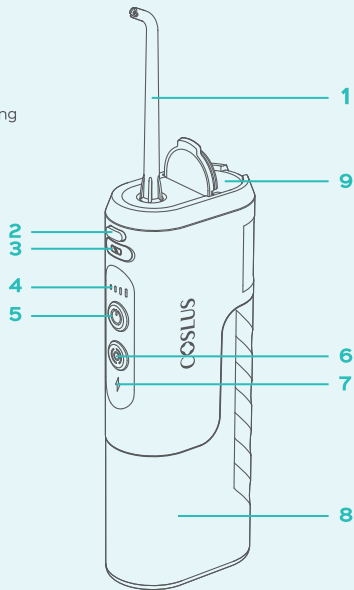
02 Packing List

EN

- ① Water flosser x 1
- ② Charging cable x 1
- ③ Storage bag x 1
- ④ Instruction manual certificate x 1
- ⑤ Nozzles x 4
(standard nozzle x 2, orthodontic nozzle x 1, tongue scraper x 1)

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03 Tips EN

Gum bleeding may occur at the beginning of use, so it is recommended to use the soft mode first to adapt. After getting used to it, choose the appropriate mode according to personal needs. Short press the mode button to switch gears. Keep the nozzle close to the gaps between teeth to remove food debris. Use with your mouth closed to prevent splashing and reduce gum discomfort.

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04 Safety Instructions EN

Danger

1. Do not use the product without understanding the operation procedures and safety instructions to avoid injury.
2. Do not disassemble or modify the product, as it may cause damage and void the warranty.
3. Do not use water above 40°C in the tank to avoid burns.
4. If the power cord or plug is damaged, replace it before use.
5. Do not throw the product into fire, heat it, or charge/use/store it in high temperatures.

Warning

1. Do use this product for purposes other than oral cleaning.
2. Do not press the nozzle forcefully against gums or teeth to avoid injury.
3. Do not immerse the product in water for long periods to prevent malfunction, fire, or electric shock.
4. Do not damage, modify, or misuse the power cord. Avoid pulling, bending, or twisting it, and do not place heavy objects it or trap it between objects. Keep the product and power cord away from heat sources.
5. Do not store in high-temperature or humid places, and avoid direct sunlight.
6. Do not charge the in high-temperature or humid places, such as bathrooms.
7. If not in use for a long time, unplug the power and charging cords to prevent fire or hazards.

Note

1. This product has no replaceable or repairable parts. For any repairs, please return it to the manufacturer.
2. Do not use without a nozzle.
3. This product is not suitable for children or infants unless supervised by a responsible adult. Keep out of their reach and do not use as a toy.
4. Consult a doctor before use if you have periodontal disease or have had oral surgery in the 2 months.
5. Do not add antiseptic mouthwash, as it can corrode the tank and cause it to stick to the body.
6. Before charging, ensure the power adapter's port is Type-C and that the charging port is free of debris and moisture.
7. This product provides deep oral cleaning and helps prevent tartar but does not replace regular tooth brushing.

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Name	Portable oral irrigator	Model	T40
Nozzle diameter	0.6mm	Pause frequency	1450±15%times/min
Battery capacity	1100mAh	Water tank capacity	180ml
Input	5V 2A (Max)	Noise	≤72dB
Charging port	Type-C	Rated power	8W
Waterproof level	IPX7	Charging time	About 4h
Size (Height*Width*Length)	164.7*38*64.9mm	Gear	4 gears
		Net weight	About 207g

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Problem	Cause	Solution
Not working	Newly purchased or not used for over 3 months	Please charge for at least 4 hours
Low water pressure	Something is blocking the nozzle, causing it to deform	Replace the nozzle
Only lasts a short time after a full charge	The battery has reached the end of its service life	Please contact customer service
	Insufficient charging time	Please charge for more than 4 hours before use
The charging indicator does not light up when charging	The power plug is not properly connected	Please connect the power properly
	There is a fault in the product's charging circuit	Please contact customer service

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Mode	Applicable People
Soft	For sensitive teeth or first-time use
Clean	Daily cleaning of food debris
Deep clean	Remove interdental debris
Pulse	Alternating strength for deep cleaning and gum massage
Battery indicator	When low, the indicator turns red. When charging, the red indicator blinks. When fully charged, the indicator turns green.

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EN Tongue scraper: clean tongue surface, remove bad breath, food residue and thick tongue fur. When cleaning tongue, please use the soft mode.

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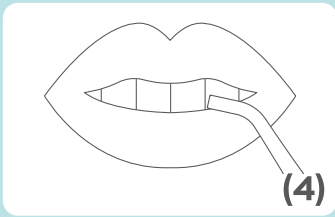
EN 4. Place the nozzle in your mouth, press the on/off button, and start flossing.

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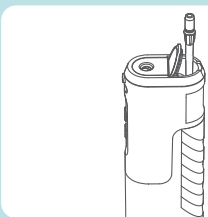
EN 5. After flossing for 1-2 min, remove the nozzle and store it by opening the fill cap.

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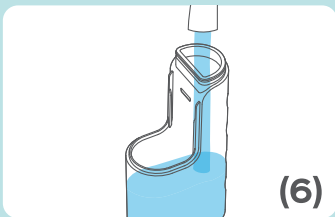
EN 6. For cleaning, detach the water tank and clean the interior.

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① **Nozzle:**
Press and hold the nozzle release button to remove the nozzle. Rinse with water and wipe dry with a cloth.

● Dentists recommend replacing the nozzle every 6 months.

● If the nozzle tip is deformed, replace it to avoid malfunction.

② **Main Unit:**

Wipe off dirt with a cloth dipped in liquid soap, then dry with a clean cloth. If dirty, you can clean with water.

● Do not use thinner, gasoline, alcohol, etc., to clean the unit.

● Do not use water hotter than 40°C.

● After cleaning with water, dry thoroughly.

③ **Water Inlet:**

Wipe off dirt to prevent leakage.

④ **Water Tank:**

Regularly remove the tank cover and clean the water tank.

● Use a brush scrub the inside, then dry with a clean cloth.

● If not used for a long time, let the tank air dry.

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