
ASSEMBLY / MONTAJE / MONTAGE /
ASSEMBLAGE / ASSEMBLAGGIO / 組立

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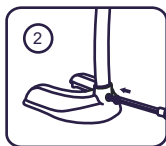K40001

K40002

The **KMINA SHOEHORN/SHOE REMOVER** is designed to facilitate the putting on and removal of footwear, especially for people with reduced mobility or difficulty bending down.

INSTRUCTIONS FOR USE

1. Ensure that the shoehorn/shoe remover is clean and free of visible damage.
2. Place the shoe on the ground, open and ready to be used.
3. Hold the handle of the shoehorn firmly. Insert the tip of the shoehorn into the shoe so that it rests on the heel of the shoe.
4. Carefully slide your foot into the shoe, using the shoehorn as a guide to avoid friction or excessive effort.
5. Remove the shoehorn once your foot is completely inside the shoe.

FEATURES

- Material: Plastic and aluminum.
- Dimensions: 4.5 x 30 x 6" (12,4 x 76 x 14,8 cm).
- Weight: 170 g (approximately 0.37 lbs).

CLEANING AND MAINTENANCE

Clean the shoehorn with a damp cloth and mild soap if necessary. Avoid using harsh chemical products that could damage the material.

WARNINGS

- Do not use the shoehorn for purposes other than those indicated.
- Regularly inspect the product to ensure it is free of damage or excessive wear.
- Use the product on a flat and stable surface to avoid falls.

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