Ergonomic Wedge Cushion

Ergonomisches Keilkissen

DE EN ES FR IT

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EN

Intended use:

The ergonomic wedge cushion from Healthfix is used to correct and optimize sitting posture. It supports the natural curvature of the back, promotes an upright posture and can help reduce back pain and discomfort. The adjustable U-shaped neckline provides pressure relief on the tailbone. Ideal for office workers, people who spend long periods of time sitting or anyone who wants to improve their seating ergonomics. Suitable for various seating furniture such as office chairs, dining chairs or in the car.

Indications:


Poor posture Back pain Coccyx problems Sitting for long periods Pelvic and hip problems Prevention of postural problems Rehabilitation after operations or injuries

Contraindications:

Lymphatic drainage disorders, sensory and critical circulatory disorders in the supplied body regions, PAD unstable fractures in the spine and pelvic area. Medical conditions for which the treating doctor or therapist advises against it.

Instructions for use:

- Place the wedge cushion on a flat seat with the matt anti-slip nubs facing down and the recess for the tailbone facing backward.
 - Make sure the wider side of the wedge cushion is at the back and the narrower side is at the front.
 - Placement: Place the cushion on a chair or seat of your choice. It should fit tightly and securely without slipping.
- After just 30 seconds of using the **Healthfix® wedge cushion**, the MEMORY-FOAM will adapt to your body shape with the help of your body heat. We still recommend that you get up regularly.

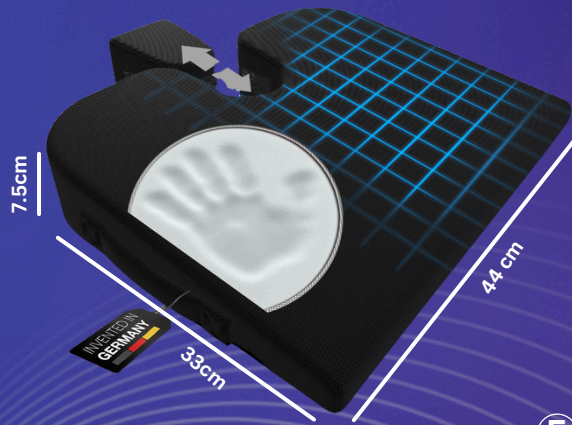Storage and transport conditions:

Constant conditions, ideally temperatures between 18 to 25 ° degrees and a humidity of 35 to 65%.

Cleaning information:

Open the zipper at the bottom and remove the foam core. The foam core should not be washed as it could lose its effectiveness.
Wash the cover in the washing machine at 40°.

Scope of delivery:

1x Healthfix® wedge cushion
1x wedge cushion cover
Instructions and brochure

Warranty and customer service

The Healthfix® wedge cushion comes with a 12-month warranty.
If you have any questions or problems, please contact our customer service at the email address provided.

Intended useful life

- 5 years

Safety instructions:

Keep the wedge cushion away from fire-inducing heat, sparks or open flames.

Only use the product as described in these instructions.

Make sure that the cushion sits firmly and does not slip.

The cushion is not a toy. Keep the cushion out of the reach of children.

Instructions for safe disposal:

Disposal as plastic waste

Further information about the application:

If you have any difficulties in use, please contact your dealer or the responsible Healthfix® customer advisor.

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All rights reserved.
Healthfix GmbH | Bentheimer Str. 66 | 48529 Nordhorn | Germany