THE DeCHOKER®
Because every second counts.™

INSTRUCTIONS FOR USE

©2022 by Dechoker LLC All Rights Reserved.

Dechoker LLC
4880 Van Gordon St. #200
Wheat Ridge, CO 80033

EC REP
MDSS GmbH
Schiffgraben 41
30175 Hannover, Germany

About The IFU

This IFU provides general instructions for the Dechoker operations. **Dechoker recommends first action Red Cross/AHA protocol. If this fails continue CPR, and use the Dechoker.** It is important that upon receipt of your device, you examine for defects. Further, it is important to read the IFU and familiarize yourself with the device and procedures on how to use the Dechoker. Observing the warnings and cautions listed within this manual will provide effective use of the device.

Requirements

The Dechoker is used within the home healthcare environment, restaurants, hospitals, and within commercial facilities. Operators of the Dechoker shall be familiar with this manual.

Indications for Use (Purposes of the Dechoker)

The Dechoker is a medical device that creates suction to clear the acute upper airway from obstruction in the relief of choking. The Dechoker device can be applied to a choking victim to generate a suction to remove the acute obstruction, and is intended for use as an aid to family members, healthcare professionals, caregivers, EMTs, restaurant workers or any average user.

Warning: Safety and effectiveness of the Dechoker has not been established for uses outside of those specified.

Contraindications

The Dechoker is contraindicated for use by patients/users with the following:

- The Dechoker should not be used without contacting 911 or emergency assistance.
- The Dechoker should not be used on infants (less than 1-years-old)

Recommended Sizes

- 01DCH01 - Toddler (1 yr to 3 yrs)
- 01DCH03 - Child (3 yrs to 12 yrs)
- 01DCH02 - Adult (12 yrs and up)

Important Symbols

The table below shows all symbols and the corresponding meaning of the symbols used within this manual and on the Dechoker labeling. Operators must pay careful attention to these symbols as they provide important safety precautions and warnings.

Label Symbol	Meaning
	Attention, Caution, Warning, Danger, Important, Note, or Refer to Accompanying Documentation
	Refer to Instruction Manual / Booklet
	Do not reuse
MD	Medical device
LOT	Batch code
	Use by date
	Manufacturer
EC REP	Authorized Representative in the European Community
	Keep dry
	Do not use if package is damaged

Table 1: Important Symbols and Descriptions

Dechoker Warranty/Returns

Dechoker LLC offers a 30 day return policy. Simply email customerservice@dechoker.com and provide your order number.

If the Dechoker is used in a life-saving event, contact life-saves@dechoker.com for a replacement device.

General Warnings

It is important to read all warnings before using the Dechoker. Reference the table below for specific warnings.

Warning	Statement
	All users must read and understand the Dechoker IFU, including indications, contraindications, operating instructions, warnings and precautions before using the device. Failure to do so may result in injury.
	Equipment is not suitable or use in the presence of flammable anesthetic mixture with air or oxygen or nitrous oxide.
	The Dechoker plastics are latex free. Skin irritation may result with the constant handling of the device for individuals who have sensitivity or allergic reaction to ABS materials.

Table 2: Warning Descriptions

Technical Support

For questions regarding the Dechoker, please contact Dechoker LLC Technical Support Team at 1-888-998-8995 in North America or the EU Authorized Representative at +44-17-04-54-49-44. Further, if a serious incident has occurred in relation to the device, please report this to the manufacturer and the competent authority of the Member State in which the user/patient is established. Please have the Dechoker Lot number (located on the Dechoker package label) available for reference when contacting Dechoker LLC.

Cleaning Specifications

The Dechoker is a patient individual signal use device. Once used on the individual to remove the acute obstruction, do not clean the device. Dispose of the device with the regulations of the country in which it was sold.

Whom to Call for Help

Always call 911 or the emergency call number of the country the Dechoker device is to be used prior to using the device on an acute airway obstructed individual during this emergency. In the event that the Dechoker is inoperable revert to other acute airway obstruction method of relief.

Caution: Do not attempt to disassemble the device. The manufacturer or an authorized representative must perform device failure evaluation.

Dechoker Setup

1. Remove the Dechoker package from the box.
2. Open the Dechoker bag and remove the Dechoker device.
3. Examine device for defects.
4. Read IFU and familiarize yourself with the Dechoker.
5. Dechoker ready to be applied to choking individual.

Dechoker recommends first action Red Cross/ AHA protocol. If this fails continue CPR, and use the Dechoker.


1. Have someone call your emergency provider. Remove Dechoker from the package and assure choking victim that help is on the way.

2. Tilt head back and lift chin for optimal access to airway.

3. Insert tongue depressor into mouth with respirator mask covering mouth and nose.

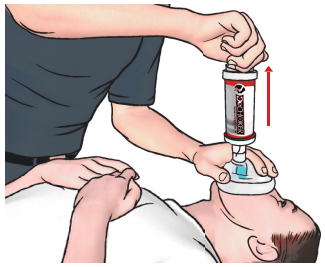4. Hold respirator mask between your thumb and fingers and apply light but firm pressure to form a seal between the face and mask.

5. Continue to apply light pressure on the mask while pulling the handle. Remove mask for ventilation for victim.


6. If victim is still choking, repeat steps 4 and 5. Multiple attempts can be performed to create additional suction if necessary.

7. Once debris has been dislodged, roll individual over on their side to allow additional debris and/or fluid out of the mouth.

Dechoker Components

1. **Respiratory Mask** - The respiratory mask covers the mouth and nose of the individual the device is being used on.
2. **Suction Barrel** - The sealed suction barrel encapsulates the piston to create the appropriate set suction level to remove the acute obstruction.
3. **Pull Handle** - The pull handle provides a simple grip to drawing back the piston.
4. **Exhaust Valve** - The exhaust valve permits the pull handle to move the piston back to its original position to allow another suction attempt in removing the acute obstruction.
5. **Enclosure (Outer Cylinder)** - The enclosure is used as a protective housing separating the user from the internal suction components of the Dechoker. Do not attempt to disassemble the device.
6. **Device Labeling** - The device has a label on the package identifying the device manufacturer, and reference the device lot number when making any inquiries.
7. **Device Quick Reference Guide** - The Dechoker has an operator's quick reference device instructions wrapped around the barrel of the device. A hard copy quick reference guide is also enclosed within the shipping box.
8. **Cross Slit Valve** - Prevents air or obstructions from going back into the patient.
9. **Tongue Depressor** - Prevents mouth cavity from collapse.

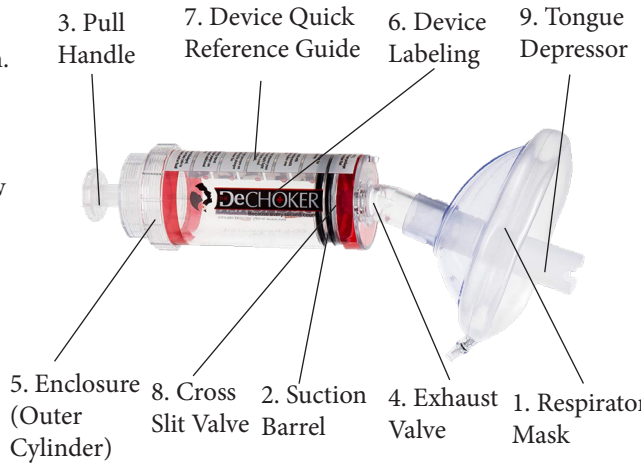THE DeCHOKER®
Because every second counts.™

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