easy@Home

# Ovulation Test

- EN -

REF EZW2-S

## INTENDED USE

Easy@Home ovulation test is a qualitative test that is designed to help predict when a woman will experience a surge in luteinizing hormone (LH). The test detects the surge in LH to identify the most fertile days of a woman's menstrual cycle.
For *in vitro* diagnostic use.
For self-testing use.

## HOW DOES IT WORK?

The LH (Luteinizing Hormone) in the urine of most women will increase dramatically in the middle of the menstrual cycle. The LH increase triggers ovulation, which is when a woman's ovary releases an egg. WHO experts state that LH testing is a reliable way to detect ovulation. The LH surge occurs about 24 to 48 hours before you actually ovulate. A positive LH test result may indicate that you are entering your most fertile period of your cycle.

The ovulation test contains antibodies which specifically react with LH. When the strip is immersed in a urine specimen capillary action carries the specimen to migrate along the membrane. When the luteinizing hormone (LH) in the urine reaches the test zone, a colored line will appear in the test zone.

If the test has been performed properly, a colored line will appear at the control zone. The control line also works as a reference line for you to find the LH peak time when the test line is as dark or darker than the control line.

## CONTENTS OF THE TEST KIT

1. One pouch containing a test strip and a desiccant. The desiccant is for storage purpose only and is not used in the test.
2. Instructions for use.

## WHAT ELSE DO YOU NEED?

1. A clean, dry, plastic or glass container to collect the urine.
2. Timer (watch, clock or smartphone)

## PRECAUTIONS

1. This kit is for external use only. Do not swallow.
2. Do not use the test kit beyond expiration date.
3. Do not use the kit if the test pouch is punctured or not fully sealed.
4. Keep out of the reach of children.
5. Discard after first use. The test strip cannot be used more than once.
6. Please follow local regulations to discard used tests.

## STORAGE AND STABILITY

1. Store at 4°C to 30°C (39°F to 86°F) in the sealed pouch up to the expiration date.
2. The test should be used immediately after opening.
3. Keep away from sunlight, moisture and heat.
4. DO NOT FREEZE


## COLLECTION AND STORAGE OF SPECIMENS

***WHEN TO COLLECT URINE FOR THE TEST?***
Any urine specimen is appropriate for ovulation testing.

***HOW TO COLLECT THE URINE SAMPLE?***
This test is designed to have urine collected in a clean container just prior to testing. Urine may be collected in any clean, dry, plastic or glass container.

## WHEN TO START TESTING?

The length of the menstrual cycle is the duration from your first menstrual bleeding day to the day before the next bleeding begins. Determine the length of your menstrual cycle before testing. Count the first day of bleeding as cycle day 1. Please refer to the chart below to determine when you should start testing. If you do not know your cycle length, start testing from cycle day 5 until the LH surge has been detected.

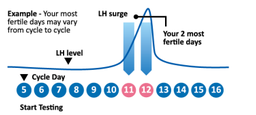
**Example:** If your cycle is normally 26 days, according to the Cycle Chart, testing should begin on day 10 of your cycle. To determine day 10, count 10 days from the first day of your menstrual bleeding (cycle day 1). For example, if your first day of menstrual bleeding was January 3, then day 10 would be January 12. You would start testing for ovulation on January 12 and continue until the LH surge has been detected.

**Sample calendar**


| | | | | | | |
|---|---|---|---|---|---|---|
| 1 | 2 | 3+ day 1 | 4 | 5 | 6 | 7 |
| 8 | 9 | 10 | 11 | 12 Giorno 10 | 13 | 14 |
| 15 | 16 | 17 | 18 | 19 | 20 | 21 |
| 22 | 23 | 24 | 25 | 26 | 27 | 28 |

3+ = The first day of menstrual bleeding
12# = The day to begin testing (Day 10)
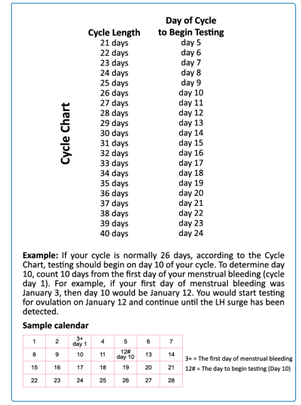
## HOW TO DO THE TEST?

1. Remove the test strip from the sealed pouch.
2. Dip the white absorbent end into the urine sample up to the black 'Max' line until you see the dye begins to rise into the testing window (Est. time: 5-10 seconds).
3. Lay the test strip flat on a clean, dry, and non-absorbent surface.
4. Read the results in 5 minutes after immersion or after any excess dye has completely cleared the testing window. Do not read after 5 minutes.

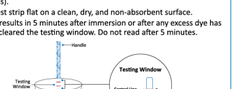
The test strip and urine sample should be at room temperature (59-86°F or 15-30°C) before performing the test to ensure accurate results.

## HOW TO READ THE RESULTS?

**Positive(LH Surge)**
If two color lines are visible, and the test line is as dark or darker than the control line, you will probably ovulate in the next 24-48 hours. If trying to get pregnant, the best time to have intercourse is after 24 but before 48 hours.

**Negative (No LH Surge)**
Only one color line appears in the control zone, or the test line appears but is lighter than the control line. This means there is no LH surge.

**Invalid**
No visible lines at all, or there is a visible line only in the test zone but not in the control zone.
Repeat with a new test strip. If test still fails, please contact us for assistance.
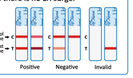
## LIMITATIONS

1. The tests results are more accurate when the test procedures are precisely followed.
2. This test may not be used as a form of birth control.
3. Consult a doctor if you experience irregular or unusually long cycles.
4. This test is a qualitative screening assay. It is not designed to determine the quantitative concentration of LH.
5. The device is only suitable for human urine test.
6. Unclean specimen or interfering substances in the urine may cause false result.
7. Some drugs for infertility treatment, such as 17 β-Estradiol, Progesterone and Estradiol Valerate, may cause false result.
8. Unclean test device may cause false result.

## PERFORMANCE CHARACTERISTICS

**SENSITIVITY AND SPECIFICITY**
The Easy@Home Ovulation Test will detect hLH in human urine at the level higher than 25 mIU/ml. The test has been calibrated by the W.H.O. Standard (2nd WHO International Standard (IS) 80/552).
No cross-reaction appears when there is hFSH (200 mIU/mL) or hTSH (250 mIU/L) in specimen.

**INTERFERING SUBSTANCE**
The following potentially interfering substances were added to LH negative and positive specimens.

| | | | |
|---|---|---|---|
| Acetaminophen | 20mg/dl | Gentisic acid | 20mg/dl |
| Acetylsalicylic acid | 20mg/dl | Glucose | 2mg/dl |
| Ascorbic acid | 20mg/dl | Hemoglobin | 20mg/dl |
| Atropine | 20mg/dl | Tetracycline | 20mg/dl |
| Caffeine | 20mg/dl | Ampicilline | 20mg/dl |

None of the substances at the concentration tested interfered in the assay.

## QUESTIONS AND ANSWERS

**Q Is it a concern if there is a big difference between the control lines of two different ovulation tests?**
No. Variations in the color of the control line will not affect the test result. It is only effective to compare the color of the test line to that of the control line of the same test on the day the test is performed.

**Q How long should I continue to perform the test?**
Continue testing for at least 5 days or until the LH surge has been detected.

**Q A pink background color and vertical streaking appeared in the result area during the testing period. Is this a concern?**
No. Each urine sample will vary in its chemical makeup, as will the humidity of the air in the testing chamber. Such variations in physical conditions can cause the vertical streaking and/or the pink-rose background color, but will not affect the test results. As long as the control line appears within five minutes, the test result is valid.

**Q Does alcohol or common medications affect the test?**
No, but you should consult your physician if you are taking any medication.

**Q Once I see a positive result, when is the best time to have intercourse?**
Ovulation is likely to occur within 24-48 hours. This is your most fertile time. Sexual intercourse is advised within this time frame.

**Q I am now using the basal body temperature method (BBT). Does this test replace BBT?**
The elevated basal body temperature primarily indicates that ovulation has already occurred. This test indicates that ovulation is about to occur within 24 to 48 hours.

**Q I have received a positive result and had intercourse during these fertile days. I have not become pregnant. What shall I do?**
Many factors can affect your ability to become pregnant. Often you may need to use the test kit for 3-4 months. You and your partner should consult your physician if pregnancy is not achieved after 3-4 months.

**Q Will the amount of liquid I drink affect the result?**
Yes. We suggest that you limit your fluid intake for about two hours before you collect your urine. Heavy intake of fluids prior to testing will dilute the hormone in your urine.

**Q Can I use first morning urine for testing?**
Unlike pregnancy testing, the FMU (first morning urine) is NOT recommended to use for ovulation testing. Luteinizing hormone(LH) levels usually start to rise in the early morning and may not be easily detectable in the first morning urine. The ideal time to collect a urine sample is from 10 a.m. to 8 p.m., although testing can safely be performed at any time.

**NOTE:** It is suggested to test twice a day when it is close to the LH surge day so that you won't miss it; the LH surge can be very short for some women. If you choose to test twice a day, you may test once around noon and another time in the early evening. We recommend collecting the urine specimen at about the same time each day.

## BIBLIOGRAPHY

1. Sievert LL, Dubois CA. "Validating signals of ovulation: Do women who think they know, really know?" Department of Anthropology, Machmer Hall, UMass Amherst, Amherst, Massachusetts. Am J Hum Biol. 2005 Apr 22;17(3):310-320.
2. Ruiz Angus J, Carballo Mondragon E, Anta Jaen E, Diaz Calvillo J, Kably Ambe A. "Usefulness of luteinizant hormone levels to determine the oocyte quality" Ginecol Obstet Mex. 2005 Jan;73:28-35

## MEANING OF SYMBOLS ON PACKAGE

## QUESTIONS OR COMMENTS?

Please contact us via email: **support@premom.com**
Easy Healthcare Corporation
360 Shore Dr., Burr Ridge, IL USA 60527
Made in China

www. premom.com

Scan the code to get the manual in more languages

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REF EZW1-S

## INTENDED USE

Easy@Home Pregnancy Test is a rapid pregnancy test, which you can easily carry out yourself. It detects the presence of human chorionic gonadotropin ("hCG"), which appears in urine very early during pregnancy.

For *in vitro* diagnostic use. For self-testing.

## HOW DOES IT WORK?

Human chorionic gonadotropin (hCG) is a hormone produced by the developing placenta shortly after conception and secreted into the urine. The pregnancy test contains antibodies that specifically react with this hormone.

When the strip is immersed in a urine specimen, capillary action carries the specimen to migrate along the membrane. When hCG in the sample reaches the Test Zone of the membrane, it will form a colored line. The absence of this colored line suggests a negative result.
To serve as a procedure control, a colored line will appear at the Control Zone if the test has been performed properly.

## CONTENTS OF THE TEST KIT

1. One pouch containing a test strip and a desiccant. The desiccant is for storage purposes only and is not used in the test.
2. Instructions for use.

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## PRECAUTIONS

1. This kit is for external use only. Do not swallow.
2. Do not use the test kit beyond expiration date.
3. Do not use the kit if the test pouch is punctured or not fully sealed.
4. Keep out of the reach of children.
5. Discard after first use. The test strip cannot be used more than once.
6. Urine specimens may be infectious; Ensure proper handling and disposal of all used reaction devices as per the local statute.

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
***WHEN TO COLLECT URINE FOR THE TEST?***
Any urine specimen is appropriate for pregnancy testing. However, the first morning urine specimen is preferred because of its highest concentration of HCG.

***HOW TO COLLECT THE URINE SAMPLE?***
This test is designed to have urine collected in a clean container just prior to testing.
Urine may be collected in any clean, dry, plastic or glass container.

## HOW TO DO THE TEST?

1. Remove the test strip from the sealed pouch.
2. Dip the white absorbent tip of the test into the urine up to the black "MAX" line and keep the tip in urine until dye rises into the results window. (Est. time: 5-10 seconds)
3. Lay the test strip flat on a clean, dry, and non-absorbent surface.
4. Read the results in 3-5 minutes after immersion. Do not read the result after 5 minutes.

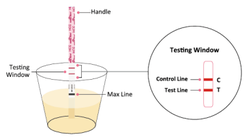
The test strip and urine have to be at room temperature 15°C-30°C (59°F-86°F) for testing.

## HOW TO READ THE RESULTS?

**Positive(Pregnant)**
Two distinct color lines appear in the Control and Test Zones. One line should be in the control zone (C) and another line should be in the test zone (T).

**Negative(Not pregnant)**
Only one color line appears in the control zone (C). No apparent Rosy Pink line in the test zone(T).

**Invalid**
No visible lines at all, or there is a visible line only in the test zone (T) but not in the control zone (C). Repeat with a new test strip. If test still fails, please contact us for assistance.

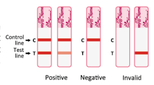***NOTE:***
***The color intensity of the test lines may vary since different stages of pregnancy have different concentrations of hCG hormone. If the color of the test lines is faint, it is recommended to repeat the test in 48 hours.***

## WHAT TO DO WHEN YOU KNOW THE RESULT?

***IF THE RESULT IS POSITIVE (TWO LINES)***
Consult your physician to confirm your pregnancy after evaluating all clinical and laboratory findings. In certain conditions, the test may give a false positive result, and your physician will also give you further information and advise you on your pregnancy.
See also "Limitations".

***IF THE RESULT IS NEGATIVE (ONE LINE)***
If the result is negative and you still suspect you may be pregnant, you should wait a few more days and repeat the test. If you continue to have questions or concerns, you should consult your physician.
See also "Limitations".

## LIMITATIONS

1. As it is with any diagnostic procedure, a confirmed pregnancy diagnosis should only be made by a physician after evaluating all clinical and laboratory findings.
2. If a urine sample is too dilute (ie, low specific gravity) it may not contain a representative level of hCG. If pregnancy is still suspected, another urine specimen should be collected 48 hours later and tested.
3. Low concentration of hCG in a very early pregnancy can give a negative result. In this case, another specimen should be obtained at least 48 hours later and tested.
4. Elevated levels of hCG can be caused by a few conditions other than pregnancy. Therefore, the presence of hCG in urine specimen should not be used to diagnose pregnancy unless these conditions have been ruled out.
5. A normal pregnancy can not be distinguished from an ectopic pregnancy based on hCG levels alone. Also, spontaneous miscarriage may cause confusion in interpreting test results.
6. Sterility treatments, based on hCG, may cause false results. Consult your physician.

## EXPECTED RESULTS

Negative results are expected in healthy non-pregnant women and healthy men. Although the hCG levels in normal early pregnant women are variable, Easy@home Pregnancy Test is capable of detecting pregnancy as early as 1 day after the first missed menses.

## PERFORMANCE CHARACTERISTICS

**Sensitivity:**
Pregnancy Test will detect urine hCG concentration of 10mIU/mL or greater.

**Specificity:**
The specificity of the Pregnancy Test was determined from cross-reactivity studies with known amounts of Luteinizing Hormone (hLH), Follicle Stimulating Hormone (hFSH), and Thyroid Stimulating Hormone (hTSH). Negative results were obtained from all tests conducted with 300 mIU/ml hLH, 300 mIU/ml hFSH and 1000 µIU/ml hTSH.

## QUESTIONS AND ANSWERS

1. **There is much difference between the control lines of two tests. Is this a concern?**
No. Variations in the color of the control line will not affect the test result.


2. **Can test results be interpreted after more than five minutes?**
No. Test results must be read in 5 minutes. Though a positive result should not change for several days, a negative result may change to a false positive within minutes after the end of the testing period, which would not be an accurate reading. It is always best to read the results in 5 minutes and then discard the test to avoid confusion.

3. **A pink background color and vertical streaking appeared in the result area during the testing period. Is this a concern?**
No. Each urine sample will vary in its chemical makeup, as will the humidity of the air in the testing chamber (room). Such variations in physical conditions can cause the vertical streaking and/or the pink-rose background color but will not affect the test results. As long as the control line appears within 5 minutes, the test is working properly.

## BIBLIOGRAPHY

1. Canfield RE, Morgan FJ, Kammerman S, Bell JJ, Agosto GM. Studies of human chorionic gonadotropin. Recent Prog Horm Res. 1971; 27:121–164.
2. Ross GT. Clinical relevance of research on the structure of human chorionic gonadotropin. American Journal of Obstetrics and Gynecology 1977; 129:795.

## MEANING OF SYMBOLS ON PACKAGE
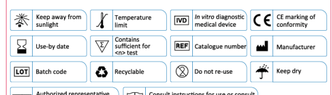
## QUESTIONS OR COMMENTS?

Please contact us via email: **support@premom.com**
Easy Healthcare Corporation
360 Shore Dr., Burr Ridge, IL USA 60527
Made in China
www. premom.com

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